

Acumed Implants in the MR Environment



Table of Contents

| | |
|---|-----------|
| Introduction..... | 5 |
| Terminology | 6 |
| Warnings & Precautions | 7 |
| IFU and Surgical Technique | 7 |
| Anatomic Radial Head Portfolio..... | 8 |
| Anatomic Radial Head System..... | 8 |
| Anatomic Radial Head Solutions..... | 10 |
| Anatomic Radial Head Solutions 2 (stems and heads only)..... | 11 |
| Intramedullary Nails Portfolio..... | 15 |
| Polarus® Humeral Rod System..... | 15 |
| Polarus® Plus Humeral Rod System | 16 |
| Polarus® 3 Solution (Locking and Proximal Locking Nails)..... | 18 |
| Fibula Rod System | 20 |
| Forearm Fracture Solutions (Radius & Ulna Rods)..... | 20 |
| Fibula Nail 2 System | 21 |
| Ulna Nail 2 System | 24 |
| Hand and Wrist Portfolio..... | 28 |
| Acu-Loc® Wrist Plating System | 28 |
| Acu-Loc® 2 Wrist Plating System..... | 31 |
| Acu-Loc® Wrist Spanning Plate | 36 |
| Forearm Fracture Solutions (plates) | 40 |
| Hand Fracture System..... | 41 |
| Modular Hand System..... | 44 |
| Osteotomy System..... | 45 |
| Total Wrist Fusion Plating System | 46 |
| Screw and Pin Portfolio..... | 49 |
| Acutrak® Headless Compression Screw System | 50 |
| Acutrak 2® Headless Compression Screw System | 52 |
| Acutrak® 3 Headless Compression Screw System | 54 |
| AcuTwist® Acutrak Compression Screw System..... | 58 |
| Biotrak® Headless Resorbable Compression Screw System..... | 58 |
| Cannulated Screw System | 59 |
| Extremity Screw System..... | 66 |
| Tension Band Pin System / Tension Band Pin System 2..... | 67 |

Table of Contents [continued]

| | |
|--|------------|
| Screw and Pin Portfolio - ExsoMed™ | 68 |
| ArcPhix™ Functional Flexion Compression Screw System | 69 |
| InFrame™ Intramedullary Micro Nail System | 69 |
| INnate™ Intramedullary Threaded Nail System | 69 |
| NanoPhix™ Cannulated Lag Screw | 70 |
| Shoulder and Elbow Portfolio | 72 |
| Elbow Plating System | 73 |
| Congruent Elbow Plating System | 81 |
| Radial Head Plating System | 85 |
| Clavicle Plating System | 86 |
| Clavicle Hook Plating System | 93 |
| Locking Clavicle Plating System | 95 |
| Acu-Sinch® Repair System | 98 |
| Dual-Trak Clavicle Screw System | 99 |
| Polarus® 3 Solution (plates) | 99 |
| Polarus® Proximal Humeral Plating System | 102 |
| Scapula Plating System | 104 |
| Foot and Ankle Portfolio | 108 |
| Ankle Plating System | 108 |
| Ankle Plating System 3 | 115 |
| Calcaneal Plating System | 119 |
| Forefoot-Midfoot Plating System | 126 |
| Small Fragment Base Set | 133 |
| Hip and Pelvic Portfolio: | 141 |
| Pelvic Plating System | 141 |
| External Fixation Systems | 145 |
| Stableloc External Fixation System | 145 |
| Small Bone Distractor | 145 |
| Small Bone Fixator | 145 |
| Scapholunate Repair System | 146 |
| Scapholunate Repair System | 146 |
| Tendon and Ligament Portfolio | 147 |
| Acu-Sinch Knotless System | 147 |
| Acu-Sinch Knotless Mini System | 147 |
| PEEK Suture Anchor System | 148 |

Table of Contents [continued]

| | |
|---|------------|
| Chest Wall Portfolio | 149 |
| RibLoc U Plus Chest Wall Plating System | 150 |
| Craniomaxillofacial (CMF) Portfolio | 151 |
| 2.0 Orthognathic | 152 |
| OSA™ | 154 |
| OsteoMed CFx™ Craniomaxillofacial Fixation System | 156 |
| OsteoMed MFx™ Mandibular Fixation System | 157 |
| OsteoMed ICON™ Facial Plating System | 174 |
| 2.0 mm Auto Drive MMF Screws | 175 |
| Neurosurgical Portfolio | 176 |
| OsteoMed Fast Flap™ Neurosurgical Fixation System | 177 |
| OsteoMed Profile Zero™ Low Profile Neuro Plating System | 180 |
| OsteoMed® Upper Extremities Portfolio | 182 |
| OsteoMed® Hand Plating System (hps™) | 184 |
| OsteoMed® hps™ Hand Fusion System | 195 |
| OsteoMed® M3-X™ Extremity Fixation System | 196 |
| OsteoMed® InstaFix™ Shape Memory Fixation System | 200 |
| OsteoMed® ExtremiLock™ Wrist Plating System | 201 |
| OsteoMed® Lower Extremities Portfolio | 207 |
| OsteoMed® Talar-Fit™ Subtalar Arthroereisis Implant | 209 |
| OsteoMed® Hemi™ Great Toe Implant | 210 |
| OsteoMed® FPS™ Foot Plating System | 211 |
| OsteoMed® ExtremiFuse™ Hammertoe Fixation System | 227 |
| OsteoMed® ExtremiLock™ Ankle Plating System | 228 |
| OsteoMed® EnCompass™ Metatarsal Resurfacing Implant | 237 |
| OsteoMed® ReFlexion™ 1st MTP Implant System | 238 |

Introduction

Magnetic Resonance Imaging (MRI) is one of many imaging modalities used by radiologists to diagnose and treat various conditions within the body. MRI is often favored for its detail and clarity when imaging soft tissues, and doing so without delivering the ionizing radiation inherent with X-ray and CT. Because of the strong, static magnetic field, and varying magnetic and electric fields during scanning, the safety of patients in the MR environment is of great importance to the healthcare professionals who perform these procedures.

Safety concerns arise when a patient has orthopaedic implants manufactured from metal or carbon fiber, or when external fixation devices are in place, because of the potential for movement or heating of the devices in the MR environment. In some cases, safety concerns may be greater due to the type, location, or configuration of the implant or the presence of multiple implants in close proximity to one another. Consequently, the use of MRI technology may not always be possible.

Purpose

This document will communicate MRI safety information regarding Acumed implants in the MR environment. While we cannot make specific MR claims for all of our implants at this time, Acumed will present the results of our testing in this document as additional information becomes available. When allowed by regulatory authorities, Acumed will use one of the internationally recognized ASTM symbols in our implant labeling.

Implant Materials

Acumed implants are commonly manufactured from commercial pure (CP) titanium (Ti) and titanium alloys, such as Ti-6Al-4V ELI. A recent FDA article¹ states:

“...[Test] results indicate that for alloys Ti-6Al-7Nb, CP titanium Grades 1 through 4, Ti-6Al-4V ELI, Ti-15Mo, Biodur 108, and Co-28Cr-6Mo ..., magnetically induced force is less than the device weight in any current clinical 1.5, 3.0, and 7.0-T MR system.”

Accordingly, many Acumed devices do not pose significant risk for magnetically induced force or torque. Refer to the device label and Instructions for Use to determine the implant material.

Terminology

ASTM International is an international standards organization that has published the following voluntary consensus standards for the testing and labeling of devices that may be used in MRI:

- ▶ F2052 Test Method for Measurement of Magnetically Induced Displacement Force on Medical Devices in the Magnetic Resonance Environment
- ▶ F2182 Test Method for Measurement of Radio Frequency Induced Heating On or Near Passive Implants During Magnetic Resonance Imaging
- ▶ F2213 Test Method for Measurement of Magnetically Induced Torque on Medical Devices in the Magnetic Resonance Environment
- ▶ F2503 Practice for Marking Medical Devices and Other Items for Safety in the Magnetic Resonance Environment

ASTM F2503 defines the three classifications for products used in MRI as MR Safe, MR Conditional, and MR Unsafe. Products in these classifications may be labeled with the following symbols:



MR Safe

MR Safe is an item that poses no known hazards resulting from exposure to any MR environment. MR Safe items are composed of materials that are electrically nonconductive, nonmetallic, and nonmagnetic. The MR Safe icon consists of the letters “MR” surrounded by a green square.



MR Conditional

MR Conditional is an item with demonstrated safety in the MR environment within defined conditions. At minimum, the conditions of the static magnetic field, the switched gradient magnetic field, and the radiofrequency fields must be disclosed. Additional conditions may be required. The MR Conditional icon consists of the letters “MR” within a yellow equilateral triangle with a thick black band around the perimeter.



MR Unsafe

MR Unsafe is an item that poses unacceptable risks to the patient, medical staff, or other persons within the MR environment. The MR Unsafe marking consists of the letters “MR” surrounded by a red circle with a diagonal red bar across the letters extending from the upper left quadrant to the lower right quadrant of the circle and oriented at 45° from the horizontal.

The symbols may appear in black and white when not practical to print in the standard colors.

Warnings & Precautions

Unless explicitly stated in this document, Acumed devices have not been tested for heating, migration, or image artifact in the magnetic resonance (MR) environment and the safety of these devices in the MR environment is unknown. Scanning a patient with untested implants may result in patient injury.

Patients with implants should seek a medical opinion before entering the MR environment.

For Devices That Have Been Tested

Warning

- ▶ Failure to follow the specific MR conditions may result in injury to the patient.
- ▶ Unconscious patients or those insensitive to heat or pain may not inform the technologist that the implant is being affected by the MR environment.
- ▶ Long scanning time or higher specific absorption rate (SAR) increases the risk of heating. Pause the scan periodically and solicit feedback from the patient.
- ▶ Unknown and untested implant configurations or the addition of implants from other manufacturers may increase the risks of heating and image artifact.

Caution

- ▶ Local transmit coils are not recommended for use on the area with implants as their effects have not been tested.
- ▶ Maximize the distance from the implant to the image area to minimize any potential artifact, when possible.
- ▶ Maximize the distance between the implant and sides of the bore, coils, and conductors to minimize the risk of heating, when possible.
- ▶ The implants have not been assessed for compatibility with MRI systems above 3 Tesla.

Important

- ▶ An increase in ventilation through the bore may help to lessen the increase in body temperature during scanning and reduce the effects of implant heating.
- ▶ Nonclinical testing does not account for the variations in MRI systems, pulse sequences, multi-channel RF coils, and local transmit coils that might be used. Some MRI systems have image artifact reduction sequences if imaging near a device is necessary.
- ▶ The use of low field MRI systems (less than 1 T) can mitigate risks for patients with implants by lowering the frequency and SAR.

IFU and Surgical Technique

Acumed offers Instructions for Use (IFU) and Surgical Techniques to promote the safe and effective use of our systems. Visit www.acumed.net for the most current versions.

Important: IFUs and Surgical Techniques may contain important safety information. They are intended for the Operating Surgeon and supporting Healthcare Professionals.

Anatomic Radial Head Portfolio

- ▶ Anatomic Radial Head System
- ▶ Anatomic Radial Head Solutions
- ▶ Anatomic Radial Head Solutions 2



Conditions for Safe Scanning in the MRI

A patient with the head or stem implants from the Anatomic Radial Head, Anatomic Radial Head Solutions, and Anatomic Radial Head Solutions 2 systems may be safely scanned at 1.5 T or 3.0 T under the following conditions. Failure to follow these conditions may result in injury to the patient.

Anatomic Radial Head Portfolio

| Parameter | Condition |
|---|---|
| Nominal values(s) of Static Magnetic Field | 1.5 T or 3 T |
| Maximum Spatial Field Gradient | 19 T/m (1900 gauss/cm) |
| RF Excitation | Circularly Polarized (CP) |
| RF Transmit Coil Type | Whole body transmit coil, Head RF transmit-receive coil |
| Maximum Whole Body SAR at 1.5 T and 3.0 T | 2.0 W/kg |
| Limits on Scan Duration 1.5 T | 2.0 W/kg whole body average SAR for 15 minutes of continuous RF (a sequence or back to back series/scan without breaks) |
| Limits on Scan Duration 3 T | 2.0 W/kg whole body average SAR for 12 minutes of continuous RF (a sequence or back to back series/scan without breaks) |
| MR Image Artifact | The presence of this implant may produce an image artifact of 71 mm. |
| If more information about a specific parameter is not included, there are no conditions associated with that parameter. | |

Anatomic Radial Head System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|---------------------|
| TR-H200L-S | N/A | 20.0 mm Head, Left |
| TR-H200R-S | N/A | 20.0 mm Head, Right |
| TR-H220L-S | N/A | 22.0 mm Head, Left |
| TR-H220R-S | N/A | 22.0 mm Head, Right |
| TR-H240L-S | N/A | 24.0 mm Head, Left |
| TR-H240R-S | N/A | 24.0 mm Head, Right |
| TR-H260L-S | N/A | 26.0 mm Head, Left |
| TR-H260R-S | N/A | 26.0 mm Head, Right |

Anatomic Radial Head Portfolio [continued]

Anatomic Radial Head System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|-----------------------|
| TR-H280L-S | N/A | 28.0 mm Head, Left |
| TR-H280R-S | N/A | 28.0 mm Head, Right |
| TR-S0600-S | N/A | 6.0 mm x 0.0 mm Stem |
| TR-S0602-S | N/A | 6.0 mm x 2.0 mm Stem |
| TR-S0604-S | N/A | 6.0 mm x 4.0 mm Stem |
| TR-S0606-S | N/A | 6.0 mm x 6.0 mm Stem |
| TR-S0608-S | N/A | 6.0 mm x 8.0 mm Stem |
| TR-S0700-S | N/A | 7.0 mm x 0.0 mm Stem |
| TR-S0702-S | N/A | 7.0 mm x 2.0 mm Stem |
| TR-S0704-S | N/A | 7.0 mm x 4.0 mm Stem |
| TR-S0706-S | N/A | 7.0 mm x 6.0 mm Stem |
| TR-S0708-S | N/A | 7.0 mm x 8.0 mm Stem |
| TR-S0800-S | N/A | 8.0 mm x 0.0 mm Stem |
| TR-S0802-S | N/A | 8.0 mm x 2.0 mm Stem |
| TR-S0804-S | N/A | 8.0 mm x 4.0 mm Stem |
| TR-S0806-S | N/A | 8.0 mm x 6.0 mm Stem |
| TR-S0808-S | N/A | 8.0 mm x 8.0 mm Stem |
| TR-S0900-S | N/A | 9.0 mm x 0.0 mm Stem |
| TR-S0902-S | N/A | 9.0 mm x 2.0 mm Stem |
| TR-S0904-S | N/A | 9.0 mm x 4.0 mm Stem |
| TR-S0906-S | N/A | 9.0 mm x 6.0 mm Stem |
| TR-S0908-S | N/A | 9.0 mm x 8.0 mm Stem |
| TR-S1000-S | N/A | 10.0 mm x 0.0 mm Stem |
| TR-S1002-S | N/A | 10.0 mm x 2.0 mm Stem |
| TR-S1004-S | N/A | 10.0 mm x 4.0 mm Stem |
| TR-S1006-S | N/A | 10.0 mm x 6.0 mm Stem |
| TR-S1008-S | N/A | 10.0 mm x 8.0 mm Stem |

Anatomic Radial Head Portfolio [continued]

Anatomic Radial Head Solutions

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|----------------------|
| TR-H200L-S | N/A | 20.0 mm Head, Left |
| TR-H200R-S | N/A | 20.0 mm Head, Right |
| TR-H220L-S | N/A | 22.0 mm Head, Left |
| TR-H220R-S | N/A | 22.0 mm Head, Right |
| TR-H240L-S | N/A | 24.0 mm Head, Left |
| TR-H240R-S | N/A | 24.0 mm Head, Right |
| TR-H260L-S | N/A | 26.0 mm Head, Left |
| TR-H260R-S | N/A | 26.0 mm Head, Right |
| TR-H280L-S | N/A | 28.0 mm Head, Left |
| TR-H280R-S | N/A | 28.0 mm Head, Right |
| TR-S0600-S | N/A | 6.0 mm x 0.0 mm Stem |
| TR-S0602-S | N/A | 6.0 mm x 2.0 mm Stem |
| TR-S0604-S | N/A | 6.0 mm x 4.0 mm Stem |
| TR-S0606-S | N/A | 6.0 mm x 6.0 mm Stem |
| TR-S0608-S | N/A | 6.0 mm x 8.0 mm Stem |
| TR-S0700-S | N/A | 7.0 mm x 0.0 mm Stem |
| TR-S0702-S | N/A | 7.0 mm x 2.0 mm Stem |
| TR-S0704-S | N/A | 7.0 mm x 4.0 mm Stem |
| TR-S0706-S | N/A | 7.0 mm x 6.0 mm Stem |
| TR-S0708-S | N/A | 7.0 mm x 8.0 mm Stem |
| TR-S0800-S | N/A | 8.0 mm x 0.0 mm Stem |
| TR-S0802-S | N/A | 8.0 mm x 2.0 mm Stem |
| TR-S0804-S | N/A | 8.0 mm x 4.0 mm Stem |
| TR-S0806-S | N/A | 8.0 mm x 6.0 mm Stem |
| TR-S0808-S | N/A | 8.0 mm x 8.0 mm Stem |
| TR-S0900-S | N/A | 9.0 mm x 0.0 mm Stem |
| TR-S0902-S | N/A | 9.0 mm x 2.0 mm Stem |
| TR-S0904-S | N/A | 9.0 mm x 4.0 mm Stem |
| TR-S0906-S | N/A | 9.0 mm x 6.0 mm Stem |

Anatomic Radial Head Portfolio [continued]

Anatomic Radial Head Solutions

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|-----------------------------|
| TR-S0908-S | N/A | 9.0 mm x 8.0 mm Stem |
| TR-S1000-S | N/A | 10.0 mm x 0.0 mm Stem |
| TR-S1002-S | N/A | 10.0 mm x 2.0 mm Stem |
| TR-S1004-S | N/A | 10.0 mm x 4.0 mm Stem |
| TR-S1006-S | N/A | 10.0 mm x 6.0 mm Stem |
| TR-S1008-S | N/A | 10.0 mm x 8.0 mm Stem |
| TR-SL06-S | N/A | 6 mm Morse Taper Long Stem |
| TR-SL08-S | N/A | 8 mm Morse Taper Long Stem |
| TR-SL10-S | N/A | 10 mm Morse Taper Long Stem |
| TR-SL12-S | N/A | 12 mm Morse Taper Long Stem |

Anatomic Radial Head Solutions 2 (stems and heads only)

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|---|
| 50-0056-S | N/A | 6.0 mm x 0.0 mm Anatomic Radial Head Partial Blast Stem |
| 50-0057-S | N/A | 6.0 mm x 2.0 mm Anatomic Radial Head Partial Blast Stem |
| 50-0058-S | N/A | 6.0 mm x 4.0 mm Anatomic Radial Head Partial Blast Stem |
| 50-0059-S | N/A | 6.0 mm x 6.0 mm Anatomic Radial Head Partial Blast Stem |
| 50-0060-S | N/A | 6.0 mm x 8.0 mm Anatomic Radial Head Partial Blast Stem |
| 50-0061-S | N/A | 7.0 mm x 0.0 mm Anatomic Radial Head Partial Blast Stem |
| 50-0062-S | N/A | 7.0 mm x 2.0 mm Anatomic Radial Head Partial Blast Stem |
| 50-0063-S | N/A | 7.0 mm x 4.0 mm Anatomic Radial Head Partial Blast Stem |
| 50-0064-S | N/A | 7.0 mm x 6.0 mm Anatomic Radial Head Partial Blast Stem |
| 50-0065-S | N/A | 7.0 mm x 8.0 mm Anatomic Radial Head Partial Blast Stem |
| 50-0066-S | N/A | 8.0 mm x 0.0 mm Anatomic Radial Head Partial Blast Stem |
| 50-0067-S | N/A | 8.0 mm x 2.0 mm Anatomic Radial Head Partial Blast Stem |
| 50-0068-S | N/A | 8.0 mm x 4.0 mm Anatomic Radial Head Partial Blast Stem |
| 50-0069-S | N/A | 8.0 mm x 6.0 mm Anatomic Radial Head Partial Blast Stem |
| 50-0070-S | N/A | 8.0 mm x 8.0 mm Anatomic Radial Head Partial Blast Stem |
| 50-0071-S | N/A | 9.0 mm x 0.0 mm Anatomic Radial Head Partial Blast Stem |

Anatomic Radial Head Portfolio [continued]

Anatomic Radial Head Solutions 2 (stems and heads only)

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|--|
| 50-0072-S | N/A | 9.0 mm x 2.0 mm Anatomic Radial Head Partial Blast Stem |
| 50-0073-S | N/A | 9.0 mm x 4.0 mm Anatomic Radial Head Partial Blast Stem |
| 50-0074-S | N/A | 9.0 mm x 6.0 mm Anatomic Radial Head Partial Blast Stem |
| 50-0075-S | N/A | 9.0 mm x 8.0 mm Anatomic Radial Head Partial Blast Stem |
| 50-0076-S | N/A | 10.0 mm x 0.0 mm Anatomic Radial Head Partial Blast Stem |
| 50-0077-S | N/A | 10.0 mm x 2.0 mm Anatomic Radial Head Partial Blast Stem |
| 50-0078-S | N/A | 10.0 mm x 4.0 mm Anatomic Radial Head Partial Blast Stem |
| 50-0079-S | N/A | 10.0 mm x 6.0 mm Anatomic Radial Head Partial Blast Stem |
| 50-0080-S | N/A | 10.0 mm x 8.0 mm Anatomic Radial Head Partial Blast Stem |
| 5001-0518L-S | N/A | Anatomic Radial Head Solutions 2 Head 18 mm, Left |
| 5001-0518R-S | N/A | Anatomic Radial Head Solutions 2 Head 18 mm, Right |
| 5001-0520L-S | N/A | Anatomic Radial Head Solutions 2 Head 20 mm, Left |
| 5001-0520R-S | N/A | Anatomic Radial Head Solutions 2 Head 20 mm, Right |
| 5001-0522L-S | N/A | Anatomic Radial Head Solutions 2 Head 22 mm, Left |
| 5001-0522R-S | N/A | Anatomic Radial Head Solutions 2 Head 22 mm, Right |
| 5001-0524L-S | N/A | Anatomic Radial Head Solutions 2 Head 24 mm, Left |
| 5001-0524R-S | N/A | Anatomic Radial Head Solutions 2 Head 24 mm, Right |
| 5001-0526L-S | N/A | Anatomic Radial Head Solutions 2 Head 26 mm, Left |
| 5001-0526R-S | N/A | Anatomic Radial Head Solutions 2 Head 26 mm, Right |
| 5001-0528L-S | N/A | Anatomic Radial Head Solutions 2 Head 28 mm, Left |
| 5001-0528R-S | N/A | Anatomic Radial Head Solutions 2 Head 28 mm, Right |
| 50-0109-S | N/A | 11.0 mm x 0.0 mm Anatomic Radial Head Partial Blast Stem |
| 50-0110-S | N/A | 11.0 mm x 2.0 mm Anatomic Radial Head Partial Blast Stem |
| 50-0111-S | N/A | 11.0 mm x 4.0 mm Anatomic Radial Head Partial Blast Stem |
| 50-0112-S | N/A | 11.0 mm x 6.0 mm Anatomic Radial Head Partial Blast Stem |
| 50-0113-S | N/A | 11.0 mm x 8.0 mm Anatomic Radial Head Partial Blast Stem |
| 50-0114-S | N/A | 12.0 mm x 0.0 mm Anatomic Radial Head Partial Blast Stem |
| 50-0115-S | N/A | 12.0 mm x 2.0 mm Anatomic Radial Head Partial Blast Stem |
| 50-0116-S | N/A | 12.0 mm x 4.0 mm Anatomic Radial Head Partial Blast Stem |

Anatomic Radial Head Portfolio [continued]

Anatomic Radial Head Solutions 2 (stems and heads only)

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|--|
| 50-0117-S | N/A | 12.0 mm x 6.0 mm Anatomic Radial Head Partial Blast Stem |
| 50-0118-S | N/A | 12.0 mm x 8.0 mm Anatomic Radial Head Partial Blast Stem |
| TR-S0600-S | N/A | 6.0 mm x 0.0 mm Stem |
| TR-S0602-S | N/A | 6.0 mm x 2.0 mm Stem |
| TR-S0604-S | N/A | 6.0 mm x 4.0 mm Stem |
| TR-S0606-S | N/A | 6.0 mm x 6.0 mm Stem |
| TR-S0608-S | N/A | 6.0 mm x 8.0 mm Stem |
| TR-S0700-S | N/A | 7.0 mm x 0.0 mm Stem |
| TR-S0702-S | N/A | 7.0 mm x 2.0 mm Stem |
| TR-S0704-S | N/A | 7.0 mm x 4.0 mm Stem |
| TR-S0706-S | N/A | 7.0 mm x 6.0 mm Stem |
| TR-S0708-S | N/A | 7.0 mm x 8.0 mm Stem |
| TR-S0800-S | N/A | 8.0 mm x 0.0 mm Stem |
| TR-S0802-S | N/A | 8.0 mm x 2.0 mm Stem |
| TR-S0804-S | N/A | 8.0 mm x 4.0 mm Stem |
| TR-S0806-S | N/A | 8.0 mm x 6.0 mm Stem |
| TR-S0808-S | N/A | 8.0 mm x 8.0 mm Stem |
| TR-S0900-S | N/A | 9.0 mm x 0.0 mm Stem |
| TR-S0902-S | N/A | 9.0 mm x 2.0 mm Stem |
| TR-S0904-S | N/A | 9.0 mm x 4.0 mm Stem |
| TR-S0906-S | N/A | 9.0 mm x 6.0 mm Stem |
| TR-S0908-S | N/A | 9.0 mm x 8.0 mm Stem |
| TR-S1000-S | N/A | 10.0 mm x 0.0 mm Stem |
| TR-S1002-S | N/A | 10.0 mm x 2.0 mm Stem |
| TR-S1004-S | N/A | 10.0 mm x 4.0 mm Stem |
| TR-S1006-S | N/A | 10.0 mm x 6.0 mm Stem |
| TR-S1008-S | N/A | 10.0 mm x 8.0 mm Stem |
| TR-S1100-S | N/A | 11.0 mm x 0.0 mm Stem |
| TR-S1102-S | N/A | 11.0 mm x 2.0 mm Stem |

Anatomic Radial Head Portfolio [continued]

Anatomic Radial Head Solutions 2 (stems and heads only)

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|-----------------------------|
| TR-S1104-S | N/A | 11.0 mm x 4.0 mm Stem |
| TR-S1106-S | N/A | 11.0 mm x 6.0 mm Stem |
| TR-S1108-S | N/A | 11.0 mm x 8.0 mm Stem |
| TR-S1200-S | N/A | 12.0 mm x 0.0 mm Stem |
| TR-S1202-S | N/A | 12.0 mm x 2.0 mm Stem |
| TR-S1204-S | N/A | 12.0 mm x 4.0 mm Stem |
| TR-S1206-S | N/A | 12.0 mm x 6.0 mm Stem |
| TR-S1208-S | N/A | 12.0 mm x 8.0 mm Stem |
| TR-SL06-S | N/A | 6 mm Morse Taper Long Stem |
| TR-SL07-S | N/A | 7 mm Morse Taper Long Stem |
| TR-SL08-S | N/A | 8 mm Morse Taper Long Stem |
| TR-SL09-S | N/A | 9 mm Morse Taper Long Stem |
| TR-SL10-S | N/A | 10 mm Morse Taper Long Stem |
| TR-SL11-S | N/A | 11 mm Morse Taper Long Stem |
| TR-SL12-S | N/A | 12 mm Morse Taper Long Stem |

Intramedullary Nails Portfolio

- ▶ Polarus® Humeral Rod System
- ▶ Polarus Plus Humeral Rod
- ▶ Polarus 3 Solution (nails)
- ▶ Fibula Rod System
- ▶ Forearm Fracture Solutions (nails)
- ▶ Fibula Nail 2 System
- ▶ Ulna Nail 2 System



Conditions for Safe Scanning in the MRI

A patient with the nail/rod, screw, washer, end cap, bushing, and set screw implants from the Polarus Humeral Rod, Polarus Plus Humeral Rod, Polarus 3 Solution, Fibula Rod, Forearm Fracture Solutions, Fibula Nail 2, and Ulna Nail 2 Systems may be safely scanned at 1.5 T or 3 T under the following conditions. Failure to follow these conditions may result in injury to the patient.

Intramedullary Nails Portfolio

| Parameter | Condition |
|---|--|
| Nominal values(s) of Static Magnetic Field | 1.5 T or 3 T |
| Maximum Spatial Field Gradient | 19 T/m (1900 gauss/cm) |
| RF Excitation | Circularly Polarized (CP) |
| RF Transmit Coil Type | Whole body transmit coil, Head RF transmit-receive coil |
| Maximum Whole Body SAR at 1.5 T and 3.0 T | 2.0 W/kg |
| Limits on Scan Duration | 2.0 W/kg whole body average SAR for 6 minutes of continuous RF (a sequence or back to back series/scan without breaks) |
| MR Image Artifact | The presence of this implant may produce an image artifact of 24 mm. |
| If more information about a specific parameter is not included, there are no conditions associated with that parameter. | |

Polarus® Humeral Rod System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|-----------------------------------|
| HCA5200-S | N/A | 5.0 mm x 20.0 mm Cancellous Screw |
| HCA5250-S | N/A | 5.0 mm x 25.0 mm Cancellous Screw |
| HCA5300-S | N/A | 5.0 mm x 30.0 mm Cancellous Screw |
| HCA5325-S | N/A | 5.0 mm x 32.5 mm Cancellous Screw |
| HCA5350-S | N/A | 5.0 mm x 35.0 mm Cancellous Screw |
| HCA5375-S | N/A | 5.0 mm x 37.5 mm Cancellous Screw |
| HCA5400-S | N/A | 5.0 mm x 40.0 mm Cancellous Screw |
| HCA5450-S | N/A | 5.0 mm x 45.0 mm Cancellous Screw |
| HCA5500-S | N/A | 5.0 mm x 50.0 mm Cancellous Screw |

Intramedullary Nails Portfolio [continued]

Polarus® Humeral Rod System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|--|
| HCA5550-S | N/A | 5.0 mm x 55.0 mm Cancellous Screw |
| HCA-L5325-S | N/A | 5.3 mm x 25 mm Cancellous Locking Screw |
| HCA-L5330-S | N/A | 5.3 mm x 30 mm Cancellous Locking Screw |
| HCA-L5335-S | N/A | 5.3 mm x 35 mm Cancellous Locking Screw |
| HCA-L5340-S | N/A | 5.3 mm x 40 mm Cancellous Locking Screw |
| HCA-L5345-S | N/A | 5.3 mm x 45 mm Cancellous Locking Screw |
| HCO3175-S | N/A | 3.5 mm x 17.5 mm Cortical Screw |
| HCO3200-S | N/A | 3.5 mm x 20.0 mm Cortical Screw |
| HCO3225-S | N/A | 3.5 mm x 22.5 mm Cortical Screw |
| HCO3250-S | N/A | 3.5 mm x 25.0 mm Cortical Screw |
| HCO3275-S | N/A | 3.5 mm x 27.5 mm Cortical Screw |
| HCO3300-S | N/A | 3.5 mm x 30.0 mm Cortical Screw |
| HCO3325-S | N/A | 3.5 mm x 32.5 mm Cortical Screw |
| HCO3350-S | N/A | 3.5 mm x 35.0 mm Cortical Screw |
| HCO3400-S | N/A | 3.5 mm x 40.0 mm Cortical Screw |
| HCO3450-S | N/A | 3.5 mm x 45.0 mm Cortical Screw |
| HCO3500-S | N/A | 3.5 mm x 50.0 mm Cortical Screw |
| HCO3550-S | N/A | 3.5 mm x 55.0 mm Cortical Screw |
| HCO3600-S | N/A | 3.5 mm x 60.0 mm Cortical Screw |
| HR-0001-S | N/A | Polarus Cap Screw Assembly |
| HR-0050-S | N/A | Polarus Cap Screw Assembly |
| HR-L1115-S | N/A | 11 mm x 150 mm Polarus Locking Humeral Rod |

Polarus® Plus Humeral Rod System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|-----------------------------------|
| HCA5200-S | N/A | 5.0 mm x 20.0 mm Cancellous Screw |
| HCA5250-S | N/A | 5.0 mm x 25.0 mm Cancellous Screw |
| HCA5300-S | N/A | 5.0 mm x 30.0 mm Cancellous Screw |
| HCA5325-S | N/A | 5.0 mm x 32.5 mm Cancellous Screw |

Intramedullary Nails Portfolio [continued]

Polarus® Plus Humeral Rod System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|--|
| HCA5350-S | N/A | 5.0 mm x 35.0 mm Cancellous Screw |
| HCA5375-S | N/A | 5.0 mm x 37.5 mm Cancellous Screw |
| HCA5400-S | N/A | 5.0 mm x 40.0 mm Cancellous Screw |
| HCA5450-S | N/A | 5.0 mm x 45.0 mm Cancellous Screw |
| HCA5500-S | N/A | 5.0 mm x 50.0 mm Cancellous Screw |
| HCA5550-S | N/A | 5.0 mm x 55.0 mm Cancellous Screw |
| HCO3175-S | N/A | 3.5 mm x 17.5 mm Cortical Screw |
| HCO3200-S | N/A | 3.5 mm x 20.0 mm Cortical Screw |
| HCO3225-S | N/A | 3.5 mm x 22.5 mm Cortical Screw |
| HCO3250-S | N/A | 3.5 mm x 25.0 mm Cortical Screw |
| HCO3275-S | N/A | 3.5 mm x 27.5 mm Cortical Screw |
| HCO3300-S | N/A | 3.5 mm x 30.0 mm Cortical Screw |
| HCO3325-S | N/A | 3.5 mm x 32.5 mm Cortical Screw |
| HCO3350-S | N/A | 3.5 mm x 35.0 mm Cortical Screw |
| HCO3400-S | N/A | 3.5 mm x 40.0 mm Cortical Screw |
| HCO3450-S | N/A | 3.5 mm x 45.0 mm Cortical Screw |
| HCO3500-S | N/A | 3.5 mm x 50.0 mm Cortical Screw |
| HCO3550-S | N/A | 3.5 mm x 55.0 mm Cortical Screw |
| HCO35WA | N/A | 3.5 mm Cortical Screw Washer |
| HCO35WA-S | N/A | 3.5 mm Cortical Screw Washer |
| HCO3600-S | N/A | 3.5 mm x 60.0 mm Cortical Screw |
| HR-0001-S | N/A | Polarus Cap Screw Assembly |
| HR-0050-S | N/A | Polarus Cap Screw Assembly |
| HR-0820-S | N/A | 8 mm x 200 mm Polarus Plus Humeral Rod |
| HR-0822-S | N/A | 8 mm x 220 mm Polarus Plus Humeral Rod |
| HR-0824-S | N/A | 8 mm x 240 mm Polarus Plus Humeral Rod |
| HR-0826-S | N/A | 8 mm x 260 mm Polarus Plus Humeral Rod |
| HR-0828-S | N/A | 8 mm x 280 mm Polarus Plus Humeral Rod |

Intramedullary Nails Portfolio [continued]

Polarus® 3 Solution (Locking and Proximal Locking Nails)

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|--|
| 4001-1015L-S | N/A | Polarus 3 Proximal Locking Nail 150 mm Left |
| 4001-1015R-S | N/A | Polarus 3 Proximal Locking Nail 150 mm Right |
| 4002-10200-S | N/A | Polarus 3 Locking Nail 200 mm |
| 4002-10220-S | N/A | Polarus 3 Locking Nail 220 mm |
| 4002-10240-S | N/A | Polarus 3 Locking Nail 240 mm |
| 4002-10260-S | N/A | Polarus 3 Locking Nail 260 mm |
| 4002-10280-S | N/A | Polarus 3 Locking Nail 280 mm |
| 4004-10000-S | N/A | Polarus 3 Cap Screw 0 mm |
| 4004-10002-S | N/A | Polarus 3 Cap Screw 2 mm |
| 4004-10004-S | N/A | Polarus 3 Cap Screw 4 mm |
| 4004-10006-S | N/A | Polarus 3 Cap Screw 6 mm |
| 3011-43018-S | 3011-43018 | 4.3 mm x 18 mm Low-Profile Hexalobe Screw |
| 3011-43020-S | 3011-43020 | 4.3 mm x 20 mm Low-Profile Hexalobe Screw |
| 3011-43022-S | 3011-43022 | 4.3 mm x 22 mm Low-Profile Hexalobe Screw |
| 3011-43024-S | 3011-43024 | 4.3 mm x 24 mm Low-Profile Hexalobe Screw |
| 3011-43026-S | 3011-43026 | 4.3 mm x 26 mm Low-Profile Hexalobe Screw |
| 3011-43028-S | 3011-43028 | 4.3 mm x 28 mm Low-Profile Hexalobe Screw |
| 3011-43030-S | 3011-43030 | 4.3 mm x 30 mm Low-Profile Hexalobe Screw |
| 3011-43032-S | 3011-43032 | 4.3 mm x 32 mm Low-Profile Hexalobe Screw |
| 3011-43034-S | 3011-43034 | 4.3 mm x 34 mm Low-Profile Hexalobe Screw |
| 3011-43036-S | 3011-43036 | 4.3 mm x 36 mm Low-Profile Hexalobe Screw |
| 3011-43038-S | 3011-43038 | 4.3 mm x 38 mm Low-Profile Hexalobe Screw |
| 3011-43040-S | 3011-43040 | 4.3 mm x 40 mm Low-Profile Hexalobe Screw |
| 3011-43042-S | 3011-43042 | 4.3 mm x 42 mm Low-Profile Hexalobe Screw |
| 3011-43044-S | 3011-43044 | 4.3 mm x 44 mm Low-Profile Hexalobe Screw |
| 3011-43046-S | 3011-43046 | 4.3 mm x 46 mm Low-Profile Hexalobe Screw |
| 3011-43048-S | 3011-43048 | 4.3 mm x 48 mm Low-Profile Hexalobe Screw |
| 3011-43050-S | 3011-43050 | 4.3 mm x 50 mm Low-Profile Hexalobe Screw |
| 3011-43052-S | 3011-43052 | 4.3 mm x 52 mm Low-Profile Hexalobe Screw |

Intramedullary Nails Portfolio [continued]

Polarus® 3 Solution (Locking and Proximal Locking Nails)

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|---|
| 3011-43054-S | 3011-43054 | 4.3 mm x 54 mm Low-Profile Hexalobe Screw |
| 3011-43056-S | 3011-43056 | 4.3 mm x 56 mm Low-Profile Hexalobe Screw |
| 3011-43058-S | 3011-43058 | 4.3 mm x 58 mm Low-Profile Hexalobe Screw |
| 3011-43060-S | 3011-43060 | 4.3 mm x 60 mm Low-Profile Hexalobe Screw |
| 3011-43062-S | 3011-43062 | 4.3 mm x 62 mm Low-Profile Hexalobe Screw |
| 3011-43064-S | 3011-43064 | 4.3 mm x 64 mm Low-Profile Hexalobe Screw |
| 3025-35018-S | 3025-35018 | 3.5 x 18 mm Nonlocking Low-Profile Hexalobe Screw |
| 3025-35020-S | 3025-35020 | 3.5 x 20 mm Nonlocking Low-Profile Hexalobe Screw |
| 3025-35022-S | 3025-35022 | 3.5 x 22 mm Nonlocking Low-Profile Hexalobe Screw |
| 3025-35024-S | 3025-35024 | 3.5 x 24 mm Nonlocking Low-Profile Hexalobe Screw |
| 3025-35026-S | 3025-35026 | 3.5 x 26 mm Nonlocking Low-Profile Hexalobe Screw |
| 3025-35028-S | 3025-35028 | 3.5 x 28 mm Nonlocking Low-Profile Hexalobe Screw |
| 3025-35030-S | 3025-35030 | 3.5 x 30 mm Nonlocking Low-Profile Hexalobe Screw |
| 3025-35032-S | 3025-35032 | 3.5 x 32 mm Nonlocking Low-Profile Hexalobe Screw |
| 3025-35034-S | 3025-35034 | 3.5 x 34 mm Nonlocking Low-Profile Hexalobe Screw |
| 3025-35036-S | 3025-35036 | 3.5 x 36 mm Nonlocking Low-Profile Hexalobe Screw |
| 3025-35038-S | 3025-35038 | 3.5 x 38 mm Nonlocking Low-Profile Hexalobe Screw |
| 3025-35040-S | 3025-35040 | 3.5 x 40 mm Nonlocking Low-Profile Hexalobe Screw |
| 3025-35042-S | 3025-35042 | 3.5 x 42 mm Nonlocking Low-Profile Hexalobe Screw |
| 3025-35044-S | 3025-35044 | 3.5 x 44 mm Nonlocking Low-Profile Hexalobe Screw |
| 3025-35046-S | 3025-35046 | 3.5 x 46 mm Nonlocking Low-Profile Hexalobe Screw |
| 3025-35048-S | 3025-35048 | 3.5 x 48 mm Nonlocking Low-Profile Hexalobe Screw |
| 3025-35050-S | 3025-35050 | 3.5 x 50 mm Nonlocking Low-Profile Hexalobe Screw |
| 3025-35052-S | 3025-35052 | 3.5 x 52 mm Nonlocking Low-Profile Hexalobe Screw |
| 3025-35054-S | 3025-35054 | 3.5 x 54 mm Nonlocking Low-Profile Hexalobe Screw |
| 3025-35056-S | 3025-35056 | 3.5 x 56 mm Nonlocking Low-Profile Hexalobe Screw |
| 3025-35058-S | 3025-35058 | 3.5 x 58 mm Nonlocking Low-Profile Hexalobe Screw |
| 3025-35060-S | 3025-35060 | 3.5 x 60 mm Nonlocking Low-Profile Hexalobe Screw |
| 3025-35062-S | 3025-35062 | 3.5 x 62 mm Nonlocking Low-Profile Hexalobe Screw |
| 3025-35064-S | 3025-35064 | 3.5 x 64 mm Nonlocking Low-Profile Hexalobe Screw |

Intramedullary Nails Portfolio [continued]

Fibula Rod System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|---------------------------------|
| 40-0026-S | N/A | 3.0 mm x 110 mm Fibula Rod |
| 40-0027-S | N/A | 3.0 mm x 145 mm Fibula Rod |
| 40-0028-S | N/A | 3.0 mm x 180 mm Fibula Rod |
| 40-0029-S | N/A | 3.6 mm x 110 mm Fibula Rod |
| 40-0030-S | N/A | 3.6 mm x 145 mm Fibula Rod |
| 40-0031-S | N/A | 3.6 mm x 180 mm Fibula Rod |
| CO-3080-S | CO-3080 | 3.5 mm x 8.0 mm Cortical Screw |
| CO-3100-S | CO-3100 | 3.5 mm x 10.0 mm Cortical Screw |
| CO-3120-S | CO-3120 | 3.5 mm x 12.0 mm Cortical Screw |
| CO-3140-S | CO-3140 | 3.5 mm x 14.0 mm Cortical Screw |
| CO-3160-S | CO-3160 | 3.5 mm x 16.0 mm Cortical Screw |
| CO-3180-S | CO-3180 | 3.5 mm x 18.0 mm Cortical Screw |
| CO-3200-S | CO-3200 | 3.5 mm x 20.0 mm Cortical Screw |
| CO-3220-S | CO-3220 | 3.5 mm x 22.0 mm Cortical Screw |
| CO-3400-S | CO-3400 | 3.5 mm x 40.0 mm Cortical Screw |
| CO-3450-S | CO-3450 | 3.5 mm x 45.0 mm Cortical Screw |
| CO-3500-S | CO-3500 | 3.5 mm x 50.0 mm Cortical Screw |
| CO-3550-S | CO-3550 | 3.5 mm x 55.0 mm Cortical Screw |
| CO-3600-S | CO-3600 | 3.5 mm x 60.0 mm Cortical Screw |
| CO-3650-S | CO-3650 | 3.5 mm x 65.0 mm Cortical Screw |

Forearm Fracture Solutions (Radius & Ulna Rods)

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|---------------------------------|
| CO-3125-S | N/A | 3.5 mm x 12.5 mm Cortical Screw |
| CO-3150-S | N/A | 3.5 mm x 15.0 mm Cortical Screw |
| CO-3175-S | N/A | 3.5 mm x 17.5 mm Cortical Screw |
| CO-3200-S | CO-3200 | 3.5 mm x 20.0 mm Cortical Screw |
| CO-3225-S | CO-3225 | 3.5 mm x 22.5 mm Cortical Screw |
| CO-3250-S | CO-3250 | 3.5 mm x 25.0 mm Cortical Screw |

Intramedullary Nails Portfolio [continued]

Forearm Fracture Solutions (Radius & Ulna Rods)

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|-----------------------------------|
| CO-3275-S | CO-3275 | 3.5 mm x 27.5 mm Cortical Screw |
| CO-3300-S | CO-3300 | 3.5 mm x 30.0 mm Cortical Screw |
| CO-3325-S | CO-3325 | 3.5 mm x 32.5 mm Cortical Screw |
| RL-3019-S | N/A | 3.0 mm x 190 mm Radius Rod, Left |
| RL-3021-S | N/A | 3.0 mm x 210 mm Radius Rod, Left |
| RL-3023-S | N/A | 3.0 mm x 230 mm Radius Rod, Left |
| RL-3619-S | N/A | 3.6 mm x 190 mm Radius Rod, Left |
| RL-3621-S | N/A | 3.6 mm x 210 mm Radius Rod, Left |
| RL-3623-S | N/A | 3.6 mm x 230 mm Radius Rod, Left |
| RR-3019-S | N/A | 3.0 mm x 190 mm Radius Rod, Right |
| RR-3021-S | N/A | 3.0 mm x 210 mm Radius Rod, Right |
| RR-3023-S | N/A | 3.0 mm x 230 mm Radius Rod, Right |
| RR-3619-S | N/A | 3.6 mm x 190 mm Radius Rod, Right |
| RR-3621-S | N/A | 3.6 mm x 210 mm Radius Rod, Right |
| RR-3623-S | N/A | 3.6 mm x 230 mm Radius Rod, Right |
| UL-3021-S | N/A | 3.0 mm x 210 mm Ulna Rod |
| UL-3023-S | N/A | 3.0 mm x 230 mm Ulna Rod |
| UL-3025-S | N/A | 3.0 mm x 250 mm Ulna Rod |
| UL-3027-S | N/A | 3.0 mm x 270 mm Ulna Rod |
| UL-3621-S | N/A | 3.6 mm x 210 mm Ulna Rod |
| UL-3623-S | N/A | 3.6 mm x 230 mm Ulna Rod |
| UL-3625-S | N/A | 3.6 mm x 250 mm Ulna Rod |
| UL-3627-S | N/A | 3.6 mm x 270 mm Ulna Rod |

Fibula Nail 2 System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|--|
| 30-0255-S | 30-0255 | 3.5 mm x 8 mm Nonlocking Hexalobe Screw |
| 30-0256-S | 30-0256 | 3.5 mm x 10 mm Nonlocking Hexalobe Screw |
| 30-0257-S | 30-0257 | 3.5 mm x 12 mm Nonlocking Hexalobe Screw |

Intramedullary Nails Portfolio [continued]

Fibula Nail 2 System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|--|
| 30-0258-S | 30-0258 | 3.5 mm x 14 mm Nonlocking Hexalobe Screw |
| 30-0259-S | 30-0259 | 3.5 mm x 16 mm Nonlocking Hexalobe Screw |
| 30-0260-S | 30-0260 | 3.5 mm x 18 mm Nonlocking Hexalobe Screw |
| 30-0261-S | 30-0261 | 3.5 mm x 20 mm Nonlocking Hexalobe Screw |
| 30-0262-S | 30-0262 | 3.5 mm x 22 mm Nonlocking Hexalobe Screw |
| 30-0263-S | 30-0263 | 3.5 mm x 24 mm Nonlocking Hexalobe Screw |
| 30-0264-S | 30-0264 | 3.5 mm x 26 mm Nonlocking Hexalobe Screw |
| 30-0265-S | 30-0265 | 3.5 mm x 28 mm Nonlocking Hexalobe Screw |
| 30-0266-S | 30-0266 | 3.5 mm x 30 mm Nonlocking Hexalobe Screw |
| 30-0267-S | 30-0267 | 3.5 mm x 32 mm Nonlocking Hexalobe Screw |
| 30-0268-S | 30-0268 | 3.5 mm x 34 mm Nonlocking Hexalobe Screw |
| 30-0269-S | 30-0269 | 3.5 mm x 36 mm Nonlocking Hexalobe Screw |
| 30-0270-S | 30-0270 | 3.5 mm x 38 mm Nonlocking Hexalobe Screw |
| 30-0271-S | 30-0271 | 3.5 mm x 40 mm Nonlocking Hexalobe Screw |
| 30-0272-S | 30-0272 | 3.5 mm x 45 mm Nonlocking Hexalobe Screw |
| 30-0273-S | 30-0273 | 3.5 mm x 50 mm Nonlocking Hexalobe Screw |
| 30-0274-S | 30-0274 | 3.5 mm x 55 mm Nonlocking Hexalobe Screw |
| 30-0275-S | 30-0275 | 3.5 mm x 60 mm Nonlocking Hexalobe Screw |
| 30-0276-S | 30-0276 | 3.5 mm x 65 mm Nonlocking Hexalobe Screw |
| N/A | 3017-25006 | Far End Locking Set Screw, 6 mm |
| N/A | 3017-25007 | Far End Locking Set Screw, 7 mm |
| N/A | 3017-25008 | Far End Locking Set Screw, 8 mm |
| N/A | 3017-25009 | Far End Locking Set Screw, 9 mm |
| N/A | 3017-25010 | Far End Locking Set Screw, 10 mm |
| N/A | 3017-25011 | Far End Locking Set Screw, 11 mm |
| N/A | 3017-25012 | Far End Locking Set Screw, 12 mm |
| N/A | 3017-25013 | Far End Locking Set Screw, 13 mm |
| N/A | 3017-25014 | Far End Locking Set Screw, 14 mm |
| N/A | 3017-25015 | Far End Locking Set Screw, 15 mm |

Intramedullary Nails Portfolio [continued]

Fibula Nail 2 System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|--|
| N/A | 3017-25016 | Far End Locking Set Screw, 16 mm |
| N/A | 3017-65006 | Far End Locking Bushing, 6 mm |
| N/A | 3017-65007 | Far End Locking Bushing, 7 mm |
| N/A | 3017-65008 | Far End Locking Bushing, 8 mm |
| N/A | 3017-65009 | Far End Locking Bushing, 9 mm |
| N/A | 3017-65010 | Far End Locking Bushing, 10 mm |
| N/A | 3017-65011 | Far End Locking Bushing, 11 mm |
| N/A | 3017-65012 | Far End Locking Bushing, 12 mm |
| N/A | 3017-65013 | Far End Locking Bushing, 13 mm |
| N/A | 3017-65014 | Far End Locking Bushing, 14 mm |
| N/A | 3017-65015 | Far End Locking Bushing, 15 mm |
| N/A | 3017-65016 | Far End Locking Bushing, 16 mm |
| N/A | 3018-47012 | 3.5 mm x 12 mm Headless Hexalobe Screw |
| N/A | 3018-47014 | 3.5 mm x 14 mm Headless Hexalobe Screw |
| N/A | 3018-47016 | 3.5 mm x 16 mm Headless Hexalobe Screw |
| N/A | 3018-47018 | 3.5 mm x 18 mm Headless Hexalobe Screw |
| N/A | 3018-47020 | 3.5 mm x 20 mm Headless Hexalobe Screw |
| N/A | 3018-47022 | 3.5 mm x 22 mm Headless Hexalobe Screw |
| N/A | 3018-47024 | 3.5 mm x 24 mm Headless Hexalobe Screw |
| N/A | 3018-47026 | 3.5 mm x 26 mm Headless Hexalobe Screw |
| N/A | 3018-47028 | 3.5 mm x 28 mm Headless Hexalobe Screw |
| N/A | 3018-47030 | 3.5 mm x 30 mm Headless Hexalobe Screw |
| N/A | 3018-47032 | 3.5 mm x 32 mm Headless Hexalobe Screw |
| N/A | 3018-47034 | 3.5 mm x 34 mm Headless Hexalobe Screw |
| N/A | 3018-47036 | 3.5 mm x 36 mm Headless Hexalobe Screw |
| N/A | 3018-47038 | 3.5 mm x 38 mm Headless Hexalobe Screw |
| N/A | 3018-47040 | 3.5 mm x 40 mm Headless Hexalobe Screw |
| N/A | 3018-47045 | 3.5 mm x 45 mm Headless Hexalobe Screw |
| N/A | 3018-47050 | 3.5 mm x 50 mm Headless Hexalobe Screw |

Intramedullary Nails Portfolio [continued]

Fibula Nail 2 System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|---|
| N/A | 3018-47055 | 3.5 mm x 55 mm Headless Hexalobe Screw |
| N/A | 3018-47060 | 3.5 mm x 60 mm Headless Hexalobe Screw |
| N/A | 3018-47065 | 3.5 mm x 65 mm Headless Hexalobe Screw |
| 4010-2611N-S | N/A | 2.6 mm x 110 mm Fibula Nail |
| 4010-2614N-S | N/A | 2.6 mm x 145 mm Fibula Nail |
| 4010-2619N-S | N/A | 2.6 mm x 190 mm Fibula Nail |
| 4010-2627N-S | N/A | 2.6 mm x 270 mm Fibula Nail |
| 4010-3011N-S | N/A | 3.0 mm x 110 mm Fibula Nail |
| 4010-3014N-S | N/A | 3.0 mm x 145 mm Fibula Nail |
| 4010-3019N-S | N/A | 3.0 mm x 190 mm Fibula Nail |
| 4010-3027N-S | N/A | 3.0 mm x 270 mm Fibula Nail |
| 4010-3611N-S | N/A | 3.6 mm x 110 mm Fibula Nail |
| 4010-3614N-S | N/A | 3.6 mm x 145 mm Fibula Nail |
| 4010-3619N-S | N/A | 3.6 mm x 190 mm Fibula Nail |
| 4010-3627N-S | N/A | 3.6 mm x 270 mm Fibula Nail |
| N/A | 4014-0600 | Fibula and Forearm Nail +0.4 mm End Cap |
| N/A | 4014-0705 | Fibula and Forearm Nail +5 mm End Cap |
| N/A | 4014-0710 | Fibula and Forearm Nail +10 mm End Cap |
| N/A | 4014-0715 | Fibula and Forearm Nail +15 mm End Cap |

Ulna Nail 2 System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|--|
| 30-0255-S | 30-0255 | 3.5 mm x 8 mm Nonlocking Hexalobe Screw |
| 30-0256-S | 30-0256 | 3.5 mm x 10 mm Nonlocking Hexalobe Screw |
| 30-0257-S | 30-0257 | 3.5 mm x 12 mm Nonlocking Hexalobe Screw |
| 30-0258-S | 30-0258 | 3.5 mm x 14 mm Nonlocking Hexalobe Screw |
| 30-0259-S | 30-0259 | 3.5 mm x 16 mm Nonlocking Hexalobe Screw |
| 30-0260-S | 30-0260 | 3.5 mm x 18 mm Nonlocking Hexalobe Screw |
| 30-0261-S | 30-0261 | 3.5 mm x 20 mm Nonlocking Hexalobe Screw |

Intramedullary Nails Portfolio [continued]

Ulna Nail 2 System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|--|
| 30-0262-S | 30-0262 | 3.5 mm x 22 mm Nonlocking Hexalobe Screw |
| 30-0263-S | 30-0263 | 3.5 mm x 24 mm Nonlocking Hexalobe Screw |
| 30-0264-S | 30-0264 | 3.5 mm x 26 mm Nonlocking Hexalobe Screw |
| 30-0265-S | 30-0265 | 3.5 mm x 28 mm Nonlocking Hexalobe Screw |
| 30-0266-S | 30-0266 | 3.5 mm x 30 mm Nonlocking Hexalobe Screw |
| 30-0267-S | 30-0267 | 3.5 mm x 32 mm Nonlocking Hexalobe Screw |
| 30-0268-S | 30-0268 | 3.5 mm x 34 mm Nonlocking Hexalobe Screw |
| 30-0269-S | 30-0269 | 3.5 mm x 36 mm Nonlocking Hexalobe Screw |
| 30-0270-S | 30-0270 | 3.5 mm x 38 mm Nonlocking Hexalobe Screw |
| 30-0271-S | 30-0271 | 3.5 mm x 40 mm Nonlocking Hexalobe Screw |
| 30-0272-S | 30-0272 | 3.5 mm x 45 mm Nonlocking Hexalobe Screw |
| 30-0273-S | 30-0273 | 3.5 mm x 50 mm Nonlocking Hexalobe Screw |
| 30-0274-S | 30-0274 | 3.5 mm x 55 mm Nonlocking Hexalobe Screw |
| 30-0275-S | 30-0275 | 3.5 mm x 60 mm Nonlocking Hexalobe Screw |
| 30-0276-S | 30-0276 | 3.5 mm x 65 mm Nonlocking Hexalobe Screw |
| N/A | 3017-25006 | Far End Locking Set Screw, 6 mm |
| N/A | 3017-25007 | Far End Locking Set Screw, 7 mm |
| N/A | 3017-25008 | Far End Locking Set Screw, 8 mm |
| N/A | 3017-25009 | Far End Locking Set Screw, 9 mm |
| N/A | 3017-25010 | Far End Locking Set Screw, 10 mm |
| N/A | 3017-25011 | Far End Locking Set Screw, 11 mm |
| N/A | 3017-25012 | Far End Locking Set Screw, 12 mm |
| N/A | 3017-25013 | Far End Locking Set Screw, 13 mm |
| N/A | 3017-25014 | Far End Locking Set Screw, 14 mm |
| N/A | 3017-25015 | Far End Locking Set Screw, 15 mm |
| N/A | 3017-25016 | Far End Locking Set Screw, 16 mm |
| N/A | 3017-65006 | Far End Locking Bushing, 6 mm |
| N/A | 3017-65007 | Far End Locking Bushing, 7 mm |
| N/A | 3017-65008 | Far End Locking Bushing, 8 mm |

Intramedullary Nails Portfolio [continued]

Ulna Nail 2 System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|--|
| N/A | 3017-65009 | Far End Locking Bushing, 9 mm |
| N/A | 3017-65010 | Far End Locking Bushing, 10 mm |
| N/A | 3017-65011 | Far End Locking Bushing, 11 mm |
| N/A | 3017-65012 | Far End Locking Bushing, 12 mm |
| N/A | 3017-65013 | Far End Locking Bushing, 13 mm |
| N/A | 3017-65014 | Far End Locking Bushing, 14 mm |
| N/A | 3017-65015 | Far End Locking Bushing, 15 mm |
| N/A | 3017-65016 | Far End Locking Bushing, 16 mm |
| N/A | 3018-47012 | 3.5 mm x 12 mm Headless Hexalobe Screw |
| N/A | 3018-47014 | 3.5 mm x 14 mm Headless Hexalobe Screw |
| N/A | 3018-47016 | 3.5 mm x 16 mm Headless Hexalobe Screw |
| N/A | 3018-47018 | 3.5 mm x 18 mm Headless Hexalobe Screw |
| N/A | 3018-47020 | 3.5 mm x 20 mm Headless Hexalobe Screw |
| N/A | 3018-47022 | 3.5 mm x 22 mm Headless Hexalobe Screw |
| N/A | 3018-47024 | 3.5 mm x 24 mm Headless Hexalobe Screw |
| N/A | 3018-47026 | 3.5 mm x 26 mm Headless Hexalobe Screw |
| N/A | 3018-47028 | 3.5 mm x 28 mm Headless Hexalobe Screw |
| N/A | 3018-47030 | 3.5 mm x 30 mm Headless Hexalobe Screw |
| N/A | 3018-47032 | 3.5 mm x 32 mm Headless Hexalobe Screw |
| N/A | 3018-47034 | 3.5 mm x 34 mm Headless Hexalobe Screw |
| N/A | 3018-47036 | 3.5 mm x 36 mm Headless Hexalobe Screw |
| N/A | 3018-47038 | 3.5 mm x 38 mm Headless Hexalobe Screw |
| N/A | 3018-47040 | 3.5 mm x 40 mm Headless Hexalobe Screw |
| N/A | 3018-47045 | 3.5 mm x 45 mm Headless Hexalobe Screw |
| N/A | 3018-47050 | 3.5 mm x 50 mm Headless Hexalobe Screw |
| N/A | 3018-47055 | 3.5 mm x 55 mm Headless Hexalobe Screw |
| N/A | 3018-47060 | 3.5 mm x 60 mm Headless Hexalobe Screw |
| N/A | 3018-47065 | 3.5 mm x 65 mm Headless Hexalobe Screw |
| 4011-3012N-S | N/A | 3.0 mm x 120 mm Short Ulna Nail |

Intramedullary Nails Portfolio [continued]

Ulna Nail 2 System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|---|
| 4011-3017N-S | N/A | 3.0 mm x 170 mm Ulna Nail |
| 4011-3019N-S | N/A | 3.0 mm x 190 mm Ulna Nail |
| 4011-3021N-S | N/A | 3.0 mm x 210 mm Ulna Nail |
| 4011-3023N-S | N/A | 3.0 mm x 230 mm Ulna Nail |
| 4011-3025N-S | N/A | 3.0 mm x 250 mm Ulna Nail |
| 4011-3027N-S | N/A | 3.0 mm x 270 mm Ulna Nail |
| 4011-3612N-S | N/A | 3.6 mm x 120 mm Short Ulna Nail |
| 4011-3617N-S | N/A | 3.6 mm x 170 mm Ulna Nail |
| 4011-3619N-S | N/A | 3.6 mm x 190 mm Ulna Nail |
| 4011-3621N-S | N/A | 3.6 mm x 210 mm Ulna Nail |
| 4011-3623N-S | N/A | 3.6 mm x 230 mm Ulna Nail |
| 4011-3625N-S | N/A | 3.6 mm x 250 mm Ulna Nail |
| 4011-3627N-S | N/A | 3.6 mm x 270 mm Ulna Nail |
| 4011-4012N-S | N/A | 4.0 mm x 120 mm Short Ulna Nail |
| 4011-4017N-S | N/A | 4.0 mm x 170 mm Ulna Nail |
| 4011-4019N-S | N/A | 4.0 mm x 190 mm Ulna Nail |
| 4011-4021N-S | N/A | 4.0 mm x 210 mm Ulna Nail |
| 4011-4023N-S | N/A | 4.0 mm x 230 mm Ulna Nail |
| 4011-4025N-S | N/A | 4.0 mm x 250 mm Ulna Nail |
| 4011-4027N-S | N/A | 4.0 mm x 270 mm Ulna Nail |
| N/A | 4014-0600 | Fibula and Forearm Nail +0.4 mm End Cap |
| N/A | 4014-0705 | Fibula and Forearm Nail +5 mm End Cap |
| N/A | 4014-0710 | Fibula and Forearm Nail +10 mm End Cap |
| N/A | 4014-0715 | Fibula and Forearm Nail +15 mm End Cap |

Hand and Wrist Portfolio

- ▶ Acu-Loc® Wrist Plating System
- ▶ Acu-Loc 2 Wrist Plating System
- ▶ Acu-Loc Wrist Spanning Plate
- ▶ Forearm Fracture Solutions (plates)
- ▶ Hand Fracture System
- ▶ Modular Hand System
- ▶ Osteotomy System
- ▶ Total Wrist Fusion Plating System



Conditions for Safe Scanning in the MRI

A patient with the peg, plate, screw, sleeve, and screw cover implants from the Acu-Loc Wrist Plating, Acu-Loc 2 Wrist Plating, Acu-Loc Wrist Spanning Plate, Anatomic Midshaft Forearm Plating, Hand Fracture, Modular Hand, Osteotomy, and Total Wrist Fusion Plating Systems may be safely scanned at 1.5 T or 3 T under the following conditions. Failure to follow these conditions may result in injury to the patient.

Hand and Wrist Portfolio

| Parameter | Condition |
|---|---|
| Nominal values(s) of Static Magnetic Field | 1.5 T or 3 T |
| Maximum Spatial Field Gradient | 30 T/m (3000 gauss/cm) |
| RF Excitation | Circularly Polarized (CP) |
| RF Transmit Coil Type | Whole body transmit coil, Head RF transmit-receive coil |
| Maximum Whole Body SAR at 1.5 T and 3.0 T | 1.0 W/kg* |
| Limits on Scan Duration | 1.0 W/kg whole body average SAR for 15 minutes of continuous RF (a sequence or back to back series/scan without breaks) with 5 minutes of cooling for an hour-long scanning session |
| MR Image Artifact | The presence of this implant may produce an image artifact of 37 mm. |
| If more information about a specific parameter is not included, there are no conditions associated with that parameter. | |

*Whole body SAR of 1.0 W/kg may limit the availability and diagnostic quality of some MR procedures

Acu-Loc® Wrist Plating System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|--|
| 70-0045-S | 70-0045 | Acu-Loc Volar Distal Ulna Plate Left Standard |
| 70-0046-S | 70-0046 | Acu-Loc Volar Distal Ulna Plate Right Standard |
| 70-0047-S | 70-0047 | Acu-Loc Volar Distal Ulna Plate Left Long |
| 70-0048-S | 70-0048 | Acu-Loc Volar Distal Ulna Plate Right Long |
| 70-0055-S | 70-0055 | Acu-Loc Dorsal Plate Standard, Left |
| 70-0056-S | 70-0056 | Acu-Loc Dorsal Plate Standard, Right |

Hand and Wrist Portfolio [continued]

Acu-Loc® Wrist Plating System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|--|
| 70-0057-S | 70-0057 | Acu-Loc Dorsal Plate Narrow Left |
| 70-0058-S | 70-0058 | Acu-Loc Dorsal Plate Narrow Right |
| 70-0063-S | 70-0063 | Acu-Loc Extra Articular Standard |
| 70-0064-S | 70-0064 | Acu-Loc Extra Articular Narrow |
| CO-3100-S | CO-3100 | 3.5 mm x 10.0 mm Cortical Screw |
| CO-3120-S | CO-3120 | 3.5 mm x 12.0 mm Cortical Screw |
| CO-3140-S | CO-3140 | 3.5 mm x 14.0 mm Cortical Screw |
| CO-3160-S | CO-3160 | 3.5 mm x 16.0 mm Cortical Screw |
| CO-3180-S | CO-3180 | 3.5 mm x 18.0 mm Cortical Screw |
| COL-3080-S | COL-3080 | 3.5 mm x 8.0 mm Locking Cortical Screw |
| COL-3100-S | COL-3100 | 3.5 mm x 10.0 mm Locking Cortical Screw |
| COL-3120-S | COL-3120 | 3.5 mm x 12.0 mm Locking Cortical Screw |
| COL-3140-S | COL-3140 | 3.5 mm x 14.0 mm Locking Cortical Screw |
| COL-3160-S | COL-3160 | 3.5 mm x 16.0 mm Locking Cortical Screw |
| COL-3180-S | COL-3180 | 3.5 mm x 18.0 mm Locking Cortical Screw |
| CO-N2308-S | CO-N2308 | 2.3 mm x 8 mm Non-Toggling Cortical Screw |
| CO-N2310-S | CO-N2310 | 2.3 mm x 10 mm Non-Toggling Cortical Screw |
| CO-N2312-S | CO-N2312 | 2.3 mm x 12 mm Non-Toggling Cortical Screw |
| CO-N2314-S | CO-N2314 | 2.3 mm x 14 mm Non-Toggling Cortical Screw |
| CO-N2316-S | CO-N2316 | 2.3 mm x 16 mm Non-Toggling Cortical Screw |
| CO-N2318-S | CO-N2318 | 2.3 mm x 18 mm Non-Toggling Cortical Screw |
| CO-N2320-S | CO-N2320 | 2.3 mm x 20 mm Non-Toggling Cortical Screw |
| CO-N2322-S | CO-N2322 | 2.3 mm x 22 mm Non-Toggling Cortical Screw |
| CO-N2324-S | CO-N2324 | 2.3 mm x 24 mm Non-Toggling Cortical Screw |
| CO-N2326-S | CO-N2326 | 2.3 mm x 26 mm Non-Toggling Cortical Screw |
| CO-N2328-S | CO-N2328 | 2.3 mm x 28 mm Non-Toggling Cortical Screw |
| CO-N2330-S | CO-N2330 | 2.3 mm x 30 mm Non-Toggling Cortical Screw |
| CO-N2332-S | CO-N2332 | 2.3 mm x 32 mm Non-Toggling Cortical Screw |
| CO-S2308-S | CO-S2308 | 2.3 mm x 8 mm Locking Cortical Peg |
| CO-S2310-S | CO-S2310 | 2.3 mm x 10 mm Locking Cortical Peg |

Hand and Wrist Portfolio [continued]

Acu-Loc® Wrist Plating System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|--|
| CO-S2312-S | CO-S2312 | 2.3 mm x 12 mm Locking Cortical Peg |
| CO-S2314-S | CO-S2314 | 2.3 mm x 14 mm Locking Cortical Peg |
| CO-S2316-S | CO-S2316 | 2.3 mm x 16 mm Locking Cortical Peg |
| CO-S2318-S | CO-S2318 | 2.3 mm x 18 mm Locking Cortical Peg |
| CO-S2320-S | CO-S2320 | 2.3 mm x 20 mm Locking Cortical Peg |
| CO-S2322-S | CO-S2322 | 2.3 mm x 22 mm Locking Cortical Peg |
| CO-S2324-S | CO-S2324 | 2.3 mm x 24 mm Locking Cortical Peg |
| CO-S2326-S | CO-S2326 | 2.3 mm x 26 mm Locking Cortical Peg |
| CO-S2328-S | CO-S2328 | 2.3 mm x 28 mm Locking Cortical Peg |
| CO-T2308-S | CO-T2308 | 2.3 mm x 8 mm Locking Cortical Screw |
| CO-T2310-S | CO-T2310 | 2.3 mm x 10 mm Locking Cortical Screw |
| CO-T2312-S | CO-T2312 | 2.3 mm x 12 mm Locking Cortical Screw |
| CO-T2314-S | CO-T2314 | 2.3 mm x 14 mm Locking Cortical Screw |
| CO-T2316-S | CO-T2316 | 2.3 mm x 16 mm Locking Cortical Screw |
| CO-T2318-S | CO-T2318 | 2.3 mm x 18 mm Locking Cortical Screw |
| CO-T2320-S | CO-T2320 | 2.3 mm x 20 mm Locking Cortical Screw |
| CO-T2322-S | CO-T2322 | 2.3 mm x 22 mm Locking Cortical Screw |
| CO-T2324-S | CO-T2324 | 2.3 mm x 24 mm Locking Cortical Screw |
| CO-T2326-S | CO-T2326 | 2.3 mm x 26 mm Locking Cortical Screw |
| CO-T2328-S | CO-T2328 | 2.3 mm x 28 mm Locking Cortical Screw |
| PL-DR30L-S | PL-DR30L | Acu-Loc Volar Distal Radius Plate, Narrow, Left |
| PL-DR30R-S | PL-DR30R | Acu-Loc Volar Distal Radius Plate, Narrow, Right |
| PL-DR50L-S | PL-DR50L | Acu-Loc Volar Distal Radius Plate, Standard, Left |
| PL-DR50R-S | PL-DR50R | Acu-Loc Volar Distal Radius Plate, Standard, Right |
| PL-DR60L-S | PL-DR60L | Acu-Loc Volar Distal Radius Plate, Long, Left |
| PL-DR60R-S | PL-DR60R | Acu-Loc Volar Distal Radius Plate, Long, Right |
| PL-DR65L-S | PL-DR65L | Acu-Loc Volar Distal Radius Plate, X-Long, Left |
| PL-DR65R-S | PL-DR65R | Acu-Loc Volar Distal Radius Plate, X-Long, Right |
| PL-DR70L-S | PL-DR70L | Acu-Loc Volar Distal Radius Plate, Wide, Left |
| PL-DR70R-S | PL-DR70R | Acu-Loc Volar Distal Radius Plate, Wide, Right |

Hand and Wrist Portfolio [continued]

Acu-Loc® 2 Wrist Plating System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|--|
| 30-0093-S | 30-0093 | Acu-Loc 2 Volar Distal Radius Plate Extension Link Screw |
| 30-0100-S | 30-0100 | Acu-Loc 2 Volar Distal Radius Hexalobe Ext Link Screw |
| 30-0218-S | N/A | 3.5 mm x 9 mm Locking Hexalobe Screw |
| 30-0219-S | N/A | 3.5 mm x 11 mm Locking Hexalobe Screw |
| 30-0220-S | N/A | 3.5 mm x 13 mm Locking Hexalobe Screw |
| 30-0221-S | N/A | 3.5 mm x 15 mm Locking Hexalobe Screw |
| 30-0222-S | N/A | 3.5 mm x 17 mm Locking Hexalobe Screw |
| 30-0223-S | N/A | 3.5 mm x 19 mm Locking Hexalobe Screw |
| 30-0224-S | N/A | 3.5 mm x 9 mm Nonlocking Hexalobe Screw |
| 30-0225-S | N/A | 3.5 mm x 11 mm Nonlocking Hexalobe Screw |
| 30-0226-S | N/A | 3.5 mm x 13 mm Nonlocking Hexalobe Screw |
| 30-0227-S | N/A | 3.5 mm x 15 mm Nonlocking Hexalobe Screw |
| 30-0228-S | N/A | 3.5 mm x 17 mm Nonlocking Hexalobe Screw |
| 30-0229-S | N/A | 3.5 mm x 19 mm Nonlocking Hexalobe Screw |
| 30-0232-S | 30-0232 | 3.5 mm x 8 mm Locking Hexalobe Screw |
| 30-0233-S | 30-0233 | 3.5 mm x 10 mm Locking Hexalobe Screw |
| 30-0234-S | 30-0234 | 3.5 mm x 12 mm Locking Hexalobe Screw |
| 30-0235-S | 30-0235 | 3.5 mm x 14 mm Locking Hexalobe Screw |
| 30-0236-S | 30-0236 | 3.5 mm x 16 mm Locking Hexalobe Screw |
| 30-0237-S | 30-0237 | 3.5 mm x 18 mm Locking Hexalobe Screw |
| 30-0256-S | 30-0256 | 3.5 mm x 10 mm Nonlocking Hexalobe Screw |
| 30-0257-S | 30-0257 | 3.5 mm x 12 mm Nonlocking Hexalobe Screw |
| 30-0258-S | 30-0258 | 3.5 mm x 14 mm Nonlocking Hexalobe Screw |
| 30-0259-S | 30-0259 | 3.5 mm x 16 mm Nonlocking Hexalobe Screw |
| 30-0260-S | 30-0260 | 3.5 mm x 18 mm Nonlocking Hexalobe Screw |
| 30-0370-S | 30-0370 | Frag-Loc® Compression Sleeve |
| 30-0371-S | 30-0371 | Frag-Loc Compression Screw |
| 30-0372-S | 30-0372 | Frag-Loc Compression Screw, Long |
| 30-2314-S | 30-2314 | 2.3 mm x 14 mm Locking Variable Angle Screw |

Hand and Wrist Portfolio [continued]

Acu-Loc® 2 Wrist Plating System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|---|
| 30-2316-S | 30-2316 | 2.3 mm x 16 mm Locking Variable Angle Screw |
| 30-2318-S | 30-2318 | 2.3 mm x 18 mm Locking Variable Angle Screw |
| 30-2320-S | 30-2320 | 2.3 mm x 20 mm Locking Variable Angle Screw |
| 30-2322-S | 30-2322 | 2.3 mm x 22 mm Locking Variable Angle Screw |
| 30-2324-S | 30-2324 | 2.3 mm x 24 mm Locking Variable Angle Screw |
| 30-2326-S | 30-2326 | 2.3 mm x 26 mm Locking Variable Angle Screw |
| 30-2328-S | 30-2328 | 2.3 mm x 28 mm Locking Variable Angle Screw |
| N/A | 3040-23008 | 2.7 x 8 mm Locking Low-Profile Hexalobe Screw |
| N/A | 3040-23009 | 2.7 x 9 mm Locking Low-Profile Hexalobe Screw |
| N/A | 3040-23010 | 2.7 x 10 mm Locking Low-Profile Hexalobe Screw |
| N/A | 3040-23011 | 2.7 x 11 mm Locking Low-Profile Hexalobe Screw |
| N/A | 3040-23012 | 2.7 x 12 mm Locking Low-Profile Hexalobe Screw |
| N/A | 3040-23013 | 2.7 x 13 mm Locking Low-Profile Hexalobe Screw |
| N/A | 3040-23014 | 2.7 x 14 mm Locking Low-Profile Hexalobe Screw |
| N/A | 3040-23015 | 2.7 x 15 mm Locking Low-Profile Hexalobe Screw |
| N/A | 3040-23016 | 2.7 x 16 mm Locking Low-Profile Hexalobe Screw |
| N/A | 3040-23017 | 2.7 x 17 mm Locking Low-Profile Hexalobe Screw |
| N/A | 3040-23018 | 2.7 x 18 mm Locking Low-Profile Hexalobe Screw |
| N/A | 3040-23019 | 2.7 x 19 mm Locking Low-Profile Hexalobe Screw |
| N/A | 3041-23009 | 2.7 x 9 mm Nonlocking Low-Profile Hexalobe Screw |
| N/A | 3041-23010 | 2.7 x 10 mm Nonlocking Low-Profile Hexalobe Screw |
| N/A | 3041-23011 | 2.7 x 11 mm Nonlocking Low-Profile Hexalobe Screw |
| N/A | 3041-23012 | 2.7 x 12 mm Nonlocking Low-Profile Hexalobe Screw |
| N/A | 3041-23013 | 2.7 x 13 mm Nonlocking Low-Profile Hexalobe Screw |
| N/A | 3041-23014 | 2.7 x 14 mm Nonlocking Low-Profile Hexalobe Screw |
| N/A | 3041-23015 | 2.7 x 15 mm Nonlocking Low-Profile Hexalobe Screw |
| N/A | 3041-23016 | 2.7 x 16 mm Nonlocking Low-Profile Hexalobe Screw |
| N/A | 3041-23017 | 2.7 x 17 mm Nonlocking Low-Profile Hexalobe Screw |
| N/A | 3041-23018 | 2.7 x 18 mm Nonlocking Low-Profile Hexalobe Screw |

Hand and Wrist Portfolio [continued]

Acu-Loc® 2 Wrist Plating System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|---|
| N/A | 3041-23019 | 2.7 x 19 mm Nonlocking Low-Profile Hexalobe Screw |
| 70-0045-S | 70-0045 | Acu-Loc Volar Distal Ulna Plate Left Standard |
| 70-0046-S | 70-0046 | Acu-Loc Volar Distal Ulna Plate Right Standard |
| 70-0047-S | 70-0047 | Acu-Loc Volar Distal Ulna Plate Left Long |
| 70-0048-S | 70-0048 | Acu-Loc Volar Distal Ulna Plate Right Long |
| 70-0055-S | 70-0055 | Acu-Loc Dorsal Plate Standard, Left |
| 70-0056-S | 70-0056 | Acu-Loc Dorsal Plate Standard, Right |
| 70-0057-S | 70-0057 | Acu-Loc Dorsal Plate Narrow Left |
| 70-0058-S | 70-0058 | Acu-Loc Dorsal Plate Narrow Right |
| 70-0063-S | 70-0063 | Acu-Loc Extra Articular Standard |
| 70-0064-S | 70-0064 | Acu-Loc Extra Articular Narrow |
| 70-0331-S | 70-0331 | Divergent Radial Styloid Plate |
| 70-0334-S | 70-0334 | Volar Lunate Suture Plate |
| 70-0335-S | 70-0335 | Dorsal Rim Buttress Plate, Right |
| 70-0336-S | 70-0336 | Dorsal Rim Buttress Plate, Left |
| 70-0337-S | 70-0337 | Dorsal Lunate Plate, Right |
| 70-0338-S | 70-0338 | Dorsal Lunate Plate, Left |
| 70-0350-S | 70-0350 | Acu-Loc 2 Volar Distal Radius Proximal Plate, Standard, Left |
| 70-0351-S | 70-0351 | Acu-Loc 2 Volar Distal Radius Proximal Plate, Standard, Right |
| 70-0352-S | 70-0352 | Acu-Loc 2 Volar Distal Radius Proximal Plate, Narrow, Left |
| 70-0353-S | 70-0353 | Acu-Loc 2 Volar Distal Radius Proximal Plate, Narrow, Right |
| 70-0354-S | 70-0354 | Acu-Loc 2 Volar Distal Radius Proximal Plate, Wide, Left |
| 70-0355-S | 70-0355 | Acu-Loc 2 Volar Distal Radius Proximal Plate, Wide, Right |
| 70-0356-S | 70-0356 | Acu-Loc 2 Volar Distal Radius Plate Standard, Left |
| 70-0357-S | 70-0357 | Acu-Loc 2 Volar Distal Radius Plate Standard, Right |
| 70-0358-S | 70-0358 | Acu-Loc 2 Volar Distal Radius Plate, Narrow, Left |
| 70-0359-S | 70-0359 | Acu-Loc 2 Volar Distal Radius Plate, Narrow, Right |
| 70-0360-S | 70-0360 | Acu-Loc 2 Volar Distal Radius Plate, Wide, Left |
| 70-0361-S | 70-0361 | Acu-Loc 2 Volar Distal Radius Plate, Wide, Right |

Hand and Wrist Portfolio [continued]

Acu-Loc® 2 Wrist Plating System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|--|
| 70-0364-S | 70-0364 | Acu-Loc 2 Volar Distal Radius Extension Plate, Neutral |
| 70-0365-S | 70-0365 | Acu-Loc 2 Volar Distal Radius Extension Plate, Long, Left |
| 70-0366-S | 70-0366 | Acu-Loc 2 Volar Distal Radius Extension Plate, Long, Right |
| 70-0368-S | 70-0368 | Acu-Loc 2 Volar Distal Radius Plate Standard Long, Left |
| 70-0369-S | 70-0369 | Acu-Loc 2 Volar Distal Radius Plate Standard Long, Right |
| 70-0370-S | 70-0370 | Acu-Loc 2 Volar Distal Radius Plate, Narrow Long, Left |
| 70-0371-S | 70-0371 | Acu-Loc 2 Volar Distal Radius Plate, Narrow Long, Right |
| 70-0372-S | 70-0372 | Acu-Loc 2 Volar Distal Radius Proximal Plate, Standard Long, Left |
| 70-0373-S | 70-0373 | Acu-Loc 2 Volar Distal Radius Proximal Plate, Standard Long, Right |
| 70-0382-S | 70-0382 | Acu-Loc 2 Volar Distal Radius Proximal Plate, Narrow Long, Left |
| 70-0383-S | 70-0383 | Acu-Loc 2 Volar Distal Radius Proximal Plate, Narrow Long, Right |
| 7006-3370N-S | N/A | Acu-Loc 2 Narrow Wrist Spanning Plate, Short |
| 7006-3390N-S | N/A | Acu-Loc 2 Narrow Wrist Spanning Plate, Long |
| 7005-08001-S | 7005-08001 | 0.8 mm Avulsion Hook Plate |
| CO-3100-S | CO-3100 | 3.5 mm x 10.0 mm Cortical Screw |
| CO-3120-S | CO-3120 | 3.5 mm x 12.0 mm Cortical Screw |
| CO-3140-S | CO-3140 | 3.5 mm x 14.0 mm Cortical Screw |
| CO-3160-S | CO-3160 | 3.5 mm x 16.0 mm Cortical Screw |
| CO-3180-S | CO-3180 | 3.5 mm x 18.0 mm Cortical Screw |
| COL-3080-S | COL-3080 | 3.5 mm x 8.0 mm Locking Cortical Screw |
| COL-3100-S | COL-3100 | 3.5 mm x 10.0 mm Locking Cortical Screw |
| COL-3120-S | COL-3120 | 3.5 mm x 12.0 mm Locking Cortical Screw |
| COL-3140-S | COL-3140 | 3.5 mm x 14.0 mm Locking Cortical Screw |
| COL-3160-S | COL-3160 | 3.5 mm x 16.0 mm Locking Cortical Screw |
| COL-3180-S | COL-3180 | 3.5 mm x 18.0 mm Locking Cortical Screw |
| CO-N2308-S | CO-N2308 | 2.3 mm x 8 mm Non-Toggling Cortical Screw |
| CO-N2310-S | CO-N2310 | 2.3 mm x 10 mm Non-Toggling Cortical Screw |
| CO-N2312-S | CO-N2312 | 2.3 mm x 12 mm Non-Toggling Cortical Screw |
| CO-N2314-S | CO-N2314 | 2.3 mm x 14 mm Non-Toggling Cortical Screw |

Hand and Wrist Portfolio [continued]

Acu-Loc® 2 Wrist Plating System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|--|
| CO-N2316-S | CO-N2316 | 2.3 mm x 16 mm Non-Toggling Cortical Screw |
| CO-N2318-S | CO-N2318 | 2.3 mm x 18 mm Non-Toggling Cortical Screw |
| CO-N2320-S | CO-N2320 | 2.3 mm x 20 mm Non-Toggling Cortical Screw |
| CO-N2322-S | CO-N2322 | 2.3 mm x 22 mm Non-Toggling Cortical Screw |
| CO-N2324-S | CO-N2324 | 2.3 mm x 24 mm Non-Toggling Cortical Screw |
| CO-N2326-S | CO-N2326 | 2.3 mm x 26 mm Non-Toggling Cortical Screw |
| CO-N2328-S | CO-N2328 | 2.3 mm x 28 mm Non-Toggling Cortical Screw |
| CO-N2330-S | CO-N2330 | 2.3 mm x 30 mm Non-Toggling Cortical Screw |
| CO-N2332-S | CO-N2332 | 2.3 mm x 32 mm Non-Toggling Cortical Screw |
| CO-N2334-S | CO-N2334 | 2.3 mm x 34 mm Non-Toggling Cortical Screw |
| CO-N2336-S | CO-N2336 | 2.3 mm x 36 mm Non-Toggling Cortical Screw |
| CO-N2338-S | CO-N2338 | 2.3 mm x 38 mm Non-Toggling Cortical Screw |
| CO-N2340-S | CO-N2340 | 2.3 mm x 40 mm Non-Toggling Cortical Screw |
| CO-N2342-S | CO-N2342 | 2.3 mm x 42 mm Non-Toggling Cortical Screw |
| CO-N2344-S | CO-N2344 | 2.3 mm x 44 mm Non-Toggling Cortical Screw |
| CO-N2346-S | CO-N2346 | 2.3 mm x 46 mm Non-Toggling Cortical Screw |
| CO-S2308-S | CO-S2308 | 2.3 mm x 8 mm Locking Cortical Peg |
| CO-S2310-S | CO-S2310 | 2.3 mm x 10 mm Locking Cortical Peg |
| CO-S2312-S | CO-S2312 | 2.3 mm x 12 mm Locking Cortical Peg |
| CO-S2314-S | CO-S2314 | 2.3 mm x 14 mm Locking Cortical Peg |
| CO-S2316-S | CO-S2316 | 2.3 mm x 16 mm Locking Cortical Peg |
| CO-S2318-S | CO-S2318 | 2.3 mm x 18 mm Locking Cortical Peg |
| CO-S2320-S | CO-S2320 | 2.3 mm x 20 mm Locking Cortical Peg |
| CO-S2322-S | CO-S2322 | 2.3 mm x 22 mm Locking Cortical Peg |
| CO-S2324-S | CO-S2324 | 2.3 mm x 24 mm Locking Cortical Peg |
| CO-S2326-S | CO-S2326 | 2.3 mm x 26 mm Locking Cortical Peg |
| CO-S2328-S | CO-S2328 | 2.3 mm x 28 mm Locking Cortical Peg |
| CO-T2308-S | CO-T2308 | 2.3 mm x 8 mm Locking Cortical Screw |
| CO-T2310-S | CO-T2310 | 2.3 mm x 10 mm Locking Cortical Screw |
| CO-T2312-S | CO-T2312 | 2.3 mm x 12 mm Locking Cortical Screw |

Hand and Wrist Portfolio [continued]

Acu-Loc® 2 Wrist Plating System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|---------------------------------------|
| CO-T2314-S | CO-T2314 | 2.3 mm x 14 mm Locking Cortical Screw |
| CO-T2316-S | CO-T2316 | 2.3 mm x 16 mm Locking Cortical Screw |
| CO-T2318-S | CO-T2318 | 2.3 mm x 18 mm Locking Cortical Screw |
| CO-T2320-S | CO-T2320 | 2.3 mm x 20 mm Locking Cortical Screw |
| CO-T2322-S | CO-T2322 | 2.3 mm x 22 mm Locking Cortical Screw |
| CO-T2324-S | CO-T2324 | 2.3 mm x 24 mm Locking Cortical Screw |
| CO-T2326-S | CO-T2326 | 2.3 mm x 26 mm Locking Cortical Screw |
| CO-T2328-S | CO-T2328 | 2.3 mm x 28 mm Locking Cortical Screw |
| CO-T2330-S | CO-T2330 | 2.3 mm x 30 mm Locking Cortical Screw |
| CO-T2332-S | CO-T2332 | 2.3 mm x 32 mm Locking Cortical Screw |
| CO-T2334-S | CO-T2334 | 2.3 mm x 34 mm Locking Cortical Screw |
| CO-T2336-S | CO-T2336 | 2.3 mm x 36 mm Locking Cortical Screw |
| CO-T2338-S | CO-T2338 | 2.3 mm x 38 mm Locking Cortical Screw |
| CO-T2340-S | CO-T2340 | 2.3 mm x 40 mm Locking Cortical Screw |
| CO-T2342-S | CO-T2342 | 2.3 mm x 42 mm Locking Cortical Screw |
| CO-T2344-S | CO-T2344 | 2.3 mm x 44 mm Locking Cortical Screw |
| CO-T2346-S | CO-T2346 | 2.3 mm x 46 mm Locking Cortical Screw |

Acu-Loc® Wrist Spanning Plate

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|--|
| 30-0232-S | 30-0232 | 3.5 mm x 8 mm Locking Hexalobe Screw |
| 30-0233-S | 30-0233 | 3.5 mm x 10 mm Locking Hexalobe Screw |
| 30-0234-S | 30-0234 | 3.5 mm x 12 mm Locking Hexalobe Screw |
| 30-0235-S | 30-0235 | 3.5 mm x 14 mm Locking Hexalobe Screw |
| 30-0236-S | 30-0236 | 3.5 mm x 16 mm Locking Hexalobe Screw |
| 30-0237-S | 30-0237 | 3.5 mm x 18 mm Locking Hexalobe Screw |
| 30-0255-S | 30-0255 | 3.5 mm x 8 mm Nonlocking Hexalobe Screw |
| 30-0256-S | 30-0256 | 3.5 mm x 10 mm Nonlocking Hexalobe Screw |
| 30-0257-S | 30-0257 | 3.5 mm x 12 mm Nonlocking Hexalobe Screw |

Hand and Wrist Portfolio [continued]

Acu-Loc® Wrist Spanning Plate

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|--|
| 30-0258-S | 30-0258 | 3.5 mm x 14 mm Nonlocking Hexalobe Screw |
| 30-0259-S | 30-0259 | 3.5 mm x 16 mm Nonlocking Hexalobe Screw |
| 30-0260-S | 30-0260 | 3.5 mm x 18 mm Nonlocking Hexalobe Screw |
| 30-0324-S | 30-0324 | 2.7 mm x 8 mm Locking Hexalobe Screw |
| 30-0325-S | 30-0325 | 2.7 mm x 10 mm Locking Hexalobe Screw |
| 30-0326-S | 30-0326 | 2.7 mm x 12 mm Locking Hexalobe Screw |
| 30-0327-S | 30-0327 | 2.7 mm x 14 mm Locking Hexalobe Screw |
| 30-0328-S | 30-0328 | 2.7 mm x 16 mm Locking Hexalobe Screw |
| 30-0329-S | 30-0329 | 2.7 mm x 18 mm Locking Hexalobe Screw |
| 30-0330-S | 30-0330 | 2.7 mm x 20 mm Locking Hexalobe Screw |
| 30-0331-S | 30-0331 | 2.7 mm x 22 mm Locking Hexalobe Screw |
| 30-0343-S | 30-0343 | 2.7 mm x 8 mm Nonlocking Hexalobe Screw |
| 30-0344-S | 30-0344 | 2.7 mm x 10 mm Nonlocking Hexalobe Screw |
| 30-0345-S | 30-0345 | 2.7 mm x 12 mm Nonlocking Hexalobe Screw |
| 30-0346-S | 30-0346 | 2.7 mm x 14 mm Nonlocking Hexalobe Screw |
| 30-0347-S | 30-0347 | 2.7 mm x 16 mm Nonlocking Hexalobe Screw |
| 30-0348-S | 30-0348 | 2.7 mm x 18 mm Nonlocking Hexalobe Screw |
| 30-0349-S | 30-0349 | 2.7 mm x 20 mm Nonlocking Hexalobe Screw |
| 30-0350-S | 30-0350 | 2.7 mm x 22 mm Nonlocking Hexalobe Screw |
| 7006-1170N-S | N/A | Acu-Loc Wrist Spanning Plate, Short |
| 7006-1190N-S | N/A | Acu-Loc Wrist Spanning Plate, Long |
| N/A | 3060-27008 | 2.7 mm x 8 mm Hexalobe Locking CoCr Screw |
| N/A | 3060-27010 | 2.7 mm x 10 mm Hexalobe Locking CoCr Screw |
| N/A | 3060-27012 | 2.7 mm x 12 mm Hexalobe Locking CoCr Screw |
| N/A | 3060-27014 | 2.7 mm x 14 mm Hexalobe Locking CoCr Screw |
| N/A | 3060-27016 | 2.7 mm x 16 mm Hexalobe Locking CoCr Screw |
| N/A | 3060-27018 | 2.7 mm x 18 mm Hexalobe Locking CoCr Screw |
| N/A | 3060-27020 | 2.7 mm x 20 mm Hexalobe Locking CoCr Screw |
| N/A | 3060-27022 | 2.7 mm x 22 mm Hexalobe Locking CoCr Screw |
| N/A | 3060-27024 | 2.7 mm x 24 mm Hexalobe Locking CoCr Screw |

Hand and Wrist Portfolio [continued]

Acu-Loc® Wrist Spanning Plate

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|--|
| N/A | 3060-27026 | 2.7 mm x 26 mm Hexalobe Locking CoCr Screw |
| N/A | 3060-27028 | 2.7 mm x 28 mm Hexalobe Locking CoCr Screw |
| N/A | 3060-27030 | 2.7 mm x 30 mm Hexalobe Locking CoCr Screw |
| N/A | 3060-27032 | 2.7 mm x 32 mm Hexalobe Locking CoCr Screw |
| N/A | 3061-27008 | 2.7 mm x 8 mm Hexalobe Non-Locking CoCr Screw |
| N/A | 3061-27010 | 2.7 mm x 10 mm Hexalobe Non-Locking CoCr Screw |
| N/A | 3061-27012 | 2.7 mm x 12 mm Hexalobe Non-Locking CoCr Screw |
| N/A | 3061-27014 | 2.7 mm x 14 mm Hexalobe Non-Locking CoCr Screw |
| N/A | 3061-27016 | 2.7 mm x 16 mm Hexalobe Non-Locking CoCr Screw |
| N/A | 3061-27018 | 2.7 mm x 18 mm Hexalobe Non-Locking CoCr Screw |
| N/A | 3061-27020 | 2.7 mm x 20 mm Hexalobe Non-Locking CoCr Screw |
| N/A | 3061-27022 | 2.7 mm x 22 mm Hexalobe Non-Locking CoCr Screw |
| N/A | 3061-27024 | 2.7 mm x 24 mm Hexalobe Non-Locking CoCr Screw |
| N/A | 3061-27026 | 2.7 mm x 26 mm Hexalobe Non-Locking CoCr Screw |
| N/A | 3061-27028 | 2.7 mm x 28 mm Hexalobe Non-Locking CoCr Screw |
| N/A | 3061-27030 | 2.7 mm x 30 mm Hexalobe Non-Locking CoCr Screw |
| N/A | 3061-27032 | 2.7 mm x 32 mm Hexalobe Non-Locking CoCr Screw |
| N/A | 3062-35008 | 3.5 mm x 8 mm Hexalobe Locking CoCr Screw |
| N/A | 3062-35010 | 3.5 mm x 10 mm Hexalobe Locking CoCr Screw |
| N/A | 3062-35012 | 3.5 mm x 12 mm Hexalobe Locking CoCr Screw |
| N/A | 3062-35014 | 3.5 mm x 14 mm Hexalobe Locking CoCr Screw |
| N/A | 3062-35016 | 3.5 mm x 16 mm Hexalobe Locking CoCr Screw |
| N/A | 3062-35018 | 3.5 mm x 18 mm Hexalobe Locking CoCr Screw |
| N/A | 3062-35020 | 3.5 mm x 20 mm Hexalobe Locking CoCr Screw |
| N/A | 3062-35022 | 3.5 mm x 22 mm Hexalobe Locking CoCr Screw |
| N/A | 3062-35024 | 3.5 mm x 24 mm Hexalobe Locking CoCr Screw |
| N/A | 3062-35026 | 3.5 mm x 26 mm Hexalobe Locking CoCr Screw |
| N/A | 3062-35028 | 3.5 mm x 28 mm Hexalobe Locking CoCr Screw |
| N/A | 3062-35030 | 3.5 mm x 30 mm Hexalobe Locking CoCr Screw |
| N/A | 3062-35032 | 3.5 mm x 32 mm Hexalobe Locking CoCr Screw |

Hand and Wrist Portfolio [continued]

Acu-Loc® Wrist Spanning Plate

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|--|
| N/A | 3062-35034 | 3.5 mm x 34 mm Hexalobe Locking CoCr Screw |
| N/A | 3062-35036 | 3.5 mm x 36 mm Hexalobe Locking CoCr Screw |
| N/A | 3062-35038 | 3.5 mm x 38 mm Hexalobe Locking CoCr Screw |
| N/A | 3062-35040 | 3.5 mm x 40 mm Hexalobe Locking CoCr Screw |
| N/A | 3062-35045 | 3.5 mm x 45 mm Hexalobe Locking CoCr Screw |
| N/A | 3062-35050 | 3.5 mm x 50 mm Hexalobe Locking CoCr Screw |
| N/A | 3062-35055 | 3.5 mm x 55 mm Hexalobe Locking CoCr Screw |
| N/A | 3062-35060 | 3.5 mm x 60 mm Hexalobe Locking CoCr Screw |
| N/A | 3063-35008 | 3.5 mm x 8 mm Hexalobe Non-Locking CoCr Screw |
| N/A | 3063-35010 | 3.5 mm x 10 mm Hexalobe Non-Locking CoCr Screw |
| N/A | 3063-35012 | 3.5 mm x 12 mm Hexalobe Non-Locking CoCr Screw |
| N/A | 3063-35014 | 3.5 mm x 14 mm Hexalobe Non-Locking CoCr Screw |
| N/A | 3063-35016 | 3.5 mm x 16 mm Hexalobe Non-Locking CoCr Screw |
| N/A | 3063-35018 | 3.5 mm x 18 mm Hexalobe Non-Locking CoCr Screw |
| N/A | 3063-35020 | 3.5 mm x 20 mm Hexalobe Non-Locking CoCr Screw |
| N/A | 3063-35022 | 3.5 mm x 22 mm Hexalobe Non-Locking CoCr Screw |
| N/A | 3063-35024 | 3.5 mm x 24 mm Hexalobe Non-Locking CoCr Screw |
| N/A | 3063-35026 | 3.5 mm x 26 mm Hexalobe Non-Locking CoCr Screw |
| N/A | 3063-35028 | 3.5 mm x 28 mm Hexalobe Non-Locking CoCr Screw |
| N/A | 3063-35030 | 3.5 mm x 30 mm Hexalobe Non-Locking CoCr Screw |
| N/A | 3063-35032 | 3.5 mm x 32 mm Hexalobe Non-Locking CoCr Screw |
| N/A | 3063-35034 | 3.5 mm x 34 mm Hexalobe Non-Locking CoCr Screw |
| N/A | 3063-35036 | 3.5 mm x 36 mm Hexalobe Non-Locking CoCr Screw |
| N/A | 3063-35038 | 3.5 mm x 38 mm Hexalobe Non-Locking CoCr Screw |
| N/A | 3063-35040 | 3.5 mm x 40 mm Hexalobe Non-Locking CoCr Screw |
| N/A | 3063-35045 | 3.5 mm x 45 mm Hexalobe Non-Locking CoCr Screw |
| N/A | 3063-35050 | 3.5 mm x 50 mm Hexalobe Non-Locking CoCr Screw |
| N/A | 3063-35055 | 3.5 mm x 55 mm Hexalobe Non-Locking CoCr Screw |
| N/A | 3063-35060 | 3.5 mm x 60 mm Hexalobe Non-Locking CoCr Screw |
| N/A | 3063-35065 | 3.5 mm x 65 mm Hexalobe Non-Locking CoCr Screw |

Hand and Wrist Portfolio [continued]

Forearm Fracture Solutions (plates)

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|--|
| 30-0232-S | 30-0232 | 3.5 mm x 8 mm Locking Hexalobe Screw |
| 30-0233-S | 30-0233 | 3.5 mm x 10 mm Locking Hexalobe Screw |
| 30-0234-S | 30-0234 | 3.5 mm x 12 mm Locking Hexalobe Screw |
| 30-0235-S | 30-0235 | 3.5 mm x 14 mm Locking Hexalobe Screw |
| 30-0236-S | 30-0236 | 3.5 mm x 16 mm Locking Hexalobe Screw |
| 30-0237-S | 30-0237 | 3.5 mm x 18 mm Locking Hexalobe Screw |
| 30-0238-S | 30-0238 | 3.5 mm x 20 mm Locking Hexalobe Screw |
| 30-0239-S | 30-0239 | 3.5 mm x 22 mm Locking Hexalobe Screw |
| 30-0256-S | 30-0256 | 3.5 mm x 10 mm Nonlocking Hexalobe Screw |
| 30-0257-S | 30-0257 | 3.5 mm x 12 mm Nonlocking Hexalobe Screw |
| 30-0258-S | 30-0258 | 3.5 mm x 14 mm Nonlocking Hexalobe Screw |
| 30-0259-S | 30-0259 | 3.5 mm x 16 mm Nonlocking Hexalobe Screw |
| 30-0260-S | 30-0260 | 3.5 mm x 18 mm Nonlocking Hexalobe Screw |
| 30-0261-S | 30-0261 | 3.5 mm x 20 mm Nonlocking Hexalobe Screw |
| 30-0262-S | 30-0262 | 3.5 mm x 22 mm Nonlocking Hexalobe Screw |
| 30-0302-S | 30-0302 | 3.0 mm x 10 mm Nonlocking Hexalobe Screw |
| 30-0303-S | 30-0303 | 3.0 mm x 12 mm Nonlocking Hexalobe Screw |
| 30-0304-S | 30-0304 | 3.0 mm x 14 mm Nonlocking Hexalobe Screw |
| 30-0305-S | 30-0305 | 3.0 mm x 16 mm Nonlocking Hexalobe Screw |
| 30-0306-S | 30-0306 | 3.0 mm x 18 mm Nonlocking Hexalobe Screw |
| 30-0307-S | 30-0307 | 3.0 mm x 20 mm Nonlocking Hexalobe Screw |
| 30-0308-S | 30-0308 | 3.0 mm x 22 mm Nonlocking Hexalobe Screw |
| 70-0066-S | 70-0066 | 6-Hole Volar Midshaft Radius Plate |
| 70-0067-S | 70-0067 | 8-Hole Volar Midshaft Radius Plate |
| 70-0068-S | 70-0068 | 10-Hole Volar Midshaft Radius Plate |
| 70-0069-S | 70-0069 | 12-Hole Volar Midshaft Radius Plate |
| 70-0070-S | 70-0070 | 6-Hole Midshaft Ulna Plate |
| 70-0071-S | 70-0071 | 8-Hole Midshaft Ulna Plate |
| 70-0072-S | 70-0072 | 10-Hole Midshaft Ulna Plate |

Hand and Wrist Portfolio [continued]

Forearm Fracture Solutions (plates)

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|--|
| 70-0073-S | 70-0073 | 12-Hole Midshaft Ulna Plate |
| 70-0074-S | 70-0074 | 6-Hole Dorsolateral Midshaft Radius Plate |
| 70-0075-S | 70-0075 | 8-Hole Dorsolateral Midshaft Radius Plate |
| 70-0076-S | 70-0076 | 10-Hole Dorsolateral Midshaft Radius Plate |
| 70-0077-S | 70-0077 | 12-Hole Dorsolateral Midshaft Radius Plate |
| 70-0463-S | N/A | 14-Hole Midshaft Ulna Plate |
| 70-0464-S | N/A | 16-Hole Midshaft Ulna Plate |
| 70-0466-S | N/A | 14-Hole Midshaft Dorsolateral Radius Plate |
| 70-0467-S | N/A | 16-Hole Midshaft Dorsolateral Radius Plate |
| 70-0469-S | N/A | 14-Hole Volar Midshaft Radius Plate |
| 70-0470-S | N/A | 16-Hole Volar Midshaft Radius Plate |
| CO-3125-S | N/A | 3.5 mm x 12.5 mm Cortical Screw |
| CO-3150-S | N/A | 3.5 mm x 15.0 mm Cortical Screw |
| CO-3175-S | N/A | 3.5 mm x 17.5 mm Cortical Screw |
| CO-3200-S | CO-3200 | 3.5 mm x 20.0 mm Cortical Screw |
| CO-3225-S | CO-3225 | 3.5 mm x 22.5 mm Cortical Screw |
| CO-3250-S | CO-3250 | 3.5 mm x 25.0 mm Cortical Screw |
| CO-3275-S | CO-3275 | 3.5 mm x 27.5 mm Cortical Screw |
| CO-3300-S | CO-3300 | 3.5 mm x 30.0 mm Cortical Screw |
| CO-3325-S | CO-3325 | 3.5 mm x 32.5 mm Cortical Screw |
| N/A | CO-T2326 | 2.3 mm x 26 mm Locking Cortical Screw |
| N/A | CO-T2328 | 2.3 mm x 28 mm Locking Cortical Screw |

Hand Fracture System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|-----------------------------------|
| 3004-15005-S | 3004-15005 | 1.5 mm x 5 mm Hexalobe MultiScrew |
| 3004-15006-S | 3004-15006 | 1.5 mm x 6 mm Hexalobe MultiScrew |
| 3004-15007-S | 3004-15007 | 1.5 mm x 7 mm Hexalobe MultiScrew |
| 3004-15008-S | 3004-15008 | 1.5 mm x 8 mm Hexalobe MultiScrew |

Hand and Wrist Portfolio [continued]

Hand Fracture System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|------------------------------------|
| 3004-15009-S | 3004-15009 | 1.5 mm x 9 mm Hexalobe MultiScrew |
| 3004-15010-S | 3004-15010 | 1.5 mm x 10 mm Hexalobe MultiScrew |
| 3004-15011-S | 3004-15011 | 1.5 mm x 11 mm Hexalobe MultiScrew |
| 3004-15012-S | 3004-15012 | 1.5 mm x 12 mm Hexalobe MultiScrew |
| 3004-15013-S | 3004-15013 | 1.5 mm x 13 mm Hexalobe MultiScrew |
| 3004-15014-S | 3004-15014 | 1.5 mm x 14 mm Hexalobe MultiScrew |
| 3004-15016-S | 3004-15016 | 1.5 mm x 16 mm Hexalobe MultiScrew |
| 3004-15018-S | 3004-15018 | 1.5 mm x 18 mm Hexalobe MultiScrew |
| 3004-15020-S | 3004-15020 | 1.5 mm x 20 mm Hexalobe MultiScrew |
| 3004-23005-S | 3004-23005 | 2.3 mm x 5 mm Hexalobe MultiScrew |
| 3004-23006-S | 3004-23006 | 2.3 mm x 6 mm Hexalobe MultiScrew |
| 3004-23007-S | 3004-23007 | 2.3 mm x 7 mm Hexalobe MultiScrew |
| 3004-23008-S | 3004-23008 | 2.3 mm x 8 mm Hexalobe MultiScrew |
| 3004-23009-S | 3004-23009 | 2.3 mm x 9 mm Hexalobe MultiScrew |
| 3004-23010-S | 3004-23010 | 2.3 mm x 10 mm Hexalobe MultiScrew |
| 3004-23011-S | 3004-23011 | 2.3 mm x 11 mm Hexalobe MultiScrew |
| 3004-23012-S | 3004-23012 | 2.3 mm x 12 mm Hexalobe MultiScrew |
| 3004-23013-S | 3004-23013 | 2.3 mm x 13 mm Hexalobe MultiScrew |
| 3004-23014-S | 3004-23014 | 2.3 mm x 14 mm Hexalobe MultiScrew |
| 3004-23016-S | 3004-23016 | 2.3 mm x 16 mm Hexalobe MultiScrew |
| 3004-23018-S | 3004-23018 | 2.3 mm x 18 mm Hexalobe MultiScrew |
| 3004-23020-S | 3004-23020 | 2.3 mm x 20 mm Hexalobe MultiScrew |
| 3012-15005-S | 3012-15005 | 1.5 mm x 5 mm Hexalobe Lag Screw |
| 3012-15006-S | 3012-15006 | 1.5 mm x 6 mm Hexalobe Lag Screw |
| 3012-15007-S | 3012-15007 | 1.5 mm x 7 mm Hexalobe Lag Screw |
| 3012-15008-S | 3012-15008 | 1.5 mm x 8 mm Hexalobe Lag Screw |
| 3012-15009-S | 3012-15009 | 1.5 mm x 9 mm Hexalobe Lag Screw |
| 3012-15010-S | 3012-15010 | 1.5 mm x 10 mm Hexalobe Lag Screw |
| 3012-15011-S | 3012-15011 | 1.5 mm x 11 mm Hexalobe Lag Screw |

Hand and Wrist Portfolio [continued]

Hand Fracture System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|------------------------------------|
| 3012-15012-S | 3012-15012 | 1.5 mm x 12 mm Hexalobe Lag Screw |
| 3012-15013-S | 3012-15013 | 1.5 mm x 13 mm Hexalobe Lag Screw |
| 3012-15014-S | 3012-15014 | 1.5 mm x 14 mm Hexalobe Lag Screw |
| 3012-15016-S | 3012-15016 | 1.5 mm x 16 mm Hexalobe Lag Screw |
| 3012-15018-S | 3012-15018 | 1.5 mm x 18 mm Hexalobe Lag Screw |
| 3012-15020-S | 3012-15020 | 1.5 mm x 20 mm Hexalobe Lag Screw |
| 3012-23005-S | 3012-23005 | 2.3 mm x 5 mm Hexalobe Lag Screw |
| 3012-23006-S | 3012-23006 | 2.3 mm x 6 mm Hexalobe Lag Screw |
| 3012-23007-S | 3012-23007 | 2.3 mm x 7 mm Hexalobe Lag Screw |
| 3012-23008-S | 3012-23008 | 2.3 mm x 8 mm Hexalobe Lag Screw |
| 3012-23009-S | 3012-23009 | 2.3 mm x 9 mm Hexalobe Lag Screw |
| 3012-23010-S | 3012-23010 | 2.3 mm x 10 mm Hexalobe Lag Screw |
| 3012-23011-S | 3012-23011 | 2.3 mm x 11 mm Hexalobe Lag Screw |
| 3012-23012-S | 3012-23012 | 2.3 mm x 12 mm Hexalobe Lag Screw |
| 3012-23013-S | 3012-23013 | 2.3 mm x 13 mm Hexalobe Lag Screw |
| 3012-23014-S | 3012-23014 | 2.3 mm x 14 mm Hexalobe Lag Screw |
| 3012-23016-S | 3012-23016 | 2.3 mm x 16 mm Hexalobe Lag Screw |
| 3012-23018-S | 3012-23018 | 2.3 mm x 18 mm Hexalobe Lag Screw |
| 3012-23020-S | 3012-23020 | 2.3 mm x 20 mm Hexalobe Lag Screw |
| 7005-08001-S | 7005-08001 | 0.8 mm Avulsion Hook Plate |
| 7005-08003-S | 7005-08003 | 0.8 mm T-Plate |
| 7005-08004-S | 7005-08004 | 0.8 mm Offset Plate |
| 7005-08006-S | 7005-08006 | 0.8 mm Compression Plate, 6-Hole |
| 7005-08007-S | 7005-08007 | 0.8 mm Curved Medial/Lateral Plate |
| 7005-08010-S | 7005-08010 | 0.8 mm Straight Plate, 10-Hole |
| 7005-13002-S | 7005-13002 | 1.3 mm T-Plate |
| 7005-13004-S | 7005-13004 | 1.3 mm Rolando Fracture Hook Plate |
| 7005-13005-S | 7005-13005 | 1.3 mm Rotational Correction Plate |
| 7005-13006-S | 7005-13006 | 1.3 mm Compression Plate, 6-Hole |

Hand and Wrist Portfolio [continued]

Hand Fracture System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|-------------------------------------|
| 7005-13010-S | 7005-13010 | 1.3 mm Straight Plate, 10-Hole |
| 7005-1303L-S | 7005-1303L | 1.3 mm Metacarpal Neck Plate, Left |
| 7005-1303R-S | 7005-1303R | 1.3 mm Metacarpal Neck Plate, Right |

Modular Hand System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|---|
| CO-F2106-S | CO-F2106 | 2.1 mm x 6.0 mm Cruciform Screw |
| CO-F2108-S | CO-F2108 | 2.1 mm x 8.0 mm Cruciform Screw |
| CO-F2110-S | CO-F2110 | 2.1 mm x 10.0 mm Cruciform Screw |
| CO-F2112-S | CO-F2112 | 2.1 mm x 12.0 mm Cruciform Screw |
| CO-F2114-S | CO-F2114 | 2.1 mm x 14.0 mm Cruciform Screw |
| CO-F2116-S | CO-F2116 | 2.1 mm x 16.0 mm Cruciform Screw |
| CO-F2706-S | CO-F2706 | 2.7 mm x 6.0 mm Cruciform Screw |
| CO-F2708-S | CO-F2708 | 2.7 mm x 8.0 mm Cruciform Screw |
| CO-F2710-S | CO-F2710 | 2.7 mm x 10.0 mm Cruciform Screw |
| CO-F2712-S | CO-F2712 | 2.7 mm x 12.0 mm Cruciform Screw |
| CO-F2714-S | CO-F2714 | 2.7 mm x 14.0 mm Cruciform Screw |
| CO-F2716-S | CO-F2716 | 2.7 mm x 16.0 mm Cruciform Screw |
| PL-MCPL-S | PL-MCPL | 1st Metacarpophalangeal Fusion Plate, Left |
| PL-MCPR-S | PL-MCPR | 1st Metacarpophalangeal Fusion Plate, Right |
| PL-WF33-S | PL-WF33 | Hub Cap® Wrist Fusion Plate, Mini, 3 corner |
| PL-WF40-S | PL-WF40 | Hub Cap Wrist Fusion Plate |
| PL-WF44-S | PL-WF44 | Hub Cap Wrist Fusion Plate, Mini, 4 corner |
| PL-WF60-S | PL-WF60 | Hub Cap Wrist Fusion Plate, Screw Cover |

Hand and Wrist Portfolio [continued]

Osteotomy System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|--|
| 30-0232-S | 30-0232 | 3.5 mm x 8 mm Locking Hexalobe Screw |
| 30-0233-S | 30-0233 | 3.5 mm x 10 mm Locking Hexalobe Screw |
| 30-0234-S | 30-0234 | 3.5 mm x 12 mm Locking Hexalobe Screw |
| 30-0235-S | 30-0235 | 3.5 mm x 14 mm Locking Hexalobe Screw |
| 30-0236-S | 30-0236 | 3.5 mm x 16 mm Locking Hexalobe Screw |
| 30-0237-S | 30-0237 | 3.5 mm x 18 mm Locking Hexalobe Screw |
| 30-0238-S | 30-0238 | 3.5 mm x 20 mm Locking Hexalobe Screw |
| 30-0255-S | 30-0255 | 3.5 mm x 8 mm Nonlocking Hexalobe Screw |
| 30-0256-S | 30-0256 | 3.5 mm x 10 mm Nonlocking Hexalobe Screw |
| 30-0257-S | 30-0257 | 3.5 mm x 12 mm Nonlocking Hexalobe Screw |
| 30-0258-S | 30-0258 | 3.5 mm x 14 mm Nonlocking Hexalobe Screw |
| 30-0259-S | 30-0259 | 3.5 mm x 16 mm Nonlocking Hexalobe Screw |
| 30-0260-S | 30-0260 | 3.5 mm x 18 mm Nonlocking Hexalobe Screw |
| 30-0261-S | 30-0261 | 3.5 mm x 20 mm Nonlocking Hexalobe Screw |
| CO-3080-S | CO-3080 | 3.5 mm x 8.0 mm Cortical Screw |
| CO-3100-S | CO-3100 | 3.5 mm x 10.0 mm Cortical Screw |
| CO-3120-S | CO-3120 | 3.5 mm x 12.0 mm Cortical Screw |
| CO-3140-S | CO-3140 | 3.5 mm x 14.0 mm Cortical Screw |
| CO-3160-S | CO-3160 | 3.5 mm x 16.0 mm Cortical Screw |
| CO-3180-S | CO-3180 | 3.5 mm x 18.0 mm Cortical Screw |
| CO-3200-S | CO-3200 | 3.5 mm x 20.0 mm Cortical Screw |
| COL-3080-S | COL-3080 | 3.5 mm x 8.0 mm Locking Cortical Screw |
| COL-3100-S | COL-3100 | 3.5 mm x 10.0 mm Locking Cortical Screw |
| COL-3120-S | COL-3120 | 3.5 mm x 12.0 mm Locking Cortical Screw |
| COL-3140-S | COL-3140 | 3.5 mm x 14.0 mm Locking Cortical Screw |
| COL-3160-S | COL-3160 | 3.5 mm x 16.0 mm Locking Cortical Screw |
| COL-3180-S | COL-3180 | 3.5 mm x 18.0 mm Locking Cortical Screw |
| COL-3200-S | COL-3200 | 3.5 mm x 20.0 mm Locking Cortical Screw |
| PL-UL06-S | PL-UL06 | 6-Hole Ulna Shortening Plate |

Hand and Wrist Portfolio [continued]

Total Wrist Fusion Plating System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|--|
| 30-0232-S | 30-0232 | 3.5 mm x 8 mm Locking Hexalobe Screw |
| 30-0233-S | 30-0233 | 3.5 mm x 10 mm Locking Hexalobe Screw |
| 30-0234-S | 30-0234 | 3.5 mm x 12 mm Locking Hexalobe Screw |
| 30-0235-S | 30-0235 | 3.5 mm x 14 mm Locking Hexalobe Screw |
| 30-0236-S | 30-0236 | 3.5 mm x 16 mm Locking Hexalobe Screw |
| 30-0237-S | 30-0237 | 3.5 mm x 18 mm Locking Hexalobe Screw |
| 30-0238-S | 30-0238 | 3.5 mm x 20 mm Locking Hexalobe Screw |
| 30-0239-S | 30-0239 | 3.5 mm x 22 mm Locking Hexalobe Screw |
| 30-0240-S | 30-0240 | 3.5 mm x 24 mm Locking Hexalobe Screw |
| 30-0241-S | 30-0241 | 3.5 mm x 26 mm Locking Hexalobe Screw |
| 30-0255-S | 30-0255 | 3.5 mm x 8 mm Nonlocking Hexalobe Screw |
| 30-0256-S | 30-0256 | 3.5 mm x 10 mm Nonlocking Hexalobe Screw |
| 30-0257-S | 30-0257 | 3.5 mm x 12 mm Nonlocking Hexalobe Screw |
| 30-0258-S | 30-0258 | 3.5 mm x 14 mm Nonlocking Hexalobe Screw |
| 30-0259-S | 30-0259 | 3.5 mm x 16 mm Nonlocking Hexalobe Screw |
| 30-0260-S | 30-0260 | 3.5 mm x 18 mm Nonlocking Hexalobe Screw |
| 30-0261-S | 30-0261 | 3.5 mm x 20 mm Nonlocking Hexalobe Screw |
| 30-0262-S | 30-0262 | 3.5 mm x 22 mm Nonlocking Hexalobe Screw |
| 30-0263-S | 30-0263 | 3.5 mm x 24 mm Nonlocking Hexalobe Screw |
| 30-0264-S | 30-0264 | 3.5 mm x 26 mm Nonlocking Hexalobe Screw |
| 3004-23005-S | 3004-23005 | 2.3 mm x 5 mm Hexalobe MultiScrew |
| 3004-23006-S | 3004-23006 | 2.3 mm x 6 mm Hexalobe MultiScrew |
| 3004-23007-S | 3004-23007 | 2.3 mm x 7 mm Hexalobe MultiScrew |
| 3004-23008-S | 3004-23008 | 2.3 mm x 8 mm Hexalobe MultiScrew |
| 3004-23009-S | 3004-23009 | 2.3 mm x 9 mm Hexalobe MultiScrew |
| 3004-23010-S | 3004-23010 | 2.3 mm x 10 mm Hexalobe MultiScrew |
| 3004-23011-S | 3004-23011 | 2.3 mm x 11 mm Hexalobe MultiScrew |
| 3004-23012-S | 3004-23012 | 2.3 mm x 12 mm Hexalobe MultiScrew |
| 3004-23013-S | 3004-23013 | 2.3 mm x 13 mm Hexalobe MultiScrew |

Hand and Wrist Portfolio [continued]

Total Wrist Fusion Plating System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|---|
| 3004-23014-S | 3004-23014 | 2.3 mm x 14 mm Hexalobe MultiScrew |
| 3004-23016-S | 3004-23016 | 2.3 mm x 16 mm Hexalobe MultiScrew |
| 3004-23018-S | 3004-23018 | 2.3 mm x 18 mm Hexalobe MultiScrew |
| 3004-23020-S | 3004-23020 | 2.3 mm x 20 mm Hexalobe MultiScrew |
| 3012-23008-S | 3012-23008 | 2.3 mm x 8 mm Hexalobe Lag Screw |
| 3012-23009-S | 3012-23009 | 2.3 mm x 9 mm Hexalobe Lag Screw |
| 3012-23010-S | 3012-23010 | 2.3 mm x 10 mm Hexalobe Lag Screw |
| 3012-23011-S | 3012-23011 | 2.3 mm x 11 mm Hexalobe Lag Screw |
| 3012-23012-S | 3012-23012 | 2.3 mm x 12 mm Hexalobe Lag Screw |
| 3012-23014-S | 3012-23014 | 2.3 mm x 14 mm Hexalobe Lag Screw |
| 3012-23016-S | 3012-23016 | 2.3 mm x 16 mm Hexalobe Lag Screw |
| 3012-23018-S | 3012-23018 | 2.3 mm x 18 mm Hexalobe Lag Screw |
| 70-0325-S | 70-0325 | Total Wrist Fusion Plate, Standard, Left |
| 70-0326-S | 70-0326 | Total Wrist Fusion Plate, Standard, Right |
| 70-0327-S | 70-0327 | Total Wrist Fusion Plate, Small, Left |
| 70-0328-S | 70-0328 | Total Wrist Fusion Plate, Small, Right |
| 70-0362-S | 70-0362 | Total Wrist Fusion Plate, Neutral |
| CO-3080-S | CO-3080 | 3.5 mm x 8.0 mm Cortical Screw |
| CO-3100-S | CO-3100 | 3.5 mm x 10.0 mm Cortical Screw |
| CO-3120-S | CO-3120 | 3.5 mm x 12.0 mm Cortical Screw |
| CO-3140-S | CO-3140 | 3.5 mm x 14.0 mm Cortical Screw |
| CO-3160-S | CO-3160 | 3.5 mm x 16.0 mm Cortical Screw |
| CO-3180-S | CO-3180 | 3.5 mm x 18.0 mm Cortical Screw |
| CO-3200-S | CO-3200 | 3.5 mm x 20.0 mm Cortical Screw |
| CO-3220-S | CO-3220 | 3.5 mm x 22.0 mm Cortical Screw |
| CO-3240-S | CO-3240 | 3.5 mm x 24.0 mm Cortical Screw |
| CO-3260-S | CO-3260 | 3.5 mm x 26.0 mm Cortical Screw |
| COL-3080-S | COL-3080 | 3.5 mm x 8.0 mm Locking Cortical Screw |
| COL-3100-S | COL-3100 | 3.5 mm x 10.0 mm Locking Cortical Screw |

Hand and Wrist Portfolio [continued]

Total Wrist Fusion Plating System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|---|
| COL-3120-S | COL-3120 | 3.5 mm x 12.0 mm Locking Cortical Screw |
| COL-3140-S | COL-3140 | 3.5 mm x 14.0 mm Locking Cortical Screw |
| COL-3160-S | COL-3160 | 3.5 mm x 16.0 mm Locking Cortical Screw |
| COL-3180-S | COL-3180 | 3.5 mm x 18.0 mm Locking Cortical Screw |
| COL-3200-S | COL-3200 | 3.5 mm x 20.0 mm Locking Cortical Screw |
| COL-3220-S | COL-3220 | 3.5 mm x 22.0 mm Locking Cortical Screw |
| COL-3240-S | COL-3240 | 3.5 mm x 24.0 mm Locking Cortical Screw |
| COL-3260-S | COL-3260 | 3.5 mm x 26.0 mm Locking Cortical Screw |

Screw and Pin Portfolio

- ▶ Acutrak® Headless Compression Screw System
- ▶ Acutrak 2® Headless Compression Screw System
- ▶ Acutrak 3 Headless Compression Screw System
- ▶ AcuTwist® Acutrak Compression Screw System
- ▶ Biotrak® Headless Resorbable Compression System
- ▶ Cannulated Screw System
- ▶ Extremity Screw System
- ▶ Tension Band Pin System/Tension Band Pin System 2



Conditions for Safe Scanning in the MRI

A patient with the screw, fusion device, helical nail, pin, and washer implants from the Acutrak Headless Compression Screw, Acutrak 2 Headless Compression Screw, Acutrak 3 Headless Compression Screw, AcuTwist Acutrak Compression Screw, Biotrak Headless Resorbable Compression, Cannulated Screw, Tension Band Pin/Tension Band Pin 2, and Extremity Screw Systems may be safely scanned at 1.5 T or 3 T under the following conditions. Failure to follow these conditions may result in injury to the patient.

Screw and Pin Portfolio

| Parameter | Condition |
|---|---|
| Nominal values(s) of Static Magnetic Field | 1.5 T or 3 T |
| Maximum Spatial Field Gradient | 30 T/m (3000 gauss/cm) |
| RF Excitation | Circularly Polarized (CP) |
| RF Transmit Coil Type | Whole body transmit coil, Head RF transmit-receive coil |
| Maximum Whole Body SAR | 1.0 W/kg* |
| Limits on Scan Duration | 1.0 W/kg whole body average SAR for 10 minutes of continuous RF (a sequence or back to back series/scan without breaks) followed by a 5 minute cooling period for a one hour session. |
| MR Image Artifact | The presence of this implant may produce an image artifact of 35 mm. |
| If more information about a specific parameter is not included, there are no conditions associated with that parameter. | |

* Whole body SAR of 1.0 W/kg may limit the availability and diagnostic quality of some MR procedures

Screw and Pin Portfolio [continued]

Acutrak® Headless Compression Screw System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|---------------------------------------|
| 30-0033-S | 30-0033 | 30 mm Acutrak Fusion Device 6 mm Nose |
| AM-0025-S | AM-0025 | 25.0 mm Acutrak 4/5 Bone Screw |
| AM-0030-S | AM-0030 | 30.0 mm Acutrak 4/5 Bone Screw |
| AM-0035-S | AM-0035 | 35.0 mm Acutrak 4/5 Bone Screw |
| AM-0040-S | AM-0040 | 40.0 mm Acutrak 4/5 Bone Screw |
| AM-0045-S | AM-0045 | 45.0 mm Acutrak 4/5 Bone Screw |
| AM-0050-S | AM-0050 | 50.0 mm Acutrak 4/5 Bone Screw |
| AP-0035-S | AP-0035 | 35.0 mm Acutrak Plus Fixation Screw |
| AP-0040-S | AP-0040 | 40.0 mm Acutrak Plus Fixation Screw |
| AP-0045-S | AP-0045 | 45.0 mm Acutrak Plus Fixation Screw |
| AP-0050-S | AP-0050 | 50.0 mm Acutrak Plus Fixation Screw |
| AP-0055-S | AP-0055 | 55.0 mm Acutrak Plus Fixation Screw |
| AP-0060-S | AP-0060 | 60.0 mm Acutrak Plus Fixation Screw |
| AP-0065-S | AP-0065 | 65.0 mm Acutrak Plus Fixation Screw |
| AP-0070-S | AP-0070 | 70.0 mm Acutrak Plus Fixation Screw |
| AP-0075-S | AP-0075 | 75.0 mm Acutrak Plus Fixation Screw |
| AP-0080-S | AP-0080 | 80.0 mm Acutrak Plus Fixation Screw |
| AP-6740-S | N/A | 40 mm Acutrak 6/7 Fixation Screw |
| AP-6745-S | N/A | 45 mm Acutrak 6/7 Fixation Screw |
| AP-6750-S | AP-6750 | 50 mm Acutrak 6/7 Fixation Screw |
| AP-6755-S | AP-6755 | 55 mm Acutrak 6/7 Fixation Screw |
| AP-6760-S | AP-6760 | 60 mm Acutrak 6/7 Fixation Screw |
| AP-6765-S | AP-6765 | 65 mm Acutrak 6/7 Fixation Screw |
| AP-6770-S | AP-6770 | 70 mm Acutrak 6/7 Fixation Screw |
| AP-6775-S | AP-6775 | 75 mm Acutrak 6/7 Fixation Screw |
| AP-6780-S | AP-6780 | 80 mm Acutrak 6/7 Fixation Screw |
| AP-6785-S | AP-6785 | 85 mm Acutrak 6/7 Fixation Screw |
| AP-6790-S | AP-6790 | 90 mm Acutrak 6/7 Fixation Screw |
| AP-6795-S | AP-6795 | 95 mm Acutrak 6/7 Fixation Screw |

Screw and Pin Portfolio [continued]

Acutrak® Headless Compression Screw System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|---|
| AP-67100-S | AP-67100 | 100 mm Acutrak 6/7 Fixation Screw |
| AT-1125-S | AT-1125 | 12.5 mm Standard Acutrak Fixation Screw |
| AT-1150-S | AT-1150 | 15.0 mm Standard Acutrak Fixation Screw |
| AT-1175-S | AT-1175 | 17.5 mm Standard Acutrak Fixation Screw |
| AT-1200-S | AT-1200 | 20.0 mm Standard Acutrak Fixation Screw |
| AT-1225-S | AT-1225 | 22.5 mm Standard Acutrak Fixation Screw |
| AT-1250-S | AT-1250 | 25.0 mm Standard Acutrak Fixation Screw |
| AT-1275-S | AT-1275 | 27.5 mm Standard Acutrak Fixation Screw |
| AT-1300-S | AT-1300 | 30.0 mm Standard Acutrak Fixation Screw |
| ATF-140-S | ATF-140 | 14.0 mm Acutrak Fusion Device |
| ATF-160-S | ATF-160 | 16.0 mm Acutrak Fusion Device |
| ATF-180-S | ATF-180 | 18.0 mm Acutrak Fusion Device |
| ATF-200-S | ATF-200 | 20.0 mm Acutrak Fusion Device |
| ATF-220-S | ATF-220 | 22.0 mm Acutrak Fusion Device |
| ATF-240-S | ATF-240 | 24.0 mm Acutrak Fusion Device |
| ATF-270-S | ATF-270 | 27.0 mm Acutrak Fusion Device |
| ATF-320-S | ATF-320 | 32.0 mm Acutrak Fusion Device |
| ATF-370-S | ATF-370 | 37.0 mm Acutrak Fusion Device |
| ATM-008-S | ATM-008 | 8.0 mm Mini Acutrak Fixation Screw |
| ATM-100-S | ATM-100 | 10.0 mm Mini Acutrak Fixation Screw |
| ATM-120-S | ATM-120 | 12.0 mm Mini Acutrak Fixation Screw |
| ATM-140-S | ATM-140 | 14.0 mm Mini Acutrak Fixation Screw |
| ATM-160-S | ATM-160 | 16.0 mm Mini Acutrak Fixation Screw |
| ATM-180-S | ATM-180 | 18.0 mm Mini Acutrak Fixation Screw |
| ATM-200-S | ATM-200 | 20.0 mm Mini Acutrak Fixation Screw |
| ATM-220-S | ATM-220 | 22.0 mm Mini Acutrak Fixation Screw |
| ATM-240-S | ATM-240 | 24.0 mm Mini Acutrak Fixation Screw |
| ATM-260-S | ATM-260 | 26.0 mm Mini Acutrak Fixation Screw |

Screw and Pin Portfolio [continued]

Acutrak 2® Headless Compression Screw System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|----------------------------------|
| 30-0021-S | 30-0021 | 25.0 mm Acutrak 2 - 5.5 mm Screw |
| 30-0023-S | 30-0023 | 30.0 mm Acutrak 2 - 5.5 mm Screw |
| 30-0025-S | 30-0025 | 35.0 mm Acutrak 2 - 5.5 mm Screw |
| 30-0027-S | 30-0027 | 40.0 mm Acutrak 2 - 5.5 mm Screw |
| 30-0029-S | 30-0029 | 45.0 mm Acutrak 2 - 5.5 mm Screw |
| 30-0031-S | 30-0031 | 50.0 mm Acutrak 2 - 5.5 mm Screw |
| 30-0084-S | 30-0084 | 55.0 mm Acutrak 2 - 5.5 mm Screw |
| 30-0085-S | 30-0085 | 60.0 mm Acutrak 2 - 5.5 mm Screw |
| 30-0620-S | 30-0620 | 20.0 mm Acutrak 2 - 4.7 mm Screw |
| 30-0622-S | 30-0622 | 22.0 mm Acutrak 2 - 4.7 mm Screw |
| 30-0624-S | 30-0624 | 24.0 mm Acutrak 2 - 4.7 mm Screw |
| 30-0626-S | 30-0626 | 26.0 mm Acutrak 2 - 4.7 mm Screw |
| 30-0628-S | 30-0628 | 28.0 mm Acutrak 2 - 4.7 mm Screw |
| 30-0630-S | 30-0630 | 30.0 mm Acutrak 2 - 4.7 mm Screw |
| 30-0635-S | 30-0635 | 35.0 mm Acutrak 2 - 4.7 mm Screw |
| 30-0640-S | 30-0640 | 40.0 mm Acutrak 2 - 4.7 mm Screw |
| 30-0645-S | 30-0645 | 45.0 mm Acutrak 2 - 4.7 mm Screw |
| 30-0650-S | 30-0650 | 50.0 mm Acutrak 2 - 4.7 mm Screw |
| 30-0740-S | 30-0740 | 40.0 mm Acutrak 2 - 7.5 mm Screw |
| 30-0745-S | 30-0745 | 45.0 mm Acutrak 2 - 7.5 mm Screw |
| 30-0750-S | 30-0750 | 50.0 mm Acutrak 2 - 7.5 mm Screw |
| 30-0755-S | 30-0755 | 55.0 mm Acutrak 2 - 7.5 mm Screw |
| 30-0760-S | 30-0760 | 60.0 mm Acutrak 2 - 7.5 mm Screw |
| 30-0765-S | 30-0765 | 65.0 mm Acutrak 2 - 7.5 mm Screw |
| 30-0770-S | 30-0770 | 70.0 mm Acutrak 2 - 7.5 mm Screw |
| 30-0775-S | 30-0775 | 75.0 mm Acutrak 2 - 7.5 mm Screw |
| 30-0780-S | 30-0780 | 80.0 mm Acutrak 2 - 7.5 mm Screw |
| 30-0785-S | 30-0785 | 85.0 mm Acutrak 2 - 7.5 mm Screw |
| 30-0790-S | 30-0790 | 90.0 mm Acutrak 2 - 7.5 mm Screw |

Screw and Pin Portfolio [continued]

Acutrak 2® Headless Compression Screw System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|-------------------------------------|
| 30-0795-S | 30-0795 | 95.0 mm Acutrak 2 - 7.5 mm Screw |
| 30-0800-S | 30-0800 | 100.0 mm Acutrak 2 - 7.5 mm Screw |
| 30-0805-S | 30-0805 | 105.0 mm Acutrak 2 - 7.5 mm Screw |
| 30-0810-S | 30-0810 | 110.0 mm Acutrak 2 - 7.5 mm Screw |
| 30-0815-S | 30-0815 | 115.0 mm Acutrak 2 - 7.5 mm Screw |
| 30-0820-S | 30-0820 | 120.0 mm Acutrak 2 - 7.5 mm Screw |
| AT2-C08-S | AT2-C08 | 8.0 mm, Micro Acutrak 2 Bone Screw |
| AT2-C09-S | AT2-C09 | 9.0 mm, Micro Acutrak 2 Bone Screw |
| AT2-C10-S | AT2-C10 | 10.0 mm, Micro Acutrak 2 Bone Screw |
| AT2-C11-S | AT2-C11 | 11.0 mm, Micro Acutrak 2 Bone Screw |
| AT2-C12-S | AT2-C12 | 12.0 mm, Micro Acutrak 2 Bone Screw |
| AT2-C13-S | AT2-C13 | 13.0 mm, Micro Acutrak 2 Bone Screw |
| AT2-C14-S | AT2-C14 | 14.0 mm, Micro Acutrak 2 Bone Screw |
| AT2-C16-S | AT2-C16 | 16.0 mm, Micro Acutrak 2 Bone Screw |
| AT2-C18-S | AT2-C18 | 18.0 mm, Micro Acutrak 2 Bone Screw |
| AT2-C20-S | AT2-C20 | 20.0 mm, Micro Acutrak 2 Bone Screw |
| AT2-C22-S | AT2-C22 | 22.0 mm, Micro Acutrak 2 Bone Screw |
| AT2-C24-S | AT2-C24 | 24.0 mm, Micro Acutrak 2 Bone Screw |
| AT2-C26-S | AT2-C26 | 26.0 mm, Micro Acutrak 2 Bone Screw |
| AT2-C28-S | AT2-C28 | 28.0 mm, Micro Acutrak 2 Bone Screw |
| AT2-C30-S | AT2-C30 | 30.0 mm, Micro Acutrak 2 Bone Screw |
| AT2-M16-S | AT2-M16 | 16.0 mm, Mini Acutrak 2 Bone Screw |
| AT2-M18-S | AT2-M18 | 18.0 mm, Mini Acutrak 2 Bone Screw |
| AT2-M20-S | AT2-M20 | 20.0 mm, Mini Acutrak 2 Bone Screw |
| AT2-M22-S | AT2-M22 | 22.0 mm, Mini Acutrak 2 Bone Screw |
| AT2-M24-S | AT2-M24 | 24.0 mm, Mini Acutrak 2 Bone Screw |
| AT2-M26-S | AT2-M26 | 26.0 mm, Mini Acutrak 2 Bone Screw |
| AT2-M28-S | AT2-M28 | 28.0 mm, Mini Acutrak 2 Bone Screw |
| AT2-M30-S | AT2-M30 | 30.0 mm, Mini Acutrak 2 Bone Screw |

Screw and Pin Portfolio [continued]

Acutrak 2® Headless Compression Screw System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|--|
| AT2-S16-S | AT2-S16 | 16.0 mm, Standard Acutrak 2 Bone Screw |
| AT2-S18-S | AT2-S18 | 18.0 mm, Standard Acutrak 2 Bone Screw |
| AT2-S20-S | AT2-S20 | 20.0 mm, Standard Acutrak 2 Bone Screw |
| AT2-S22-S | AT2-S22 | 22.0 mm, Standard Acutrak 2 Bone Screw |
| AT2-S24-S | AT2-S24 | 24.0 mm, Standard Acutrak 2 Bone Screw |
| AT2-S26-S | AT2-S26 | 26.0 mm, Standard Acutrak 2 Bone Screw |
| AT2-S28-S | AT2-S28 | 28.0 mm, Standard Acutrak 2 Bone Screw |
| AT2-S30-S | AT2-S30 | 30.0 mm, Standard Acutrak 2 Bone Screw |
| AT2-S32-S | AT2-S32 | 32.0 mm, Standard Acutrak 2 Bone Screw |
| AT2-S34-S | AT2-S34 | 34.0 mm, Standard Acutrak 2 Bone Screw |

Acutrak® 3 Headless Compression Screw System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|------------------------------------|
| 3050-20008-S | 3050-20008 | 8.0 mm, Nano Acutrak 3 Bone Screw |
| 3050-20009-S | 3050-20009 | 9.0 mm, Nano Acutrak 3 Bone Screw |
| 3050-20010-S | 3050-20010 | 10.0 mm, Nano Acutrak 3 Bone Screw |
| 3050-20011-S | 3050-20011 | 11.0 mm, Nano Acutrak 3 Bone Screw |
| 3050-20012-S | 3050-20012 | 12.0 mm, Nano Acutrak 3 Bone Screw |
| 3050-20013-S | 3050-20013 | 13.0 mm, Nano Acutrak 3 Bone Screw |
| 3050-20014-S | 3050-20014 | 14.0 mm, Nano Acutrak 3 Bone Screw |
| 3050-20016-S | 3050-20016 | 16.0 mm, Nano Acutrak 3 Bone Screw |
| 3050-20018-S | 3050-20018 | 18.0 mm, Nano Acutrak 3 Bone Screw |
| 3050-20020-S | 3050-20020 | 20.0 mm, Nano Acutrak 3 Bone Screw |
| 3050-20022-S | 3050-20022 | 22.0 mm, Nano Acutrak 3 Bone Screw |
| 3050-20024-S | 3050-20024 | 24.0 mm, Nano Acutrak 3 Bone Screw |
| 3050-20026-S | 3050-20026 | 26.0 mm, Nano Acutrak 3 Bone Screw |
| 3050-20028-S | 3050-20028 | 28.0 mm, Nano Acutrak 3 Bone Screw |
| 3050-20030-S | 3050-20030 | 30.0 mm, Nano Acutrak 3 Bone Screw |
| 3050-20032-S | 3050-20032 | 32.0 mm, Nano Acutrak 3 Bone Screw |

Screw and Pin Portfolio [continued]

Acutrak® 3 Headless Compression Screw System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|-------------------------------------|
| 3050-20034-S | 3050-20034 | 34.0 mm, Nano Acutrak 3 Bone Screw |
| 3050-20036-S | 3050-20036 | 36.0 mm, Nano Acutrak 3 Bone Screw |
| 3050-20038-S | 3050-20038 | 38.0 mm, Nano Acutrak 3 Bone Screw |
| 3050-20040-S | 3050-20040 | 40.0 mm, Nano Acutrak 3 Bone Screw |
| 3051-25008-S | 3051-25008 | 8.0 mm, Micro Acutrak 3 Bone Screw |
| 3051-25009-S | 3051-25009 | 9.0 mm, Micro Acutrak 3 Bone Screw |
| 3051-25010-S | 3051-25010 | 10.0 mm, Micro Acutrak 3 Bone Screw |
| 3051-25011-S | 3051-25011 | 11.0 mm, Micro Acutrak 3 Bone Screw |
| 3051-25012-S | 3051-25012 | 12.0 mm, Micro Acutrak 3 Bone Screw |
| 3051-25013-S | 3051-25013 | 13.0 mm, Micro Acutrak 3 Bone Screw |
| 3051-25014-S | 3051-25014 | 14.0 mm, Micro Acutrak 3 Bone Screw |
| 3051-25016-S | 3051-25016 | 16.0 mm, Micro Acutrak 3 Bone Screw |
| 3051-25018-S | 3051-25018 | 18.0 mm, Micro Acutrak 3 Bone Screw |
| 3051-25020-S | 3051-25020 | 20.0 mm, Micro Acutrak 3 Bone Screw |
| 3051-25022-S | 3051-25022 | 22.0 mm, Micro Acutrak 3 Bone Screw |
| 3051-25024-S | 3051-25024 | 24.0 mm, Micro Acutrak 3 Bone Screw |
| 3051-25026-S | 3051-25026 | 26.0 mm, Micro Acutrak 3 Bone Screw |
| 3051-25028-S | 3051-25028 | 28.0 mm, Micro Acutrak 3 Bone Screw |
| 3051-25030-S | 3051-25030 | 30.0 mm, Micro Acutrak 3 Bone Screw |
| 3051-25032-S | 3051-25032 | 32.0 mm, Micro Acutrak 3 Bone Screw |
| 3051-25034-S | 3051-25034 | 34.0 mm, Micro Acutrak 3 Bone Screw |
| 3051-25036-S | 3051-25036 | 36.0 mm, Micro Acutrak 3 Bone Screw |
| 3051-25038-S | 3051-25038 | 38.0 mm, Micro Acutrak 3 Bone Screw |
| 3051-25040-S | 3051-25040 | 40.0 mm, Micro Acutrak 3 Bone Screw |
| 3051-25042-S | 3051-25042 | 42.0 mm, Micro Acutrak 3 Bone Screw |
| 3051-25044-S | 3051-25044 | 44.0 mm, Micro Acutrak 3 Bone Screw |
| 3051-25046-S | 3051-25046 | 46.0 mm, Micro Acutrak 3 Bone Screw |
| 3051-25048-S | 3051-25048 | 48.0 mm, Micro Acutrak 3 Bone Screw |
| 3051-25050-S | 3051-25050 | 50.0 mm, Micro Acutrak 3 Bone Screw |
| 3052-35012-S | 3052-35012 | 12.0 mm, Mini Acutrak 3 Bone Screw |

Screw and Pin Portfolio [continued]

Acutrak® 3 Headless Compression Screw System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|--|
| 3052-35014-S | 3052-35014 | 14.0 mm, Mini Acutrak 3 Bone Screw |
| 3052-35016-S | 3052-35016 | 16.0 mm, Mini Acutrak 3 Bone Screw |
| 3052-35018-S | 3052-35018 | 18.0 mm, Mini Acutrak 3 Bone Screw |
| 3052-35020-S | 3052-35020 | 20.0 mm, Mini Acutrak 3 Bone Screw |
| 3052-35022-S | 3052-35022 | 22.0 mm, Mini Acutrak 3 Bone Screw |
| 3052-35024-S | 3052-35024 | 24.0 mm, Mini Acutrak 3 Bone Screw |
| 3052-35026-S | 3052-35026 | 26.0 mm, Mini Acutrak 3 Bone Screw |
| 3052-35028-S | 3052-35028 | 28.0 mm, Mini Acutrak 3 Bone Screw |
| 3052-35030-S | 3052-35030 | 30.0 mm, Mini Acutrak 3 Bone Screw |
| 3052-35032-S | 3052-35032 | 32.0 mm, Mini Acutrak 3 Bone Screw |
| 3052-35034-S | 3052-35034 | 34.0 mm, Mini Acutrak 3 Bone Screw |
| 3052-35036-S | 3052-35036 | 36.0 mm, Mini Acutrak 3 Bone Screw |
| 3052-35038-S | 3052-35038 | 38.0 mm, Mini Acutrak 3 Bone Screw |
| 3052-35040-S | 3052-35040 | 40.0 mm, Mini Acutrak 3 Bone Screw |
| 3052-35042-S | 3052-35042 | 42.0 mm, Mini Acutrak 3 Bone Screw |
| 3052-35044-S | 3052-35044 | 44.0 mm, Mini Acutrak 3 Bone Screw |
| 3052-35046-S | 3052-35046 | 46.0 mm, Mini Acutrak 3 Bone Screw |
| 3052-35048-S | 3052-35048 | 48.0 mm, Mini Acutrak 3 Bone Screw |
| 3052-35050-S | 3052-35050 | 50.0 mm, Mini Acutrak 3 Bone Screw |
| 3052-35052-S | 3052-35052 | 52.0 mm, Mini Acutrak 3 Bone Screw |
| 3052-35054-S | 3052-35054 | 54.0 mm, Mini Acutrak 3 Bone Screw |
| 3052-35056-S | 3052-35056 | 56.0 mm, Mini Acutrak 3 Bone Screw |
| 3052-35058-S | 3052-35058 | 58.0 mm, Mini Acutrak 3 Bone Screw |
| 3052-35060-S | 3052-35060 | 60.0 mm, Mini Acutrak 3 Bone Screw |
| 3053-40016-S | 3053-40016 | 16.0 mm, Standard Acutrak 3 Bone Screw |
| 3053-40018-S | 3053-40018 | 18.0 mm, Standard Acutrak 3 Bone Screw |
| 3053-40020-S | 3053-40020 | 20.0 mm, Standard Acutrak 3 Bone Screw |
| 3053-40022-S | 3053-40022 | 22.0 mm, Standard Acutrak 3 Bone Screw |
| 3053-40024-S | 3053-40024 | 24.0 mm, Standard Acutrak 3 Bone Screw |

Screw and Pin Portfolio [continued]

Acutrak® 3 Headless Compression Screw System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|--|
| 3053-40026-S | 3053-40026 | 26.0 mm, Standard Acutrak 3 Bone Screw |
| 3053-40028-S | 3053-40028 | 28.0 mm, Standard Acutrak 3 Bone Screw |
| 3053-40030-S | 3053-40030 | 30.0 mm, Standard Acutrak 3 Bone Screw |
| 3053-40032-S | 3053-40032 | 32.0 mm, Standard Acutrak 3 Bone Screw |
| 3053-40034-S | 3053-40034 | 34.0 mm, Standard Acutrak 3 Bone Screw |
| 3053-40036-S | 3053-40036 | 36.0 mm, Standard Acutrak 3 Bone Screw |
| 3053-40038-S | 3053-40038 | 38.0 mm, Standard Acutrak 3 Bone Screw |
| 3053-40040-S | 3053-40040 | 40.0 mm, Standard Acutrak 3 Bone Screw |
| 3053-40042-S | 3053-40042 | 42.0 mm, Standard Acutrak 3 Bone Screw |
| 3053-40044-S | 3053-40044 | 44.0 mm, Standard Acutrak 3 Bone Screw |
| 3053-40046-S | 3053-40046 | 46.0 mm, Standard Acutrak 3 Bone Screw |
| 3053-40048-S | 3053-40048 | 48.0 mm, Standard Acutrak 3 Bone Screw |
| 3053-40050-S | 3053-40050 | 50.0 mm, Standard Acutrak 3 Bone Screw |
| 3053-40052-S | 3053-40052 | 52.0 mm, Standard Acutrak 3 Bone Screw |
| 3053-40054-S | 3053-40054 | 54.0 mm, Standard Acutrak 3 Bone Screw |
| 3053-40056-S | 3053-40056 | 56.0 mm, Standard Acutrak 3 Bone Screw |
| 3053-40058-S | 3053-40058 | 58.0 mm, Standard Acutrak 3 Bone Screw |
| 3053-40060-S | 3053-40060 | 60.0 mm, Standard Acutrak 3 Bone Screw |

Screw and Pin Portfolio [continued]

AcuTwist® Acutrak Compression Screw System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|--|
| AI-0010-S | AI-0010 | 10 mm AcuTwist Acutrak Compression Screw |
| AI-0012-S | AI-0012 | 12 mm AcuTwist Acutrak Compression Screw |
| AI-0014-S | AI-0014 | 14 mm AcuTwist Acutrak Compression Screw |
| AI-0016-S | AI-0016 | 16 mm AcuTwist Acutrak Compression Screw |
| AI-0018-S | AI-0018 | 18 mm AcuTwist Acutrak Compression Screw |
| AI-0020-S | AI-0020 | 20 mm AcuTwist Acutrak Compression Screw |
| AI-0022-S | AI-0022 | 22 mm AcuTwist Acutrak Compression Screw |
| AI-0024-S | AI-0024 | 24 mm AcuTwist Acutrak Compression Screw |
| AI-0026-S | AI-0026 | 26 mm AcuTwist Acutrak Compression Screw |
| AI-0028-S | AI-0028 | 28 mm AcuTwist Acutrak Compression Screw |
| AI-0030-S | AI-0030 | 30 mm AcuTwist Acutrak Compression Screw |

Biotrak® Headless Resorbable Compression Screw System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|----------------------------|
| 30170007-S | N/A | Biotrak 16 mm Screw |
| 30170008-S | N/A | Biotrak 18 mm Screw |
| 30170009-S | N/A | Biotrak 20 mm Screw |
| 30170010-S | N/A | Biotrak 22 mm Screw |
| 30170011-S | N/A | Biotrak 24 mm Screw |
| 30170057-S | N/A | Biotrak Mini 16 mm Screw |
| 30170058-S | N/A | Biotrak Mini 18 mm Screw |
| 30170059-S | N/A | Biotrak Mini 20 mm Screw |
| 30170060-S | N/A | Biotrak Mini 22 mm Screw |
| 30170061-S | N/A | Biotrak Mini 24 mm Screw |
| 30170320-S | N/A | Biotrak 2.0 mm x 20 mm Pin |
| 30170330-S | N/A | Biotrak 2.0 mm x 30 mm Pin |
| 30170340-S | N/A | Biotrak 2.0 mm x 40 mm Pin |
| 55-0001-S | 55-0001 | Biotrak 20 mm Helical Nail |
| 55-0002-S | 55-0002 | Biotrak 30 mm Helical Nail |
| 55-0003-S | 55-0003 | Biotrak 40 mm Helical Nail |

Screw and Pin Portfolio [continued]

Cannulated Screw System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|---|
| N/A | 3005-40010 | 4.0 mm x 10 mm Cannulated Screw, Short Thread |
| N/A | 3005-40012 | 4.0 mm x 12 mm Cannulated Screw, Short Thread |
| N/A | 3005-40014 | 4.0 mm x 14 mm Cannulated Screw, Short Thread |
| N/A | 3005-40016 | 4.0 mm x 16 mm Cannulated Screw, Short Thread |
| N/A | 3005-40018 | 4.0 mm x 18 mm Cannulated Screw, Short Thread |
| N/A | 3005-40020 | 4.0 mm x 20 mm Cannulated Screw, Short Thread |
| N/A | 3005-40022 | 4.0 mm x 22 mm Cannulated Screw, Short Thread |
| N/A | 3005-40024 | 4.0 mm x 24 mm Cannulated Screw, Short Thread |
| N/A | 3005-40026 | 4.0 mm x 26 mm Cannulated Screw, Short Thread |
| N/A | 3005-40028 | 4.0 mm x 28 mm Cannulated Screw, Short Thread |
| N/A | 3005-40030 | 4.0 mm x 30 mm Cannulated Screw, Short Thread |
| N/A | 3005-40032 | 4.0 mm x 32 mm Cannulated Screw, Short Thread |
| N/A | 3005-40034 | 4.0 mm x 34 mm Cannulated Screw, Short Thread |
| N/A | 3005-40036 | 4.0 mm x 36 mm Cannulated Screw, Short Thread |
| N/A | 3005-40038 | 4.0 mm x 38 mm Cannulated Screw, Short Thread |
| N/A | 3005-40040 | 4.0 mm x 40 mm Cannulated Screw, Short Thread |
| N/A | 3005-40042 | 4.0 mm x 42 mm Cannulated Screw, Short Thread |
| N/A | 3005-40044 | 4.0 mm x 44 mm Cannulated Screw, Short Thread |
| N/A | 3005-40046 | 4.0 mm x 46 mm Cannulated Screw, Short Thread |
| N/A | 3005-40048 | 4.0 mm x 48 mm Cannulated Screw, Short Thread |
| N/A | 3005-40050 | 4.0 mm x 50 mm Cannulated Screw, Short Thread |
| N/A | 3005-40052 | 4.0 mm x 52 mm Cannulated Screw, Short Thread |
| N/A | 3005-40054 | 4.0 mm x 54 mm Cannulated Screw, Short Thread |
| N/A | 3005-40056 | 4.0 mm x 56 mm Cannulated Screw, Short Thread |
| N/A | 3005-40058 | 4.0 mm x 58 mm Cannulated Screw, Short Thread |
| N/A | 3005-40060 | 4.0 mm x 60 mm Cannulated Screw, Short Thread |
| N/A | 3005-40064 | 4.0 mm x 64 mm Cannulated Screw, Short Thread |
| N/A | 3005-40068 | 4.0 mm x 68 mm Cannulated Screw, Short Thread |
| N/A | 3005-40072 | 4.0 mm x 72 mm Cannulated Screw, Short Thread |

Screw and Pin Portfolio [continued]

Cannulated Screw System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|--|
| N/A | 3005-65030 | 6.5 mm x 30 mm Cannulated Screw, 16 mm Thread |
| N/A | 3005-65035 | 6.5 mm x 35 mm Cannulated Screw, 16 mm Thread |
| N/A | 3005-65040 | 6.5 mm x 40 mm Cannulated Screw, 16 mm Thread |
| N/A | 3005-65045 | 6.5 mm x 45 mm Cannulated Screw, 16 mm Thread |
| N/A | 3005-65050 | 6.5 mm x 50 mm Cannulated Screw, 16 mm Thread |
| N/A | 3005-65055 | 6.5 mm x 55 mm Cannulated Screw, 16 mm Thread |
| N/A | 3005-65060 | 6.5 mm x 60 mm Cannulated Screw, 16 mm Thread |
| N/A | 3005-65065 | 6.5 mm x 65 mm Cannulated Screw, 16 mm Thread |
| N/A | 3005-65070 | 6.5 mm x 70 mm Cannulated Screw, 16 mm Thread |
| N/A | 3005-65075 | 6.5 mm x 75 mm Cannulated Screw, 16 mm Thread |
| N/A | 3005-65080 | 6.5 mm x 80 mm Cannulated Screw, 16 mm Thread |
| N/A | 3005-65085 | 6.5 mm x 85 mm Cannulated Screw, 16 mm Thread |
| N/A | 3005-65090 | 6.5 mm x 90 mm Cannulated Screw, 16 mm Thread |
| N/A | 3005-65095 | 6.5 mm x 95 mm Cannulated Screw, 16 mm Thread |
| N/A | 3005-65100 | 6.5 mm x 100 mm Cannulated Screw, 16 mm Thread |
| N/A | 3005-65105 | 6.5 mm x 105 mm Cannulated Screw, 16 mm Thread |
| N/A | 3005-65110 | 6.5 mm x 110 mm Cannulated Screw, 16 mm Thread |
| N/A | 3005-65115 | 6.5 mm x 115 mm Cannulated Screw, 16 mm Thread |
| N/A | 3005-65120 | 6.5 mm x 120 mm Cannulated Screw, 16 mm Thread |
| N/A | 3005-65125 | 6.5 mm x 125 mm Cannulated Screw, 16 mm Thread |
| N/A | 3005-65130 | 6.5 mm x 130 mm Cannulated Screw, 16 mm Thread |
| N/A | 3005-65135 | 6.5 mm x 135 mm Cannulated Screw, 16 mm Thread |
| N/A | 3005-65140 | 6.5 mm x 140 mm Cannulated Screw, 16 mm Thread |
| N/A | 3005-65145 | 6.5 mm x 145 mm Cannulated Screw, 16 mm Thread |
| N/A | 3005-65150 | 6.5 mm x 150 mm Cannulated Screw, 16 mm Thread |
| N/A | 3005-73030 | 7.3 mm x 30 mm Cannulated Screw, 16 mm Thread |
| N/A | 3005-73035 | 7.3 mm x 35 mm Cannulated Screw, 16 mm Thread |
| N/A | 3005-73040 | 7.3 mm x 40 mm Cannulated Screw, 16 mm Thread |
| N/A | 3005-73045 | 7.3 mm x 45 mm Cannulated Screw, 16 mm Thread |

Screw and Pin Portfolio [continued]

Cannulated Screw System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|--|
| N/A | 3005-73050 | 7.3 mm x 50 mm Cannulated Screw, 16 mm Thread |
| N/A | 3005-73055 | 7.3 mm x 55 mm Cannulated Screw, 16 mm Thread |
| N/A | 3005-73060 | 7.3 mm x 60 mm Cannulated Screw, 16 mm Thread |
| N/A | 3005-73065 | 7.3 mm x 65 mm Cannulated Screw, 16 mm Thread |
| N/A | 3005-73070 | 7.3 mm x 70 mm Cannulated Screw, 16 mm Thread |
| N/A | 3005-73075 | 7.3 mm x 75 mm Cannulated Screw, 16 mm Thread |
| N/A | 3005-73080 | 7.3 mm x 80 mm Cannulated Screw, 16 mm Thread |
| N/A | 3005-73085 | 7.3 mm x 85 mm Cannulated Screw, 16 mm Thread |
| N/A | 3005-73090 | 7.3 mm x 90 mm Cannulated Screw, 16 mm Thread |
| N/A | 3005-73095 | 7.3 mm x 95 mm Cannulated Screw, 16 mm Thread |
| N/A | 3005-73100 | 7.3 mm x 100 mm Cannulated Screw, 16 mm Thread |
| N/A | 3005-73105 | 7.3 mm x 105 mm Cannulated Screw, 16 mm Thread |
| N/A | 3005-73110 | 7.3 mm x 110 mm Cannulated Screw, 16 mm Thread |
| N/A | 3005-73115 | 7.3 mm x 115 mm Cannulated Screw, 16 mm Thread |
| N/A | 3005-73120 | 7.3 mm x 120 mm Cannulated Screw, 16 mm Thread |
| N/A | 3005-73125 | 7.3 mm x 125 mm Cannulated Screw, 16 mm Thread |
| N/A | 3005-73130 | 7.3 mm x 130 mm Cannulated Screw, 16 mm Thread |
| N/A | 3005-73135 | 7.3 mm x 135 mm Cannulated Screw, 16 mm Thread |
| N/A | 3005-73140 | 7.3 mm x 140 mm Cannulated Screw, 16 mm Thread |
| N/A | 3005-73145 | 7.3 mm x 145 mm Cannulated Screw, 16 mm Thread |
| N/A | 3005-73150 | 7.3 mm x 150 mm Cannulated Screw, 16 mm Thread |
| N/A | 3006-40016 | 4.0 mm x 16 mm Cannulated Screw, Long Thread |
| N/A | 3006-40018 | 4.0 mm x 18 mm Cannulated Screw, Long Thread |
| N/A | 3006-40020 | 4.0 mm x 20 mm Cannulated Screw, Long Thread |
| N/A | 3006-40022 | 4.0 mm x 22 mm Cannulated Screw, Long Thread |
| N/A | 3006-40024 | 4.0 mm x 24 mm Cannulated Screw, Long Thread |
| N/A | 3006-40026 | 4.0 mm x 26 mm Cannulated Screw, Long Thread |
| N/A | 3006-40028 | 4.0 mm x 28 mm Cannulated Screw, Long Thread |
| N/A | 3006-40030 | 4.0 mm x 30 mm Cannulated Screw, Long Thread |

Screw and Pin Portfolio [continued]

Cannulated Screw System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|---|
| N/A | 3006-40032 | 4.0 mm x 32 mm Cannulated Screw, Long Thread |
| N/A | 3006-40034 | 4.0 mm x 34 mm Cannulated Screw, Long Thread |
| 3006-40036-S | 3006-40036 | 4.0 mm x 36 mm Cannulated Screw, Long Thread |
| N/A | 3006-40038 | 4.0 mm x 38 mm Cannulated Screw, Long Thread |
| N/A | 3006-40040 | 4.0 mm x 40 mm Cannulated Screw, Long Thread |
| 3006-40042-S | 3006-40042 | 4.0 mm x 42 mm Cannulated Screw, Long Thread |
| N/A | 3006-40044 | 4.0 mm x 44 mm Cannulated Screw, Long Thread |
| N/A | 3006-40046 | 4.0 mm x 46 mm Cannulated Screw, Long Thread |
| 3006-40048-S | 3006-40048 | 4.0 mm x 48 mm Cannulated Screw, Long Thread |
| N/A | 3006-40050 | 4.0 mm x 50 mm Cannulated Screw, Long Thread |
| N/A | 3006-40052 | 4.0 mm x 52 mm Cannulated Screw, Long Thread |
| N/A | 3006-40054 | 4.0 mm x 54 mm Cannulated Screw, Long Thread |
| N/A | 3006-40056 | 4.0 mm x 56 mm Cannulated Screw, Long Thread |
| N/A | 3006-40058 | 4.0 mm x 58 mm Cannulated Screw, Long Thread |
| N/A | 3006-40060 | 4.0 mm x 60 mm Cannulated Screw, Long Thread |
| N/A | 3006-40064 | 4.0 mm x 64 mm Cannulated Screw, Long Thread |
| N/A | 3006-40068 | 4.0 mm x 68 mm Cannulated Screw, Long Thread |
| N/A | 3006-40072 | 4.0 mm x 72 mm Cannulated Screw, Long Thread |
| N/A | 3006-65045 | 6.5 mm x 45 mm Cannulated Screw, 32 mm Thread |
| N/A | 3006-65050 | 6.5 mm x 50 mm Cannulated Screw, 32 mm Thread |
| N/A | 3006-65055 | 6.5 mm x 55 mm Cannulated Screw, 32 mm Thread |
| N/A | 3006-65060 | 6.5 mm x 60 mm Cannulated Screw, 32 mm Thread |
| N/A | 3006-65065 | 6.5 mm x 65 mm Cannulated Screw, 32 mm Thread |
| N/A | 3006-65070 | 6.5 mm x 70 mm Cannulated Screw, 32 mm Thread |
| N/A | 3006-65075 | 6.5 mm x 75 mm Cannulated Screw, 32 mm Thread |
| N/A | 3006-65080 | 6.5 mm x 80 mm Cannulated Screw, 32 mm Thread |
| N/A | 3006-65085 | 6.5 mm x 85 mm Cannulated Screw, 32 mm Thread |
| N/A | 3006-65090 | 6.5 mm x 90 mm Cannulated Screw, 32 mm Thread |
| N/A | 3006-65095 | 6.5 mm x 95 mm Cannulated Screw, 32 mm Thread |

Screw and Pin Portfolio [continued]

Cannulated Screw System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|--|
| N/A | 3006-65100 | 6.5 mm x 100 mm Cannulated Screw, 32 mm Thread |
| N/A | 3006-65105 | 6.5 mm x 105 mm Cannulated Screw, 32 mm Thread |
| N/A | 3006-65110 | 6.5 mm x 110 mm Cannulated Screw, 32 mm Thread |
| N/A | 3006-65115 | 6.5 mm x 115 mm Cannulated Screw, 32 mm Thread |
| N/A | 3006-65120 | 6.5 mm x 120 mm Cannulated Screw, 32 mm Thread |
| N/A | 3006-65125 | 6.5 mm x 125 mm Cannulated Screw, 32 mm Thread |
| N/A | 3006-65130 | 6.5 mm x 130 mm Cannulated Screw, 32 mm Thread |
| N/A | 3006-65135 | 6.5 mm x 135 mm Cannulated Screw, 32 mm Thread |
| N/A | 3006-65140 | 6.5 mm x 140 mm Cannulated Screw, 32 mm Thread |
| N/A | 3006-65145 | 6.5 mm x 145 mm Cannulated Screw, 32 mm Thread |
| N/A | 3006-65150 | 6.5 mm x 150 mm Cannulated Screw, 32 mm Thread |
| N/A | 3006-73045 | 7.3 mm x 45 mm Cannulated Screw, 32 mm Thread |
| N/A | 3006-73050 | 7.3 mm x 50 mm Cannulated Screw, 32 mm Thread |
| N/A | 3006-73055 | 7.3 mm x 55 mm Cannulated Screw, 32 mm Thread |
| N/A | 3006-73060 | 7.3 mm x 60 mm Cannulated Screw, 32 mm Thread |
| N/A | 3006-73065 | 7.3 mm x 65 mm Cannulated Screw, 32 mm Thread |
| N/A | 3006-73070 | 7.3 mm x 70 mm Cannulated Screw, 32 mm Thread |
| N/A | 3006-73075 | 7.3 mm x 75 mm Cannulated Screw, 32 mm Thread |
| N/A | 3006-73080 | 7.3 mm x 80 mm Cannulated Screw, 32 mm Thread |
| N/A | 3006-73085 | 7.3 mm x 85 mm Cannulated Screw, 32 mm Thread |
| N/A | 3006-73090 | 7.3 mm x 90 mm Cannulated Screw, 32 mm Thread |
| N/A | 3006-73095 | 7.3 mm x 95 mm Cannulated Screw, 32 mm Thread |
| N/A | 3006-73100 | 7.3 mm x 100 mm Cannulated Screw, 32 mm Thread |
| N/A | 3006-73105 | 7.3 mm x 105 mm Cannulated Screw, 32 mm Thread |
| N/A | 3006-73110 | 7.3 mm x 110 mm Cannulated Screw, 32 mm Thread |
| N/A | 3006-73115 | 7.3 mm x 115 mm Cannulated Screw, 32 mm Thread |
| N/A | 3006-73120 | 7.3 mm x 120 mm Cannulated Screw, 32 mm Thread |
| N/A | 3006-73125 | 7.3 mm x 125 mm Cannulated Screw, 32 mm Thread |
| N/A | 3006-73130 | 7.3 mm x 130 mm Cannulated Screw, 32 mm Thread |

Screw and Pin Portfolio [continued]

Cannulated Screw System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|--|
| N/A | 3006-73135 | 7.3 mm x 135 mm Cannulated Screw, 32 mm Thread |
| N/A | 3006-73140 | 7.3 mm x 140 mm Cannulated Screw, 32 mm Thread |
| N/A | 3006-73145 | 7.3 mm x 145 mm Cannulated Screw, 32 mm Thread |
| N/A | 3006-73150 | 7.3 mm x 150 mm Cannulated Screw, 32 mm Thread |
| N/A | 3007-65030 | 6.5 mm x 30 mm Cannulated Screw, Full Thread |
| N/A | 3007-65035 | 6.5 mm x 35 mm Cannulated Screw, Full Thread |
| N/A | 3007-65040 | 6.5 mm x 40 mm Cannulated Screw, Full Thread |
| N/A | 3007-65045 | 6.5 mm x 45 mm Cannulated Screw, Full Thread |
| N/A | 3007-65050 | 6.5 mm x 50 mm Cannulated Screw, Full Thread |
| N/A | 3007-65055 | 6.5 mm x 55 mm Cannulated Screw, Full Thread |
| N/A | 3007-65060 | 6.5 mm x 60 mm Cannulated Screw, Full Thread |
| N/A | 3007-65065 | 6.5 mm x 65 mm Cannulated Screw, Full Thread |
| N/A | 3007-65070 | 6.5 mm x 70 mm Cannulated Screw, Full Thread |
| N/A | 3007-65075 | 6.5 mm x 75 mm Cannulated Screw, Full Thread |
| N/A | 3007-65080 | 6.5 mm x 80 mm Cannulated Screw, Full Thread |
| N/A | 3007-65085 | 6.5 mm x 85 mm Cannulated Screw, Full Thread |
| N/A | 3007-65090 | 6.5 mm x 90 mm Cannulated Screw, Full Thread |
| N/A | 3007-65095 | 6.5 mm x 95 mm Cannulated Screw, Full Thread |
| N/A | 3007-65100 | 6.5 mm x 100 mm Cannulated Screw, Full Thread |
| N/A | 3007-65105 | 6.5 mm x 105 mm Cannulated Screw, Full Thread |
| N/A | 3007-65110 | 6.5 mm x 110 mm Cannulated Screw, Full Thread |
| N/A | 3007-65115 | 6.5 mm x 115 mm Cannulated Screw, Full Thread |
| N/A | 3007-65120 | 6.5 mm x 120 mm Cannulated Screw, Full Thread |
| N/A | 3007-65125 | 6.5 mm x 125 mm Cannulated Screw, Full Thread |
| N/A | 3007-65130 | 6.5 mm x 130 mm Cannulated Screw, Full Thread |
| N/A | 3007-65135 | 6.5 mm x 135 mm Cannulated Screw, Full Thread |
| N/A | 3007-65140 | 6.5 mm x 140 mm Cannulated Screw, Full Thread |
| N/A | 3007-65145 | 6.5 mm x 145 mm Cannulated Screw, Full Thread |
| N/A | 3007-65150 | 6.5 mm x 150 mm Cannulated Screw, Full Thread |

Screw and Pin Portfolio [continued]

Cannulated Screw System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|--|
| N/A | 3007-73030 | 7.3 mm x 30 mm Cannulated Screw, Full Thread |
| N/A | 3007-73035 | 7.3 mm x 35 mm Cannulated Screw, Full Thread |
| N/A | 3007-73040 | 7.3 mm x 40 mm Cannulated Screw, Full Thread |
| N/A | 3007-73045 | 7.3 mm x 45 mm Cannulated Screw, Full Thread |
| N/A | 3007-73050 | 7.3 mm x 50 mm Cannulated Screw, Full Thread |
| N/A | 3007-73055 | 7.3 mm x 55 mm Cannulated Screw, Full Thread |
| N/A | 3007-73060 | 7.3 mm x 60 mm Cannulated Screw, Full Thread |
| N/A | 3007-73065 | 7.3 mm x 65 mm Cannulated Screw, Full Thread |
| N/A | 3007-73070 | 7.3 mm x 70 mm Cannulated Screw, Full Thread |
| N/A | 3007-73075 | 7.3 mm x 75 mm Cannulated Screw, Full Thread |
| N/A | 3007-73080 | 7.3 mm x 80 mm Cannulated Screw, Full Thread |
| N/A | 3007-73085 | 7.3 mm x 85 mm Cannulated Screw, Full Thread |
| N/A | 3007-73090 | 7.3 mm x 90 mm Cannulated Screw, Full Thread |
| N/A | 3007-73095 | 7.3 mm x 95 mm Cannulated Screw, Full Thread |
| N/A | 3007-73100 | 7.3 mm x 100 mm Cannulated Screw, Full Thread |
| N/A | 3007-73105 | 7.3 mm x 105 mm Cannulated Screw, Full Thread |
| N/A | 3007-73110 | 7.3 mm x 110 mm Cannulated Screw, Full Thread |
| N/A | 3007-73115 | 7.3 mm x 115 mm Cannulated Screw, Full Thread |
| N/A | 3007-73120 | 7.3 mm x 120 mm Cannulated Screw, Full Thread |
| N/A | 3007-73125 | 7.3 mm x 125 mm Cannulated Screw, Full Thread |
| N/A | 3007-73130 | 7.3 mm x 130 mm Cannulated Screw, Full Thread |
| N/A | 3007-73135 | 7.3 mm x 135 mm Cannulated Screw, Full Thread |
| N/A | 3007-73140 | 7.3 mm x 140 mm Cannulated Screw, Full Thread |
| N/A | 3007-73145 | 7.3 mm x 145 mm Cannulated Screw, Full Thread |
| N/A | 3007-73150 | 7.3 mm x 150 mm Cannulated Screw, Full Thread |
| 7003-07036-S | 7003-07036 | Cannulated Screw Washer 7.0 mm OD x 3.6 mm ID |
| N/A | 7003-13067 | Cannulated Screw Washer 13.0 mm OD x 6.7 mm ID |

Screw and Pin Portfolio [continued]

Extremity Screw System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|---|
| 27-0120-S | 27-0120 | 2.7 mm x 12.0 mm Cannulated Extremity Screw |
| 27-0140-S | 27-0140 | 2.7 mm x 14.0 mm Cannulated Extremity Screw |
| 27-0160-S | 27-0160 | 2.7 mm x 16.0 mm Cannulated Extremity Screw |
| 27-0180-S | 27-0180 | 2.7 mm x 18.0 mm Cannulated Extremity Screw |
| 27-0200-S | 27-0200 | 2.7 mm x 20.0 mm Cannulated Extremity Screw |
| 27-0220-S | 27-0220 | 2.7 mm x 22.0 mm Cannulated Extremity Screw |
| 27-0240-S | 27-0240 | 2.7 mm x 24.0 mm Cannulated Extremity Screw |
| FA-0140-S | FA-0140 | 3.5 mm x 14.0 mm Cannulated Extremity Screw |
| FA-0160-S | FA-0160 | 3.5 mm x 16.0 mm Cannulated Extremity Screw |
| FA-0180-S | FA-0180 | 3.5 mm x 18.0 mm Cannulated Extremity Screw |
| FA-0200-S | FA-0200 | 3.5 mm x 20.0 mm Cannulated Extremity Screw |
| FA-0220-S | FA-0220 | 3.5 mm x 22.0 mm Cannulated Extremity Screw |
| FA-0240-S | FA-0240 | 3.5 mm x 24.0 mm Cannulated Extremity Screw |
| FA-0260-S | FA-0260 | 3.5 mm x 26.0 mm Cannulated Extremity Screw |
| FA-0280-S | FA-0280 | 3.5 mm x 28.0 mm Cannulated Extremity Screw |
| FA-0300-S | FA-0300 | 3.5 mm x 30.0 mm Cannulated Extremity Screw |
| FA-0320-S | FA-0320 | 3.5 mm x 32.0 mm Cannulated Extremity Screw |
| FA-0340-S | FA-0340 | 3.5 mm x 34.0 mm Cannulated Extremity Screw |
| FA-0360-S | FA-0360 | 3.5 mm x 36.0 mm Cannulated Extremity Screw |
| FA-0380-S | FA-0380 | 3.5 mm x 38.0 mm Cannulated Extremity Screw |
| FA-0400-S | FA-0400 | 3.5 mm x 40.0 mm Cannulated Extremity Screw |
| FA-0420-S | FA-0420 | 3.5 mm x 42.0 mm Cannulated Extremity Screw |
| FA-4014-S | FA-4014 | 4.0 mm x 14.0 mm Cannulated Extremity Screw |
| FA-4016-S | FA-4016 | 4.0 mm x 16.0 mm Cannulated Extremity Screw |
| FA-4018-S | FA-4018 | 4.0 mm x 18.0 mm Cannulated Extremity Screw |
| FA-4020-S | FA-4020 | 4.0 mm x 20.0 mm Cannulated Extremity Screw |
| FA-4022-S | FA-4022 | 4.0 mm x 22.0 mm Cannulated Extremity Screw |
| FA-4024-S | FA-4024 | 4.0 mm x 24.0 mm Cannulated Extremity Screw |
| FA-4026-S | FA-4026 | 4.0 mm x 26.0 mm Cannulated Extremity Screw |

Screw and Pin Portfolio [continued]

Extremity Screw System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|---|
| FA-4028-S | FA-4028 | 4.0 mm x 28.0 mm Cannulated Extremity Screw |
| FA-4030-S | FA-4030 | 4.0 mm x 30.0 mm Cannulated Extremity Screw |
| FA-4032-S | FA-4032 | 4.0 mm x 32.0 mm Cannulated Extremity Screw |
| FA-4034-S | FA-4034 | 4.0 mm x 34.0 mm Cannulated Extremity Screw |
| FA-4036-S | FA-4036 | 4.0 mm x 36.0 mm Cannulated Extremity Screw |
| FA-4038-S | FA-4038 | 4.0 mm x 38.0 mm Cannulated Extremity Screw |
| FA-4040-S | FA-4040 | 4.0 mm x 40.0 mm Cannulated Extremity Screw |
| FA-4042-S | FA-4042 | 4.0 mm x 42.0 mm Cannulated Extremity Screw |

Tension Band Pin System / Tension Band Pin System 2

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|-------------------------------|
| N/A | 30-0097 | 50.0 mm Tension Band Pin |
| N/A | 30-0098 | 70.0 mm Tension Band Pin |
| N/A | 30-0099 | 90.0 mm Tension Band Pin |
| N/A | 30-0166 | 2 mm x 50 mm Tension Band Pin |
| N/A | 30-0167 | 2 mm x 70 mm Tension Band Pin |
| N/A | 30-0168 | 2 mm x 90 mm Tension Band Pin |

Screw and Pin Portfolio - ExsoMed™

- ▶ ArcPhix™ Functional Flexion Compression Screw System
- ▶ InFrame™ Intramedullary Micro Nail System
- ▶ INnate™ Intramedullary Threaded Nail System
- ▶ NanoPhix™ Cannulated Lag Screw



Conditions for Safe Scanning in the MRI

A patient with the Acumed ExsoMed product line may be safely scanned at 1.5 T or 3 T under the following conditions. Failure to follow these conditions may result in injury to the patient.

Acumed ExsoMed™ Product Line

| Parameter | Condition |
|---|---|
| Nominal values(s) of Static Magnetic Field | 1.5 T or 3 T |
| Maximum Spatial Field Gradient | 30 T/m (3000 gauss/cm) |
| RF Excitation | Circularly Polarized (CP) |
| RF Transmit Coil Type | There are no transmit coil restrictions |
| Operating Mode | Normal operating mode |
| Maximum Whole Body SAR | 2.0 W/kg (Normal operating mode) |
| Maximum Head SAR | 3.2 W/kg (Normal operating mode) |
| Limits on Scan Duration | 2.0 W/kg whole body average SAR for 60 minutes of continuous RF (a sequence or back to back series/scan without breaks) |
| MR Image Artifact | The presence of this implant may produce an image artifact of 32 mm. |
| If more information about a specific parameter is not included, there are no conditions associated with that parameter. | |

Screw and Pin Portfolio - ExsoMed™ [continued]

ArcPhix™ Functional Flexion Compression Screw System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|-------------------------------|
| EXARC903028 | N/A | ArcPhix System 3.0 mm x 28 mm |

InFrame™ Intramedullary Micro Nail System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|--------------------------------|
| EXINF922012 | N/A | InFrame Implant 2.0 mm x 12 mm |
| EXINF922014 | N/A | InFrame Implant 2.0 mm x 14 mm |
| EXINF922016 | N/A | InFrame Implant 2.0 mm x 16 mm |
| EXINF922018 | N/A | InFrame Implant 2.0 mm x 18 mm |
| EXINF922020 | N/A | InFrame Implant 2.0 mm x 20 mm |
| EXINF922022 | N/A | InFrame Implant 2.0 mm x 22 mm |
| EXINF922024 | N/A | InFrame Implant 2.0 mm x 24 mm |
| EXINF922026 | N/A | InFrame Implant 2.0 mm x 26 mm |
| EXINF922028 | N/A | InFrame Implant 2.0 mm x 28 mm |
| EXINF922030 | N/A | InFrame Implant 2.0 mm x 30 mm |
| EXINF922032 | N/A | InFrame Implant 2.0 mm x 32 mm |
| EXINF922034 | N/A | InFrame Implant 2.0 mm x 34 mm |
| EXINF922036 | N/A | InFrame Implant 2.0 mm x 36 mm |
| EXINF922038 | N/A | InFrame Implant 2.0 mm x 38 mm |
| EXINF922040 | N/A | InFrame Implant 2.0 mm x 40 mm |
| EXINF922042 | N/A | InFrame Implant 2.0 mm x 42 mm |
| EXINF922044 | N/A | InFrame Implant 2.0 mm x 44 mm |
| EXINF922046 | N/A | InFrame Implant 2.0 mm x 46 mm |
| EXINF922048 | N/A | InFrame Implant 2.0 mm x 48 mm |

INnate™ Intramedullary Threaded Nail System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|-------------------------------|
| EXINN923625 | N/A | INnate Implant 3.6 mm x 25 mm |
| EXINN923630 | N/A | INnate Implant 3.6 mm x 30 mm |

Screw and Pin Portfolio - ExsoMed™ [continued]

INnate™ Intramedullary Threaded Nail System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|-------------------------------|
| EXINN923635 | N/A | INnate Implant 3.6 mm x 35 mm |
| EXINN923640 | N/A | INnate Implant 3.6 mm x 40 mm |
| EXINN923645 | N/A | INnate Implant 3.6 mm x 45 mm |
| EXINN923650 | N/A | INnate Implant 3.6 mm x 50 mm |
| EXINN923655 | N/A | INnate Implant 3.6 mm x 55 mm |
| EXINN924535 | N/A | INnate Implant 4.5 mm x 35 mm |
| EXINN924540 | N/A | INnate Implant 4.5 mm x 40 mm |
| EXINN924545 | N/A | INnate Implant 4.5 mm x 45 mm |
| EXINN924550 | N/A | INnate Implant 4.5 mm x 50 mm |
| EXINN924555 | N/A | INnate Implant 4.5 mm x 55 mm |
| EXINN924565 | N/A | INnate Implant 4.5 mm x 65 mm |
| EXINN924575 | N/A | INnate Implant 4.5 mm x 75 mm |

NanoPhix™ Cannulated Lag Screw

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|---------------------------------|
| EXINN921506 | N/A | NanoPhix Implant 1.5 mm x 6 mm |
| EXINN921507 | N/A | NanoPhix Implant 1.5 mm x 7 mm |
| EXINN921508 | N/A | NanoPhix Implant 1.5 mm x 8 mm |
| EXINN921509 | N/A | NanoPhix Implant 1.5 mm x 9 mm |
| EXINN921510 | N/A | NanoPhix Implant 1.5 mm x 10 mm |
| EXINN921511 | N/A | NanoPhix Implant 1.5 mm x 11 mm |
| EXINN921512 | N/A | NanoPhix Implant 1.5 mm x 12 mm |
| EXINN921513 | N/A | NanoPhix Implant 1.5 mm x 13 mm |
| EXINN921514 | N/A | NanoPhix Implant 1.5 mm x 14 mm |
| EXINN921515 | N/A | NanoPhix Implant 1.5 mm x 15 mm |
| EXINN921516 | N/A | NanoPhix Implant 1.5 mm x 16 mm |
| EXINN921517 | N/A | NanoPhix Implant 1.5 mm x 17 mm |
| EXINN921518 | N/A | NanoPhix Implant 1.5 mm x 18 mm |
| EXINN921519 | N/A | NanoPhix Implant 1.5 mm x 19 mm |

Screw and Pin Portfolio - ExsoMed™ [continued]

NanoPhix™ Cannulated Lag Screw

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|---------------------------------|
| EXINN921520 | N/A | NanoPhix Implant 1.5 mm x 20 mm |
| EXINN921521 | N/A | NanoPhix Implant 1.5 mm x 21 mm |
| EXINN921522 | N/A | NanoPhix Implant 1.5 mm x 22 mm |
| EXINN921523 | N/A | NanoPhix Implant 1.5 mm x 23 mm |
| EXINN921524 | N/A | NanoPhix Implant 1.5 mm x 24 mm |
| EXINN921525 | N/A | NanoPhix Implant 1.5 mm x 25 mm |
| EXINN921526 | N/A | NanoPhix Implant 1.5 mm x 26 mm |

Shoulder and Elbow Portfolio

- ▶ Elbow Plating System
- ▶ Congruent Elbow Plating System
- ▶ Radial Head Plating System
- ▶ Clavicle Plating System
- ▶ Clavicle Hook Plating System
- ▶ Locking Clavicle Plating System
- ▶ Acu-Sinch® Repair System
- ▶ Dual-Trak Clavicle Screw System
- ▶ Polarus® 3 Solution (plates)
- ▶ Polarus Proximal Humeral Plating System
- ▶ Scapula Plating System



Conditions for Safe Scanning in the MRI

A patient with the screw, plate, pin, suture anchor, suture retainer, suture, and washer implants from the Elbow Plating, Congruent Elbow Plating, Radial Head Plating, Clavicle Plating, Locking Clavicle Plating, Clavicle Hook Plating, Acu-Sinch Repair, Dual-Trak Clavicle Screw, Polarus 3 Solution, Polarus Proximal Humeral Plating, and Scapula Plating Systems may be safely scanned at 1.5 T or 3 T under the following conditions. Failure to follow these conditions may result in injury to the patient.

Shoulder and Elbow Portfolio

| Parameter | Condition |
|---|---|
| Nominal values(s) of Static Magnetic Field | 1.5 T or 3 T |
| Maximum Spatial Field Gradient | 30 T/m (3000 gauss/cm) |
| RF Excitation | Circularly Polarized (CP) |
| RF Transmit Coil Type | Whole body transmit coil, Head RF transmit-receive coil |
| Maximum Whole Body SAR | 1.0 W/kg* |
| Limits on Scan Duration | 1.0 W/kg whole body average SAR for 20 minutes of continuous RF (a sequence or back to back series/scan without breaks) followed by a 3 minute cooling period for a one hour session. |
| MR Image Artifact | The presence of this implant may produce an image artifact of 37 mm. |
| If more information about a specific parameter is not included, there are no conditions associated with that parameter. | |

* Whole body SAR of 1.0 W/kg may limit the availability and diagnostic quality of some MR procedures

Shoulder and Elbow Portfolio [continued]

Elbow Plating System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|--|
| N/A | 30-0098 | 70.0 mm Tension Band Pin |
| N/A | 30-0099 | 90.0 mm Tension Band Pin |
| 30-0232-S | 30-0232 | 3.5 mm x 8 mm Locking Hexalobe Screw |
| 30-0233-S | 30-0233 | 3.5 mm x 10 mm Locking Hexalobe Screw |
| 30-0234-S | 30-0234 | 3.5 mm x 12 mm Locking Hexalobe Screw |
| 30-0235-S | 30-0235 | 3.5 mm x 14 mm Locking Hexalobe Screw |
| 30-0236-S | 30-0236 | 3.5 mm x 16 mm Locking Hexalobe Screw |
| 30-0237-S | 30-0237 | 3.5 mm x 18 mm Locking Hexalobe Screw |
| 30-0238-S | 30-0238 | 3.5 mm x 20 mm Locking Hexalobe Screw |
| 30-0239-S | 30-0239 | 3.5 mm x 22 mm Locking Hexalobe Screw |
| 30-0240-S | 30-0240 | 3.5 mm x 24 mm Locking Hexalobe Screw |
| 30-0241-S | 30-0241 | 3.5 mm x 26 mm Locking Hexalobe Screw |
| 30-0242-S | 30-0242 | 3.5 mm x 28 mm Locking Hexalobe Screw |
| 30-0243-S | 30-0243 | 3.5 mm x 30 mm Locking Hexalobe Screw |
| 30-0244-S | 30-0244 | 3.5 mm x 32 mm Locking Hexalobe Screw |
| 30-0245-S | 30-0245 | 3.5 mm x 34 mm Locking Hexalobe Screw |
| 30-0246-S | 30-0246 | 3.5 mm x 36 mm Locking Hexalobe Screw |
| 30-0247-S | 30-0247 | 3.5 mm x 38 mm Locking Hexalobe Screw |
| 30-0248-S | 30-0248 | 3.5 mm x 40 mm Locking Hexalobe Screw |
| 30-0249-S | 30-0249 | 3.5 mm x 45 mm Locking Hexalobe Screw |
| 30-0250-S | 30-0250 | 3.5 mm x 50 mm Locking Hexalobe Screw |
| 30-0251-S | 30-0251 | 3.5 mm x 55 mm Locking Hexalobe Screw |
| 30-0252-S | 30-0252 | 3.5 mm x 60 mm Locking Hexalobe Screw |
| 30-0255-S | 30-0255 | 3.5 mm x 8 mm Nonlocking Hexalobe Screw |
| 30-0256-S | 30-0256 | 3.5 mm x 10 mm Nonlocking Hexalobe Screw |
| 30-0257-S | 30-0257 | 3.5 mm x 12 mm Nonlocking Hexalobe Screw |
| 30-0258-S | 30-0258 | 3.5 mm x 14 mm Nonlocking Hexalobe Screw |
| 30-0259-S | 30-0259 | 3.5 mm x 16 mm Nonlocking Hexalobe Screw |
| 30-0260-S | 30-0260 | 3.5 mm x 18 mm Nonlocking Hexalobe Screw |

Shoulder and Elbow Portfolio [continued]

Elbow Plating System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|--|
| 30-0261-S | 30-0261 | 3.5 mm x 20 mm Nonlocking Hexalobe Screw |
| 30-0262-S | 30-0262 | 3.5 mm x 22 mm Nonlocking Hexalobe Screw |
| 30-0263-S | 30-0263 | 3.5 mm x 24 mm Nonlocking Hexalobe Screw |
| 30-0264-S | 30-0264 | 3.5 mm x 26 mm Nonlocking Hexalobe Screw |
| 30-0265-S | 30-0265 | 3.5 mm x 28 mm Nonlocking Hexalobe Screw |
| 30-0266-S | 30-0266 | 3.5 mm x 30 mm Nonlocking Hexalobe Screw |
| 30-0267-S | 30-0267 | 3.5 mm x 32 mm Nonlocking Hexalobe Screw |
| 30-0268-S | 30-0268 | 3.5 mm x 34 mm Nonlocking Hexalobe Screw |
| 30-0269-S | 30-0269 | 3.5 mm x 36 mm Nonlocking Hexalobe Screw |
| 30-0270-S | 30-0270 | 3.5 mm x 38 mm Nonlocking Hexalobe Screw |
| 30-0271-S | 30-0271 | 3.5 mm x 40 mm Nonlocking Hexalobe Screw |
| 30-0272-S | 30-0272 | 3.5 mm x 45 mm Nonlocking Hexalobe Screw |
| 30-0273-S | 30-0273 | 3.5 mm x 50 mm Nonlocking Hexalobe Screw |
| 30-0274-S | 30-0274 | 3.5 mm x 55 mm Nonlocking Hexalobe Screw |
| 30-0275-S | 30-0275 | 3.5 mm x 60 mm Nonlocking Hexalobe Screw |
| 30-0276-S | 30-0276 | 3.5 mm x 65 mm Nonlocking Hexalobe Screw |
| 30-0278-S | 30-0278 | 3.0 mm x 8 mm Locking Hexalobe Screw |
| 30-0279-S | 30-0279 | 3.0 mm x 10 mm Locking Hexalobe Screw |
| 30-0280-S | 30-0280 | 3.0 mm x 12 mm Locking Hexalobe Screw |
| 30-0281-S | 30-0281 | 3.0 mm x 14 mm Locking Hexalobe Screw |
| 30-0282-S | 30-0282 | 3.0 mm x 16 mm Locking Hexalobe Screw |
| 30-0283-S | 30-0283 | 3.0 mm x 18 mm Locking Hexalobe Screw |
| 30-0284-S | 30-0284 | 3.0 mm x 20 mm Locking Hexalobe Screw |
| 30-0285-S | 30-0285 | 3.0 mm x 22 mm Locking Hexalobe Screw |
| 30-0286-S | 30-0286 | 3.0 mm x 24 mm Locking Hexalobe Screw |
| 30-0287-S | 30-0287 | 3.0 mm x 26 mm Locking Hexalobe Screw |
| 30-0288-S | 30-0288 | 3.0 mm x 28 mm Locking Hexalobe Screw |
| 30-0289-S | 30-0289 | 3.0 mm x 30 mm Locking Hexalobe Screw |
| 30-0290-S | 30-0290 | 3.0 mm x 32 mm Locking Hexalobe Screw |

Shoulder and Elbow Portfolio [continued]

Elbow Plating System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|--|
| 30-0291-S | 30-0291 | 3.0 mm x 34 mm Locking Hexalobe Screw |
| 30-0292-S | 30-0292 | 3.0 mm x 36 mm Locking Hexalobe Screw |
| 30-0293-S | 30-0293 | 3.0 mm x 38 mm Locking Hexalobe Screw |
| 30-0294-S | 30-0294 | 3.0 mm x 40 mm Locking Hexalobe Screw |
| 30-0295-S | 30-0295 | 3.0 mm x 45 mm Locking Hexalobe Screw |
| 30-0296-S | 30-0296 | 3.0 mm x 50 mm Locking Hexalobe Screw |
| 30-0297-S | 30-0297 | 3.0 mm x 55 mm Locking Hexalobe Screw |
| 30-0298-S | 30-0298 | 3.0 mm x 60 mm Locking Hexalobe Screw |
| 30-0301-S | 30-0301 | 3.0 mm x 8 mm Nonlocking Hexalobe Screw |
| 30-0302-S | 30-0302 | 3.0 mm x 10 mm Nonlocking Hexalobe Screw |
| 30-0303-S | 30-0303 | 3.0 mm x 12 mm Nonlocking Hexalobe Screw |
| 30-0304-S | 30-0304 | 3.0 mm x 14 mm Nonlocking Hexalobe Screw |
| 30-0305-S | 30-0305 | 3.0 mm x 16 mm Nonlocking Hexalobe Screw |
| 30-0306-S | 30-0306 | 3.0 mm x 18 mm Nonlocking Hexalobe Screw |
| 30-0307-S | 30-0307 | 3.0 mm x 20 mm Nonlocking Hexalobe Screw |
| 30-0308-S | 30-0308 | 3.0 mm x 22 mm Nonlocking Hexalobe Screw |
| 30-0309-S | 30-0309 | 3.0 mm x 24 mm Nonlocking Hexalobe Screw |
| 30-0310-S | 30-0310 | 3.0 mm x 26 mm Nonlocking Hexalobe Screw |
| 30-0311-S | 30-0311 | 3.0 mm x 28 mm Nonlocking Hexalobe Screw |
| 30-0312-S | 30-0312 | 3.0 mm x 30 mm Nonlocking Hexalobe Screw |
| 30-0313-S | 30-0313 | 3.0 mm x 32 mm Nonlocking Hexalobe Screw |
| 30-0314-S | 30-0314 | 3.0 mm x 34 mm Nonlocking Hexalobe Screw |
| 30-0315-S | 30-0315 | 3.0 mm x 36 mm Nonlocking Hexalobe Screw |
| 30-0316-S | 30-0316 | 3.0 mm x 38 mm Nonlocking Hexalobe Screw |
| 30-0317-S | 30-0317 | 3.0 mm x 40 mm Nonlocking Hexalobe Screw |
| 30-0318-S | 30-0318 | 3.0 mm x 45 mm Nonlocking Hexalobe Screw |
| 30-0319-S | 30-0319 | 3.0 mm x 50 mm Nonlocking Hexalobe Screw |
| 30-0320-S | 30-0320 | 3.0 mm x 55 mm Nonlocking Hexalobe Screw |
| 30-0321-S | 30-0321 | 3.0 mm x 60 mm Nonlocking Hexalobe Screw |

Shoulder and Elbow Portfolio [continued]

Elbow Plating System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|--|
| 30-0322-S | 30-0322 | 3.0 mm x 65 mm Nonlocking Hexalobe Screw |
| 30-0324-S | 30-0324 | 2.7 mm x 8 mm Locking Hexalobe Screw |
| 30-0325-S | 30-0325 | 2.7 mm x 10 mm Locking Hexalobe Screw |
| 30-0326-S | 30-0326 | 2.7 mm x 12 mm Locking Hexalobe Screw |
| 30-0327-S | 30-0327 | 2.7 mm x 14 mm Locking Hexalobe Screw |
| 30-0328-S | 30-0328 | 2.7 mm x 16 mm Locking Hexalobe Screw |
| 30-0329-S | 30-0329 | 2.7 mm x 18 mm Locking Hexalobe Screw |
| 30-0330-S | 30-0330 | 2.7 mm x 20 mm Locking Hexalobe Screw |
| 30-0331-S | 30-0331 | 2.7 mm x 22 mm Locking Hexalobe Screw |
| 30-0332-S | 30-0332 | 2.7 mm x 24 mm Locking Hexalobe Screw |
| 30-0333-S | 30-0333 | 2.7 mm x 26 mm Locking Hexalobe Screw |
| 30-0334-S | 30-0334 | 2.7 mm x 28 mm Locking Hexalobe Screw |
| 30-0335-S | 30-0335 | 2.7 mm x 30 mm Locking Hexalobe Screw |
| 30-0336-S | 30-0336 | 2.7 mm x 32 mm Locking Hexalobe Screw |
| 30-0343-S | 30-0343 | 2.7 mm x 8 mm Nonlocking Hexalobe Screw |
| 30-0344-S | 30-0344 | 2.7 mm x 10 mm Nonlocking Hexalobe Screw |
| 30-0345-S | 30-0345 | 2.7 mm x 12 mm Nonlocking Hexalobe Screw |
| 30-0346-S | 30-0346 | 2.7 mm x 14 mm Nonlocking Hexalobe Screw |
| 30-0347-S | 30-0347 | 2.7 mm x 16 mm Nonlocking Hexalobe Screw |
| 30-0348-S | 30-0348 | 2.7 mm x 18 mm Nonlocking Hexalobe Screw |
| 30-0349-S | 30-0349 | 2.7 mm x 20 mm Nonlocking Hexalobe Screw |
| 30-0350-S | 30-0350 | 2.7 mm x 22 mm Nonlocking Hexalobe Screw |
| 30-0351-S | 30-0351 | 2.7 mm x 24 mm Nonlocking Hexalobe Screw |
| 30-0352-S | 30-0352 | 2.7 mm x 26 mm Nonlocking Hexalobe Screw |
| 30-0353-S | 30-0353 | 2.7 mm x 28 mm Nonlocking Hexalobe Screw |
| 30-0354-S | 30-0354 | 2.7 mm x 30 mm Nonlocking Hexalobe Screw |
| 30-0355-S | 30-0355 | 2.7 mm x 32 mm Nonlocking Hexalobe Screw |
| N/A | 3060-27008 | 2.7 mm x 8 mm Hexalobe Locking CoCr Screw |
| N/A | 3060-27010 | 2.7 mm x 10 mm Hexalobe Locking CoCr Screw |

Shoulder and Elbow Portfolio [continued]

Elbow Plating System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|--|
| N/A | 3060-27012 | 2.7 mm x 12 mm Hexalobe Locking CoCr Screw |
| N/A | 3060-27014 | 2.7 mm x 14 mm Hexalobe Locking CoCr Screw |
| N/A | 3060-27016 | 2.7 mm x 16 mm Hexalobe Locking CoCr Screw |
| N/A | 3060-27018 | 2.7 mm x 18 mm Hexalobe Locking CoCr Screw |
| N/A | 3060-27020 | 2.7 mm x 20 mm Hexalobe Locking CoCr Screw |
| N/A | 3060-27022 | 2.7 mm x 22 mm Hexalobe Locking CoCr Screw |
| N/A | 3060-27024 | 2.7 mm x 24 mm Hexalobe Locking CoCr Screw |
| N/A | 3060-27026 | 2.7 mm x 26 mm Hexalobe Locking CoCr Screw |
| N/A | 3060-27028 | 2.7 mm x 28 mm Hexalobe Locking CoCr Screw |
| N/A | 3060-27030 | 2.7 mm x 30 mm Hexalobe Locking CoCr Screw |
| N/A | 3060-27032 | 2.7 mm x 32 mm Hexalobe Locking CoCr Screw |
| N/A | 3061-27008 | 2.7 mm x 8 mm Hexalobe Non-Locking CoCr Screw |
| N/A | 3061-27010 | 2.7 mm x 10 mm Hexalobe Non-Locking CoCr Screw |
| N/A | 3061-27012 | 2.7 mm x 12 mm Hexalobe Non-Locking CoCr Screw |
| N/A | 3061-27014 | 2.7 mm x 14 mm Hexalobe Non-Locking CoCr Screw |
| N/A | 3061-27016 | 2.7 mm x 16 mm Hexalobe Non-Locking CoCr Screw |
| N/A | 3061-27018 | 2.7 mm x 18 mm Hexalobe Non-Locking CoCr Screw |
| N/A | 3061-27020 | 2.7 mm x 20 mm Hexalobe Non-Locking CoCr Screw |
| N/A | 3061-27022 | 2.7 mm x 22 mm Hexalobe Non-Locking CoCr Screw |
| N/A | 3061-27024 | 2.7 mm x 24 mm Hexalobe Non-Locking CoCr Screw |
| N/A | 3061-27026 | 2.7 mm x 26 mm Hexalobe Non-Locking CoCr Screw |
| N/A | 3061-27028 | 2.7 mm x 28 mm Hexalobe Non-Locking CoCr Screw |
| N/A | 3061-27030 | 2.7 mm x 30 mm Hexalobe Non-Locking CoCr Screw |
| N/A | 3061-27032 | 2.7 mm x 32 mm Hexalobe Non-Locking CoCr Screw |
| N/A | 3062-35008 | 3.5 mm x 8 mm Hexalobe Locking CoCr Screw |
| N/A | 3062-35010 | 3.5 mm x 10 mm Hexalobe Locking CoCr Screw |
| N/A | 3062-35012 | 3.5 mm x 12 mm Hexalobe Locking CoCr Screw |
| N/A | 3062-35014 | 3.5 mm x 14 mm Hexalobe Locking CoCr Screw |

Shoulder and Elbow Portfolio [continued]

Elbow Plating System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|--|
| N/A | 3062-35016 | 3.5 mm x 16 mm Hexalobe Locking CoCr Screw |
| N/A | 3062-35018 | 3.5 mm x 18 mm Hexalobe Locking CoCr Screw |
| N/A | 3062-35020 | 3.5 mm x 20 mm Hexalobe Locking CoCr Screw |
| N/A | 3062-35022 | 3.5 mm x 22 mm Hexalobe Locking CoCr Screw |
| N/A | 3062-35024 | 3.5 mm x 24 mm Hexalobe Locking CoCr Screw |
| N/A | 3062-35026 | 3.5 mm x 26 mm Hexalobe Locking CoCr Screw |
| N/A | 3062-35028 | 3.5 mm x 28 mm Hexalobe Locking CoCr Screw |
| N/A | 3062-35030 | 3.5 mm x 30 mm Hexalobe Locking CoCr Screw |
| N/A | 3062-35032 | 3.5 mm x 32 mm Hexalobe Locking CoCr Screw |
| N/A | 3062-35034 | 3.5 mm x 34 mm Hexalobe Locking CoCr Screw |
| N/A | 3062-35036 | 3.5 mm x 36 mm Hexalobe Locking CoCr Screw |
| N/A | 3062-35038 | 3.5 mm x 38 mm Hexalobe Locking CoCr Screw |
| N/A | 3062-35040 | 3.5 mm x 40 mm Hexalobe Locking CoCr Screw |
| N/A | 3062-35045 | 3.5 mm x 45 mm Hexalobe Locking CoCr Screw |
| N/A | 3062-35050 | 3.5 mm x 50 mm Hexalobe Locking CoCr Screw |
| N/A | 3062-35055 | 3.5 mm x 55 mm Hexalobe Locking CoCr Screw |
| N/A | 3062-35060 | 3.5 mm x 60 mm Hexalobe Locking CoCr Screw |
| N/A | 3063-35008 | 3.5 mm x 8 mm Hexalobe Non-Locking CoCr Screw |
| N/A | 3063-35010 | 3.5 mm x 10 mm Hexalobe Non-Locking CoCr Screw |
| N/A | 3063-35012 | 3.5 mm x 12 mm Hexalobe Non-Locking CoCr Screw |
| N/A | 3063-35014 | 3.5 mm x 14 mm Hexalobe Non-Locking CoCr Screw |
| N/A | 3063-35016 | 3.5 mm x 16 mm Hexalobe Non-Locking CoCr Screw |
| N/A | 3063-35018 | 3.5 mm x 18 mm Hexalobe Non-Locking CoCr Screw |
| N/A | 3063-35020 | 3.5 mm x 20 mm Hexalobe Non-Locking CoCr Screw |
| N/A | 3063-35022 | 3.5 mm x 22 mm Hexalobe Non-Locking CoCr Screw |
| N/A | 3063-35024 | 3.5 mm x 24 mm Hexalobe Non-Locking CoCr Screw |
| N/A | 3063-35026 | 3.5 mm x 26 mm Hexalobe Non-Locking CoCr Screw |
| N/A | 3063-35028 | 3.5 mm x 28 mm Hexalobe Non-Locking CoCr Screw |
| N/A | 3063-35030 | 3.5 mm x 30 mm Hexalobe Non-Locking CoCr Screw |

Shoulder and Elbow Portfolio [continued]

Elbow Plating System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|---|
| N/A | 3063-35032 | 3.5 mm x 32 mm Hexalobe Non-Locking CoCr Screw |
| N/A | 3063-35034 | 3.5 mm x 34 mm Hexalobe Non-Locking CoCr Screw |
| N/A | 3063-35036 | 3.5 mm x 36 mm Hexalobe Non-Locking CoCr Screw |
| N/A | 3063-35038 | 3.5 mm x 38 mm Hexalobe Non-Locking CoCr Screw |
| N/A | 3063-35040 | 3.5 mm x 40 mm Hexalobe Non-Locking CoCr Screw |
| N/A | 3063-35045 | 3.5 mm x 45 mm Hexalobe Non-Locking CoCr Screw |
| N/A | 3063-35050 | 3.5 mm x 50 mm Hexalobe Non-Locking CoCr Screw |
| N/A | 3063-35055 | 3.5 mm x 55 mm Hexalobe Non-Locking CoCr Screw |
| N/A | 3063-35060 | 3.5 mm x 60 mm Hexalobe Non-Locking CoCr Screw |
| N/A | 3063-35065 | 3.5 mm x 65 mm Hexalobe Non-Locking CoCr Screw |
| 70-0302-S | 70-0302 | Olecranon Plate, Standard 3-Hole, Left |
| 70-0303-S | 70-0303 | Olecranon Plate, Standard 3-Hole, Right |
| 70-0304-S | 70-0304 | Olecranon Plate, Standard 5-Hole, Left |
| 70-0305-S | 70-0305 | Olecranon Plate, Standard 5-Hole, Right |
| 70-0306-S | 70-0306 | Olecranon Plate, Standard 7-Hole, Left |
| 70-0307-S | 70-0307 | Olecranon Plate, Standard 7-Hole, Right |
| 70-0308-S | 70-0308 | Olecranon Plate, Standard 11-Hole, Left |
| 70-0309-S | 70-0309 | Olecranon Plate, Standard 11-Hole, Right |
| 70-0310-S | 70-0310 | Olecranon Plate, Standard 15-Hole, Left |
| 70-0311-S | 70-0311 | Olecranon Plate, Standard 15-Hole, Right |
| 70-0312-S | 70-0312 | Olecranon Plate, Extended 5-Hole, Left |
| 70-0313-S | 70-0313 | Olecranon Plate, Extended 5-Hole, Right |
| 70-0314-S | 70-0314 | Olecranon Plate, Extended 9-Hole, Left |
| 70-0315-S | 70-0315 | Olecranon Plate, Extended 9-Hole, Right |
| 70-0316-S | 70-0316 | Olecranon Plate, Narrow 5-Hole, Left |
| 70-0317-S | 70-0317 | Olecranon Plate, Narrow 5-Hole, Right |
| 70-0374-S | 70-0374 | Posterolateral Distal Humerus Plate 5-Hole, Left |
| 70-0375-S | 70-0375 | Posterolateral Distal Humerus Plate 5-Hole, Right |
| 70-0376-S | 70-0376 | Posterolateral Distal Humerus Plate 7-Hole, Left |

Shoulder and Elbow Portfolio [continued]

Elbow Plating System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|--|
| 70-0377-S | 70-0377 | Posterolateral Distal Humerus Plate 7-Hole, Right |
| 70-0378-S | 70-0378 | Posterolateral Distal Humerus Plate 11-Hole, Left |
| 70-0379-S | 70-0379 | Posterolateral Distal Humerus Plate 11-Hole, Right |
| 70-0380-S | 70-0380 | Posterolateral Distal Humerus Plate 15-Hole, Left |
| 70-0381-S | 70-0381 | Posterolateral Distal Humerus Plate 15-Hole, Right |
| 70-0413-S | 70-0413 | Coronoid Plate, Standard, Left |
| 70-0414-S | 70-0414 | Coronoid Plate, Standard, Right |
| 70-0415-S | 70-0415 | Coronoid Plate, Sm, Left |
| 70-0416-S | 70-0416 | Coronoid Plate, Sm, Right |
| PL-LEL10L-S | PL-LEL10L | Left 10-Hole Locking Lateral Plate |
| PL-LEL10R-S | PL-LEL10R | Right 10-Hole Locking Lateral Plate |
| PL-LEL14L-S | PL-LEL14L | Left 14-Hole Locking Lateral Plate |
| PL-LEL14R-S | PL-LEL14R | Right 14-Hole Locking Lateral Plate |
| PL-LEL20L-S | PL-LEL20L | Left 20-Hole Locking Lateral Plate |
| PL-LEL20R-S | PL-LEL20R | Right 20-Hole Locking Lateral Plate |
| PL-LEL6L-S | PL-LEL6L | Left 6-Hole Locking Lateral Plate |
| PL-LEL6R-S | PL-LEL6R | Right 6-Hole Locking Lateral Plate |
| PL-LEM12-S | PL-LEM12 | 12-Hole Locking Medial Plate |
| PL-LEM16-S | PL-LEM16 | 16-Hole Locking Medial Plate |
| PL-LEM7-S | PL-LEM7 | 7-Hole Locking Medial Plate |
| PL-LEM8-S | PL-LEM8 | 8-Hole Locking Medial Plate |
| PL-LEM9L-S | PL-LEM9L | Long 9-Hole Locking Medial Plate |
| PL-LEM9S-S | PL-LEM9S | Short 9-Hole Locking Medial Plate |

Shoulder and Elbow Portfolio [continued]

Congruent Elbow Plating System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|-----------------------------------|
| CA-4120-S | CA-4120 | 4.0 mm x 12.0 mm Cancellous Screw |
| CA-4140-S | CA-4140 | 4.0 mm x 14.0 mm Cancellous Screw |
| CA-4160-S | CA-4160 | 4.0 mm x 16.0 mm Cancellous Screw |
| CA-4180-S | CA-4180 | 4.0 mm x 18.0 mm Cancellous Screw |
| CA-4200-S | CA-4200 | 4.0 mm x 20.0 mm Cancellous Screw |
| CA-4220-S | CA-4220 | 4.0 mm x 22.0 mm Cancellous Screw |
| CA-4240-S | CA-4240 | 4.0 mm x 24.0 mm Cancellous Screw |
| CA-4260-S | CA-4260 | 4.0 mm x 26.0 mm Cancellous Screw |
| CA-4280-S | CA-4280 | 4.0 mm x 28.0 mm Cancellous Screw |
| CA-4300-S | CA-4300 | 4.0 mm x 30.0 mm Cancellous Screw |
| CA-4350-S | CA-4350 | 4.0 mm x 35.0 mm Cancellous Screw |
| CA-4400-S | CA-4400 | 4.0 mm x 40.0 mm Cancellous Screw |
| CA-4450-S | CA-4450 | 4.0 mm x 45.0 mm Cancellous Screw |
| CA-4500-S | CA-4500 | 4.0 mm x 50.0 mm Cancellous Screw |
| CA-4550-S | CA-4550 | 4.0 mm x 55.0 mm Cancellous Screw |
| CA-4600-S | CA-4600 | 4.0 mm x 60.0 mm Cancellous Screw |
| CO-2712-S | CO-2712 | 2.7 mm x 12.0 mm Cortical Screw |
| CO-2714-S | CO-2714 | 2.7 mm x 14.0 mm Cortical Screw |
| CO-2716-S | CO-2716 | 2.7 mm x 16.0 mm Cortical Screw |
| CO-2718-S | CO-2718 | 2.7 mm x 18.0 mm Cortical Screw |
| CO-2720-S | CO-2720 | 2.7 mm x 20.0 mm Cortical Screw |
| CO-2722-S | CO-2722 | 2.7 mm x 22.0 mm Cortical Screw |
| CO-2724-S | CO-2724 | 2.7 mm x 24.0 mm Cortical Screw |
| CO-2726-S | CO-2726 | 2.7 mm x 26.0 mm Cortical Screw |
| CO-2728-S | CO-2728 | 2.7 mm x 28.0 mm Cortical Screw |
| CO-2730-S | CO-2730 | 2.7 mm x 30.0 mm Cortical Screw |
| CO-2732-S | CO-2732 | 2.7 mm x 32.0 mm Cortical Screw |
| CO-2734-S | CO-2734 | 2.7 mm x 34.0 mm Cortical Screw |
| CO-2736-S | CO-2736 | 2.7 mm x 36.0 mm Cortical Screw |
| CO-2738-S | CO-2738 | 2.7 mm x 38.0 mm Cortical Screw |

Shoulder and Elbow Portfolio [continued]

Congruent Elbow Plating System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|---|
| CO-2740-S | CO-2740 | 2.7 mm x 40.0 mm Cortical Screw |
| CO-2745-S | CO-2745 | 2.7 mm x 45.0 mm Cortical Screw |
| CO-2750-S | CO-2750 | 2.7 mm x 50.0 mm Cortical Screw |
| CO-2755-S | CO-2755 | 2.7 mm x 55.0 mm Cortical Screw |
| CO-2760-S | CO-2760 | 2.7 mm x 60.0 mm Cortical Screw |
| CO-2765-S | CO-2765 | 2.7 mm x 65.0 mm Cortical Screw |
| CO-3120-S | CO-3120 | 3.5 mm x 12.0 mm Cortical Screw |
| CO-3140-S | CO-3140 | 3.5 mm x 14.0 mm Cortical Screw |
| CO-3160-S | CO-3160 | 3.5 mm x 16.0 mm Cortical Screw |
| CO-3180-S | CO-3180 | 3.5 mm x 18.0 mm Cortical Screw |
| CO-3200-S | CO-3200 | 3.5 mm x 20.0 mm Cortical Screw |
| CO-3220-S | CO-3220 | 3.5 mm x 22.0 mm Cortical Screw |
| CO-3240-S | CO-3240 | 3.5 mm x 24.0 mm Cortical Screw |
| CO-3260-S | CO-3260 | 3.5 mm x 26.0 mm Cortical Screw |
| CO-3280-S | CO-3280 | 3.5 mm x 28.0 mm Cortical Screw |
| CO-3300-S | CO-3300 | 3.5 mm x 30.0 mm Cortical Screw |
| CO-3320-S | CO-3320 | 3.5 mm x 32.0 mm Cortical Screw |
| CO-3340-S | CO-3340 | 3.5 mm x 34.0 mm Cortical Screw |
| CO-3360-S | CO-3360 | 3.5 mm x 36.0 mm Cortical Screw |
| CO-3380-S | CO-3380 | 3.5 mm x 38.0 mm Cortical Screw |
| CO-3400-S | CO-3400 | 3.5 mm x 40.0 mm Cortical Screw |
| CO-3450-S | CO-3450 | 3.5 mm x 45.0 mm Cortical Screw |
| CO-3500-S | CO-3500 | 3.5 mm x 50.0 mm Cortical Screw |
| CO-3550-S | CO-3550 | 3.5 mm x 55.0 mm Cortical Screw |
| CO-3600-S | CO-3600 | 3.5 mm x 60.0 mm Cortical Screw |
| CO-3650-S | CO-3650 | 3.5 mm x 65.0 mm Cortical Screw |
| COL-2120-S | COL-2120 | 2.7 mm x 12.0 mm Locking Cortical Screw |
| COL-2140-S | COL-2140 | 2.7 mm x 14.0 mm Locking Cortical Screw |
| COL-2160-S | COL-2160 | 2.7 mm x 16.0 mm Locking Cortical Screw |
| COL-2180-S | COL-2180 | 2.7 mm x 18.0 mm Locking Cortical Screw |

Shoulder and Elbow Portfolio [continued]

Congruent Elbow Plating System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|---|
| COL-2200-S | COL-2200 | 2.7 mm x 20.0 mm Locking Cortical Screw |
| COL-2220-S | COL-2220 | 2.7 mm x 22.0 mm Locking Cortical Screw |
| COL-3120-S | COL-3120 | 3.5 mm x 12.0 mm Locking Cortical Screw |
| COL-3140-S | COL-3140 | 3.5 mm x 14.0 mm Locking Cortical Screw |
| COL-3160-S | COL-3160 | 3.5 mm x 16.0 mm Locking Cortical Screw |
| COL-3180-S | COL-3180 | 3.5 mm x 18.0 mm Locking Cortical Screw |
| COL-3200-S | COL-3200 | 3.5 mm x 20.0 mm Locking Cortical Screw |
| COL-3220-S | COL-3220 | 3.5 mm x 22.0 mm Locking Cortical Screw |
| COL-3240-S | COL-3240 | 3.5 mm x 24.0 mm Locking Cortical Screw |
| COL-3260-S | COL-3260 | 3.5 mm x 26.0 mm Locking Cortical Screw |
| COL-3280-S | COL-3280 | 3.5 mm x 28.0 mm Locking Cortical Screw |
| COL-3300-S | COL-3300 | 3.5 mm x 30.0 mm Locking Cortical Screw |
| COL-3320-S | COL-3320 | 3.5 mm x 32.0 mm Locking Cortical Screw |
| COL-3340-S | COL-3340 | 3.5 mm x 34.0 mm Locking Cortical Screw |
| COL-3360-S | COL-3360 | 3.5 mm x 36.0 mm Locking Cortical Screw |
| COL-3380-S | COL-3380 | 3.5 mm x 38.0 mm Locking Cortical Screw |
| COL-3400-S | COL-3400 | 3.5 mm x 40.0 mm Locking Cortical Screw |
| COL-3450-S | COL-3450 | 3.5 mm x 45.0 mm Locking Cortical Screw |
| COL-3500-S | COL-3500 | 3.5 mm x 50.0 mm Locking Cortical Screw |
| FA-CO2736-S | FA-CO2736 | 2.7 mm x 36 mm Tap-Loc® Screw |
| FA-CO2738-S | FA-CO2738 | 2.7 mm x 38 mm Tap-Loc Screw |
| FA-CO2740-S | FA-CO2740 | 2.7 mm x 40 mm Tap-Loc Screw |
| FA-CO2745-S | FA-CO2745 | 2.7 mm x 45 mm Tap-Loc Screw |
| FA-CO2750-S | FA-CO2750 | 2.7 mm x 50 mm Tap-Loc Screw |
| FA-CO2755-S | FA-CO2755 | 2.7 mm x 55 mm Tap-Loc Screw |
| FA-CO2760-S | FA-CO2760 | 2.7 mm x 60 mm Tap-Loc Screw |
| FA-CO3536-S | FA-CO3536 | 3.5 mm x 36 mm Tap-Loc Screw |
| FA-CO3538-S | FA-CO3538 | 3.5 mm x 38 mm Tap-Loc Screw |
| FA-CO3540-S | FA-CO3540 | 3.5 mm x 40 mm Tap-Loc Screw |
| FA-CO3545-S | FA-CO3545 | 3.5 mm x 45 mm Tap-Loc Screw |

Shoulder and Elbow Portfolio [continued]

Congruent Elbow Plating System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|--|
| FA-CO3550-S | FA-CO3550 | 3.5 mm x 50 mm Tap-Loc Screw |
| FA-CO3555-S | FA-CO3555 | 3.5 mm x 55 mm Tap-Loc Screw |
| FA-CO3560-S | FA-CO3560 | 3.5 mm x 60 mm Tap-Loc Screw |
| PL-ELCLL-S | PL-ELCLL | Coronoid Plate, Left, Extended |
| PL-ELCLR-S | PL-ELCLR | Coronoid Plate, Right, Extended |
| PL-ELCOL-S | PL-ELCOL | Coronoid Plate, Left |
| PL-ELCOR-S | PL-ELCOR | Coronoid Plate, Right |
| PL-ELPO-S | PL-ELPO | Elbow, Posterior Plate |
| PL-LEL10L-S | PL-LEL10L | Left 10-Hole Locking Lateral Plate |
| PL-LEL10R-S | PL-LEL10R | Right 10-Hole Locking Lateral Plate |
| PL-LEL14L-S | PL-LEL14L | Left 14-Hole Locking Lateral Plate |
| PL-LEL14R-S | PL-LEL14R | Right 14-Hole Locking Lateral Plate |
| PL-LEL20L-S | PL-LEL20L | Left 20-Hole Locking Lateral Plate |
| PL-LEL20R-S | PL-LEL20R | Right 20-Hole Locking Lateral Plate |
| PL-LEL6L-S | PL-LEL6L | Left 6-Hole Locking Lateral Plate |
| PL-LEL6R-S | PL-LEL6R | Right 6-Hole Locking Lateral Plate |
| PL-LEM12-S | PL-LEM12 | 12-Hole Locking Medial Plate |
| PL-LEM16-S | PL-LEM16 | 16-Hole Locking Medial Plate |
| PL-LEM7-S | PL-LEM7 | 7-Hole Locking Medial Plate |
| PL-LEM8-S | PL-LEM8 | 8-Hole Locking Medial Plate |
| PL-LEM9L-S | PL-LEM9L | Long 9-Hole Locking Medial Plate |
| PL-LEM9S-S | PL-LEM9S | Short 9-Hole Locking Medial Plate |
| PL-LEO11-S | PL-LEO11 | 11-Hole Locking Olecranon Plate |
| PL-LEO13-S | PL-LEO13 | 13-Hole Locking Olecranon Plate |
| PL-LEO13E-S | PL-LEO13E | Extended 13-Hole Locking Olecranon Plate |
| PL-LEO17L-S | PL-LEO17L | Left 17-Hole Locking Olecranon Plate |
| PL-LEO17R-S | PL-LEO17R | Right 17-Hole Locking Olecranon Plate |
| PL-LEO9-S | PL-LEO9 | 9-Hole Locking Olecranon Plate |

Shoulder and Elbow Portfolio [continued]

Radial Head Plating System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|---|
| 70-0097-S | 70-0097 | Locking Radial Head Plate 3-Hole Standard Curvature |
| 70-0098-S | 70-0098 | Locking Radial Head Plate 5-Hole Standard Curvature |
| 70-0099-S | 70-0099 | Locking Radial Head Plate 3-Hole Small Curvature |
| 70-0100-S | 70-0100 | Locking Radial Head Plate 5-Hole Small Curvature |
| CO-N2312-S | CO-N2312 | 2.3 mm x 12 mm Non-Toggling Cortical Screw |
| CO-N2314-S | CO-N2314 | 2.3 mm x 14 mm Non-Toggling Cortical Screw |
| CO-N2316-S | CO-N2316 | 2.3 mm x 16 mm Non-Toggling Cortical Screw |
| CO-N2318-S | CO-N2318 | 2.3 mm x 18 mm Non-Toggling Cortical Screw |
| CO-N2320-S | CO-N2320 | 2.3 mm x 20 mm Non-Toggling Cortical Screw |
| CO-N2322-S | CO-N2322 | 2.3 mm x 22 mm Non-Toggling Cortical Screw |
| CO-N2324-S | CO-N2324 | 2.3 mm x 24 mm Non-Toggling Cortical Screw |
| CO-N2326-S | CO-N2326 | 2.3 mm x 26 mm Non-Toggling Cortical Screw |
| CO-N2328-S | CO-N2328 | 2.3 mm x 28 mm Non-Toggling Cortical Screw |
| CO-T2312-S | CO-T2312 | 2.3 mm x 12 mm Locking Cortical Screw |
| CO-T2314-S | CO-T2314 | 2.3 mm x 14 mm Locking Cortical Screw |
| CO-T2316-S | CO-T2316 | 2.3 mm x 16 mm Locking Cortical Screw |
| CO-T2318-S | CO-T2318 | 2.3 mm x 18 mm Locking Cortical Screw |
| CO-T2320-S | CO-T2320 | 2.3 mm x 20 mm Locking Cortical Screw |
| CO-T2322-S | CO-T2322 | 2.3 mm x 22 mm Locking Cortical Screw |
| CO-T2324-S | CO-T2324 | 2.3 mm x 24 mm Locking Cortical Screw |
| CO-T2326-S | CO-T2326 | 2.3 mm x 26 mm Locking Cortical Screw |
| CO-T2328-S | CO-T2328 | 2.3 mm x 28 mm Locking Cortical Screw |

Shoulder and Elbow Portfolio [continued]

Clavicle Plating System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|--|
| 30-0232-S | 30-0232 | 3.5 mm x 8 mm Locking Hexalobe Screw |
| 30-0233-S | 30-0233 | 3.5 mm x 10 mm Locking Hexalobe Screw |
| 30-0234-S | 30-0234 | 3.5 mm x 12 mm Locking Hexalobe Screw |
| 30-0235-S | 30-0235 | 3.5 mm x 14 mm Locking Hexalobe Screw |
| 30-0236-S | 30-0236 | 3.5 mm x 16 mm Locking Hexalobe Screw |
| 30-0237-S | 30-0237 | 3.5 mm x 18 mm Locking Hexalobe Screw |
| 30-0238-S | 30-0238 | 3.5 mm x 20 mm Locking Hexalobe Screw |
| 30-0239-S | 30-0239 | 3.5 mm x 22 mm Locking Hexalobe Screw |
| 30-0240-S | 30-0240 | 3.5 mm x 24 mm Locking Hexalobe Screw |
| 30-0241-S | 30-0241 | 3.5 mm x 26 mm Locking Hexalobe Screw |
| 30-0255-S | 30-0255 | 3.5 mm x 8 mm Nonlocking Hexalobe Screw |
| 30-0256-S | 30-0256 | 3.5 mm x 10 mm Nonlocking Hexalobe Screw |
| 30-0257-S | 30-0257 | 3.5 mm x 12 mm Nonlocking Hexalobe Screw |
| 30-0258-S | 30-0258 | 3.5 mm x 14 mm Nonlocking Hexalobe Screw |
| 30-0259-S | 30-0259 | 3.5 mm x 16 mm Nonlocking Hexalobe Screw |
| 30-0260-S | 30-0260 | 3.5 mm x 18 mm Nonlocking Hexalobe Screw |
| 30-0261-S | 30-0261 | 3.5 mm x 20 mm Nonlocking Hexalobe Screw |
| 30-0262-S | 30-0262 | 3.5 mm x 22 mm Nonlocking Hexalobe Screw |
| 30-0263-S | 30-0263 | 3.5 mm x 24 mm Nonlocking Hexalobe Screw |
| 30-0264-S | 30-0264 | 3.5 mm x 26 mm Nonlocking Hexalobe Screw |
| 30-0278-S | 30-0278 | 3.0 mm x 8 mm Locking Hexalobe Screw |
| 30-0279-S | 30-0279 | 3.0 mm x 10 mm Locking Hexalobe Screw |
| 30-0280-S | 30-0280 | 3.0 mm x 12 mm Locking Hexalobe Screw |
| 30-0281-S | 30-0281 | 3.0 mm x 14 mm Locking Hexalobe Screw |
| 30-0282-S | 30-0282 | 3.0 mm x 16 mm Locking Hexalobe Screw |
| 30-0283-S | 30-0283 | 3.0 mm x 18 mm Locking Hexalobe Screw |
| 30-0284-S | 30-0284 | 3.0 mm x 20 mm Locking Hexalobe Screw |
| 30-0285-S | 30-0285 | 3.0 mm x 22 mm Locking Hexalobe Screw |
| 30-0286-S | 30-0286 | 3.0 mm x 24 mm Locking Hexalobe Screw |

Shoulder and Elbow Portfolio [continued]

Clavicle Plating System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|---|
| 30-0287-S | 30-0287 | 3.0 mm x 26 mm Locking Hexalobe Screw |
| 30-0301-S | 30-0301 | 3.0 mm x 8 mm Nonlocking Hexalobe Screw |
| 30-0302-S | 30-0302 | 3.0 mm x 10 mm Nonlocking Hexalobe Screw |
| 30-0303-S | 30-0303 | 3.0 mm x 12 mm Nonlocking Hexalobe Screw |
| 30-0304-S | 30-0304 | 3.0 mm x 14 mm Nonlocking Hexalobe Screw |
| 30-0305-S | 30-0305 | 3.0 mm x 16 mm Nonlocking Hexalobe Screw |
| 30-0306-S | 30-0306 | 3.0 mm x 18 mm Nonlocking Hexalobe Screw |
| 30-0307-S | 30-0307 | 3.0 mm x 20 mm Nonlocking Hexalobe Screw |
| 30-0308-S | 30-0308 | 3.0 mm x 22 mm Nonlocking Hexalobe Screw |
| 30-0309-S | 30-0309 | 3.0 mm x 24 mm Nonlocking Hexalobe Screw |
| 30-0310-S | 30-0310 | 3.0 mm x 26 mm Nonlocking Hexalobe Screw |
| 70-0111-S | 70-0111 | Distal Clavicle Plate 3.5 mm 12-Hole, Right |
| 70-0112-S | 70-0112 | Distal Clavicle Plate 3.5 mm 12-Hole, Left |
| 70-0116-S | 70-0116 | Distal Clavicle Plate 3.5 mm 9-Hole, Right |
| 70-0117-S | 70-0117 | Distal Clavicle Plate 3.5 mm 9-Hole, Left |
| 70-0118-S | 70-0118 | 8-Hole Lateral Anterior Clavicle Plate |
| 70-0119-S | 70-0119 | 8-Hole Medial Anterior Clavicle Plate |
| 70-0120-S | 70-0120 | 6-Hole Medial Anterior Clavicle Plate |
| 70-0121-S | 70-0121 | 10-Hole Anterior Clavicle Plate |
| 70-0122-S | 70-0122 | 6-Hole Lateral Anterior Clavicle Plate |
| 70-0123-S | 70-0123 | Distal Clavicle Plate 2.3 mm 16-Hole, Right |
| 70-0124-S | 70-0124 | Distal Clavicle Plate 2.3 mm 16-Hole, Left |
| 70-0125-S | 70-0125 | Distal Clavicle Plate 2.3 mm 13-Hole, Right |
| 70-0126-S | 70-0126 | Distal Clavicle Plate 2.3 mm 13-Hole, Left |
| 70-0286-S | 70-0286 | Low-Profile Clavicle Plate, 8-Hole Straight Left |
| 70-0287-S | 70-0287 | Low-Profile Clavicle Plate, 8-Hole Straight Right |
| 70-0288-S | 70-0288 | Low-Profile Clavicle Plate, 8-Hole Large Left |
| 70-0289-S | 70-0289 | Low-Profile Clavicle Plate, 8-Hole Large Right |
| 70-0290-S | 70-0290 | Low-Profile Clavicle Plate, 8-Hole Medium Left |

Shoulder and Elbow Portfolio [continued]

Clavicle Plating System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|--|
| 70-0291-S | 70-0291 | Low-Profile Clavicle Plate, 8-Hole Medium Right |
| 70-0292-S | 70-0292 | Low-Profile Clavicle Plate, 8-Hole Small Left |
| 70-0293-S | 70-0293 | Low-Profile Clavicle Plate, 8-Hole Small Right |
| 70-0294-S | 70-0294 | Low-Profile Clavicle Plate, 10-Hole Left |
| 70-0295-S | 70-0295 | Low-Profile Clavicle Plate, 10-Hole Right |
| 70-0296-S | 70-0296 | Narrow Profile Clavicle Plate, 6-Hole Left |
| 70-0297-S | 70-0297 | Narrow Profile Clavicle Plate, 6-Hole Right |
| 70-0298-S | 70-0298 | Narrow Profile Clavicle Plate, 8-Hole Straight Left |
| 70-0299-S | 70-0299 | Narrow Profile Clavicle Plate, 8-Hole Straight Right |
| 70-0300-S | 70-0300 | Narrow Profile Clavicle Plate, 8-Hole Large Left |
| 70-0301-S | 70-0301 | Narrow Profile Clavicle Plate, 8-Hole Large Right |
| 70-0319-S | 70-0319 | Low-Profile Clavicle J-Plate, 8-Hole Left |
| 70-0320-S | 70-0320 | Low-Profile Clavicle J-Plate, 8-Hole Right |
| 7002-0416L-S | 7002-0416L | Distal Clavicle Plate 3.5 mm 16-Hole Left |
| 7002-0416R-S | 7002-0416R | Distal Clavicle Plate 3.5 mm 16-Hole Right |
| CA-4120-S | CA-4120 | 4.0 mm x 12.0 mm Cancellous Screw |
| CA-4140-S | CA-4140 | 4.0 mm x 14.0 mm Cancellous Screw |
| CA-4160-S | CA-4160 | 4.0 mm x 16.0 mm Cancellous Screw |
| CA-4180-S | CA-4180 | 4.0 mm x 18.0 mm Cancellous Screw |
| CA-4200-S | CA-4200 | 4.0 mm x 20.0 mm Cancellous Screw |
| CA-4220-S | CA-4220 | 4.0 mm x 22.0 mm Cancellous Screw |
| CA-4240-S | CA-4240 | 4.0 mm x 24.0 mm Cancellous Screw |
| CA-4260-S | CA-4260 | 4.0 mm x 26.0 mm Cancellous Screw |
| CA-4280-S | CA-4280 | 4.0 mm x 28.0 mm Cancellous Screw |
| CA-4300-S | CA-4300 | 4.0 mm x 30.0 mm Cancellous Screw |
| CA-4350-S | CA-4350 | 4.0 mm x 35.0 mm Cancellous Screw |
| CA-4400-S | CA-4400 | 4.0 mm x 40.0 mm Cancellous Screw |
| CA-4450-S | CA-4450 | 4.0 mm x 45.0 mm Cancellous Screw |
| CA-4500-S | CA-4500 | 4.0 mm x 50.0 mm Cancellous Screw |
| CA-4550-S | CA-4550 | 4.0 mm x 55.0 mm Cancellous Screw |

Shoulder and Elbow Portfolio [continued]

Clavicle Plating System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|-----------------------------------|
| CA-4600-S | CA-4600 | 4.0 mm x 60.0 mm Cancellous Screw |
| CO-2708-S | CO-2708 | 2.7 mm x 8.0 mm Cortical Screw |
| CO-2710-S | CO-2710 | 2.7 mm x 10.0 mm Cortical Screw |
| CO-2712-S | CO-2712 | 2.7 mm x 12.0 mm Cortical Screw |
| CO-2714-S | CO-2714 | 2.7 mm x 14.0 mm Cortical Screw |
| CO-2716-S | CO-2716 | 2.7 mm x 16.0 mm Cortical Screw |
| CO-2718-S | CO-2718 | 2.7 mm x 18.0 mm Cortical Screw |
| CO-2720-S | CO-2720 | 2.7 mm x 20.0 mm Cortical Screw |
| CO-2722-S | CO-2722 | 2.7 mm x 22.0 mm Cortical Screw |
| CO-2724-S | CO-2724 | 2.7 mm x 24.0 mm Cortical Screw |
| CO-2726-S | CO-2726 | 2.7 mm x 26.0 mm Cortical Screw |
| CO-2728-S | CO-2728 | 2.7 mm x 28.0 mm Cortical Screw |
| CO-2730-S | CO-2730 | 2.7 mm x 30.0 mm Cortical Screw |
| CO-2732-S | CO-2732 | 2.7 mm x 32.0 mm Cortical Screw |
| CO-2734-S | CO-2734 | 2.7 mm x 34.0 mm Cortical Screw |
| CO-2736-S | CO-2736 | 2.7 mm x 36.0 mm Cortical Screw |
| CO-2738-S | CO-2738 | 2.7 mm x 38.0 mm Cortical Screw |
| CO-2740-S | CO-2740 | 2.7 mm x 40.0 mm Cortical Screw |
| CO-2745-S | CO-2745 | 2.7 mm x 45.0 mm Cortical Screw |
| CO-2750-S | CO-2750 | 2.7 mm x 50.0 mm Cortical Screw |
| CO-2755-S | CO-2755 | 2.7 mm x 55.0 mm Cortical Screw |
| CO-2760-S | CO-2760 | 2.7 mm x 60.0 mm Cortical Screw |
| CO-2765-S | CO-2765 | 2.7 mm x 65.0 mm Cortical Screw |
| CO-3060-S | CO-3060 | 3.5 mm x 6.0 mm Cortical Screw |
| CO-3080-S | CO-3080 | 3.5 mm x 8.0 mm Cortical Screw |
| CO-3100-S | CO-3100 | 3.5 mm x 10.0 mm Cortical Screw |
| CO-3120-S | CO-3120 | 3.5 mm x 12.0 mm Cortical Screw |
| CO-3140-S | CO-3140 | 3.5 mm x 14.0 mm Cortical Screw |
| CO-3160-S | CO-3160 | 3.5 mm x 16.0 mm Cortical Screw |

Shoulder and Elbow Portfolio [continued]

Clavicle Plating System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|--|
| CO-3180-S | CO-3180 | 3.5 mm x 18.0 mm Cortical Screw |
| CO-3200-S | CO-3200 | 3.5 mm x 20.0 mm Cortical Screw |
| CO-3220-S | CO-3220 | 3.5 mm x 22.0 mm Cortical Screw |
| CO-3240-S | CO-3240 | 3.5 mm x 24.0 mm Cortical Screw |
| CO-3260-S | CO-3260 | 3.5 mm x 26.0 mm Cortical Screw |
| CO-3280-S | CO-3280 | 3.5 mm x 28.0 mm Cortical Screw |
| CO-3300-S | CO-3300 | 3.5 mm x 30.0 mm Cortical Screw |
| CO-3320-S | CO-3320 | 3.5 mm x 32.0 mm Cortical Screw |
| CO-3340-S | CO-3340 | 3.5 mm x 34.0 mm Cortical Screw |
| CO-3360-S | CO-3360 | 3.5 mm x 36.0 mm Cortical Screw |
| CO-3380-S | CO-3380 | 3.5 mm x 38.0 mm Cortical Screw |
| CO-3400-S | CO-3400 | 3.5 mm x 40.0 mm Cortical Screw |
| CO-3450-S | CO-3450 | 3.5 mm x 45.0 mm Cortical Screw |
| CO-3500-S | CO-3500 | 3.5 mm x 50.0 mm Cortical Screw |
| CO-3550-S | CO-3550 | 3.5 mm x 55.0 mm Cortical Screw |
| CO-3600-S | CO-3600 | 3.5 mm x 60.0 mm Cortical Screw |
| CO-3650-S | CO-3650 | 3.5 mm x 65.0 mm Cortical Screw |
| CO-N2308-S | CO-N2308 | 2.3 mm x 8 mm Non-Toggling Cortical Screw |
| CO-N2310-S | CO-N2310 | 2.3 mm x 10 mm Non-Toggling Cortical Screw |
| CO-N2312-S | CO-N2312 | 2.3 mm x 12 mm Non-Toggling Cortical Screw |
| CO-N2314-S | CO-N2314 | 2.3 mm x 14 mm Non-Toggling Cortical Screw |
| CO-N2316-S | CO-N2316 | 2.3 mm x 16 mm Non-Toggling Cortical Screw |
| CO-N2318-S | CO-N2318 | 2.3 mm x 18 mm Non-Toggling Cortical Screw |
| CO-N2320-S | CO-N2320 | 2.3 mm x 20 mm Non-Toggling Cortical Screw |
| CO-N2322-S | CO-N2322 | 2.3 mm x 22 mm Non-Toggling Cortical Screw |
| CO-N2324-S | CO-N2324 | 2.3 mm x 24 mm Non-Toggling Cortical Screw |
| CO-N2326-S | CO-N2326 | 2.3 mm x 26 mm Non-Toggling Cortical Screw |
| CO-T2308-S | CO-T2308 | 2.3 mm x 8 mm Locking Cortical Screw |
| CO-T2310-S | CO-T2310 | 2.3 mm x 10 mm Locking Cortical Screw |

Shoulder and Elbow Portfolio [continued]

Clavicle Plating System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|---|
| CO-T2312-S | CO-T2312 | 2.3 mm x 12 mm Locking Cortical Screw |
| CO-T2314-S | CO-T2314 | 2.3 mm x 14 mm Locking Cortical Screw |
| CO-T2316-S | CO-T2316 | 2.3 mm x 16 mm Locking Cortical Screw |
| CO-T2318-S | CO-T2318 | 2.3 mm x 18 mm Locking Cortical Screw |
| CO-T2320-S | CO-T2320 | 2.3 mm x 20 mm Locking Cortical Screw |
| CO-T2322-S | CO-T2322 | 2.3 mm x 22 mm Locking Cortical Screw |
| CO-T2324-S | CO-T2324 | 2.3 mm x 24 mm Locking Cortical Screw |
| CO-T2326-S | CO-T2326 | 2.3 mm x 26 mm Locking Cortical Screw |
| COL-2080-S | COL-2080 | 2.7 mm x 8.0 mm Locking Cortical Screw |
| COL-2100-S | COL-2100 | 2.7 mm x 10.0 mm Locking Cortical Screw |
| COL-2120-S | COL-2120 | 2.7 mm x 12.0 mm Locking Cortical Screw |
| COL-2140-S | COL-2140 | 2.7 mm x 14.0 mm Locking Cortical Screw |
| COL-2160-S | COL-2160 | 2.7 mm x 16.0 mm Locking Cortical Screw |
| COL-2180-S | COL-2180 | 2.7 mm x 18.0 mm Locking Cortical Screw |
| COL-2200-S | COL-2200 | 2.7 mm x 20.0 mm Locking Cortical Screw |
| COL-2220-S | COL-2220 | 2.7 mm x 22.0 mm Locking Cortical Screw |
| COL-2240-S | COL-2240 | 2.7 mm x 24.0 mm Locking Cortical Screw |
| COL-2260-S | COL-2260 | 2.7 mm x 26.0 mm Locking Cortical Screw |
| COL-2280-S | COL-2280 | 2.7 mm x 28.0 mm Locking Cortical Screw |
| COL-2300-S | COL-2300 | 2.7 mm x 30.0 mm Locking Cortical Screw |
| COL-2320-S | COL-2320 | 2.7 mm x 32.0 mm Locking Cortical Screw |
| COL-2340-S | COL-2340 | 2.7 mm x 34.0 mm Locking Cortical Screw |
| COL-2360-S | COL-2360 | 2.7 mm x 36.0 mm Locking Cortical Screw |
| COL-2380-S | COL-2380 | 2.7 mm x 38.0 mm Locking Cortical Screw |
| COL-2400-S | COL-2400 | 2.7 mm x 40.0 mm Locking Cortical Screw |
| COL-2450-S | COL-2450 | 2.7 mm x 45.0 mm Locking Cortical Screw |
| COL-2500-S | COL-2500 | 2.7 mm x 50.0 mm Locking Cortical Screw |
| COL-2550-S | COL-2550 | 2.7 mm x 55.0 mm Locking Cortical Screw |
| COL-2600-S | COL-2600 | 2.7 mm x 60.0 mm Locking Cortical Screw |

Shoulder and Elbow Portfolio [continued]

Clavicle Plating System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|---|
| COL-2650-S | COL-2650 | 2.7 mm x 65.0 mm Locking Cortical Screw |
| COL-3060-S | COL-3060 | 3.5 mm x 6.0 mm Locking Cortical Screw |
| COL-3080-S | COL-3080 | 3.5 mm x 8.0 mm Locking Cortical Screw |
| COL-3100-S | COL-3100 | 3.5 mm x 10.0 mm Locking Cortical Screw |
| COL-3120-S | COL-3120 | 3.5 mm x 12.0 mm Locking Cortical Screw |
| COL-3140-S | COL-3140 | 3.5 mm x 14.0 mm Locking Cortical Screw |
| COL-3160-S | COL-3160 | 3.5 mm x 16.0 mm Locking Cortical Screw |
| COL-3180-S | COL-3180 | 3.5 mm x 18.0 mm Locking Cortical Screw |
| COL-3200-S | COL-3200 | 3.5 mm x 20.0 mm Locking Cortical Screw |
| COL-3220-S | COL-3220 | 3.5 mm x 22.0 mm Locking Cortical Screw |
| COL-3240-S | COL-3240 | 3.5 mm x 24.0 mm Locking Cortical Screw |
| COL-3260-S | COL-3260 | 3.5 mm x 26.0 mm Locking Cortical Screw |
| COL-3280-S | COL-3280 | 3.5 mm x 28.0 mm Locking Cortical Screw |
| COL-3300-S | COL-3300 | 3.5 mm x 30.0 mm Locking Cortical Screw |
| COL-3320-S | COL-3320 | 3.5 mm x 32.0 mm Locking Cortical Screw |
| COL-3340-S | COL-3340 | 3.5 mm x 34.0 mm Locking Cortical Screw |
| COL-3360-S | COL-3360 | 3.5 mm x 36.0 mm Locking Cortical Screw |
| COL-3380-S | COL-3380 | 3.5 mm x 38.0 mm Locking Cortical Screw |
| COL-3400-S | COL-3400 | 3.5 mm x 40.0 mm Locking Cortical Screw |
| COL-3450-S | COL-3450 | 3.5 mm x 45.0 mm Locking Cortical Screw |
| COL-3500-S | COL-3500 | 3.5 mm x 50.0 mm Locking Cortical Screw |
| COL-3550-S | COL-3550 | 3.5 mm x 55.0 mm Locking Cortical Screw |
| COL-3600-S | COL-3600 | 3.5 mm x 60.0 mm Locking Cortical Screw |
| COL-3650-S | COL-3650 | 3.5 mm x 65.0 mm Locking Cortical Screw |

Shoulder and Elbow Portfolio [continued]

Clavicle Hook Plating System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|--|
| 7014-0512L-S | 7014-0512L | Clavicle Hook Plate, 5-Hole 12 mm, Left |
| 7014-0512R-S | 7014-0512R | Clavicle Hook Plate, 5-Hole 12 mm, Right |
| 7014-0516L-S | 7014-0516L | Clavicle Hook Plate, 5-Hole 16 mm, Left |
| 7014-0516R-S | 7014-0516R | Clavicle Hook Plate, 5-Hole 16 mm, Right |
| 7014-0520L-S | 7014-0520L | Clavicle Hook Plate, 5-Hole 20 mm, Left |
| 7014-0520R-S | 7014-0520R | Clavicle Hook Plate, 5-Hole 20 mm, Right |
| 7014-0612L-S | 7014-0612L | Clavicle Hook Plate, 6-Hole 12 mm, Left |
| 7014-0612R-S | 7014-0612R | Clavicle Hook Plate, 6-Hole 12 mm, Right |
| 7014-0616L-S | 7014-0616L | Clavicle Hook Plate, 6-Hole 16 mm, Left |
| 7014-0616R-S | 7014-0616R | Clavicle Hook Plate, 6-Hole 16 mm, Right |
| 7014-0620L-S | 7014-0620L | Clavicle Hook Plate, 6-Hole 20 mm, Left |
| 7014-0620R-S | 7014-0620R | Clavicle Hook Plate, 6-Hole 20 mm, Right |
| 7014-0712L-S | 7014-0712L | Clavicle Hook Plate, 7-Hole 12 mm, Left |
| 7014-0712R-S | 7014-0712R | Clavicle Hook Plate, 7-Hole 12 mm, Right |
| 7014-0716L-S | 7014-0716L | Clavicle Hook Plate, 7-Hole 16 mm, Left |
| 7014-0716R-S | 7014-0716R | Clavicle Hook Plate, 7-Hole 16 mm, Right |
| 7014-0720L-S | 7014-0720L | Clavicle Hook Plate, 7-Hole 20 mm, Left |
| 7014-0720R-S | 7014-0720R | Clavicle Hook Plate, 7-Hole 20 mm, Right |
| 7014-0912L-S | 7014-0912L | Clavicle Hook Plate, 9-Hole 12 mm, Left |
| 7014-0912R-S | 7014-0912R | Clavicle Hook Plate, 9-Hole 12 mm, Right |
| 7014-0916L-S | 7014-0916L | Clavicle Hook Plate, 9-Hole 16 mm, Left |
| 7014-0916R-S | 7014-0916R | Clavicle Hook Plate, 9-Hole 16 mm, Right |
| 7014-0920L-S | 7014-0920L | Clavicle Hook Plate, 9-Hole 20 mm, Left |
| 7014-0920R-S | 7014-0920R | Clavicle Hook Plate, 9-Hole 20 mm, Right |
| 30-0232-S | 30-0232 | 3.5 mm x 8 mm Locking Hexalobe Screw |
| 30-0233-S | 30-0233 | 3.5 mm x 10 mm Locking Hexalobe Screw |
| 30-0234-S | 30-0234 | 3.5 mm x 12 mm Locking Hexalobe Screw |
| 30-0235-S | 30-0235 | 3.5 mm x 14 mm Locking Hexalobe Screw |
| 30-0236-S | 30-0236 | 3.5 mm x 16 mm Locking Hexalobe Screw |

Shoulder and Elbow Portfolio [continued]

Clavicle Hook Plating System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|---|
| 30-0237-S | 30-0237 | 3.5 mm x 18 mm Locking Hexalobe Screw |
| 30-0238-S | 30-0238 | 3.5 mm x 20 mm Locking Hexalobe Screw |
| 30-0239-S | 30-0239 | 3.5 mm x 22 mm Locking Hexalobe Screw |
| 30-0240-S | 30-0240 | 3.5 mm x 24 mm Locking Hexalobe Screw |
| 30-0241-S | 30-0241 | 3.5 mm x 26 mm Locking Hexalobe Screw |
| 30-0255-S | 30-0255 | 3.5 mm x 8 mm Non-Locking Hexalobe Screw |
| 30-0256-S | 30-0256 | 3.5 mm x 10 mm Non-Locking Hexalobe Screw |
| 30-0257-S | 30-0257 | 3.5 mm x 12 mm Non-Locking Hexalobe Screw |
| 30-0258-S | 30-0258 | 3.5 mm x 14 mm Non-Locking Hexalobe Screw |
| 30-0259-S | 30-0259 | 3.5 mm x 16 mm Non-Locking Hexalobe Screw |
| 30-0260-S | 30-0260 | 3.5 mm x 18 mm Non-Locking Hexalobe Screw |
| 30-0261-S | 30-0261 | 3.5 mm x 20 mm Non-Locking Hexalobe Screw |
| 30-0262-S | 30-0262 | 3.5 mm x 22 mm Non-Locking Hexalobe Screw |
| 30-0263-S | 30-0263 | 3.5 mm x 24 mm Non-Locking Hexalobe Screw |
| 30-0264-S | 30-0264 | 3.5 mm x 26 mm Non-Locking Hexalobe Screw |
| CO-3080-S | CO-3080 | 3.5 mm x 8.0 mm Cortical Screw |
| CO-3100-S | CO-3100 | 3.5 mm x 10.0 mm Cortical Screw |
| CO-3120-S | CO-3120 | 3.5 mm x 12.0 mm Cortical Screw |
| CO-3140-S | CO-3140 | 3.5 mm x 14.0 mm Cortical Screw |
| CO-3160-S | CO-3160 | 3.5 mm x 16.0 mm Cortical Screw |
| CO-3180-S | CO-3180 | 3.5 mm x 18.0 mm Cortical Screw |
| CO-3200-S | CO-3200 | 3.5 mm x 20.0 mm Cortical Screw |
| CO-3220-S | CO-3220 | 3.5 mm x 22.0 mm Cortical Screw |
| CO-3240-S | CO-3240 | 3.5 mm x 24.0 mm Cortical Screw |
| CO-3260-S | CO-3260 | 3.5 mm x 26.0 mm Cortical Screw |
| COL-3080-S | COL-3080 | 3.5 mm x 8.0 mm Locking Cortical Screw |
| COL-3100-S | COL-3100 | 3.5 mm x 10.0 mm Locking Cortical Screw |
| COL-3120-S | COL-3120 | 3.5 mm x 12.0 mm Locking Cortical Screw |
| COL-3140-S | COL-3140 | 3.5 mm x 14.0 mm Locking Cortical Screw |

Shoulder and Elbow Portfolio [continued]

Clavicle Hook Plating System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|---|
| COL-3160-S | COL-3160 | 3.5 mm x 16.0 mm Locking Cortical Screw |
| COL-3180-S | COL-3180 | 3.5 mm x 18.0 mm Locking Cortical Screw |
| COL-3200-S | COL-3200 | 3.5 mm x 20.0 mm Locking Cortical Screw |
| COL-3220-S | COL-3220 | 3.5 mm x 22.0 mm Locking Cortical Screw |
| COL-3240-S | COL-3240 | 3.5 mm x 24.0 mm Locking Cortical Screw |
| COL-3260-S | COL-3260 | 3.5 mm x 26.0 mm Locking Cortical Screw |

Locking Clavicle Plating System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|---|
| 30-0232-S | 30-0232 | 3.5 mm x 8 mm Locking Hexalobe Screw |
| 30-0233-S | 30-0233 | 3.5 mm x 10 mm Locking Hexalobe Screw |
| 30-0234-S | 30-0234 | 3.5 mm x 12 mm Locking Hexalobe Screw |
| 30-0235-S | 30-0235 | 3.5 mm x 14 mm Locking Hexalobe Screw |
| 30-0236-S | 30-0236 | 3.5 mm x 16 mm Locking Hexalobe Screw |
| 30-0237-S | 30-0237 | 3.5 mm x 18 mm Locking Hexalobe Screw |
| 30-0238-S | 30-0238 | 3.5 mm x 20 mm Locking Hexalobe Screw |
| 30-0239-S | 30-0239 | 3.5 mm x 22 mm Locking Hexalobe Screw |
| 30-0240-S | 30-0240 | 3.5 mm x 24 mm Locking Hexalobe Screw |
| 30-0241-S | 30-0241 | 3.5 mm x 26 mm Locking Hexalobe Screw |
| 30-0255-S | 30-0255 | 3.5 mm x 8 mm Nonlocking Hexalobe Screw |
| 30-0232-S | 30-0232 | 3.5 mm x 8 mm Locking Hexalobe Screw |
| 30-0233-S | 30-0233 | 3.5 mm x 10 mm Locking Hexalobe Screw |
| 30-0234-S | 30-0234 | 3.5 mm x 12 mm Locking Hexalobe Screw |
| 30-0235-S | 30-0235 | 3.5 mm x 14 mm Locking Hexalobe Screw |
| 30-0236-S | 30-0236 | 3.5 mm x 16 mm Locking Hexalobe Screw |
| 30-0237-S | 30-0237 | 3.5 mm x 18 mm Locking Hexalobe Screw |
| 30-0238-S | 30-0238 | 3.5 mm x 20 mm Locking Hexalobe Screw |
| 30-0239-S | 30-0239 | 3.5 mm x 22 mm Locking Hexalobe Screw |
| 30-0240-S | 30-0240 | 3.5 mm x 24 mm Locking Hexalobe Screw |

Shoulder and Elbow Portfolio [continued]

Locking Clavicle Plating System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|--|
| 30-0241-S | 30-0241 | 3.5 mm x 26 mm Locking Hexalobe Screw |
| 30-0255-S | 30-0255 | 3.5 mm x 8 mm Nonlocking Hexalobe Screw |
| 30-0256-S | 30-0256 | 3.5 mm x 10 mm Nonlocking Hexalobe Screw |
| 30-0257-S | 30-0257 | 3.5 mm x 12 mm Nonlocking Hexalobe Screw |
| 30-0258-S | 30-0258 | 3.5 mm x 14 mm Nonlocking Hexalobe Screw |
| 30-0259-S | 30-0259 | 3.5 mm x 16 mm Nonlocking Hexalobe Screw |
| 30-0260-S | 30-0260 | 3.5 mm x 18 mm Nonlocking Hexalobe Screw |
| 30-0261-S | 30-0261 | 3.5 mm x 20 mm Nonlocking Hexalobe Screw |
| 30-0262-S | 30-0262 | 3.5 mm x 22 mm Nonlocking Hexalobe Screw |
| 30-0263-S | 30-0263 | 3.5 mm x 24 mm Nonlocking Hexalobe Screw |
| 30-0264-S | 30-0264 | 3.5 mm x 26 mm Nonlocking Hexalobe Screw |
| 30-0301-S | 30-0301 | 3.0 mm x 8 mm Nonlocking Hexalobe Screw |
| 30-0302-S | 30-0302 | 3.0 mm x 10 mm Nonlocking Hexalobe Screw |
| 30-0303-S | 30-0303 | 3.0 mm x 12 mm Nonlocking Hexalobe Screw |
| 30-0304-S | 30-0304 | 3.0 mm x 14 mm Nonlocking Hexalobe Screw |
| 30-0305-S | 30-0305 | 3.0 mm x 16 mm Nonlocking Hexalobe Screw |
| 30-0306-S | 30-0306 | 3.0 mm x 18 mm Nonlocking Hexalobe Screw |
| 30-0307-S | 30-0307 | 3.0 mm x 20 mm Nonlocking Hexalobe Screw |
| 30-0308-S | 30-0308 | 3.0 mm x 22 mm Nonlocking Hexalobe Screw |
| 30-0309-S | 30-0309 | 3.0 mm x 24 mm Nonlocking Hexalobe Screw |
| 30-0310-S | 30-0310 | 3.0 mm x 26 mm Nonlocking Hexalobe Screw |
| 70-0095-S | 70-0095 | Locking Clavicle Plate 8-Hole Straight Right |
| 70-0096-S | 70-0096 | Locking Clavicle Plate 8-Hole Straight Left |
| CO-2708-S | CO-2708 | 2.7 mm x 8.0 mm Cortical Screw |
| CO-2710-S | CO-2710 | 2.7 mm x 10.0 mm Cortical Screw |
| CO-2712-S | CO-2712 | 2.7 mm x 12.0 mm Cortical Screw |
| CO-2714-S | CO-2714 | 2.7 mm x 14.0 mm Cortical Screw |
| CO-2716-S | CO-2716 | 2.7 mm x 16.0 mm Cortical Screw |
| CO-2718-S | CO-2718 | 2.7 mm x 18.0 mm Cortical Screw |

Shoulder and Elbow Portfolio [continued]

Locking Clavicle Plating System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|---|
| CO-2720-S | CO-2720 | 2.7 mm x 20.0 mm Cortical Screw |
| CO-2722-S | CO-2722 | 2.7 mm x 22.0 mm Cortical Screw |
| CO-2724-S | CO-2724 | 2.7 mm x 24.0 mm Cortical Screw |
| CO-2726-S | CO-2726 | 2.7 mm x 26.0 mm Cortical Screw |
| CO-2728-S | CO-2728 | 2.7 mm x 28.0 mm Cortical Screw |
| CO-2730-S | CO-2730 | 2.7 mm x 30.0 mm Cortical Screw |
| CO-2732-S | CO-2732 | 2.7 mm x 32.0 mm Cortical Screw |
| CO-2734-S | CO-2734 | 2.7 mm x 34.0 mm Cortical Screw |
| CO-2736-S | CO-2736 | 2.7 mm x 36.0 mm Cortical Screw |
| CO-2738-S | CO-2738 | 2.7 mm x 38.0 mm Cortical Screw |
| CO-2740-S | CO-2740 | 2.7 mm x 40.0 mm Cortical Screw |
| CO-2745-S | CO-2745 | 2.7 mm x 45.0 mm Cortical Screw |
| CO-2750-S | CO-2750 | 2.7 mm x 50.0 mm Cortical Screw |
| CO-2755-S | CO-2755 | 2.7 mm x 55.0 mm Cortical Screw |
| CO-2760-S | CO-2760 | 2.7 mm x 60.0 mm Cortical Screw |
| CO-2765-S | CO-2765 | 2.7 mm x 65.0 mm Cortical Screw |
| CO-3080-S | CO-3080 | 3.5 mm x 8.0 mm Cortical Screw |
| CO-3100-S | CO-3100 | 3.5 mm x 10.0 mm Cortical Screw |
| CO-3120-S | CO-3120 | 3.5 mm x 12.0 mm Cortical Screw |
| CO-3140-S | CO-3140 | 3.5 mm x 14.0 mm Cortical Screw |
| CO-3160-S | CO-3160 | 3.5 mm x 16.0 mm Cortical Screw |
| CO-3180-S | CO-3180 | 3.5 mm x 18.0 mm Cortical Screw |
| CO-3200-S | CO-3200 | 3.5 mm x 20.0 mm Cortical Screw |
| CO-3220-S | CO-3220 | 3.5 mm x 22.0 mm Cortical Screw |
| CO-3240-S | CO-3240 | 3.5 mm x 24.0 mm Cortical Screw |
| CO-3260-S | CO-3260 | 3.5 mm x 26.0 mm Cortical Screw |
| COL-3080-S | COL-3080 | 3.5 mm x 8.0 mm Locking Cortical Screw |
| COL-3100-S | COL-3100 | 3.5 mm x 10.0 mm Locking Cortical Screw |
| COL-3120-S | COL-3120 | 3.5 mm x 12.0 mm Locking Cortical Screw |

Shoulder and Elbow Portfolio [continued]

Locking Clavicle Plating System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|---|
| COL-3140-S | COL-3140 | 3.5 mm x 14.0 mm Locking Cortical Screw |
| COL-3160-S | COL-3160 | 3.5 mm x 16.0 mm Locking Cortical Screw |
| COL-3180-S | COL-3180 | 3.5 mm x 18.0 mm Locking Cortical Screw |
| COL-3200-S | COL-3200 | 3.5 mm x 20.0 mm Locking Cortical Screw |
| COL-3220-S | COL-3220 | 3.5 mm x 22.0 mm Locking Cortical Screw |
| COL-3240-S | COL-3240 | 3.5 mm x 24.0 mm Locking Cortical Screw |
| COL-3260-S | COL-3260 | 3.5 mm x 26.0 mm Locking Cortical Screw |
| PL-CL10LL-S | PL-CL10LL | Locking Clavicle Plate, 10-Hole Large Left |
| PL-CL10LR-S | PL-CL10LR | Locking Clavicle Plate, 10-Hole Large Right |
| PL-CL6SL-S | PL-CL6SL | Locking Clavicle Plate, 6-Hole Small Left |
| PL-CL6SR-S | PL-CL6SR | Locking Clavicle Plate, 6-Hole Small Right |
| PL-CL8JL-S | PL-CL8JL | Locking Clavicle J-Plate, 8-Hole, Left |
| PL-CL8JR-S | PL-CL8JR | Locking Clavicle J-Plate, 8-Hole, Right |
| PL-CL8LL-S | PL-CL8LL | Locking Clavicle Plate, 8-Hole Large Left |
| PL-CL8LR-S | PL-CL8LR | Locking Clavicle Plate, 8-Hole Large Right |
| PL-CL8ML-S | PL-CL8ML | Locking Clavicle Plate, 8-Hole Medium Left |
| PL-CL8MR-S | PL-CL8MR | Locking Clavicle Plate, 8-Hole Medium Right |
| PL-CL8SL-S | PL-CL8SL | Locking Clavicle Plate, 8-Hole Small Left |
| PL-CL8SR-S | PL-CL8SR | Locking Clavicle Plate, 8-Hole Small Right |
| PL-CL9JL-S | PL-CL9JL | Locking Clavicle J-Plate, 9-Hole, Left |
| PL-CL9JR-S | PL-CL9JR | Locking Clavicle J-Plate, 9-Hole, Right |

Acu-Sinch® Repair System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|------------------------------|
| N/A | 55-0004 | Acu-Sinch Anchor |
| N/A | 55-0014 | Suture, FlexBraid®, #5 x 28" |
| N/A | 55-0005 | Acu-Sinch Suture Retainer |

Shoulder and Elbow Portfolio [continued]

Dual-Trak Clavicle Screw System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|--|
| 40-0136-S | 40-0136 | 3.0 mm x 80 mm Dual-Trak Clavicle Screw |
| 40-0137-S | 40-0137 | 3.0 mm x 90 mm Dual-Trak Clavicle Screw |
| 40-0138-S | 40-0138 | 3.0 mm x 100 mm Dual-Trak Clavicle Screw |
| 40-0139-S | 40-0139 | 3.0 mm x 110 mm Dual-Trak Clavicle Screw |
| 40-0140-S | 40-0140 | 3.0 mm x 120 mm Dual-Trak Clavicle Screw |
| 40-0143-S | 40-0143 | 3.8 mm x 80 mm Dual-Trak Clavicle Screw |
| 40-0144-S | 40-0144 | 3.8 mm x 90 mm Dual-Trak Clavicle Screw |
| 40-0145-S | 40-0145 | 3.8 mm x 100 mm Dual-Trak Clavicle Screw |
| 40-0146-S | 40-0146 | 3.8 mm x 110 mm Dual-Trak Clavicle Screw |
| 40-0147-S | 40-0147 | 3.8 mm x 120 mm Dual-Trak Clavicle Screw |

Polarus® 3 Solution (plates)

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|---|
| 3011-43018-S | 3011-43018 | 4.3 mm x 18 mm Low-Profile Hexalobe Screw |
| 3011-43020-S | 3011-43020 | 4.3 mm x 20 mm Low-Profile Hexalobe Screw |
| 3011-43022-S | 3011-43022 | 4.3 mm x 22 mm Low-Profile Hexalobe Screw |
| 3011-43024-S | 3011-43024 | 4.3 mm x 24 mm Low-Profile Hexalobe Screw |
| 3011-43026-S | 3011-43026 | 4.3 mm x 26 mm Low-Profile Hexalobe Screw |
| 3011-43028-S | 3011-43028 | 4.3 mm x 28 mm Low-Profile Hexalobe Screw |
| 3011-43030-S | 3011-43030 | 4.3 mm x 30 mm Low-Profile Hexalobe Screw |
| 3011-43032-S | 3011-43032 | 4.3 mm x 32 mm Low-Profile Hexalobe Screw |
| 3011-43034-S | 3011-43034 | 4.3 mm x 34 mm Low-Profile Hexalobe Screw |
| 3011-43036-S | 3011-43036 | 4.3 mm x 36 mm Low-Profile Hexalobe Screw |
| 3011-43038-S | 3011-43038 | 4.3 mm x 38 mm Low-Profile Hexalobe Screw |
| 3011-43040-S | 3011-43040 | 4.3 mm x 40 mm Low-Profile Hexalobe Screw |
| 3011-43042-S | 3011-43042 | 4.3 mm x 42 mm Low-Profile Hexalobe Screw |
| 3011-43044-S | 3011-43044 | 4.3 mm x 44 mm Low-Profile Hexalobe Screw |
| 3011-43046-S | 3011-43046 | 4.3 mm x 46 mm Low-Profile Hexalobe Screw |

Shoulder and Elbow Portfolio [continued]

Polarus® 3 Solution (plates)

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|---|
| 3011-43048-S | 3011-43048 | 4.3 mm x 48 mm Low-Profile Hexalobe Screw |
| 3011-43050-S | 3011-43050 | 4.3 mm x 50 mm Low-Profile Hexalobe Screw |
| 3011-43052-S | 3011-43052 | 4.3 mm x 52 mm Low-Profile Hexalobe Screw |
| 3011-43054-S | 3011-43054 | 4.3 mm x 54 mm Low-Profile Hexalobe Screw |
| 3011-43056-S | 3011-43056 | 4.3 mm x 56 mm Low-Profile Hexalobe Screw |
| 3011-43058-S | 3011-43058 | 4.3 mm x 58 mm Low-Profile Hexalobe Screw |
| 3011-43060-S | 3011-43060 | 4.3 mm x 60 mm Low-Profile Hexalobe Screw |
| 3011-43062-S | 3011-43062 | 4.3 mm x 62 mm Low-Profile Hexalobe Screw |
| 3011-43064-S | 3011-43064 | 4.3 mm x 64 mm Low-Profile Hexalobe Screw |
| 3025-35018-S | 3025-35018 | 3.5 x 18 mm Nonlocking Low-Profile Hexalobe Screw |
| 3025-35020-S | 3025-35020 | 3.5 x 20 mm Nonlocking Low-Profile Hexalobe Screw |
| 3025-35022-S | 3025-35022 | 3.5 x 22 mm Nonlocking Low-Profile Hexalobe Screw |
| 3025-35024-S | 3025-35024 | 3.5 x 24 mm Nonlocking Low-Profile Hexalobe Screw |
| 3025-35026-S | 3025-35026 | 3.5 x 26 mm Nonlocking Low-Profile Hexalobe Screw |
| 3025-35028-S | 3025-35028 | 3.5 x 28 mm Nonlocking Low-Profile Hexalobe Screw |
| 3025-35030-S | 3025-35030 | 3.5 x 30 mm Nonlocking Low-Profile Hexalobe Screw |
| 3025-35032-S | 3025-35032 | 3.5 x 32 mm Nonlocking Low-Profile Hexalobe Screw |
| 3025-35034-S | 3025-35034 | 3.5 x 34 mm Nonlocking Low-Profile Hexalobe Screw |
| 3025-35036-S | 3025-35036 | 3.5 x 36 mm Nonlocking Low-Profile Hexalobe Screw |
| 3025-35038-S | 3025-35038 | 3.5 x 38 mm Nonlocking Low-Profile Hexalobe Screw |
| 3025-35040-S | 3025-35040 | 3.5 x 40 mm Nonlocking Low-Profile Hexalobe Screw |
| 3025-35042-S | 3025-35042 | 3.5 x 42 mm Nonlocking Low-Profile Hexalobe Screw |
| 3025-35044-S | 3025-35044 | 3.5 x 44 mm Nonlocking Low-Profile Hexalobe Screw |
| 3025-35046-S | 3025-35046 | 3.5 x 46 mm Nonlocking Low-Profile Hexalobe Screw |
| 3025-35048-S | 3025-35048 | 3.5 x 48 mm Nonlocking Low-Profile Hexalobe Screw |
| 3025-35050-S | 3025-35050 | 3.5 x 50 mm Nonlocking Low-Profile Hexalobe Screw |
| 3025-35052-S | 3025-35052 | 3.5 x 52 mm Nonlocking Low-Profile Hexalobe Screw |
| 3025-35054-S | 3025-35054 | 3.5 x 54 mm Nonlocking Low-Profile Hexalobe Screw |
| 3025-35056-S | 3025-35056 | 3.5 x 56 mm Nonlocking Low-Profile Hexalobe Screw |

Shoulder and Elbow Portfolio [continued]

Polarus® 3 Solution (plates)

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|---|
| 3025-35058-S | 3025-35058 | 3.5 x 58 mm Nonlocking Low-Profile Hexalobe Screw |
| 3025-35060-S | 3025-35060 | 3.5 x 60 mm Nonlocking Low-Profile Hexalobe Screw |
| 3025-35062-S | 3025-35062 | 3.5 x 62 mm Nonlocking Low-Profile Hexalobe Screw |
| 3025-35064-S | 3025-35064 | 3.5 x 64 mm Nonlocking Low-Profile Hexalobe Screw |
| 7001-0104L-S | 7001-0104L | Polarus 3 Standard Plate 4-Hole Left |
| 7001-0104R-S | 7001-0104R | Polarus 3 Standard Plate 4-Hole Right |
| 7001-0106L-S | 7001-0106L | Polarus 3 Standard Plate 6-Hole Left |
| 7001-0106R-S | 7001-0106R | Polarus 3 Standard Plate 6-Hole Right |
| 7001-0110L-S | 7001-0110L | Polarus 3 Standard Plate 10-Hole Left |
| 7001-0110R-S | 7001-0110R | Polarus 3 Standard Plate 10-Hole Right |
| N/A | 7001-0114L-S | Polarus 3 Standard Plate 14-Hole Left |
| N/A | 7001-0114R-S | Polarus 3 Standard Plate 14-Hole Right |
| N/A | 7001-0118L-S | Polarus 3 Standard Plate 18-Hole Left |
| N/A | 7001-0118R-S | Polarus 3 Standard Plate 18-Hole Right |
| N/A | 7001-0122L-S | Polarus 3 Standard Plate 22-Hole Left |
| N/A | 7001-0122R-S | Polarus 3 Standard Plate 22-Hole Right |
| 7001-0204L-S | 7001-0204L | Polarus 3 Posterior Plate 4-Hole Left |
| 7001-0204R-S | 7001-0204R | Polarus 3 Posterior Plate 4-Hole Right |
| 7001-0206L-S | 7001-0206L | Polarus 3 Posterior Plate 6-Hole Left |
| 7001-0206R-S | 7001-0206R | Polarus 3 Posterior Plate 6-Hole Right |
| N/A | 7001-03001-S | Polarus 3 8 mm Washer, Locking |

Shoulder and Elbow Portfolio [continued]

Polarus® Proximal Humeral Plating System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|---|
| 30-0426-S | 30-0426 | 5.7 x 26 mm Locking Screw |
| 30-0428-S | 30-0428 | 5.7 x 28 mm Locking Screw |
| 30-0430-S | 30-0430 | 5.7 x 30 mm Locking Screw |
| 30-0432-S | 30-0432 | 5.7 x 32 mm Locking Screw |
| 30-0434-S | 30-0434 | 5.7 x 34 mm Locking Screw |
| 30-0436-S | 30-0436 | 5.7 x 36 mm Locking Screw |
| 30-0438-S | 30-0438 | 5.7 x 38 mm Locking Screw |
| 30-0440-S | 30-0440 | 5.7 x 40 mm Locking Screw |
| 30-0442-S | 30-0442 | 5.7 x 42 mm Locking Screw |
| 30-0444-S | 30-0444 | 5.7 x 44 mm Locking Screw |
| 30-0446-S | 30-0446 | 5.7 x 46 mm Locking Screw |
| 30-0448-S | 30-0448 | 5.7 x 48 mm Locking Screw |
| 30-0450-S | 30-0450 | 5.7 x 50 mm Locking Screw |
| 30-0452-S | 30-0452 | 5.7 x 52 mm Locking Screw |
| 30-0454-S | 30-0454 | 5.7 x 54 mm Locking Screw |
| CA-PHB25-S | CA-PHB25 | 4.5 mm x 25 mm Locking Buttress Screw |
| CA-PHB275-S | CA-PHB275 | 4.5 mm x 27.5 mm Locking Buttress Screw |
| CA-PHB30-S | CA-PHB30 | 4.5 mm x 30 mm Locking Buttress Screw |
| CA-PHB325-S | CA-PHB325 | 4.5 mm x 32.5 mm Locking Buttress Screw |
| CA-PHB35-S | CA-PHB35 | 4.5 mm x 35 mm Locking Buttress Screw |
| CA-PHB375-S | CA-PHB375 | 4.5 mm x 37.5 mm Locking Buttress Screw |
| CA-PHB40-S | CA-PHB40 | 4.5 mm x 40 mm Locking Buttress Screw |
| CA-PHB45-S | CA-PHB45 | 4.5 mm x 45 mm Locking Buttress Screw |
| CA-PHB50-S | CA-PHB50 | 4.5 mm x 50 mm Locking Buttress Screw |
| CA-PHB55-S | CA-PHB55 | 4.5 mm x 55 mm Locking Buttress Screw |
| CO-3200-S | CO-3200 | 3.5 mm x 20.0 mm Cortical Screw |
| CO-3225-S | CO-3225 | 3.5 mm x 22.5 mm Cortical Screw |
| CO-3250-S | CO-3250 | 3.5 mm x 25.0 mm Cortical Screw |
| CO-3275-S | CO-3275 | 3.5 mm x 27.5 mm Cortical Screw |

Shoulder and Elbow Portfolio [continued]

Polarus® Proximal Humeral Plating System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|---|
| CO-3300-S | CO-3300 | 3.5 mm x 30.0 mm Cortical Screw |
| CO-3325-S | CO-3325 | 3.5 mm x 32.5 mm Cortical Screw |
| CO-3350-S | CO-3350 | 3.5 mm x 35.0 mm Cortical Screw |
| CO-3400-S | CO-3400 | 3.5 mm x 40.0 mm Cortical Screw |
| CO-3450-S | CO-3450 | 3.5 mm x 45.0 mm Cortical Screw |
| CO-3500-S | CO-3500 | 3.5 mm x 50.0 mm Cortical Screw |
| CO-3550-S | CO-3550 | 3.5 mm x 55.0 mm Cortical Screw |
| COL-3200-S | COL-3200 | 3.5 mm x 20.0 mm Locking Cortical Screw |
| COL-3225-S | COL-3225 | 3.5 mm x 22.5 mm Locking Cortical Screw |
| COL-3250-S | COL-3250 | 3.5 mm x 25.0 mm Locking Cortical Screw |
| COL-3275-S | COL-3275 | 3.5 mm x 27.5 mm Locking Cortical Screw |
| COL-3300-S | COL-3300 | 3.5 mm x 30.0 mm Locking Cortical Screw |
| N/A | HCA-5125 | 5.0 mm x 25.0 mm Cancellous Screw |
| N/A | HCA-5127 | 5.0 mm x 27.5 mm Cancellous Screw |
| N/A | HCA-5130 | 5.0 mm x 30.0 mm Cancellous Screw |
| N/A | HCA-5132 | 5.0 mm x 32.5 mm Cancellous Screw |
| N/A | HCA-5135 | 5.0 mm x 35.0 mm Cancellous Screw |
| N/A | HCA-5137 | 5.0 mm x 37.5 mm Cancellous Screw |
| N/A | HCA-5140 | 5.0 mm x 40.0 mm Cancellous Screw |
| N/A | HCA-5145 | 5.0 mm x 45.0 mm Cancellous Screw |
| N/A | HCA-5150 | 5.0 mm x 50.0 mm Cancellous Screw |
| N/A | HCA-5155 | 5.0 mm x 55.0 mm Cancellous Screw |
| PL-PHGL-S | PL-PHGL | Proximal Humeral Plate Large, Left |
| PL-PHGR-S | PL-PHGR | Proximal Humeral Plate Large, Right |
| PL-PHSL-S | PL-PHSL | Proximal Humeral Plate, Small, Left |
| PL-PHSR-S | PL-PHSR | Proximal Humeral Plate, Small, Right |
| PL-PHXGL-S | PL-PHXGL | Proximal Humeral Plate, Extra Long, Left |
| PL-PHXGR-S | PL-PHXGR | Proximal Humeral Plate, Extra Long, Right |

Shoulder and Elbow Portfolio [continued]

Scapula Plating System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|-----------------------------------|
| CA-4120-S | CA-4120 | 4.0 mm x 12.0 mm Cancellous Screw |
| CA-4140-S | CA-4140 | 4.0 mm x 14.0 mm Cancellous Screw |
| CA-4160-S | CA-4160 | 4.0 mm x 16.0 mm Cancellous Screw |
| CA-4180-S | CA-4180 | 4.0 mm x 18.0 mm Cancellous Screw |
| CA-4200-S | CA-4200 | 4.0 mm x 20.0 mm Cancellous Screw |
| CA-4220-S | CA-4220 | 4.0 mm x 22.0 mm Cancellous Screw |
| CA-4240-S | CA-4240 | 4.0 mm x 24.0 mm Cancellous Screw |
| CA-4260-S | CA-4260 | 4.0 mm x 26.0 mm Cancellous Screw |
| CA-4280-S | CA-4280 | 4.0 mm x 28.0 mm Cancellous Screw |
| CA-4300-S | CA-4300 | 4.0 mm x 30.0 mm Cancellous Screw |
| CA-4350-S | CA-4350 | 4.0 mm x 35.0 mm Cancellous Screw |
| CA-4400-S | CA-4400 | 4.0 mm x 40.0 mm Cancellous Screw |
| CA-4450-S | CA-4450 | 4.0 mm x 45.0 mm Cancellous Screw |
| CA-4500-S | CA-4500 | 4.0 mm x 50.0 mm Cancellous Screw |
| CA-4550-S | CA-4550 | 4.0 mm x 55.0 mm Cancellous Screw |
| CA-4600-S | CA-4600 | 4.0 mm x 60.0 mm Cancellous Screw |
| CO-2708-S | CO-2708 | 2.7 mm x 8.0 mm Cortical Screw |
| CO-2710-S | CO-2710 | 2.7 mm x 10.0 mm Cortical Screw |
| CO-2712-S | CO-2712 | 2.7 mm x 12.0 mm Cortical Screw |
| CO-2714-S | CO-2714 | 2.7 mm x 14.0 mm Cortical Screw |
| CO-2716-S | CO-2716 | 2.7 mm x 16.0 mm Cortical Screw |
| CO-2718-S | CO-2718 | 2.7 mm x 18.0 mm Cortical Screw |
| CO-2720-S | CO-2720 | 2.7 mm x 20.0 mm Cortical Screw |
| CO-2722-S | CO-2722 | 2.7 mm x 22.0 mm Cortical Screw |
| CO-2724-S | CO-2724 | 2.7 mm x 24.0 mm Cortical Screw |
| CO-2726-S | CO-2726 | 2.7 mm x 26.0 mm Cortical Screw |
| CO-2728-S | CO-2728 | 2.7 mm x 28.0 mm Cortical Screw |
| CO-2730-S | CO-2730 | 2.7 mm x 30.0 mm Cortical Screw |
| CO-2732-S | CO-2732 | 2.7 mm x 32.0 mm Cortical Screw |
| CO-2734-S | CO-2734 | 2.7 mm x 34.0 mm Cortical Screw |

Shoulder and Elbow Portfolio [continued]

Scapula Plating System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|---------------------------------|
| CO-2736-S | CO-2736 | 2.7 mm x 36.0 mm Cortical Screw |
| CO-2738-S | CO-2738 | 2.7 mm x 38.0 mm Cortical Screw |
| CO-2740-S | CO-2740 | 2.7 mm x 40.0 mm Cortical Screw |
| CO-2745-S | CO-2745 | 2.7 mm x 45.0 mm Cortical Screw |
| CO-2750-S | CO-2750 | 2.7 mm x 50.0 mm Cortical Screw |
| CO-2755-S | CO-2755 | 2.7 mm x 55.0 mm Cortical Screw |
| CO-2760-S | CO-2760 | 2.7 mm x 60.0 mm Cortical Screw |
| CO-2765-S | CO-2765 | 2.7 mm x 65.0 mm Cortical Screw |
| CO-3060-S | CO-3060 | 3.5 mm x 6.0 mm Cortical Screw |
| CO-3080-S | CO-3080 | 3.5 mm x 8.0 mm Cortical Screw |
| CO-3100-S | CO-3100 | 3.5 mm x 10.0 mm Cortical Screw |
| CO-3120-S | CO-3120 | 3.5 mm x 12.0 mm Cortical Screw |
| CO-3140-S | CO-3140 | 3.5 mm x 14.0 mm Cortical Screw |
| CO-3160-S | CO-3160 | 3.5 mm x 16.0 mm Cortical Screw |
| CO-3180-S | CO-3180 | 3.5 mm x 18.0 mm Cortical Screw |
| CO-3200-S | CO-3200 | 3.5 mm x 20.0 mm Cortical Screw |
| CO-3220-S | CO-3220 | 3.5 mm x 22.0 mm Cortical Screw |
| CO-3240-S | CO-3240 | 3.5 mm x 24.0 mm Cortical Screw |
| CO-3260-S | CO-3260 | 3.5 mm x 26.0 mm Cortical Screw |
| CO-3280-S | CO-3280 | 3.5 mm x 28.0 mm Cortical Screw |
| CO-3300-S | CO-3300 | 3.5 mm x 30.0 mm Cortical Screw |
| CO-3320-S | CO-3320 | 3.5 mm x 32.0 mm Cortical Screw |
| CO-3340-S | CO-3340 | 3.5 mm x 34.0 mm Cortical Screw |
| CO-3360-S | CO-3360 | 3.5 mm x 36.0 mm Cortical Screw |
| CO-3380-S | CO-3380 | 3.5 mm x 38.0 mm Cortical Screw |
| CO-3400-S | CO-3400 | 3.5 mm x 40.0 mm Cortical Screw |
| CO-3450-S | CO-3450 | 3.5 mm x 45.0 mm Cortical Screw |
| CO-3500-S | CO-3500 | 3.5 mm x 50.0 mm Cortical Screw |
| CO-3550-S | CO-3550 | 3.5 mm x 55.0 mm Cortical Screw |
| CO-3600-S | CO-3600 | 3.5 mm x 60.0 mm Cortical Screw |

Shoulder and Elbow Portfolio [continued]

Scapula Plating System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|---|
| CO-3650-S | CO-3650 | 3.5 mm x 65.0 mm Cortical Screw |
| COL-2080-S | COL-2080 | 2.7 mm x 8.0 mm Locking Cortical Screw |
| COL-2100-S | COL-2100 | 2.7 mm x 10.0 mm Locking Cortical Screw |
| COL-2120-S | COL-2120 | 2.7 mm x 12.0 mm Locking Cortical Screw |
| COL-2140-S | COL-2140 | 2.7 mm x 14.0 mm Locking Cortical Screw |
| COL-2160-S | COL-2160 | 2.7 mm x 16.0 mm Locking Cortical Screw |
| COL-2180-S | COL-2180 | 2.7 mm x 18.0 mm Locking Cortical Screw |
| COL-2200-S | COL-2200 | 2.7 mm x 20.0 mm Locking Cortical Screw |
| COL-2220-S | COL-2220 | 2.7 mm x 22.0 mm Locking Cortical Screw |
| COL-2240-S | COL-2240 | 2.7 mm x 24.0 mm Locking Cortical Screw |
| COL-2260-S | COL-2260 | 2.7 mm x 26.0 mm Locking Cortical Screw |
| COL-2280-S | COL-2280 | 2.7 mm x 28.0 mm Locking Cortical Screw |
| COL-2300-S | COL-2300 | 2.7 mm x 30.0 mm Locking Cortical Screw |
| COL-2320-S | COL-2320 | 2.7 mm x 32.0 mm Locking Cortical Screw |
| COL-2340-S | COL-2340 | 2.7 mm x 34.0 mm Locking Cortical Screw |
| COL-2360-S | COL-2360 | 2.7 mm x 36.0 mm Locking Cortical Screw |
| COL-2380-S | COL-2380 | 2.7 mm x 38.0 mm Locking Cortical Screw |
| COL-2400-S | COL-2400 | 2.7 mm x 40.0 mm Locking Cortical Screw |
| COL-2450-S | COL-2450 | 2.7 mm x 45.0 mm Locking Cortical Screw |
| COL-2500-S | COL-2500 | 2.7 mm x 50.0 mm Locking Cortical Screw |
| COL-2550-S | COL-2550 | 2.7 mm x 55.0 mm Locking Cortical Screw |
| COL-2600-S | COL-2600 | 2.7 mm x 60.0 mm Locking Cortical Screw |
| COL-2650-S | COL-2650 | 2.7 mm x 65.0 mm Locking Cortical Screw |
| COL-3060-S | COL-3060 | 3.5 mm x 6.0 mm Locking Cortical Screw |
| COL-3080-S | COL-3080 | 3.5 mm x 8.0 mm Locking Cortical Screw |
| COL-3100-S | COL-3100 | 3.5 mm x 10.0 mm Locking Cortical Screw |
| COL-3120-S | COL-3120 | 3.5 mm x 12.0 mm Locking Cortical Screw |
| COL-3140-S | COL-3140 | 3.5 mm x 14.0 mm Locking Cortical Screw |
| COL-3160-S | COL-3160 | 3.5 mm x 16.0 mm Locking Cortical Screw |

Shoulder and Elbow Portfolio [continued]

Scapula Plating System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|---|
| COL-3180-S | COL-3180 | 3.5 mm x 18.0 mm Locking Cortical Screw |
| COL-3200-S | COL-3200 | 3.5 mm x 20.0 mm Locking Cortical Screw |
| COL-3220-S | COL-3220 | 3.5 mm x 22.0 mm Locking Cortical Screw |
| COL-3240-S | COL-3240 | 3.5 mm x 24.0 mm Locking Cortical Screw |
| COL-3260-S | COL-3260 | 3.5 mm x 26.0 mm Locking Cortical Screw |
| COL-3280-S | COL-3280 | 3.5 mm x 28.0 mm Locking Cortical Screw |
| COL-3300-S | COL-3300 | 3.5 mm x 30.0 mm Locking Cortical Screw |
| COL-3320-S | COL-3320 | 3.5 mm x 32.0 mm Locking Cortical Screw |
| COL-3340-S | COL-3340 | 3.5 mm x 34.0 mm Locking Cortical Screw |
| COL-3360-S | COL-3360 | 3.5 mm x 36.0 mm Locking Cortical Screw |
| COL-3380-S | COL-3380 | 3.5 mm x 38.0 mm Locking Cortical Screw |
| COL-3400-S | COL-3400 | 3.5 mm x 40.0 mm Locking Cortical Screw |
| COL-3450-S | COL-3450 | 3.5 mm x 45.0 mm Locking Cortical Screw |
| COL-3500-S | COL-3500 | 3.5 mm x 50.0 mm Locking Cortical Screw |
| COL-3550-S | COL-3550 | 3.5 mm x 55.0 mm Locking Cortical Screw |
| COL-3600-S | COL-3600 | 3.5 mm x 60.0 mm Locking Cortical Screw |
| COL-3650-S | COL-3650 | 3.5 mm x 65.0 mm Locking Cortical Screw |
| PL-SA06L-S | PL-SA06L | 6-Hole Scapula Acromion Plate, Left |
| PL-SA06R-S | PL-SA06R | 6-Hole Scapula Acromion Plate, Right |
| PL-SA07L-S | PL-SA07L | 7-Hole Scapula Acromion Plate, Left |
| PL-SA07R-S | PL-SA07R | 7-Hole Scapula Acromion Plate, Right |
| PL-SG04L-S | PL-SG04L | 4-Hole Scapula Glenoid Plate, Left |
| PL-SG04R-S | PL-SG04R | 4-Hole Scapula Glenoid Plate, Right |
| PL-SLB10L-S | PL-SLB10L | 10-Hole Scapula Lateral Border Plate, Left |
| PL-SLB10R-S | PL-SLB10R | 10-Hole Scapula Lateral Border Plate, Right |
| PL-SMB09L-S | PL-SMB09L | 9-Hole Scapula Medial Border Plate, Left |
| PL-SMB09R-S | PL-SMB09R | 9-Hole Scapula Medial Border Plate, Right |
| PL-SMB13L-S | PL-SMB13L | 13-Hole Scapula Medial Border Plate, Left |
| PL-SMB13R-S | PL-SMB13R | 13-Hole Scapula Medial Border Plate, Right |

Foot and Ankle Portfolio

- ▶ Ankle Plating System
- ▶ Ankle Plating System 3
- ▶ Calcaneal Plating System
- ▶ Forefoot-Midfoot Plating System
- ▶ Small Fragment Base Set



Conditions for Safe Scanning in the MRI

A patient with the plate, screw, pin, peg, and washer implants from the Ankle Plating, Ankle Plating 3, Calcaneal Plating, Forefoot-Midfoot Plating, and Small Fragment Base Set Systems may be safely scanned at 1.5 T or 3 T under the following conditions. Failure to follow these conditions may result in injury to the patient.

Foot and Ankle Portfolio

| Parameter | Condition |
|---|---|
| Nominal values(s) of Static Magnetic Field | 1.5 T or 3 T |
| Maximum Spatial Field Gradient | 19 T/m (1900 gauss/cm) |
| RF Excitation | Circularly Polarized (CP) |
| RF Transmit Coil Type | Whole body transmit coil, Head RF transmit-receive coil |
| Maximum Whole Body SAR | 2.0 W/kg or 1.0 W/kg* |
| Limits on Scan Duration | 2.0 W/kg whole body average SAR for 5 minutes of continuous RF (a sequence or back to back series/scan without breaks) / 1.0 W/kg whole body average SAR for 15 minutes of continuous RF (a sequence or back to back series/scan without breaks) |
| MR Image Artifact | The presence of this implant may produce an image artifact of 42 mm. |
| If more information about a specific parameter is not included, there are no conditions associated with that parameter. | |

*Whole body SAR of 1.0 W/kg may limit the availability and diagnostic quality of some MR procedures

Ankle Plating System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|---------------------------------------|
| 30-0232-S | 30-0232 | 3.5 mm x 8 mm Locking Hexalobe Screw |
| 30-0233-S | 30-0233 | 3.5 mm x 10 mm Locking Hexalobe Screw |
| 30-0234-S | 30-0234 | 3.5 mm x 12 mm Locking Hexalobe Screw |
| 30-0235-S | 30-0235 | 3.5 mm x 14 mm Locking Hexalobe Screw |
| 30-0236-S | 30-0236 | 3.5 mm x 16 mm Locking Hexalobe Screw |
| 30-0237-S | 30-0237 | 3.5 mm x 18 mm Locking Hexalobe Screw |
| 30-0238-S | 30-0238 | 3.5 mm x 20 mm Locking Hexalobe Screw |

Foot and Ankle Portfolio [continued]

Ankle Plating System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|--|
| 30-0239-S | 30-0239 | 3.5 mm x 22 mm Locking Hexalobe Screw |
| 30-0240-S | 30-0240 | 3.5 mm x 24 mm Locking Hexalobe Screw |
| 30-0241-S | 30-0241 | 3.5 mm x 26 mm Locking Hexalobe Screw |
| 30-0242-S | 30-0242 | 3.5 mm x 28 mm Locking Hexalobe Screw |
| 30-0243-S | 30-0243 | 3.5 mm x 30 mm Locking Hexalobe Screw |
| 30-0244-S | 30-0244 | 3.5 mm x 32 mm Locking Hexalobe Screw |
| 30-0245-S | 30-0245 | 3.5 mm x 34 mm Locking Hexalobe Screw |
| 30-0246-S | 30-0246 | 3.5 mm x 36 mm Locking Hexalobe Screw |
| 30-0247-S | 30-0247 | 3.5 mm x 38 mm Locking Hexalobe Screw |
| 30-0248-S | 30-0248 | 3.5 mm x 40 mm Locking Hexalobe Screw |
| 30-0249-S | 30-0249 | 3.5 mm x 45 mm Locking Hexalobe Screw |
| 30-0250-S | 30-0250 | 3.5 mm x 50 mm Locking Hexalobe Screw |
| 30-0251-S | 30-0251 | 3.5 mm x 55 mm Locking Hexalobe Screw |
| 30-0252-S | 30-0252 | 3.5 mm x 60 mm Locking Hexalobe Screw |
| 30-0255-S | 30-0255 | 3.5 mm x 8 mm Nonlocking Hexalobe Screw |
| 30-0256-S | 30-0256 | 3.5 mm x 10 mm Nonlocking Hexalobe Screw |
| 30-0257-S | 30-0257 | 3.5 mm x 12 mm Nonlocking Hexalobe Screw |
| 30-0258-S | 30-0258 | 3.5 mm x 14 mm Nonlocking Hexalobe Screw |
| 30-0259-S | 30-0259 | 3.5 mm x 16 mm Nonlocking Hexalobe Screw |
| 30-0260-S | 30-0260 | 3.5 mm x 18 mm Nonlocking Hexalobe Screw |
| 30-0261-S | 30-0261 | 3.5 mm x 20 mm Nonlocking Hexalobe Screw |
| 30-0262-S | 30-0262 | 3.5 mm x 22 mm Nonlocking Hexalobe Screw |
| 30-0263-S | 30-0263 | 3.5 mm x 24 mm Nonlocking Hexalobe Screw |
| 30-0264-S | 30-0264 | 3.5 mm x 26 mm Nonlocking Hexalobe Screw |
| 30-0265-S | 30-0265 | 3.5 mm x 28 mm Nonlocking Hexalobe Screw |
| 30-0266-S | 30-0266 | 3.5 mm x 30 mm Nonlocking Hexalobe Screw |
| 30-0267-S | 30-0267 | 3.5 mm x 32 mm Nonlocking Hexalobe Screw |
| 30-0268-S | 30-0268 | 3.5 mm x 34 mm Nonlocking Hexalobe Screw |
| 30-0269-S | 30-0269 | 3.5 mm x 36 mm Nonlocking Hexalobe Screw |

Foot and Ankle Portfolio [continued]

Ankle Plating System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|--|
| 30-0270-S | 30-0270 | 3.5 mm x 38 mm Nonlocking Hexalobe Screw |
| 30-0271-S | 30-0271 | 3.5 mm x 40 mm Nonlocking Hexalobe Screw |
| 30-0272-S | 30-0272 | 3.5 mm x 45 mm Nonlocking Hexalobe Screw |
| 30-0273-S | 30-0273 | 3.5 mm x 50 mm Nonlocking Hexalobe Screw |
| 30-0274-S | 30-0274 | 3.5 mm x 55 mm Nonlocking Hexalobe Screw |
| 30-0275-S | 30-0275 | 3.5 mm x 60 mm Nonlocking Hexalobe Screw |
| 30-0278-S | 30-0278 | 3.0 mm x 8 mm Locking Hexalobe Screw |
| 30-0279-S | 30-0279 | 3.0 mm x 10 mm Locking Hexalobe Screw |
| 30-0280-S | 30-0280 | 3.0 mm x 12 mm Locking Hexalobe Screw |
| 30-0281-S | 30-0281 | 3.0 mm x 14 mm Locking Hexalobe Screw |
| 30-0282-S | 30-0282 | 3.0 mm x 16 mm Locking Hexalobe Screw |
| 30-0283-S | 30-0283 | 3.0 mm x 18 mm Locking Hexalobe Screw |
| 30-0284-S | 30-0284 | 3.0 mm x 20 mm Locking Hexalobe Screw |
| 30-0285-S | 30-0285 | 3.0 mm x 22 mm Locking Hexalobe Screw |
| 30-0286-S | 30-0286 | 3.0 mm x 24 mm Locking Hexalobe Screw |
| 30-0287-S | 30-0287 | 3.0 mm x 26 mm Locking Hexalobe Screw |
| 30-0288-S | 30-0288 | 3.0 mm x 28 mm Locking Hexalobe Screw |
| 30-0289-S | 30-0289 | 3.0 mm x 30 mm Locking Hexalobe Screw |
| 30-0290-S | 30-0290 | 3.0 mm x 32 mm Locking Hexalobe Screw |
| 30-0291-S | 30-0291 | 3.0 mm x 34 mm Locking Hexalobe Screw |
| 30-0292-S | 30-0292 | 3.0 mm x 36 mm Locking Hexalobe Screw |
| 30-0293-S | 30-0293 | 3.0 mm x 38 mm Locking Hexalobe Screw |
| 30-0294-S | 30-0294 | 3.0 mm x 40 mm Locking Hexalobe Screw |
| 30-0295-S | 30-0295 | 3.0 mm x 45 mm Locking Hexalobe Screw |
| 30-0296-S | 30-0296 | 3.0 mm x 50 mm Locking Hexalobe Screw |
| 30-0297-S | 30-0297 | 3.0 mm x 55 mm Locking Hexalobe Screw |
| 30-0301-S | 30-0301 | 3.0 mm x 8 mm Nonlocking Hexalobe Screw |
| 30-0302-S | 30-0302 | 3.0 mm x 10 mm Nonlocking Hexalobe Screw |
| 30-0303-S | 30-0303 | 3.0 mm x 12 mm Nonlocking Hexalobe Screw |

Foot and Ankle Portfolio [continued]

Ankle Plating System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|--|
| 30-0304-S | 30-0304 | 3.0 mm x 14 mm Nonlocking Hexalobe Screw |
| 30-0305-S | 30-0305 | 3.0 mm x 16 mm Nonlocking Hexalobe Screw |
| 30-0306-S | 30-0306 | 3.0 mm x 18 mm Nonlocking Hexalobe Screw |
| 30-0307-S | 30-0307 | 3.0 mm x 20 mm Nonlocking Hexalobe Screw |
| 30-0308-S | 30-0308 | 3.0 mm x 22 mm Nonlocking Hexalobe Screw |
| 30-0309-S | 30-0309 | 3.0 mm x 24 mm Nonlocking Hexalobe Screw |
| 30-0310-S | 30-0310 | 3.0 mm x 26 mm Nonlocking Hexalobe Screw |
| 30-0311-S | 30-0311 | 3.0 mm x 28 mm Nonlocking Hexalobe Screw |
| 30-0312-S | 30-0312 | 3.0 mm x 30 mm Nonlocking Hexalobe Screw |
| 30-0313-S | 30-0313 | 3.0 mm x 32 mm Nonlocking Hexalobe Screw |
| 30-0314-S | 30-0314 | 3.0 mm x 34 mm Nonlocking Hexalobe Screw |
| 30-0315-S | 30-0315 | 3.0 mm x 36 mm Nonlocking Hexalobe Screw |
| 30-0316-S | 30-0316 | 3.0 mm x 38 mm Nonlocking Hexalobe Screw |
| 30-0317-S | 30-0317 | 3.0 mm x 40 mm Nonlocking Hexalobe Screw |
| 30-0318-S | 30-0318 | 3.0 mm x 45 mm Nonlocking Hexalobe Screw |
| 30-0319-S | 30-0319 | 3.0 mm x 50 mm Nonlocking Hexalobe Screw |
| 30-0320-S | 30-0320 | 3.0 mm x 55 mm Nonlocking Hexalobe Screw |
| 70-0145-S | 70-0145 | Low-profile Locking Lateral Fibula Plate 5-Hole |
| 70-0147-S | 70-0147 | Low-profile Locking Lateral Fibula Plate 7-Hole |
| 70-0149-S | 70-0149 | Low-profile Locking Lateral Fibula Plate 9-Hole |
| 70-0151-S | 70-0151 | Low-profile Locking Lateral Fibula Plate 11-Hole |
| 70-0153-S | 70-0153 | Low-profile Locking Lateral Fibula Plate 13-Hole |
| 70-0169-S | 70-0169 | Locking Lateral Fibula Plate 9-Hole |
| 70-0171-S | 70-0171 | Locking Lateral Fibula Plate 11-Hole |
| 70-0173-S | 70-0173 | Locking Lateral Fibula Plate 13-Hole |
| 70-0227-S | 70-0227 | Low-profile Locking Medial Tibia Plate 7-Hole |
| 70-0229-S | 70-0229 | Low-profile Locking Medial Tibia Plate 9-Hole |
| 70-0245-S | 70-0245 | Low-profile Locking Anterior Tibia Plate 5-Hole |
| 70-0247-S | 70-0247 | Low-profile Locking Anterior Tibia Plate 7-Hole |

Foot and Ankle Portfolio [continued]

Ankle Plating System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|-----------------------------------|
| CA-4120-S | CA-4120 | 4.0 mm x 12.0 mm Cancellous Screw |
| CA-4140-S | CA-4140 | 4.0 mm x 14.0 mm Cancellous Screw |
| CA-4160-S | CA-4160 | 4.0 mm x 16.0 mm Cancellous Screw |
| CA-4180-S | CA-4180 | 4.0 mm x 18.0 mm Cancellous Screw |
| CA-4200-S | CA-4200 | 4.0 mm x 20.0 mm Cancellous Screw |
| CA-4220-S | CA-4220 | 4.0 mm x 22.0 mm Cancellous Screw |
| CA-4240-S | CA-4240 | 4.0 mm x 24.0 mm Cancellous Screw |
| CA-4260-S | CA-4260 | 4.0 mm x 26.0 mm Cancellous Screw |
| CA-4280-S | CA-4280 | 4.0 mm x 28.0 mm Cancellous Screw |
| CA-4300-S | CA-4300 | 4.0 mm x 30.0 mm Cancellous Screw |
| CA-4350-S | CA-4350 | 4.0 mm x 35.0 mm Cancellous Screw |
| CA-4400-S | CA-4400 | 4.0 mm x 40.0 mm Cancellous Screw |
| CA-4450-S | CA-4450 | 4.0 mm x 45.0 mm Cancellous Screw |
| CA-4500-S | CA-4500 | 4.0 mm x 50.0 mm Cancellous Screw |
| CA-4550-S | CA-4550 | 4.0 mm x 55.0 mm Cancellous Screw |
| CA-4600-S | CA-4600 | 4.0 mm x 60.0 mm Cancellous Screw |
| CO-2708-S | CO-2708 | 2.7 mm x 8.0 mm Cortical Screw |
| CO-2710-S | CO-2710 | 2.7 mm x 10.0 mm Cortical Screw |
| CO-2712-S | CO-2712 | 2.7 mm x 12.0 mm Cortical Screw |
| CO-2714-S | CO-2714 | 2.7 mm x 14.0 mm Cortical Screw |
| CO-2716-S | CO-2716 | 2.7 mm x 16.0 mm Cortical Screw |
| CO-2718-S | CO-2718 | 2.7 mm x 18.0 mm Cortical Screw |
| CO-2720-S | CO-2720 | 2.7 mm x 20.0 mm Cortical Screw |
| CO-2722-S | CO-2722 | 2.7 mm x 22.0 mm Cortical Screw |
| CO-2724-S | CO-2724 | 2.7 mm x 24.0 mm Cortical Screw |
| CO-2726-S | CO-2726 | 2.7 mm x 26.0 mm Cortical Screw |
| CO-2728-S | CO-2728 | 2.7 mm x 28.0 mm Cortical Screw |
| CO-2730-S | CO-2730 | 2.7 mm x 30.0 mm Cortical Screw |
| CO-2732-S | CO-2732 | 2.7 mm x 32.0 mm Cortical Screw |

Foot and Ankle Portfolio [continued]

Ankle Plating System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|--|
| CO-2734-S | CO-2734 | 2.7 mm x 34.0 mm Cortical Screw |
| CO-2736-S | CO-2736 | 2.7 mm x 36.0 mm Cortical Screw |
| CO-2738-S | CO-2738 | 2.7 mm x 38.0 mm Cortical Screw |
| CO-2740-S | CO-2740 | 2.7 mm x 40.0 mm Cortical Screw |
| CO-2745-S | CO-2745 | 2.7 mm x 45.0 mm Cortical Screw |
| CO-2750-S | CO-2750 | 2.7 mm x 50.0 mm Cortical Screw |
| CO-2755-S | CO-2755 | 2.7 mm x 55.0 mm Cortical Screw |
| CO-3080-S | CO-3080 | 3.5 mm x 8.0 mm Cortical Screw |
| CO-3100-S | CO-3100 | 3.5 mm x 10.0 mm Cortical Screw |
| CO-3120-S | CO-3120 | 3.5 mm x 12.0 mm Cortical Screw |
| CO-3140-S | CO-3140 | 3.5 mm x 14.0 mm Cortical Screw |
| CO-3160-S | CO-3160 | 3.5 mm x 16.0 mm Cortical Screw |
| CO-3180-S | CO-3180 | 3.5 mm x 18.0 mm Cortical Screw |
| CO-3200-S | CO-3200 | 3.5 mm x 20.0 mm Cortical Screw |
| CO-3220-S | CO-3220 | 3.5 mm x 22.0 mm Cortical Screw |
| CO-3240-S | CO-3240 | 3.5 mm x 24.0 mm Cortical Screw |
| CO-3260-S | CO-3260 | 3.5 mm x 26.0 mm Cortical Screw |
| CO-3280-S | CO-3280 | 3.5 mm x 28.0 mm Cortical Screw |
| CO-3300-S | CO-3300 | 3.5 mm x 30.0 mm Cortical Screw |
| CO-3320-S | CO-3320 | 3.5 mm x 32.0 mm Cortical Screw |
| CO-3340-S | CO-3340 | 3.5 mm x 34.0 mm Cortical Screw |
| CO-3360-S | CO-3360 | 3.5 mm x 36.0 mm Cortical Screw |
| CO-3380-S | CO-3380 | 3.5 mm x 38.0 mm Cortical Screw |
| CO-3400-S | CO-3400 | 3.5 mm x 40.0 mm Cortical Screw |
| CO-3450-S | CO-3450 | 3.5 mm x 45.0 mm Cortical Screw |
| CO-3500-S | CO-3500 | 3.5 mm x 50.0 mm Cortical Screw |
| CO-3550-S | CO-3550 | 3.5 mm x 55.0 mm Cortical Screw |
| CO-3600-S | CO-3600 | 3.5 mm x 60.0 mm Cortical Screw |
| COL-2080-S | COL-2080 | 2.7 mm x 8.0 mm Locking Cortical Screw |

Foot and Ankle Portfolio [continued]

Ankle Plating System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|---|
| COL-2100-S | COL-2100 | 2.7 mm x 10.0 mm Locking Cortical Screw |
| COL-2120-S | COL-2120 | 2.7 mm x 12.0 mm Locking Cortical Screw |
| COL-2140-S | COL-2140 | 2.7 mm x 14.0 mm Locking Cortical Screw |
| COL-2160-S | COL-2160 | 2.7 mm x 16.0 mm Locking Cortical Screw |
| COL-2180-S | COL-2180 | 2.7 mm x 18.0 mm Locking Cortical Screw |
| COL-2200-S | COL-2200 | 2.7 mm x 20.0 mm Locking Cortical Screw |
| COL-2220-S | COL-2220 | 2.7 mm x 22.0 mm Locking Cortical Screw |
| COL-2240-S | COL-2240 | 2.7 mm x 24.0 mm Locking Cortical Screw |
| COL-2260-S | COL-2260 | 2.7 mm x 26.0 mm Locking Cortical Screw |
| COL-2280-S | COL-2280 | 2.7 mm x 28.0 mm Locking Cortical Screw |
| COL-2300-S | COL-2300 | 2.7 mm x 30.0 mm Locking Cortical Screw |
| COL-2320-S | COL-2320 | 2.7 mm x 32.0 mm Locking Cortical Screw |
| COL-2340-S | COL-2340 | 2.7 mm x 34.0 mm Locking Cortical Screw |
| COL-2360-S | COL-2360 | 2.7 mm x 36.0 mm Locking Cortical Screw |
| COL-2380-S | COL-2380 | 2.7 mm x 38.0 mm Locking Cortical Screw |
| COL-2400-S | COL-2400 | 2.7 mm x 40.0 mm Locking Cortical Screw |
| COL-2450-S | COL-2450 | 2.7 mm x 45.0 mm Locking Cortical Screw |
| COL-2500-S | COL-2500 | 2.7 mm x 50.0 mm Locking Cortical Screw |
| COL-2550-S | COL-2550 | 2.7 mm x 55.0 mm Locking Cortical Screw |
| COL-3080-S | COL-3080 | 3.5 mm x 8.0 mm Locking Cortical Screw |
| COL-3100-S | COL-3100 | 3.5 mm x 10.0 mm Locking Cortical Screw |
| COL-3120-S | COL-3120 | 3.5 mm x 12.0 mm Locking Cortical Screw |
| COL-3140-S | COL-3140 | 3.5 mm x 14.0 mm Locking Cortical Screw |
| COL-3160-S | COL-3160 | 3.5 mm x 16.0 mm Locking Cortical Screw |
| COL-3180-S | COL-3180 | 3.5 mm x 18.0 mm Locking Cortical Screw |
| COL-3200-S | COL-3200 | 3.5 mm x 20.0 mm Locking Cortical Screw |
| COL-3220-S | COL-3220 | 3.5 mm x 22.0 mm Locking Cortical Screw |
| COL-3240-S | COL-3240 | 3.5 mm x 24.0 mm Locking Cortical Screw |
| COL-3260-S | COL-3260 | 3.5 mm x 26.0 mm Locking Cortical Screw |

Foot and Ankle Portfolio [continued]

Ankle Plating System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|---|
| COL-3280-S | COL-3280 | 3.5 mm x 28.0 mm Locking Cortical Screw |
| COL-3300-S | COL-3300 | 3.5 mm x 30.0 mm Locking Cortical Screw |
| COL-3320-S | COL-3320 | 3.5 mm x 32.0 mm Locking Cortical Screw |
| COL-3340-S | COL-3340 | 3.5 mm x 34.0 mm Locking Cortical Screw |
| COL-3360-S | COL-3360 | 3.5 mm x 36.0 mm Locking Cortical Screw |
| COL-3380-S | COL-3380 | 3.5 mm x 38.0 mm Locking Cortical Screw |
| COL-3400-S | COL-3400 | 3.5 mm x 40.0 mm Locking Cortical Screw |
| COL-3450-S | COL-3450 | 3.5 mm x 45.0 mm Locking Cortical Screw |
| COL-3500-S | COL-3500 | 3.5 mm x 50.0 mm Locking Cortical Screw |
| COL-3550-S | COL-3550 | 3.5 mm x 55.0 mm Locking Cortical Screw |
| COL-3600-S | COL-3600 | 3.5 mm x 60.0 mm Locking Cortical Screw |

Ankle Plating System 3

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|---|
| N/A | 30-0097 | 50.0 mm Tension Band Pin |
| N/A | 30-0098 | 70.0 mm Tension Band Pin |
| N/A | 30-0099 | 90.0 mm Tension Band Pin |
| N/A | 30-0166 | 2 mm x 50 mm Tension Band Pin |
| N/A | 30-0167 | 2 mm x 70 mm Tension Band Pin |
| N/A | 30-0168 | 2 mm x 90 mm Tension Band Pin |
| N/A | 3005-40010 | 4.0 mm x 10 mm Cannulated Screw, Short Thread |
| N/A | 3005-40012 | 4.0 mm x 12 mm Cannulated Screw, Short Thread |
| N/A | 3005-40014 | 4.0 mm x 14 mm Cannulated Screw, Short Thread |
| N/A | 3005-40016 | 4.0 mm x 16 mm Cannulated Screw, Short Thread |
| N/A | 3005-40018 | 4.0 mm x 18 mm Cannulated Screw, Short Thread |
| N/A | 3005-40020 | 4.0 mm x 20 mm Cannulated Screw, Short Thread |
| N/A | 3005-40022 | 4.0 mm x 22 mm Cannulated Screw, Short Thread |
| N/A | 3005-40024 | 4.0 mm x 24 mm Cannulated Screw, Short Thread |
| N/A | 3005-40026 | 4.0 mm x 26 mm Cannulated Screw, Short Thread |

Foot and Ankle Portfolio [continued]

Ankle Plating System 3

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|---|
| N/A | 3005-40028 | 4.0 mm x 28 mm Cannulated Screw, Short Thread |
| N/A | 3005-40030 | 4.0 mm x 30 mm Cannulated Screw, Short Thread |
| N/A | 3005-40032 | 4.0 mm x 32 mm Cannulated Screw, Short Thread |
| N/A | 3005-40034 | 4.0 mm x 34 mm Cannulated Screw, Short Thread |
| N/A | 3005-40036 | 4.0 mm x 36 mm Cannulated Screw, Short Thread |
| N/A | 3005-40038 | 4.0 mm x 38 mm Cannulated Screw, Short Thread |
| N/A | 3005-40040 | 4.0 mm x 40 mm Cannulated Screw, Short Thread |
| N/A | 3005-40042 | 4.0 mm x 42 mm Cannulated Screw, Short Thread |
| N/A | 3005-40044 | 4.0 mm x 44 mm Cannulated Screw, Short Thread |
| N/A | 3005-40046 | 4.0 mm x 46 mm Cannulated Screw, Short Thread |
| N/A | 3005-40048 | 4.0 mm x 48 mm Cannulated Screw, Short Thread |
| N/A | 3005-40050 | 4.0 mm x 50 mm Cannulated Screw, Short Thread |
| N/A | 3005-40052 | 4.0 mm x 52 mm Cannulated Screw, Short Thread |
| N/A | 3005-40054 | 4.0 mm x 54 mm Cannulated Screw, Short Thread |
| N/A | 3005-40056 | 4.0 mm x 56 mm Cannulated Screw, Short Thread |
| N/A | 3005-40058 | 4.0 mm x 58 mm Cannulated Screw, Short Thread |
| N/A | 3005-40060 | 4.0 mm x 60 mm Cannulated Screw, Short Thread |
| N/A | 3005-40064 | 4.0 mm x 64 mm Cannulated Screw, Short Thread |
| N/A | 3005-40068 | 4.0 mm x 68 mm Cannulated Screw, Short Thread |
| N/A | 3005-40072 | 4.0 mm x 72 mm Cannulated Screw, Short Thread |
| N/A | 3006-40016 | 4.0 mm x 16 mm Cannulated Screw, Long Thread |
| N/A | 3006-40018 | 4.0 mm x 18 mm Cannulated Screw, Long Thread |
| N/A | 3006-40020 | 4.0 mm x 20 mm Cannulated Screw, Long Thread |
| N/A | 3006-40022 | 4.0 mm x 22 mm Cannulated Screw, Long Thread |
| N/A | 3006-40024 | 4.0 mm x 24 mm Cannulated Screw, Long Thread |
| N/A | 3006-40026 | 4.0 mm x 26 mm Cannulated Screw, Long Thread |
| N/A | 3006-40028 | 4.0 mm x 28 mm Cannulated Screw, Long Thread |
| N/A | 3006-40030 | 4.0 mm x 30 mm Cannulated Screw, Long Thread |
| N/A | 3006-40032 | 4.0 mm x 32 mm Cannulated Screw, Long Thread |

Foot and Ankle Portfolio [continued]

Ankle Plating System 3

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|---|
| N/A | 3006-40034 | 4.0 mm x 34 mm Cannulated Screw, Long Thread |
| 3006-40036-S | 3006-40036 | 4.0 mm x 36 mm Cannulated Screw, Long Thread |
| N/A | 3006-40038 | 4.0 mm x 38 mm Cannulated Screw, Long Thread |
| N/A | 3006-40040 | 4.0 mm x 40 mm Cannulated Screw, Long Thread |
| 3006-40042-S | 3006-40042 | 4.0 mm x 42 mm Cannulated Screw, Long Thread |
| N/A | 3006-40044 | 4.0 mm x 44 mm Cannulated Screw, Long Thread |
| N/A | 3006-40046 | 4.0 mm x 46 mm Cannulated Screw, Long Thread |
| 3006-40048-S | 3006-40048 | 4.0 mm x 48 mm Cannulated Screw, Long Thread |
| N/A | 3006-40050 | 4.0 mm x 50 mm Cannulated Screw, Long Thread |
| N/A | 3006-40052 | 4.0 mm x 52 mm Cannulated Screw, Long Thread |
| N/A | 3006-40054 | 4.0 mm x 54 mm Cannulated Screw, Long Thread |
| N/A | 3006-40056 | 4.0 mm x 56 mm Cannulated Screw, Long Thread |
| N/A | 3006-40058 | 4.0 mm x 58 mm Cannulated Screw, Long Thread |
| N/A | 3006-40060 | 4.0 mm x 60 mm Cannulated Screw, Long Thread |
| N/A | 3006-40064 | 4.0 mm x 64 mm Cannulated Screw, Long Thread |
| N/A | 3006-40068 | 4.0 mm x 68 mm Cannulated Screw, Long Thread |
| N/A | 3006-40072 | 4.0 mm x 72 mm Cannulated Screw, Long Thread |
| 3014-23026-S | 3014-23026 | 2.3 mm x 26 mm Locking Cortical Peg |
| 3014-23036-S | 3014-23036 | 2.3 mm x 36 mm Locking Cortical Peg |
| 3014-23046-S | 3014-23046 | 2.3 mm x 46 mm Locking Cortical Peg |
| 7003-07036-S | 7003-07036 | Cannulated Screw Washer 7.0 mm OD x 3.6 mm ID |
| 7007-0104L-S | 7007-0104L | Lateral Fibula Plate 4-Hole Left |
| 7007-0104R-S | 7007-0104R | Lateral Fibula Plate 4-Hole Right |
| 7007-0105L-S | 7007-0105L | Lateral Fibula Plate 5-Hole Left |
| 7007-0105R-S | 7007-0105R | Lateral Fibula Plate 5-Hole Right |
| 7007-0106L-S | 7007-0106L | Lateral Fibula Plate 6-Hole Left |
| 7007-0106R-S | 7007-0106R | Lateral Fibula Plate 6-Hole Right |
| 7007-0107L-S | 7007-0107L | Lateral Fibula Plate 7-Hole Left |
| 7007-0107R-S | 7007-0107R | Lateral Fibula Plate 7-Hole Right |

Foot and Ankle Portfolio [continued]

Ankle Plating System 3

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|--|
| 7007-0109L-S | 7007-0109L | Lateral Fibula Plate 9-Hole Left |
| 7007-0109R-S | 7007-0109R | Lateral Fibula Plate 9-Hole Right |
| 7007-0111L-S | N/A | Lateral Fibula Plate 11-Hole Left |
| 7007-0111R-S | N/A | Lateral Fibula Plate 11-Hole Right |
| 7007-0113L-S | N/A | Lateral Fibula Plate 13-Hole Left |
| 7007-0113R-S | N/A | Lateral Fibula Plate 13-Hole Right |
| 7007-0203L-S | 7007-0203L | Posterolateral Fibula Plate 3-Hole Left |
| 7007-0203R-S | 7007-0203R | Posterolateral Fibula Plate 3-Hole Right |
| 7007-0204L-S | 7007-0204L | Posterolateral Fibula Plate 4-Hole Left |
| 7007-0204R-S | 7007-0204R | Posterolateral Fibula Plate 4-Hole Right |
| 7007-0205L-S | 7007-0205L | Posterolateral Fibula Plate 5-Hole Left |
| 7007-0205R-S | 7007-0205R | Posterolateral Fibula Plate 5-Hole Right |
| 7007-0206L-S | 7007-0206L | Posterolateral Fibula Plate 6-Hole Left |
| 7007-0206R-S | 7007-0206R | Posterolateral Fibula Plate 6-Hole Right |
| 7007-0207L-S | 7007-0207L | Posterolateral Fibula Plate 7-Hole Left |
| 7007-0207R-S | 7007-0207R | Posterolateral Fibula Plate 7-Hole Right |
| 7007-0303L-S | 7007-0303L | Posterolateral Distal Tibia Plate 3-Hole Left |
| 7007-0303R-S | 7007-0303R | Posterolateral Distal Tibia Plate 3-Hole Right |
| 7007-0304L-S | 7007-0304L | Posterolateral Distal Tibia Plate 4-Hole Left |
| 7007-0304R-S | 7007-0304R | Posterolateral Distal Tibia Plate 4-Hole Right |
| 7007-0403L-S | 7007-0403L | Posteromedial Distal Tibia Plate 3-Hole Left |
| 7007-0403R-S | 7007-0403R | Posteromedial Distal Tibia Plate 3-Hole Right |
| 7007-0504-S | 7007-0504 | Medial Anti-Glide Plate 4-Hole |
| 7007-0602-S | 7007-0602 | Hook Plate 2-Hole |
| 7007-0603-S | 7007-0603 | Hook Plate 3-Hole |
| 7007-0702-S | 7007-0702 | Locking Peg Hook Plate 2-Hole |
| 7007-0703-S | 7007-0703 | Locking Peg Hook Plate 3-Hole |
| AI-0010-S | AI-0010 | 10 mm AcuTwist® Acutrak Compression Screw |
| AI-0012-S | AI-0012 | 12 mm AcuTwist® Acutrak Compression Screw |

Foot and Ankle Portfolio [continued]

Ankle Plating System 3

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|--|
| AI-0014-S | AI-0014 | 14 mm AcuTwist Acutrak Compression Screw |
| AI-0016-S | AI-0016 | 16 mm AcuTwist Acutrak Compression Screw |
| AI-0018-S | AI-0018 | 18 mm AcuTwist Acutrak Compression Screw |
| AI-0020-S | AI-0020 | 20 mm AcuTwist Acutrak Compression Screw |
| AI-0022-S | AI-0022 | 22 mm AcuTwist Acutrak Compression Screw |
| AI-0024-S | AI-0024 | 24 mm AcuTwist Acutrak Compression Screw |
| AI-0026-S | AI-0026 | 26 mm AcuTwist Acutrak Compression Screw |
| AI-0028-S | AI-0028 | 28 mm AcuTwist Acutrak Compression Screw |
| AI-0030-S | AI-0030 | 30 mm AcuTwist Acutrak Compression Screw |

Calcaneal Plating System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|---------------------------------------|
| 30-0232-S | 30-0232 | 3.5 mm x 8 mm Locking Hexalobe Screw |
| 30-0233-S | 30-0233 | 3.5 mm x 10 mm Locking Hexalobe Screw |
| 30-0234-S | 30-0234 | 3.5 mm x 12 mm Locking Hexalobe Screw |
| 30-0235-S | 30-0235 | 3.5 mm x 14 mm Locking Hexalobe Screw |
| 30-0236-S | 30-0236 | 3.5 mm x 16 mm Locking Hexalobe Screw |
| 30-0237-S | 30-0237 | 3.5 mm x 18 mm Locking Hexalobe Screw |
| 30-0238-S | 30-0238 | 3.5 mm x 20 mm Locking Hexalobe Screw |
| 30-0239-S | 30-0239 | 3.5 mm x 22 mm Locking Hexalobe Screw |
| 30-0240-S | 30-0240 | 3.5 mm x 24 mm Locking Hexalobe Screw |
| 30-0241-S | 30-0241 | 3.5 mm x 26 mm Locking Hexalobe Screw |
| 30-0242-S | 30-0242 | 3.5 mm x 28 mm Locking Hexalobe Screw |
| 30-0243-S | 30-0243 | 3.5 mm x 30 mm Locking Hexalobe Screw |
| 30-0244-S | 30-0244 | 3.5 mm x 32 mm Locking Hexalobe Screw |
| 30-0245-S | 30-0245 | 3.5 mm x 34 mm Locking Hexalobe Screw |
| 30-0246-S | 30-0246 | 3.5 mm x 36 mm Locking Hexalobe Screw |
| 30-0247-S | 30-0247 | 3.5 mm x 38 mm Locking Hexalobe Screw |
| 30-0248-S | 30-0248 | 3.5 mm x 40 mm Locking Hexalobe Screw |

Foot and Ankle Portfolio [continued]

Calcaneal Plating System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|--|
| 30-0249-S | 30-0249 | 3.5 mm x 45 mm Locking Hexalobe Screw |
| 30-0250-S | 30-0250 | 3.5 mm x 50 mm Locking Hexalobe Screw |
| 30-0251-S | 30-0251 | 3.5 mm x 55 mm Locking Hexalobe Screw |
| 30-0252-S | 30-0252 | 3.5 mm x 60 mm Locking Hexalobe Screw |
| 30-0255-S | 30-0255 | 3.5 mm x 8 mm Nonlocking Hexalobe Screw |
| 30-0256-S | 30-0256 | 3.5 mm x 10 mm Nonlocking Hexalobe Screw |
| 30-0257-S | 30-0257 | 3.5 mm x 12 mm Nonlocking Hexalobe Screw |
| 30-0258-S | 30-0258 | 3.5 mm x 14 mm Nonlocking Hexalobe Screw |
| 30-0259-S | 30-0259 | 3.5 mm x 16 mm Nonlocking Hexalobe Screw |
| 30-0260-S | 30-0260 | 3.5 mm x 18 mm Nonlocking Hexalobe Screw |
| 30-0261-S | 30-0261 | 3.5 mm x 20 mm Nonlocking Hexalobe Screw |
| 30-0262-S | 30-0262 | 3.5 mm x 22 mm Nonlocking Hexalobe Screw |
| 30-0263-S | 30-0263 | 3.5 mm x 24 mm Nonlocking Hexalobe Screw |
| 30-0264-S | 30-0264 | 3.5 mm x 26 mm Nonlocking Hexalobe Screw |
| 30-0265-S | 30-0265 | 3.5 mm x 28 mm Nonlocking Hexalobe Screw |
| 30-0266-S | 30-0266 | 3.5 mm x 30 mm Nonlocking Hexalobe Screw |
| 30-0267-S | 30-0267 | 3.5 mm x 32 mm Nonlocking Hexalobe Screw |
| 30-0268-S | 30-0268 | 3.5 mm x 34 mm Nonlocking Hexalobe Screw |
| 30-0269-S | 30-0269 | 3.5 mm x 36 mm Nonlocking Hexalobe Screw |
| 30-0270-S | 30-0270 | 3.5 mm x 38 mm Nonlocking Hexalobe Screw |
| 30-0271-S | 30-0271 | 3.5 mm x 40 mm Nonlocking Hexalobe Screw |
| 30-0272-S | 30-0272 | 3.5 mm x 45 mm Nonlocking Hexalobe Screw |
| 30-0273-S | 30-0273 | 3.5 mm x 50 mm Nonlocking Hexalobe Screw |
| 30-0274-S | 30-0274 | 3.5 mm x 55 mm Nonlocking Hexalobe Screw |
| 30-0275-S | 30-0275 | 3.5 mm x 60 mm Nonlocking Hexalobe Screw |
| 30-0278-S | 30-0278 | 3.0 mm x 8 mm Locking Hexalobe Screw |
| 30-0279-S | 30-0279 | 3.0 mm x 10 mm Locking Hexalobe Screw |
| 30-0280-S | 30-0280 | 3.0 mm x 12 mm Locking Hexalobe Screw |
| 30-0281-S | 30-0281 | 3.0 mm x 14 mm Locking Hexalobe Screw |

Foot and Ankle Portfolio [continued]

Calcaneal Plating System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|--|
| 30-0282-S | 30-0282 | 3.0 mm x 16 mm Locking Hexalobe Screw |
| 30-0283-S | 30-0283 | 3.0 mm x 18 mm Locking Hexalobe Screw |
| 30-0284-S | 30-0284 | 3.0 mm x 20 mm Locking Hexalobe Screw |
| 30-0285-S | 30-0285 | 3.0 mm x 22 mm Locking Hexalobe Screw |
| 30-0286-S | 30-0286 | 3.0 mm x 24 mm Locking Hexalobe Screw |
| 30-0287-S | 30-0287 | 3.0 mm x 26 mm Locking Hexalobe Screw |
| 30-0288-S | 30-0288 | 3.0 mm x 28 mm Locking Hexalobe Screw |
| 30-0289-S | 30-0289 | 3.0 mm x 30 mm Locking Hexalobe Screw |
| 30-0290-S | 30-0290 | 3.0 mm x 32 mm Locking Hexalobe Screw |
| 30-0291-S | 30-0291 | 3.0 mm x 34 mm Locking Hexalobe Screw |
| 30-0292-S | 30-0292 | 3.0 mm x 36 mm Locking Hexalobe Screw |
| 30-0293-S | 30-0293 | 3.0 mm x 38 mm Locking Hexalobe Screw |
| 30-0294-S | 30-0294 | 3.0 mm x 40 mm Locking Hexalobe Screw |
| 30-0295-S | 30-0295 | 3.0 mm x 45 mm Locking Hexalobe Screw |
| 30-0296-S | 30-0296 | 3.0 mm x 50 mm Locking Hexalobe Screw |
| 30-0297-S | 30-0297 | 3.0 mm x 55 mm Locking Hexalobe Screw |
| 30-0301-S | 30-0301 | 3.0 mm x 8 mm Nonlocking Hexalobe Screw |
| 30-0302-S | 30-0302 | 3.0 mm x 10 mm Nonlocking Hexalobe Screw |
| 30-0303-S | 30-0303 | 3.0 mm x 12 mm Nonlocking Hexalobe Screw |
| 30-0304-S | 30-0304 | 3.0 mm x 14 mm Nonlocking Hexalobe Screw |
| 30-0305-S | 30-0305 | 3.0 mm x 16 mm Nonlocking Hexalobe Screw |
| 30-0306-S | 30-0306 | 3.0 mm x 18 mm Nonlocking Hexalobe Screw |
| 30-0307-S | 30-0307 | 3.0 mm x 20 mm Nonlocking Hexalobe Screw |
| 30-0308-S | 30-0308 | 3.0 mm x 22 mm Nonlocking Hexalobe Screw |
| 30-0309-S | 30-0309 | 3.0 mm x 24 mm Nonlocking Hexalobe Screw |
| 30-0310-S | 30-0310 | 3.0 mm x 26 mm Nonlocking Hexalobe Screw |
| 30-0311-S | 30-0311 | 3.0 mm x 28 mm Nonlocking Hexalobe Screw |
| 30-0312-S | 30-0312 | 3.0 mm x 30 mm Nonlocking Hexalobe Screw |
| 30-0313-S | 30-0313 | 3.0 mm x 32 mm Nonlocking Hexalobe Screw |

Foot and Ankle Portfolio [continued]

Calcaneal Plating System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|---|
| 30-0314-S | 30-0314 | 3.0 mm x 34 mm Nonlocking Hexalobe Screw |
| 30-0315-S | 30-0315 | 3.0 mm x 36 mm Nonlocking Hexalobe Screw |
| 30-0316-S | 30-0316 | 3.0 mm x 38 mm Nonlocking Hexalobe Screw |
| 30-0317-S | 30-0317 | 3.0 mm x 40 mm Nonlocking Hexalobe Screw |
| 30-0318-S | 30-0318 | 3.0 mm x 45 mm Nonlocking Hexalobe Screw |
| 30-0319-S | 30-0319 | 3.0 mm x 50 mm Nonlocking Hexalobe Screw |
| 30-0320-S | 30-0320 | 3.0 mm x 55 mm Nonlocking Hexalobe Screw |
| N/A | 35-0002 | 4.0 mm Steinmann Pin, Smooth |
| N/A | 35-0003 | 5.0 mm Steinmann Pin, Smooth |
| 70-0022-S | 70-0022 | Locking Calcaneal Plate, Small Left |
| 70-0023-S | 70-0023 | Locking Calcaneal Plate, Small Right |
| 70-0024-S | 70-0024 | Locking Calcaneal Plate, Medium Left |
| 70-0025-S | 70-0025 | Locking Calcaneal Plate, Medium Right |
| 70-0026-S | 70-0026 | Locking Calcaneal Plate, Large Left |
| 70-0027-S | 70-0027 | Locking Calcaneal Plate, Large Right |
| 70-0386-S | 70-0386 | Anterior Process Calcaneal Plate, Medium, Left |
| 70-0387-S | 70-0387 | Anterior Process Calcaneal Plate, Medium, Right |
| 70-0388-S | 70-0388 | Anterior Process Calcaneal Plate, Large, Left |
| 70-0389-S | 70-0389 | Anterior Process Calcaneal Plate, Large, Right |
| 70-0394-S | 70-0394 | Posterior Tuberosity Calcaneal Plate, 5-Hole, Left |
| 70-0395-S | 70-0395 | Posterior Tuberosity Calcaneal Plate, 5-Hole, Right |
| 70-0396-S | 70-0396 | Posterior Tuberosity Calcaneal Plate, 6-Hole, Left |
| 70-0397-S | 70-0397 | Posterior Tuberosity Calcaneal Plate, 6-Hole, Right |
| 70-0400-S | 70-0400 | Combo Calcaneal Plate, Medium, 8-Hole, Left |
| 70-0401-S | 70-0401 | Combo Calcaneal Plate, Medium, 8-Hole, Right |
| 70-0404-S | 70-0404 | Combo Calcaneal Plate, Large, 8-Hole, Left |
| 70-0405-S | 70-0405 | Combo Calcaneal Plate, Large, 8-Hole, Right |
| 70-0406-S | 70-0406 | Combo Calcaneal Plate, Large, 9-Hole, Left |
| 70-0407-S | 70-0407 | Combo Calcaneal Plate, Large, 9-Hole, Right |

Foot and Ankle Portfolio [continued]

Calcaneal Plating System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|-----------------------------------|
| CA-4120-S | CA-4120 | 4.0 mm x 12.0 mm Cancellous Screw |
| CA-4140-S | CA-4140 | 4.0 mm x 14.0 mm Cancellous Screw |
| CA-4160-S | CA-4160 | 4.0 mm x 16.0 mm Cancellous Screw |
| CA-4180-S | CA-4180 | 4.0 mm x 18.0 mm Cancellous Screw |
| CA-4200-S | CA-4200 | 4.0 mm x 20.0 mm Cancellous Screw |
| CA-4220-S | CA-4220 | 4.0 mm x 22.0 mm Cancellous Screw |
| CA-4240-S | CA-4240 | 4.0 mm x 24.0 mm Cancellous Screw |
| CA-4260-S | CA-4260 | 4.0 mm x 26.0 mm Cancellous Screw |
| CA-4280-S | CA-4280 | 4.0 mm x 28.0 mm Cancellous Screw |
| CA-4300-S | CA-4300 | 4.0 mm x 30.0 mm Cancellous Screw |
| CA-4350-S | CA-4350 | 4.0 mm x 35.0 mm Cancellous Screw |
| CA-4400-S | CA-4400 | 4.0 mm x 40.0 mm Cancellous Screw |
| CA-4450-S | CA-4450 | 4.0 mm x 45.0 mm Cancellous Screw |
| CA-4500-S | CA-4500 | 4.0 mm x 50.0 mm Cancellous Screw |
| CA-4550-S | CA-4550 | 4.0 mm x 55.0 mm Cancellous Screw |
| CA-4600-S | CA-4600 | 4.0 mm x 60.0 mm Cancellous Screw |
| CO-2708-S | CO-2708 | 2.7 mm x 8.0 mm Cortical Screw |
| CO-2710-S | CO-2710 | 2.7 mm x 10.0 mm Cortical Screw |
| CO-2712-S | CO-2712 | 2.7 mm x 12.0 mm Cortical Screw |
| CO-2714-S | CO-2714 | 2.7 mm x 14.0 mm Cortical Screw |
| CO-2716-S | CO-2716 | 2.7 mm x 16.0 mm Cortical Screw |
| CO-2718-S | CO-2718 | 2.7 mm x 18.0 mm Cortical Screw |
| CO-2720-S | CO-2720 | 2.7 mm x 20.0 mm Cortical Screw |
| CO-2722-S | CO-2722 | 2.7 mm x 22.0 mm Cortical Screw |
| CO-2724-S | CO-2724 | 2.7 mm x 24.0 mm Cortical Screw |
| CO-2726-S | CO-2726 | 2.7 mm x 26.0 mm Cortical Screw |
| CO-2728-S | CO-2728 | 2.7 mm x 28.0 mm Cortical Screw |
| CO-2730-S | CO-2730 | 2.7 mm x 30.0 mm Cortical Screw |
| CO-2732-S | CO-2732 | 2.7 mm x 32.0 mm Cortical Screw |

Foot and Ankle Portfolio [continued]

Calcaneal Plating System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|--|
| CO-2734-S | CO-2734 | 2.7 mm x 34.0 mm Cortical Screw |
| CO-2736-S | CO-2736 | 2.7 mm x 36.0 mm Cortical Screw |
| CO-2738-S | CO-2738 | 2.7 mm x 38.0 mm Cortical Screw |
| CO-2740-S | CO-2740 | 2.7 mm x 40.0 mm Cortical Screw |
| CO-2745-S | CO-2745 | 2.7 mm x 45.0 mm Cortical Screw |
| CO-2750-S | CO-2750 | 2.7 mm x 50.0 mm Cortical Screw |
| CO-2755-S | CO-2755 | 2.7 mm x 55.0 mm Cortical Screw |
| CO-3080-S | CO-3080 | 3.5 mm x 8.0 mm Cortical Screw |
| CO-3100-S | CO-3100 | 3.5 mm x 10.0 mm Cortical Screw |
| CO-3120-S | CO-3120 | 3.5 mm x 12.0 mm Cortical Screw |
| CO-3140-S | CO-3140 | 3.5 mm x 14.0 mm Cortical Screw |
| CO-3160-S | CO-3160 | 3.5 mm x 16.0 mm Cortical Screw |
| CO-3180-S | CO-3180 | 3.5 mm x 18.0 mm Cortical Screw |
| CO-3200-S | CO-3200 | 3.5 mm x 20.0 mm Cortical Screw |
| CO-3220-S | CO-3220 | 3.5 mm x 22.0 mm Cortical Screw |
| CO-3240-S | CO-3240 | 3.5 mm x 24.0 mm Cortical Screw |
| CO-3260-S | CO-3260 | 3.5 mm x 26.0 mm Cortical Screw |
| CO-3280-S | CO-3280 | 3.5 mm x 28.0 mm Cortical Screw |
| CO-3300-S | CO-3300 | 3.5 mm x 30.0 mm Cortical Screw |
| CO-3320-S | CO-3320 | 3.5 mm x 32.0 mm Cortical Screw |
| CO-3340-S | CO-3340 | 3.5 mm x 34.0 mm Cortical Screw |
| CO-3360-S | CO-3360 | 3.5 mm x 36.0 mm Cortical Screw |
| CO-3380-S | CO-3380 | 3.5 mm x 38.0 mm Cortical Screw |
| CO-3400-S | CO-3400 | 3.5 mm x 40.0 mm Cortical Screw |
| CO-3450-S | CO-3450 | 3.5 mm x 45.0 mm Cortical Screw |
| CO-3500-S | CO-3500 | 3.5 mm x 50.0 mm Cortical Screw |
| CO-3550-S | CO-3550 | 3.5 mm x 55.0 mm Cortical Screw |
| CO-3600-S | CO-3600 | 3.5 mm x 60.0 mm Cortical Screw |
| COL-2080-S | COL-2080 | 2.7 mm x 8.0 mm Locking Cortical Screw |

Foot and Ankle Portfolio [continued]

Calcaneal Plating System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|---|
| COL-2100-S | COL-2100 | 2.7 mm x 10.0 mm Locking Cortical Screw |
| COL-2120-S | COL-2120 | 2.7 mm x 12.0 mm Locking Cortical Screw |
| COL-2140-S | COL-2140 | 2.7 mm x 14.0 mm Locking Cortical Screw |
| COL-2160-S | COL-2160 | 2.7 mm x 16.0 mm Locking Cortical Screw |
| COL-2180-S | COL-2180 | 2.7 mm x 18.0 mm Locking Cortical Screw |
| COL-2200-S | COL-2200 | 2.7 mm x 20.0 mm Locking Cortical Screw |
| COL-2220-S | COL-2220 | 2.7 mm x 22.0 mm Locking Cortical Screw |
| COL-2240-S | COL-2240 | 2.7 mm x 24.0 mm Locking Cortical Screw |
| COL-2260-S | COL-2260 | 2.7 mm x 26.0 mm Locking Cortical Screw |
| COL-2280-S | COL-2280 | 2.7 mm x 28.0 mm Locking Cortical Screw |
| COL-2300-S | COL-2300 | 2.7 mm x 30.0 mm Locking Cortical Screw |
| COL-2320-S | COL-2320 | 2.7 mm x 32.0 mm Locking Cortical Screw |
| COL-2340-S | COL-2340 | 2.7 mm x 34.0 mm Locking Cortical Screw |
| COL-2360-S | COL-2360 | 2.7 mm x 36.0 mm Locking Cortical Screw |
| COL-2380-S | COL-2380 | 2.7 mm x 38.0 mm Locking Cortical Screw |
| COL-2400-S | COL-2400 | 2.7 mm x 40.0 mm Locking Cortical Screw |
| COL-2450-S | COL-2450 | 2.7 mm x 45.0 mm Locking Cortical Screw |
| COL-2500-S | COL-2500 | 2.7 mm x 50.0 mm Locking Cortical Screw |
| COL-2550-S | COL-2550 | 2.7 mm x 55.0 mm Locking Cortical Screw |
| COL-3080-S | COL-3080 | 3.5 mm x 8.0 mm Locking Cortical Screw |
| COL-3100-S | COL-3100 | 3.5 mm x 10.0 mm Locking Cortical Screw |
| COL-3120-S | COL-3120 | 3.5 mm x 12.0 mm Locking Cortical Screw |
| COL-3140-S | COL-3140 | 3.5 mm x 14.0 mm Locking Cortical Screw |
| COL-3160-S | COL-3160 | 3.5 mm x 16.0 mm Locking Cortical Screw |
| COL-3180-S | COL-3180 | 3.5 mm x 18.0 mm Locking Cortical Screw |
| COL-3200-S | COL-3200 | 3.5 mm x 20.0 mm Locking Cortical Screw |
| COL-3220-S | COL-3220 | 3.5 mm x 22.0 mm Locking Cortical Screw |
| COL-3240-S | COL-3240 | 3.5 mm x 24.0 mm Locking Cortical Screw |
| COL-3260-S | COL-3260 | 3.5 mm x 26.0 mm Locking Cortical Screw |

Foot and Ankle Portfolio [continued]

Calcaneal Plating System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|---|
| COL-3280-S | COL-3280 | 3.5 mm x 28.0 mm Locking Cortical Screw |
| COL-3300-S | COL-3300 | 3.5 mm x 30.0 mm Locking Cortical Screw |
| COL-3320-S | COL-3320 | 3.5 mm x 32.0 mm Locking Cortical Screw |
| COL-3340-S | COL-3340 | 3.5 mm x 34.0 mm Locking Cortical Screw |
| COL-3360-S | COL-3360 | 3.5 mm x 36.0 mm Locking Cortical Screw |
| COL-3380-S | COL-3380 | 3.5 mm x 38.0 mm Locking Cortical Screw |
| COL-3400-S | COL-3400 | 3.5 mm x 40.0 mm Locking Cortical Screw |
| COL-3450-S | COL-3450 | 3.5 mm x 45.0 mm Locking Cortical Screw |
| COL-3500-S | COL-3500 | 3.5 mm x 50.0 mm Locking Cortical Screw |
| COL-3550-S | COL-3550 | 3.5 mm x 55.0 mm Locking Cortical Screw |
| COL-3600-S | COL-3600 | 3.5 mm x 60.0 mm Locking Cortical Screw |

Forefoot-Midfoot Plating System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|---------------------------------------|
| N/A | 30-0232 | 3.5 mm x 8 mm Locking Hexalobe Screw |
| N/A | 30-0233 | 3.5 mm x 10 mm Locking Hexalobe Screw |
| N/A | 30-0234 | 3.5 mm x 12 mm Locking Hexalobe Screw |
| N/A | 30-0235 | 3.5 mm x 14 mm Locking Hexalobe Screw |
| N/A | 30-0236 | 3.5 mm x 16 mm Locking Hexalobe Screw |
| N/A | 30-0237 | 3.5 mm x 18 mm Locking Hexalobe Screw |
| N/A | 30-0238 | 3.5 mm x 20 mm Locking Hexalobe Screw |
| N/A | 30-0239 | 3.5 mm x 22 mm Locking Hexalobe Screw |
| N/A | 30-0240 | 3.5 mm x 24 mm Locking Hexalobe Screw |
| N/A | 30-0241 | 3.5 mm x 26 mm Locking Hexalobe Screw |
| N/A | 30-0242 | 3.5 mm x 28 mm Locking Hexalobe Screw |
| N/A | 30-0243 | 3.5 mm x 30 mm Locking Hexalobe Screw |
| N/A | 30-0244 | 3.5 mm x 32 mm Locking Hexalobe Screw |
| N/A | 30-0245 | 3.5 mm x 34 mm Locking Hexalobe Screw |

Foot and Ankle Portfolio [continued]

Forefoot-Midfoot Plating System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|--|
| N/A | 30-0246 | 3.5 mm x 36 mm Locking Hexalobe Screw |
| N/A | 30-0247 | 3.5 mm x 38 mm Locking Hexalobe Screw |
| N/A | 30-0248 | 3.5 mm x 40 mm Locking Hexalobe Screw |
| N/A | 30-0249 | 3.5 mm x 45 mm Locking Hexalobe Screw |
| N/A | 30-0250 | 3.5 mm x 50 mm Locking Hexalobe Screw |
| N/A | 30-0251 | 3.5 mm x 55 mm Locking Hexalobe Screw |
| N/A | 30-0252 | 3.5 mm x 60 mm Locking Hexalobe Screw |
| 30-0255-S | 30-0255 | 3.5 mm x 8 mm Nonlocking Hexalobe Screw |
| 30-0256-S | 30-0256 | 3.5 mm x 10 mm Nonlocking Hexalobe Screw |
| 30-0257-S | 30-0257 | 3.5 mm x 12 mm Nonlocking Hexalobe Screw |
| 30-0258-S | 30-0258 | 3.5 mm x 14 mm Nonlocking Hexalobe Screw |
| 30-0259-S | 30-0259 | 3.5 mm x 16 mm Nonlocking Hexalobe Screw |
| 30-0260-S | 30-0260 | 3.5 mm x 18 mm Nonlocking Hexalobe Screw |
| 30-0261-S | 30-0261 | 3.5 mm x 20 mm Nonlocking Hexalobe Screw |
| 30-0262-S | 30-0262 | 3.5 mm x 22 mm Nonlocking Hexalobe Screw |
| 30-0263-S | 30-0263 | 3.5 mm x 24 mm Nonlocking Hexalobe Screw |
| 30-0264-S | 30-0264 | 3.5 mm x 26 mm Nonlocking Hexalobe Screw |
| 30-0265-S | 30-0265 | 3.5 mm x 28 mm Nonlocking Hexalobe Screw |
| 30-0266-S | 30-0266 | 3.5 mm x 30 mm Nonlocking Hexalobe Screw |
| 30-0267-S | 30-0267 | 3.5 mm x 32 mm Nonlocking Hexalobe Screw |
| 30-0268-S | 30-0268 | 3.5 mm x 34 mm Nonlocking Hexalobe Screw |
| 30-0269-S | 30-0269 | 3.5 mm x 36 mm Nonlocking Hexalobe Screw |
| 30-0270-S | 30-0270 | 3.5 mm x 38 mm Nonlocking Hexalobe Screw |
| 30-0271-S | 30-0271 | 3.5 mm x 40 mm Nonlocking Hexalobe Screw |
| 30-0272-S | 30-0272 | 3.5 mm x 45 mm Nonlocking Hexalobe Screw |
| 30-0273-S | 30-0273 | 3.5 mm x 50 mm Nonlocking Hexalobe Screw |
| 30-0274-S | 30-0274 | 3.5 mm x 55 mm Nonlocking Hexalobe Screw |
| 30-0275-S | 30-0275 | 3.5 mm x 60 mm Nonlocking Hexalobe Screw |
| 30-0278-S | 30-0278 | 3.0 mm x 8 mm Locking Hexalobe Screw |

Foot and Ankle Portfolio [continued]

Forefoot-Midfoot Plating System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|--|
| 30-0279-S | 30-0279 | 3.0 mm x 10 mm Locking Hexalobe Screw |
| 30-0280-S | 30-0280 | 3.0 mm x 12 mm Locking Hexalobe Screw |
| 30-0281-S | 30-0281 | 3.0 mm x 14 mm Locking Hexalobe Screw |
| 30-0282-S | 30-0282 | 3.0 mm x 16 mm Locking Hexalobe Screw |
| 30-0283-S | 30-0283 | 3.0 mm x 18 mm Locking Hexalobe Screw |
| 30-0284-S | 30-0284 | 3.0 mm x 20 mm Locking Hexalobe Screw |
| 30-0285-S | 30-0285 | 3.0 mm x 22 mm Locking Hexalobe Screw |
| 30-0286-S | 30-0286 | 3.0 mm x 24 mm Locking Hexalobe Screw |
| 30-0287-S | 30-0287 | 3.0 mm x 26 mm Locking Hexalobe Screw |
| 30-0288-S | 30-0288 | 3.0 mm x 28 mm Locking Hexalobe Screw |
| 30-0289-S | 30-0289 | 3.0 mm x 30 mm Locking Hexalobe Screw |
| 30-0290-S | 30-0290 | 3.0 mm x 32 mm Locking Hexalobe Screw |
| 30-0291-S | 30-0291 | 3.0 mm x 34 mm Locking Hexalobe Screw |
| 30-0292-S | 30-0292 | 3.0 mm x 36 mm Locking Hexalobe Screw |
| 30-0293-S | 30-0293 | 3.0 mm x 38 mm Locking Hexalobe Screw |
| 30-0294-S | 30-0294 | 3.0 mm x 40 mm Locking Hexalobe Screw |
| 30-0295-S | 30-0295 | 3.0 mm x 45 mm Locking Hexalobe Screw |
| 30-0296-S | 30-0296 | 3.0 mm x 50 mm Locking Hexalobe Screw |
| 30-0297-S | 30-0297 | 3.0 mm x 55 mm Locking Hexalobe Screw |
| 30-0301-S | 30-0301 | 3.0 mm x 8 mm Nonlocking Hexalobe Screw |
| 30-0302-S | 30-0302 | 3.0 mm x 10 mm Nonlocking Hexalobe Screw |
| 30-0303-S | 30-0303 | 3.0 mm x 12 mm Nonlocking Hexalobe Screw |
| 30-0304-S | 30-0304 | 3.0 mm x 14 mm Nonlocking Hexalobe Screw |
| 30-0305-S | 30-0305 | 3.0 mm x 16 mm Nonlocking Hexalobe Screw |
| 30-0306-S | 30-0306 | 3.0 mm x 18 mm Nonlocking Hexalobe Screw |
| 30-0307-S | 30-0307 | 3.0 mm x 20 mm Nonlocking Hexalobe Screw |
| 30-0308-S | 30-0308 | 3.0 mm x 22 mm Nonlocking Hexalobe Screw |
| 30-0309-S | 30-0309 | 3.0 mm x 24 mm Nonlocking Hexalobe Screw |
| 30-0310-S | 30-0310 | 3.0 mm x 26 mm Nonlocking Hexalobe Screw |

Foot and Ankle Portfolio [continued]

Forefoot-Midfoot Plating System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|--|
| 30-0311-S | 30-0311 | 3.0 mm x 28 mm Nonlocking Hexalobe Screw |
| 30-0312-S | 30-0312 | 3.0 mm x 30 mm Nonlocking Hexalobe Screw |
| 30-0313-S | 30-0313 | 3.0 mm x 32 mm Nonlocking Hexalobe Screw |
| 30-0314-S | 30-0314 | 3.0 mm x 34 mm Nonlocking Hexalobe Screw |
| 30-0315-S | 30-0315 | 3.0 mm x 36 mm Nonlocking Hexalobe Screw |
| 30-0316-S | 30-0316 | 3.0 mm x 38 mm Nonlocking Hexalobe Screw |
| 30-0317-S | 30-0317 | 3.0 mm x 40 mm Nonlocking Hexalobe Screw |
| 30-0318-S | 30-0318 | 3.0 mm x 45 mm Nonlocking Hexalobe Screw |
| 30-0319-S | 30-0319 | 3.0 mm x 50 mm Nonlocking Hexalobe Screw |
| 30-0320-S | 30-0320 | 3.0 mm x 55 mm Nonlocking Hexalobe Screw |
| 70-0007-S | 70-0007 | 5-Hole Locking 1st Tarsometatarsal Plate |
| 70-0008-S | 70-0008 | 4-Hole Locking 1st Tarsometatarsal Plate |
| 70-0009-S | 70-0009 | 4-Hole Locking 2nd & 3rd Tarsometatarsal Plate |
| 70-0010-S | 70-0010 | Locking Proximal Metatarsophalangeal Wedge Plate, Left |
| 70-0011-S | 70-0011 | Locking Proximal Metatarsophalangeal Wedge Plate, Right |
| 70-0012-S | 70-0012 | Locking Dorsal Metatarsophalangeal Fusion Plate, Left |
| 70-0013-S | 70-0013 | Locking Dorsal Metatarsophalangeal Fusion Plate, Right |
| 70-0018-S | 70-0018 | Locking Metatarsophalangeal Revision Plate, Left |
| 70-0019-S | 70-0019 | Locking Metatarsophalangeal Revision Plate, Right |
| 70-0036-S | 70-0036 | Locking Metatarsophalangeal/Metatarsophalangeal Joint Combo Plate, Left |
| 70-0037-S | 70-0037 | Locking Metatarsophalangeal/Metatarsophalangeal Joint Combo Plate, Right |
| 70-0049-S | 70-0049 | 8-Hole Locking 1st Tarsometatarsal Plate |
| 70-0050-S | 70-0050 | 7-Hole Locking 1st Tarsometatarsal Plate |
| 70-0051-S | 70-0051 | 7-Hole Locking 2nd & 3rd Tarsometatarsal Plate |
| 70-0109-S | 70-0109 | Locking Metatarsophalangeal/Metatarsophalangeal Joint Combo Plate 4 Degree Left |
| 70-0110-S | 70-0110 | Locking Metatarsophalangeal/Metatarsophalangeal Joint Combo Plate 4 Degree Right |
| 70-0141-S | 70-0141 | Locking Dorsal Metatarsophalangeal Plate 4 Degree Left |
| 70-0142-S | 70-0142 | Locking Dorsal Metatarsophalangeal Plate 4 Degree Right |

Foot and Ankle Portfolio [continued]

Forefoot-Midfoot Plating System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|--|
| 70-0143-S | 70-0143 | Locking Metatarsophalangeal Revision Plate 4 Degree Left |
| 70-0144-S | 70-0144 | Locking Metatarsophalangeal Revision Plate 4 Degree Right |
| 70-0323-S | 70-0323 | Locking Proximal Metatarsal Wedge Plate, 0 Degree, Left |
| 70-0324-S | 70-0324 | Locking Proximal Metatarsal Wedge Plate, 0 Degree, Right |
| 70-0329-S | 70-0329 | Locking Metatarsophalangeal Joint Petite Plate 4 Degree, Left |
| 70-0330-S | 70-0330 | Locking Metatarsophalangeal Joint Petite Plate 4 Degree, Right |
| CA-4120-S | CA-4120 | 4.0 mm x 12.0 mm Cancellous Screw |
| CA-4140-S | CA-4140 | 4.0 mm x 14.0 mm Cancellous Screw |
| CA-4160-S | CA-4160 | 4.0 mm x 16.0 mm Cancellous Screw |
| CA-4180-S | CA-4180 | 4.0 mm x 18.0 mm Cancellous Screw |
| CA-4200-S | CA-4200 | 4.0 mm x 20.0 mm Cancellous Screw |
| CA-4220-S | CA-4220 | 4.0 mm x 22.0 mm Cancellous Screw |
| CA-4240-S | CA-4240 | 4.0 mm x 24.0 mm Cancellous Screw |
| CA-4260-S | CA-4260 | 4.0 mm x 26.0 mm Cancellous Screw |
| CA-4280-S | CA-4280 | 4.0 mm x 28.0 mm Cancellous Screw |
| CA-4300-S | CA-4300 | 4.0 mm x 30.0 mm Cancellous Screw |
| CA-4350-S | CA-4350 | 4.0 mm x 35.0 mm Cancellous Screw |
| CA-4400-S | CA-4400 | 4.0 mm x 40.0 mm Cancellous Screw |
| CA-4450-S | CA-4450 | 4.0 mm x 45.0 mm Cancellous Screw |
| CA-4500-S | CA-4500 | 4.0 mm x 50.0 mm Cancellous Screw |
| CA-4550-S | CA-4550 | 4.0 mm x 55.0 mm Cancellous Screw |
| CA-4600-S | CA-4600 | 4.0 mm x 60.0 mm Cancellous Screw |
| CO-2708-S | CO-2708 | 2.7 mm x 8.0 mm Cortical Screw |
| CO-2710-S | CO-2710 | 2.7 mm x 10.0 mm Cortical Screw |
| CO-2712-S | CO-2712 | 2.7 mm x 12.0 mm Cortical Screw |
| CO-2714-S | CO-2714 | 2.7 mm x 14.0 mm Cortical Screw |
| CO-2716-S | CO-2716 | 2.7 mm x 16.0 mm Cortical Screw |
| CO-2718-S | CO-2718 | 2.7 mm x 18.0 mm Cortical Screw |
| CO-2720-S | CO-2720 | 2.7 mm x 20.0 mm Cortical Screw |

Foot and Ankle Portfolio [continued]

Forefoot-Midfoot Plating System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|---------------------------------|
| CO-2722-S | CO-2722 | 2.7 mm x 22.0 mm Cortical Screw |
| CO-2724-S | CO-2724 | 2.7 mm x 24.0 mm Cortical Screw |
| CO-2726-S | CO-2726 | 2.7 mm x 26.0 mm Cortical Screw |
| CO-2728-S | CO-2728 | 2.7 mm x 28.0 mm Cortical Screw |
| CO-2730-S | CO-2730 | 2.7 mm x 30.0 mm Cortical Screw |
| CO-2732-S | CO-2732 | 2.7 mm x 32.0 mm Cortical Screw |
| CO-2734-S | CO-2734 | 2.7 mm x 34.0 mm Cortical Screw |
| CO-2736-S | CO-2736 | 2.7 mm x 36.0 mm Cortical Screw |
| CO-2738-S | CO-2738 | 2.7 mm x 38.0 mm Cortical Screw |
| CO-2740-S | CO-2740 | 2.7 mm x 40.0 mm Cortical Screw |
| CO-2745-S | CO-2745 | 2.7 mm x 45.0 mm Cortical Screw |
| CO-2750-S | CO-2750 | 2.7 mm x 50.0 mm Cortical Screw |
| CO-2755-S | CO-2755 | 2.7 mm x 55.0 mm Cortical Screw |
| CO-3080-S | CO-3080 | 3.5 mm x 8.0 mm Cortical Screw |
| CO-3100-S | CO-3100 | 3.5 mm x 10.0 mm Cortical Screw |
| CO-3120-S | CO-3120 | 3.5 mm x 12.0 mm Cortical Screw |
| CO-3140-S | CO-3140 | 3.5 mm x 14.0 mm Cortical Screw |
| CO-3160-S | CO-3160 | 3.5 mm x 16.0 mm Cortical Screw |
| CO-3180-S | CO-3180 | 3.5 mm x 18.0 mm Cortical Screw |
| CO-3200-S | CO-3200 | 3.5 mm x 20.0 mm Cortical Screw |
| CO-3220-S | CO-3220 | 3.5 mm x 22.0 mm Cortical Screw |
| CO-3240-S | CO-3240 | 3.5 mm x 24.0 mm Cortical Screw |
| CO-3260-S | CO-3260 | 3.5 mm x 26.0 mm Cortical Screw |
| CO-3280-S | CO-3280 | 3.5 mm x 28.0 mm Cortical Screw |
| CO-3300-S | CO-3300 | 3.5 mm x 30.0 mm Cortical Screw |
| CO-3320-S | CO-3320 | 3.5 mm x 32.0 mm Cortical Screw |
| CO-3340-S | CO-3340 | 3.5 mm x 34.0 mm Cortical Screw |
| CO-3360-S | CO-3360 | 3.5 mm x 36.0 mm Cortical Screw |
| CO-3380-S | CO-3380 | 3.5 mm x 38.0 mm Cortical Screw |

Foot and Ankle Portfolio [continued]

Forefoot-Midfoot Plating System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|---|
| CO-3400-S | CO-3400 | 3.5 mm x 40.0 mm Cortical Screw |
| CO-3450-S | CO-3450 | 3.5 mm x 45.0 mm Cortical Screw |
| CO-3500-S | CO-3500 | 3.5 mm x 50.0 mm Cortical Screw |
| CO-3550-S | CO-3550 | 3.5 mm x 55.0 mm Cortical Screw |
| CO-3600-S | CO-3600 | 3.5 mm x 60.0 mm Cortical Screw |
| COL-2080-S | COL-2080 | 2.7 mm x 8.0 mm Locking Cortical Screw |
| COL-2100-S | COL-2100 | 2.7 mm x 10.0 mm Locking Cortical Screw |
| COL-2120-S | COL-2120 | 2.7 mm x 12.0 mm Locking Cortical Screw |
| COL-2140-S | COL-2140 | 2.7 mm x 14.0 mm Locking Cortical Screw |
| COL-2160-S | COL-2160 | 2.7 mm x 16.0 mm Locking Cortical Screw |
| COL-2180-S | COL-2180 | 2.7 mm x 18.0 mm Locking Cortical Screw |
| COL-2200-S | COL-2200 | 2.7 mm x 20.0 mm Locking Cortical Screw |
| COL-2220-S | COL-2220 | 2.7 mm x 22.0 mm Locking Cortical Screw |
| COL-2240-S | COL-2240 | 2.7 mm x 24.0 mm Locking Cortical Screw |
| COL-2260-S | COL-2260 | 2.7 mm x 26.0 mm Locking Cortical Screw |
| COL-2280-S | COL-2280 | 2.7 mm x 28.0 mm Locking Cortical Screw |
| COL-2300-S | COL-2300 | 2.7 mm x 30.0 mm Locking Cortical Screw |
| COL-2320-S | COL-2320 | 2.7 mm x 32.0 mm Locking Cortical Screw |
| COL-2340-S | COL-2340 | 2.7 mm x 34.0 mm Locking Cortical Screw |
| COL-2360-S | COL-2360 | 2.7 mm x 36.0 mm Locking Cortical Screw |
| COL-2380-S | COL-2380 | 2.7 mm x 38.0 mm Locking Cortical Screw |
| COL-2400-S | COL-2400 | 2.7 mm x 40.0 mm Locking Cortical Screw |
| COL-2450-S | COL-2450 | 2.7 mm x 45.0 mm Locking Cortical Screw |
| COL-2500-S | COL-2500 | 2.7 mm x 50.0 mm Locking Cortical Screw |
| COL-2550-S | COL-2550 | 2.7 mm x 55.0 mm Locking Cortical Screw |
| COL-3080-S | COL-3080 | 3.5 mm x 8.0 mm Locking Cortical Screw |
| COL-3100-S | COL-3100 | 3.5 mm x 10.0 mm Locking Cortical Screw |
| COL-3120-S | COL-3120 | 3.5 mm x 12.0 mm Locking Cortical Screw |
| COL-3140-S | COL-3140 | 3.5 mm x 14.0 mm Locking Cortical Screw |

Foot and Ankle Portfolio [continued]

Forefoot-Midfoot Plating System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|---|
| COL-3160-S | COL-3160 | 3.5 mm x 16.0 mm Locking Cortical Screw |
| COL-3180-S | COL-3180 | 3.5 mm x 18.0 mm Locking Cortical Screw |
| COL-3200-S | COL-3200 | 3.5 mm x 20.0 mm Locking Cortical Screw |
| COL-3220-S | COL-3220 | 3.5 mm x 22.0 mm Locking Cortical Screw |
| COL-3240-S | COL-3240 | 3.5 mm x 24.0 mm Locking Cortical Screw |
| COL-3260-S | COL-3260 | 3.5 mm x 26.0 mm Locking Cortical Screw |
| COL-3280-S | COL-3280 | 3.5 mm x 28.0 mm Locking Cortical Screw |
| COL-3300-S | COL-3300 | 3.5 mm x 30.0 mm Locking Cortical Screw |
| COL-3320-S | COL-3320 | 3.5 mm x 32.0 mm Locking Cortical Screw |
| COL-3340-S | COL-3340 | 3.5 mm x 34.0 mm Locking Cortical Screw |
| COL-3360-S | COL-3360 | 3.5 mm x 36.0 mm Locking Cortical Screw |
| COL-3380-S | COL-3380 | 3.5 mm x 38.0 mm Locking Cortical Screw |
| COL-3400-S | COL-3400 | 3.5 mm x 40.0 mm Locking Cortical Screw |
| COL-3450-S | COL-3450 | 3.5 mm x 45.0 mm Locking Cortical Screw |
| COL-3500-S | COL-3500 | 3.5 mm x 50.0 mm Locking Cortical Screw |
| COL-3550-S | COL-3550 | 3.5 mm x 55.0 mm Locking Cortical Screw |
| COL-3600-S | COL-3600 | 3.5 mm x 60.0 mm Locking Cortical Screw |

Small Fragment Base Set

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|---------------------------------------|
| N/A | 30-0098 | 70.0 mm Tension Band Pin |
| N/A | 30-0099 | 90.0 mm Tension Band Pin |
| 30-0232-S | 30-0232 | 3.5 mm x 8 mm Locking Hexalobe Screw |
| 30-0233-S | 30-0233 | 3.5 mm x 10 mm Locking Hexalobe Screw |
| 30-0234-S | 30-0234 | 3.5 mm x 12 mm Locking Hexalobe Screw |
| 30-0235-S | 30-0235 | 3.5 mm x 14 mm Locking Hexalobe Screw |
| 30-0236-S | 30-0236 | 3.5 mm x 16 mm Locking Hexalobe Screw |
| 30-0237-S | 30-0237 | 3.5 mm x 18 mm Locking Hexalobe Screw |

Foot and Ankle Portfolio [continued]

Small Fragment Base Set

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|--|
| 30-0238-S | 30-0238 | 3.5 mm x 20 mm Locking Hexalobe Screw |
| 30-0239-S | 30-0239 | 3.5 mm x 22 mm Locking Hexalobe Screw |
| 30-0240-S | 30-0240 | 3.5 mm x 24 mm Locking Hexalobe Screw |
| 30-0241-S | 30-0241 | 3.5 mm x 26 mm Locking Hexalobe Screw |
| 30-0242-S | 30-0242 | 3.5 mm x 28 mm Locking Hexalobe Screw |
| 30-0243-S | 30-0243 | 3.5 mm x 30 mm Locking Hexalobe Screw |
| 30-0244-S | 30-0244 | 3.5 mm x 32 mm Locking Hexalobe Screw |
| 30-0245-S | 30-0245 | 3.5 mm x 34 mm Locking Hexalobe Screw |
| 30-0246-S | 30-0246 | 3.5 mm x 36 mm Locking Hexalobe Screw |
| 30-0247-S | 30-0247 | 3.5 mm x 38 mm Locking Hexalobe Screw |
| 30-0248-S | 30-0248 | 3.5 mm x 40 mm Locking Hexalobe Screw |
| 30-0250-S | 30-0250 | 3.5 mm x 50 mm Locking Hexalobe Screw |
| 30-0251-S | 30-0251 | 3.5 mm x 55 mm Locking Hexalobe Screw |
| 30-0252-S | 30-0252 | 3.5 mm x 60 mm Locking Hexalobe Screw |
| 30-0253-S | 30-0253 | 3.5 mm x 65 mm Locking Hexalobe Screw |
| 30-0255-S | 30-0255 | 3.5 mm x 8 mm Nonlocking Hexalobe Screw |
| 30-0256-S | 30-0256 | 3.5 mm x 10 mm Nonlocking Hexalobe Screw |
| 30-0257-S | 30-0257 | 3.5 mm x 12 mm Nonlocking Hexalobe Screw |
| 30-0258-S | 30-0258 | 3.5 mm x 14 mm Nonlocking Hexalobe Screw |
| 30-0259-S | 30-0259 | 3.5 mm x 16 mm Nonlocking Hexalobe Screw |
| 30-0260-S | 30-0260 | 3.5 mm x 18 mm Nonlocking Hexalobe Screw |
| 30-0261-S | 30-0261 | 3.5 mm x 20 mm Nonlocking Hexalobe Screw |
| 30-0262-S | 30-0262 | 3.5 mm x 22 mm Nonlocking Hexalobe Screw |
| 30-0263-S | 30-0263 | 3.5 mm x 24 mm Nonlocking Hexalobe Screw |
| 30-0264-S | 30-0264 | 3.5 mm x 26 mm Nonlocking Hexalobe Screw |
| 30-0265-S | 30-0265 | 3.5 mm x 28 mm Nonlocking Hexalobe Screw |
| 30-0266-S | 30-0266 | 3.5 mm x 30 mm Nonlocking Hexalobe Screw |
| 30-0267-S | 30-0267 | 3.5 mm x 32 mm Nonlocking Hexalobe Screw |
| 30-0268-S | 30-0268 | 3.5 mm x 34 mm Nonlocking Hexalobe Screw |

Foot and Ankle Portfolio [continued]

Small Fragment Base Set

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|--|
| 30-0269-S | 30-0269 | 3.5 mm x 36 mm Nonlocking Hexalobe Screw |
| 30-0270-S | 30-0270 | 3.5 mm x 38 mm Nonlocking Hexalobe Screw |
| 30-0271-S | 30-0271 | 3.5 mm x 40 mm Nonlocking Hexalobe Screw |
| 30-0273-S | 30-0273 | 3.5 mm x 50 mm Nonlocking Hexalobe Screw |
| 30-0274-S | 30-0274 | 3.5 mm x 55 mm Nonlocking Hexalobe Screw |
| 30-0275-S | 30-0275 | 3.5 mm x 60 mm Nonlocking Hexalobe Screw |
| 30-0276-S | 30-0276 | 3.5 mm x 65 mm Nonlocking Hexalobe Screw |
| 30-0324-S | 30-0324 | 2.7 mm x 8 mm Locking Hexalobe Screw |
| 30-0325-S | 30-0325 | 2.7 mm x 10 mm Locking Hexalobe Screw |
| 30-0326-S | 30-0326 | 2.7 mm x 12 mm Locking Hexalobe Screw |
| 30-0327-S | 30-0327 | 2.7 mm x 14 mm Locking Hexalobe Screw |
| 30-0328-S | 30-0328 | 2.7 mm x 16 mm Locking Hexalobe Screw |
| 30-0329-S | 30-0329 | 2.7 mm x 18 mm Locking Hexalobe Screw |
| 30-0330-S | 30-0330 | 2.7 mm x 20 mm Locking Hexalobe Screw |
| 30-0331-S | 30-0331 | 2.7 mm x 22 mm Locking Hexalobe Screw |
| 30-0332-S | 30-0332 | 2.7 mm x 24 mm Locking Hexalobe Screw |
| 30-0333-S | 30-0333 | 2.7 mm x 26 mm Locking Hexalobe Screw |
| 30-0334-S | 30-0334 | 2.7 mm x 28 mm Locking Hexalobe Screw |
| 30-0335-S | 30-0335 | 2.7 mm x 30 mm Locking Hexalobe Screw |
| 30-0336-S | 30-0336 | 2.7 mm x 32 mm Locking Hexalobe Screw |
| 30-0337-S | 30-0337 | 2.7 mm x 34 mm Locking Hexalobe Screw |
| 30-0338-S | 30-0338 | 2.7 mm x 36 mm Locking Hexalobe Screw |
| 30-0339-S | 30-0339 | 2.7 mm x 38 mm Locking Hexalobe Screw |
| 30-0340-S | 30-0340 | 2.7 mm x 40 mm Locking Hexalobe Screw |
| 30-0343-S | 30-0343 | 2.7 mm x 8 mm Nonlocking Hexalobe Screw |
| 30-0344-S | 30-0344 | 2.7 mm x 10 mm Nonlocking Hexalobe Screw |
| 30-0345-S | 30-0345 | 2.7 mm x 12 mm Nonlocking Hexalobe Screw |
| 30-0346-S | 30-0346 | 2.7 mm x 14 mm Nonlocking Hexalobe Screw |
| 30-0347-S | 30-0347 | 2.7 mm x 16 mm Nonlocking Hexalobe Screw |

Foot and Ankle Portfolio [continued]

Small Fragment Base Set

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|--|
| 30-0348-S | 30-0348 | 2.7 mm x 18 mm Nonlocking Hexalobe Screw |
| 30-0349-S | 30-0349 | 2.7 mm x 20 mm Nonlocking Hexalobe Screw |
| 30-0350-S | 30-0350 | 2.7 mm x 22 mm Nonlocking Hexalobe Screw |
| 30-0351-S | 30-0351 | 2.7 mm x 24 mm Nonlocking Hexalobe Screw |
| 30-0352-S | 30-0352 | 2.7 mm x 26 mm Nonlocking Hexalobe Screw |
| 30-0353-S | 30-0353 | 2.7 mm x 28 mm Nonlocking Hexalobe Screw |
| 30-0354-S | 30-0354 | 2.7 mm x 30 mm Nonlocking Hexalobe Screw |
| 30-0355-S | 30-0355 | 2.7 mm x 32 mm Nonlocking Hexalobe Screw |
| 30-0356-S | 30-0356 | 2.7 mm x 34 mm Nonlocking Hexalobe Screw |
| 30-0357-S | 30-0357 | 2.7 mm x 36 mm Nonlocking Hexalobe Screw |
| 30-0358-S | 30-0358 | 2.7 mm x 38 mm Nonlocking Hexalobe Screw |
| 30-0359-S | 30-0359 | 2.7 mm x 40 mm Nonlocking Hexalobe Screw |
| 30-0361-S | 30-0361 | 2.7 mm x 50 mm Nonlocking Hexalobe Screw |
| 30-0362-S | 30-0362 | 2.7 mm x 55 mm Nonlocking Hexalobe Screw |
| 30-0363-S | 30-0363 | 2.7 mm x 60 mm Nonlocking Hexalobe Screw |
| 3013-27010-S | 3013-27010 | 2.7 mm x 10 mm Hexalobe Variable Angle Screw |
| 3013-27012-S | 3013-27012 | 2.7 mm x 12 mm Hexalobe Variable Angle Screw |
| 3013-27014-S | 3013-27014 | 2.7 mm x 14 mm Hexalobe Variable Angle Screw |
| 3013-27016-S | 3013-27016 | 2.7 mm x 16 mm Hexalobe Variable Angle Screw |
| 3013-27018-S | 3013-27018 | 2.7 mm x 18 mm Hexalobe Variable Angle Screw |
| 3013-27020-S | 3013-27020 | 2.7 mm x 20 mm Hexalobe Variable Angle Screw |
| 3013-27022-S | 3013-27022 | 2.7 mm x 22 mm Hexalobe Variable Angle Screw |
| 3013-27024-S | 3013-27024 | 2.7 mm x 24 mm Hexalobe Variable Angle Screw |
| 3013-27026-S | 3013-27026 | 2.7 mm x 26 mm Hexalobe Variable Angle Screw |
| 3013-27028-S | 3013-27028 | 2.7 mm x 28 mm Hexalobe Variable Angle Screw |
| 3013-27030-S | 3013-27030 | 2.7 mm x 30 mm Hexalobe Variable Angle Screw |
| 3013-27032-S | 3013-27032 | 2.7 mm x 32 mm Hexalobe Variable Angle Screw |
| 3013-27034-S | 3013-27034 | 2.7 mm x 34 mm Hexalobe Variable Angle Screw |
| 3013-27036-S | 3013-27036 | 2.7 mm x 36 mm Hexalobe Variable Angle Screw |

Foot and Ankle Portfolio [continued]

Small Fragment Base Set

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|--|
| 3013-27038-S | 3013-27038 | 2.7 mm x 38 mm Hexalobe Variable Angle Screw |
| 3013-27040-S | 3013-27040 | 2.7 mm x 40 mm Hexalobe Variable Angle Screw |
| 3013-27042-S | 3013-27042 | 2.7 mm x 42 mm Hexalobe Variable Angle Screw |
| 3013-27044-S | 3013-27044 | 2.7 mm x 44 mm Hexalobe Variable Angle Screw |
| 3013-27046-S | 3013-27046 | 2.7 mm x 46 mm Hexalobe Variable Angle Screw |
| 3013-27048-S | 3013-27048 | 2.7 mm x 48 mm Hexalobe Variable Angle Screw |
| 3013-27050-S | 3013-27050 | 2.7 mm x 50 mm Hexalobe Variable Angle Screw |
| 3013-27055-S | 3013-27055 | 2.7 mm x 55 mm Hexalobe Variable Angle Screw |
| 3013-27060-S | 3013-27060 | 2.7 mm x 60 mm Hexalobe Variable Angle Screw |
| 3013-35010-S | 3013-35010 | 3.5 mm x 10 mm Hexalobe Variable Angle Screw |
| 3013-35012-S | 3013-35012 | 3.5 mm x 12 mm Hexalobe Variable Angle Screw |
| 3013-35014-S | 3013-35014 | 3.5 mm x 14 mm Hexalobe Variable Angle Screw |
| 3013-35016-S | 3013-35016 | 3.5 mm x 16 mm Hexalobe Variable Angle Screw |
| 3013-35018-S | 3013-35018 | 3.5 mm x 18 mm Hexalobe Variable Angle Screw |
| 3013-35020-S | 3013-35020 | 3.5 mm x 20 mm Hexalobe Variable Angle Screw |
| 3013-35022-S | 3013-35022 | 3.5 mm x 22 mm Hexalobe Variable Angle Screw |
| 3013-35024-S | 3013-35024 | 3.5 mm x 24 mm Hexalobe Variable Angle Screw |
| 3013-35026-S | 3013-35026 | 3.5 mm x 26 mm Hexalobe Variable Angle Screw |
| 3013-35028-S | 3013-35028 | 3.5 mm x 28 mm Hexalobe Variable Angle Screw |
| 3013-35030-S | 3013-35030 | 3.5 mm x 30 mm Hexalobe Variable Angle Screw |
| 3013-35032-S | 3013-35032 | 3.5 mm x 32 mm Hexalobe Variable Angle Screw |
| 3013-35034-S | 3013-35034 | 3.5 mm x 34 mm Hexalobe Variable Angle Screw |
| 3013-35036-S | 3013-35036 | 3.5 mm x 36 mm Hexalobe Variable Angle Screw |
| 3013-35038-S | 3013-35038 | 3.5 mm x 38 mm Hexalobe Variable Angle Screw |
| 3013-35040-S | 3013-35040 | 3.5 mm x 40 mm Hexalobe Variable Angle Screw |
| 3013-35042-S | 3013-35042 | 3.5 mm x 42 mm Hexalobe Variable Angle Screw |
| 3013-35044-S | 3013-35044 | 3.5 mm x 44 mm Hexalobe Variable Angle Screw |
| 3013-35046-S | 3013-35046 | 3.5 mm x 46 mm Hexalobe Variable Angle Screw |
| 3013-35048-S | 3013-35048 | 3.5 mm x 48 mm Hexalobe Variable Angle Screw |

Foot and Ankle Portfolio [continued]

Small Fragment Base Set

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|---|
| 3013-35050-S | 3013-35050 | 3.5 mm x 50 mm Hexalobe Variable Angle Screw |
| 3013-35055-S | 3013-35055 | 3.5 mm x 55 mm Hexalobe Variable Angle Screw |
| 3013-35060-S | 3013-35060 | 3.5 mm x 60 mm Hexalobe Variable Angle Screw |
| 3013-35065-S | 3013-35065 | 3.5 mm x 65 mm Hexalobe Variable Angle Screw |
| 3015-40010-S | 3015-40010 | 4.0 mm x 10 mm Cancellous Hexalobe Screw |
| 3015-40012-S | 3015-40012 | 4.0 mm x 12 mm Cancellous Hexalobe Screw |
| 3015-40014-S | 3015-40014 | 4.0 mm x 14 mm Cancellous Hexalobe Screw |
| 3015-40016-S | 3015-40016 | 4.0 mm x 16 mm Cancellous Hexalobe Screw |
| 3015-40018-S | 3015-40018 | 4.0 mm x 18 mm Cancellous Hexalobe Screw |
| 3015-40020-S | 3015-40020 | 4.0 mm x 20 mm Cancellous Hexalobe Screw |
| 3015-40022-S | 3015-40022 | 4.0 mm x 22 mm Cancellous Hexalobe Screw |
| 3015-40024-S | 3015-40024 | 4.0 mm x 24 mm Cancellous Hexalobe Screw |
| 3015-40026-S | 3015-40026 | 4.0 mm x 26 mm Cancellous Hexalobe Screw |
| 3015-40028-S | 3015-40028 | 4.0 mm x 28 mm Cancellous Hexalobe Screw |
| 3015-40030-S | 3015-40030 | 4.0 mm x 30 mm Cancellous Hexalobe Screw |
| 3015-40035-S | 3015-40035 | 4.0 mm x 35 mm Cancellous Hexalobe Screw |
| 3015-40040-S | 3015-40040 | 4.0 mm x 40 mm Cancellous Hexalobe Screw |
| 3015-40045-S | 3015-40045 | 4.0 mm x 45 mm Cancellous Hexalobe Screw |
| 3015-40050-S | 3015-40050 | 4.0 mm x 50 mm Cancellous Hexalobe Screw |
| 3015-40055-S | 3015-40055 | 4.0 mm x 55 mm Cancellous Hexalobe Screw |
| 3015-40060-S | 3015-40060 | 4.0 mm x 60 mm Cancellous Hexalobe Screw |
| 3016-40012-S | 3016-40012 | 4.0 mm x 12 mm Partially Threaded Cancellous Hexalobe Screw |
| 3016-40014-S | 3016-40014 | 4.0 mm x 14 mm Partially Threaded Cancellous Hexalobe Screw |
| 3016-40016-S | 3016-40016 | 4.0 mm x 16 mm Partially Threaded Cancellous Hexalobe Screw |
| 3016-40018-S | 3016-40018 | 4.0 mm x 18 mm Partially Threaded Cancellous Hexalobe Screw |
| 3016-40020-S | 3016-40020 | 4.0 mm x 20 mm Partially Threaded Cancellous Hexalobe Screw |
| 3016-40022-S | 3016-40022 | 4.0 mm x 22 mm Partially Threaded Cancellous Hexalobe Screw |

Foot and Ankle Portfolio [continued]

Small Fragment Base Set

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|---|
| 3016-40024-S | 3016-40024 | 4.0 mm x 24 mm Partially Threaded Cancellous Hexalobe Screw |
| 3016-40026-S | 3016-40026 | 4.0 mm x 26 mm Partially Threaded Cancellous Hexalobe Screw |
| 3016-40028-S | 3016-40028 | 4.0 mm x 28 mm Partially Threaded Cancellous Hexalobe Screw |
| 3016-40030-S | 3016-40030 | 4.0 mm x 30 mm Partially Threaded Cancellous Hexalobe Screw |
| 3016-40035-S | 3016-40035 | 4.0 mm x 35 mm Partially Threaded Cancellous Hexalobe Screw |
| 3016-40040-S | 3016-40040 | 4.0 mm x 40 mm Partially Threaded Cancellous Hexalobe Screw |
| 3016-40045-S | 3016-40045 | 4.0 mm x 45 mm Partially Threaded Cancellous Hexalobe Screw |
| 3016-40050-S | 3016-40050 | 4.0 mm x 50 mm Partially Threaded Cancellous Hexalobe Screw |
| 3016-40055-S | 3016-40055 | 4.0 mm x 55 mm Partially Threaded Cancellous Hexalobe Screw |
| 3016-40060-S | 3016-40060 | 4.0 mm x 60 mm Partially Threaded Cancellous Hexalobe Screw |
| 30-2742-S | 30-2742 | 2.7 mm x 42 mm Nonlocking Hexalobe Screw |
| 30-2744-S | 30-2744 | 2.7 mm x 44 mm Nonlocking Hexalobe Screw |
| 30-2746-S | 30-2746 | 2.7 mm x 46 mm Nonlocking Hexalobe Screw |
| 30-2748-S | 30-2748 | 2.7 mm x 48 mm Nonlocking Hexalobe Screw |
| 30-2842-S | 30-2842 | 2.7 mm x 42 mm Locking Hexalobe Screw |
| 30-2844-S | 30-2844 | 2.7 mm x 44 mm Locking Hexalobe Screw |
| 30-2846-S | 30-2846 | 2.7 mm x 46 mm Locking Hexalobe Screw |
| 30-2848-S | 30-2848 | 2.7 mm x 48 mm Locking Hexalobe Screw |
| 30-2850-S | 30-2850 | 2.7 mm x 50 mm Locking Hexalobe Screw |
| 30-2855-S | 30-2855 | 2.7 mm x 55 mm Locking Hexalobe Screw |
| 30-2860-S | 30-2860 | 2.7 mm x 60 mm Locking Hexalobe Screw |
| 30-3542-S | 30-3542 | 3.5 mm x 42 mm Nonlocking Hexalobe Screw |
| 30-3544-S | 30-3544 | 3.5 mm x 44 mm Nonlocking Hexalobe Screw |
| 30-3546-S | 30-3546 | 3.5 mm x 46 mm Nonlocking Hexalobe Screw |
| 30-3548-S | 30-3548 | 3.5 mm x 48 mm Nonlocking Hexalobe Screw |
| 30-3642-S | 30-3642 | 3.5 mm x 42 mm Locking Hexalobe Screw |

Foot and Ankle Portfolio [continued]

Small Fragment Base Set

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|---|
| 30-3644-S | 30-3644 | 3.5 mm x 44 mm Locking Hexalobe Screw |
| 30-3646-S | 30-3646 | 3.5 mm x 46 mm Locking Hexalobe Screw |
| 30-3648-S | 30-3648 | 3.5 mm x 48 mm Locking Hexalobe Screw |
| 7003-07036-S | 7003-07036 | Cannulated Screw Washer 7.0 mm OD x 3.6 mm ID |
| 7008-0103-S | 7008-0103 | One-Third Tubular Plate 3-Hole |
| 7008-0104-S | 7008-0104 | One-Third Tubular Plate 4-Hole |
| 7008-0105-S | 7008-0105 | One-Third Tubular Plate 5-Hole |
| 7008-0106-S | 7008-0106 | One-Third Tubular Plate 6-Hole |
| 7008-0107-S | 7008-0107 | One-Third Tubular Plate 7-Hole |
| 7008-0108-S | 7008-0108 | One-Third Tubular Plate 8-Hole |
| 7008-0110-S | 7008-0110 | One-Third Tubular Plate 10-Hole |
| 7008-0112-S | 7008-0112 | One-Third Tubular Plate 12-Hole |
| 7010-0106N-S | 7010-0106N | Fragment Plate 2.7 mm |
| 7010-0107L-S | 7010-0107L | L Fragment Plate 2.7 mm Left |
| 7010-0107R-S | 7010-0107R | L Fragment Plate 2.7 mm Right |
| 7010-0108N-S | 7010-0108N | T Fragment Plate 2.7 mm |
| AI-0020-S | AI-0020 | 20 mm AcuTwist® Acutrak Compression Screw |
| AI-0026-S | AI-0026 | 26 mm AcuTwist Acutrak Compression Screw |
| AI-0030-S | AI-0030 | 30 mm AcuTwist Acutrak Compression Screw |

Hip and Pelvic Portfolio:

► Pelvic Plating System



Conditions for Safe Scanning in the MRI

A patient with the plate and screw implants from the Pelvic Plating System may be safely scanned at 1.5 T or 3 T under the following conditions. Failure to follow these conditions may result in injury to the patient.

Hip and Pelvic Portfolio

| Parameter | Condition |
|---|---|
| Nominal values(s) of Static Magnetic Field | 1.5 T or 3 T |
| Maximum Spatial Field Gradient | 19 T/m (1900 gauss/cm) |
| RF Excitation | Circularly Polarized (CP) |
| RF Transmit Coil Type | Whole body transmit coil, Head RF transmit-receive coil. |
| Maximum Whole Body SAR [W/kg] | 2.0 W/kg |
| Limits on Scan Duration | 2.0 W/kg whole body average SAR for 13 minutes of continuous RF (a sequence or back to back series/scan without breaks) |
| MR Image Artifact | The presence of this implant may produce an image artifact of 28 mm. |
| If more information about a specific parameter is not included, there are no conditions associated with that parameter. | |

Pelvic Plating System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|--|
| 30-0256-S | 30-0256 | 3.5 mm x 10 mm Nonlocking Hexalobe Screw |
| 30-0257-S | 30-0257 | 3.5 mm x 12 mm Nonlocking Hexalobe Screw |
| 30-0258-S | 30-0258 | 3.5 mm x 14 mm Nonlocking Hexalobe Screw |
| 30-0259-S | 30-0259 | 3.5 mm x 16 mm Nonlocking Hexalobe Screw |
| 30-0260-S | 30-0260 | 3.5 mm x 18 mm Nonlocking Hexalobe Screw |
| 30-0261-S | 30-0261 | 3.5 mm x 20 mm Nonlocking Hexalobe Screw |
| 30-0262-S | 30-0262 | 3.5 mm x 22 mm Nonlocking Hexalobe Screw |
| 30-0263-S | 30-0263 | 3.5 mm x 24 mm Nonlocking Hexalobe Screw |
| 30-0264-S | 30-0264 | 3.5 mm x 26 mm Nonlocking Hexalobe Screw |
| 30-0265-S | 30-0265 | 3.5 mm x 28 mm Nonlocking Hexalobe Screw |
| 30-0266-S | 30-0266 | 3.5 mm x 30 mm Nonlocking Hexalobe Screw |

Hip and Pelvic Portfolio [continued]

Pelvic Plating System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|---|
| 30-0267-S | 30-0267 | 3.5 mm x 32 mm Nonlocking Hexalobe Screw |
| 30-0268-S | 30-0268 | 3.5 mm x 34 mm Nonlocking Hexalobe Screw |
| 30-0269-S | 30-0269 | 3.5 mm x 36 mm Nonlocking Hexalobe Screw |
| 30-0270-S | 30-0270 | 3.5 mm x 38 mm Nonlocking Hexalobe Screw |
| 30-0271-S | 30-0271 | 3.5 mm x 40 mm Nonlocking Hexalobe Screw |
| 30-0272-S | 30-0272 | 3.5 mm x 45 mm Nonlocking Hexalobe Screw |
| 30-0273-S | 30-0273 | 3.5 mm x 50 mm Nonlocking Hexalobe Screw |
| 30-0274-S | 30-0274 | 3.5 mm x 55 mm Nonlocking Hexalobe Screw |
| 30-0275-S | 30-0275 | 3.5 mm x 60 mm Nonlocking Hexalobe Screw |
| 30-0276-S | 30-0276 | 3.5 mm x 65 mm Nonlocking Hexalobe Screw |
| 30-0344-S | 30-0344 | 2.7 mm x 10 mm Nonlocking Hexalobe Screw |
| 30-0345-S | 30-0345 | 2.7 mm x 12 mm Nonlocking Hexalobe Screw |
| 30-0346-S | 30-0346 | 2.7 mm x 14 mm Nonlocking Hexalobe Screw |
| 30-0347-S | 30-0347 | 2.7 mm x 16 mm Nonlocking Hexalobe Screw |
| 30-0348-S | 30-0348 | 2.7 mm x 18 mm Nonlocking Hexalobe Screw |
| 30-0349-S | 30-0349 | 2.7 mm x 20 mm Nonlocking Hexalobe Screw |
| 30-0350-S | 30-0350 | 2.7 mm x 22 mm Nonlocking Hexalobe Screw |
| 30-0351-S | 30-0351 | 2.7 mm x 24 mm Nonlocking Hexalobe Screw |
| 30-0352-S | 30-0352 | 2.7 mm x 26 mm Nonlocking Hexalobe Screw |
| 30-0353-S | 30-0353 | 2.7 mm x 28 mm Nonlocking Hexalobe Screw |
| 30-0354-S | 30-0354 | 2.7 mm x 30 mm Nonlocking Hexalobe Screw |
| 30-0355-S | 30-0355 | 2.7 mm x 32 mm Nonlocking Hexalobe Screw |
| 30-0877-S | 30-0877 | 3.5 mm x 70 mm Nonlocking Hexalobe Screw |
| 30-0878-S | 30-0878 | 3.5 mm x 75 mm Nonlocking Hexalobe Screw |
| 30-0879-S | 30-0879 | 3.5 mm x 80 mm Nonlocking Hexalobe Screw |
| 30-0880-S | 30-0880 | 3.5 mm x 85 mm Nonlocking Hexalobe Screw |
| 30-0881-S | 30-0881 | 3.5 mm x 90 mm Nonlocking Hexalobe Screw |
| 30-0882-S | 30-0882 | 3.5 mm x 95 mm Nonlocking Hexalobe Screw |
| 30-0883-S | 30-0883 | 3.5 mm x 100 mm Nonlocking Hexalobe Screw |
| 30-0884-S | 30-0884 | 3.5 mm x 105 mm Nonlocking Hexalobe Screw |

Hip and Pelvic Portfolio [continued]

Pelvic Plating System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|---|
| 30-0885-S | 30-0885 | 3.5 mm x 110 mm Nonlocking Hexalobe Screw |
| 30-0886-S | 30-0886 | 3.5 mm x 115 mm Nonlocking Hexalobe Screw |
| 30-0887-S | 30-0887 | 3.5 mm x 120 mm Nonlocking Hexalobe Screw |
| 30-0888-S | 30-0888 | 3.5 mm x 125 mm Nonlocking Hexalobe Screw |
| 30-0889-S | 30-0889 | 3.5 mm x 130 mm Nonlocking Hexalobe Screw |
| 30-0890-S | 30-0890 | 3.5 mm x 135 mm Nonlocking Hexalobe Screw |
| 30-0891-S | 30-0891 | 3.5 mm x 140 mm Nonlocking Hexalobe Screw |
| 30-0892-S | 30-0892 | 3.5 mm x 145 mm Nonlocking Hexalobe Screw |
| 30-0893-S | 30-0893 | 3.5 mm x 150 mm Nonlocking Hexalobe Screw |
| 30-0894-S | 30-0894 | 4.3 mm x 50 mm Hexalobe Column Screw |
| 30-0895-S | 30-0895 | 4.3 mm x 55 mm Hexalobe Column Screw |
| 30-0896-S | 30-0896 | 4.3 mm x 60 mm Hexalobe Column Screw |
| 30-0897-S | 30-0897 | 4.3 mm x 65 mm Hexalobe Column Screw |
| 30-0898-S | 30-0898 | 4.3 mm x 70 mm Hexalobe Column Screw |
| 30-0899-S | 30-0899 | 4.3 mm x 75 mm Hexalobe Column Screw |
| 30-0900-S | 30-0900 | 4.3 mm x 80 mm Hexalobe Column Screw |
| 30-0901-S | 30-0901 | 4.3 mm x 85 mm Hexalobe Column Screw |
| 30-0902-S | 30-0902 | 4.3 mm x 90 mm Hexalobe Column Screw |
| 30-0903-S | 30-0903 | 4.3 mm x 95 mm Hexalobe Column Screw |
| 30-0904-S | 30-0904 | 4.3 mm x 100 mm Hexalobe Column Screw |
| 30-0905-S | 30-0905 | 4.3 mm x 105 mm Hexalobe Column Screw |
| 30-0906-S | 30-0906 | 4.3 mm x 110 mm Hexalobe Column Screw |
| 30-0907-S | 30-0907 | 4.3 mm x 115 mm Hexalobe Column Screw |
| 30-0908-S | 30-0908 | 4.3 mm x 120 mm Hexalobe Column Screw |
| 30-0909-S | 30-0909 | 4.3 mm x 125 mm Hexalobe Column Screw |
| 30-0910-S | 30-0910 | 4.3 mm x 130 mm Hexalobe Column Screw |
| 30-0911-S | 30-0911 | 4.3 mm x 135 mm Hexalobe Column Screw |
| 30-0912-S | 30-0912 | 4.3 mm x 140 mm Hexalobe Column Screw |
| 30-0913-S | 30-0913 | 4.3 mm x 145 mm Hexalobe Column Screw |
| 30-0914-S | 30-0914 | 4.3 mm x 150 mm Hexalobe Column Screw |

Hip and Pelvic Portfolio [continued]

Pelvic Plating System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|--|
| 70-0426-S | 70-0426 | Posterior Wall Acetabular Plate |
| 70-0427-S | 70-0427 | Posterior Wall Acetabular Fragment Plate, Left |
| 70-0428-S | 70-0428 | Posterior Wall Acetabular Fragment Plate, Right |
| 70-0429-S | 70-0429 | Acetabular Spring Plate, 2-Hole |
| 70-0430-S | 70-0430 | Acetabular Spring Plate, 3-Hole |
| 70-0431-S | 70-0431 | Anterior Brim Plate, 12-Hole, Left |
| 70-0432-S | 70-0432 | Anterior Brim Plate, 12-Hole, Right |
| 70-0433-S | 70-0433 | Anterior Brim Plate, 14-Hole, Left |
| 70-0434-S | 70-0434 | Anterior Brim Plate, 14-Hole, Right |
| 70-0435-S | 70-0435 | Quadrilateral Surface Plate, Left |
| 70-0436-S | 70-0436 | Quadrilateral Surface Plate, Right |
| 70-0437-S | 70-0437 | Intrapelvic Plate, 5-Hole, Left |
| 70-0438-S | 70-0438 | Intrapelvic Plate, 5-Hole, Right |
| 70-0439-S | 70-0439 | Intrapelvic Plate, 9-Hole, Left |
| 70-0440-S | 70-0440 | Intrapelvic Plate, 9-Hole, Right |
| 70-0441-S | 70-0441 | 3.5 mm Reconstruction Plate, 3-Hole |
| 70-0442-S | 70-0442 | 3.5 mm Reconstruction Plate, 4-Hole |
| 70-0443-S | 70-0443 | 3.5 mm Reconstruction Plate, 6-Hole |
| 70-0444-S | 70-0444 | 3.5 mm Reconstruction Plate, 8-Hole |
| 70-0445-S | 70-0445 | 3.5 mm Reconstruction Plate, 10-Hole |
| 70-0446-S | 70-0446 | 3.5 mm Reconstruction Plate, 12-Hole |
| 70-0447-S | 70-0447 | 3.5 mm Reconstruction Plate, 14-Hole |
| 70-0448-S | 70-0448 | 3.5 mm Reconstruction Plate, 16-Hole |
| 70-0449-S | 70-0449 | 3.5 mm Interlocking Recon Plate, 11-Hole |
| 70-0450-S | 70-0450 | Pubic Symphysis Plate, 4-Hole |
| 70-0451-S | 70-0451 | Pubic Symphysis Plate, 6-Hole |
| 70-0452-S | 70-0452 | Superior Sacroiliac Plate, 4-Hole |
| 70-0458-S | 70-0458 | Curved Posterior Wall Acetabular Fragment Plate, Left |
| 70-0459-S | 70-0459 | Curved Posterior Wall Acetabular Fragment Plate, Right |

External Fixation Systems


- ▶ Stableloc External Fixation System
- ▶ Small Bone Distractor
- ▶ Small Bone Fixator

The Small Bone External Fixation System (Small Bone Distractor and Small Bone Fixator) and the Stableloc External Fixation System (shown below) are MR Unsafe.



MR Unsafe

A patient with the screw and other implants from the Stableloc External Fixation System and Small Bone External Fixation System (Small Bone Distractor and Small Bone Fixator) should not receive an MR exam. Failure to follow these conditions may result in injury to the patient.

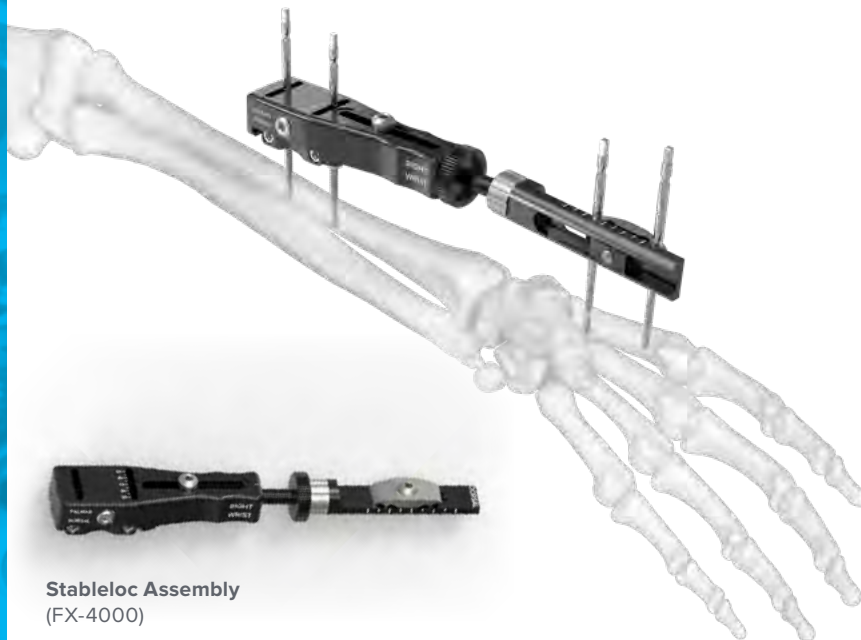


Stableloc Assembly
(FX-4000)



**3.2 mm x 30 mm
Self-Drilling Pin**
(FX-4004-S)



**2.7 mm x 25 mm
Self-Drilling Pin**
(FX-4054-S)

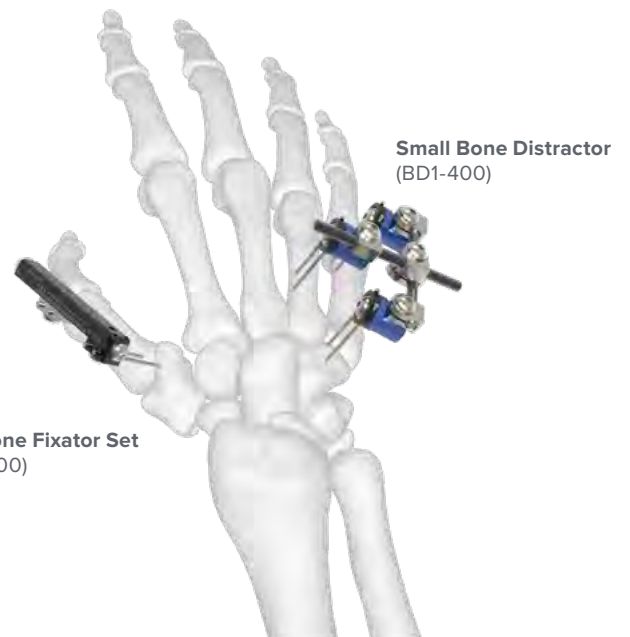

Small Bone Distractor
(BD1-400)

Small Bone Fixator Set
(SBF-0000)

Scapholunate Repair System



Conditions for Safe Scanning in the MRI

A patient with the screw implants from the Scapholunate Repair System may be safely scanned at 1.5 T or 3 T under the following conditions. Failure to follow these conditions may result in injury to the patient.

Scapholunate Repair System

| Parameter | Condition |
|---|--|
| Nominal values(s) of Static Magnetic Field | 1.5 T or 3 T |
| Maximum Spatial Field Gradient | 19 T/m |
| RF Excitation | Circularly Polarized (CP) |
| RF Transmit Coil Type | Whole body transmit coil |
| Maximum Whole Body SAR at 1.5 T and 3.0 T [W/kg] | 2.0 W/kg |
| Limits on Scan Duration 1.5 T | N/A |
| Limits on Scan Duration 3 T | N/A |
| MR Image Artifact | The presence of this implant may produce an image artifact length of 46 mm \pm 5 mm including the screw and a maximal artifact width of 10 mm \pm 3mm surrounding the screw. |
| If more information about a specific parameter is not included, there are no conditions associated with that parameter. | |

Scapholunate Repair System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|---------------------------------------|
| N/A | 55-0029 | Transrotational Screw, 4.5 mm x 22 mm |
| N/A | 55-0030 | Transrotational Screw, 4.5 mm x 23 mm |
| N/A | 55-0031 | Transrotational Screw, 4.5 mm x 24 mm |
| N/A | 55-0032 | Transrotational Screw, 4.5 mm x 25 mm |
| N/A | 55-0033 | Transrotational Screw, 4.5 mm x 26 mm |

Tendon and Ligament Portfolio

- ▶ Acu-Sinch® Knotless System
- ▶ Acu-Sinch® Knotless Mini System
- ▶ PEEK Suture Anchor System



Conditions for Safe Scanning in the MRI

A patient with the button and washer implants from the Acu-Sinch Knotless System and Acu-Sinch Knotless Mini System may be safely scanned at 1.5 T or 3 T under the following conditions. Failure to follow these conditions may result in injury to the patient.

Tendon and Ligament Portfolio

| Parameter | Condition |
|---|---|
| Nominal values(s) of Static Magnetic Field | 1.5 T or 3 T |
| Maximum Spatial Field Gradient | 20 T/m (2000 gauss/cm) |
| RF Excitation | Circularly Polarized (CP) |
| RF Transmit Coil Type | Whole body transmit coil, Head RF transmit-receive coil |
| Maximum Whole Body SAR [W/kg] | 2.0 W/kg |
| Limits on Scan Duration | 2.0 W/kg whole body average SAR for 60 minutes of continuous RF (a sequence or back to back series/scan without breaks) |
| MR Image Artifact | The presence of this implant may produce an image artifact of 18 mm. |
| If more information about a specific parameter is not included, there are no conditions associated with that parameter. | |

Acu-Sinch Knotless System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|--|
| 46-0023-S | N/A | Acu-Sinch Knotless System w/ Inserter 3.5 mm Kit |
| 46-0024-S | N/A | Acu-Sinch Knotless System w/ Inserter 3.5 mm |
| 5502-0205-S | N/A | Acu-Sinch Washer 3.5 mm Kit |

Acu-Sinch Knotless Mini System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|---|
| 46-0040-S | N/A | Acu-Sinch Knotless Mini System w/ Inserter 2.3 mm |

Tendon and Ligament Portfolio [continued]



MRI Compatibility Claim

PEEK is a thermoplastic polymer that is nonconductive, nonmetallic, nonmagnetic with no known hazards resulting from exposure of this device to any magnetic resonance environment, hence the device is considered MR Safe.

PEEK Suture Anchor System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|--|
| 46-0008-S | N/A | 4.5 mm PEEK Suture Anchor Kit |
| 46-0009-S | N/A | 5.5 mm PEEK Suture Anchor Kit |
| 46-0012-S | N/A | 4.5 mm PEEK Suture Anchor w/o Drill, Tap |
| 46-0013-S | N/A | 5.5 mm PEEK Suture Anchor w/o Drill, Tap |

Chest Wall Portfolio

► RibLoc® U Plus Chest Wall Plating System

RibLoc U Plus wall plating system **when used for indications on the sternum** is MR unsafe. Otherwise, MR conditional.



Conditions for Safe Scanning in the MRI

A patient with the screw and plate implants from the RibLoc U Plus Chest Wall Plating System implanted on the rib maybe safely scanned under the following conditions. Failure to follow these conditions may result in injury to the patient.

Chest Wall Portfolio

| Parameter | Condition |
|--|---|
| Nominal values(s) of Static Magnetic Field | 1.5 T or 3 T |
| Maximum Spatial Field Gradient | 40 T/m (4000 gauss/cm) |
| RF Excitation | Circularly Polarized (CP) |
| RF Transmit Coil Type | Whole body transmit coil, Head RF transmit-receive coil |
| Maximum Whole Body SAR [W/kg] | 2.0 W/kg |
| Limits on Scan Duration | 2.0 W/kg whole body average SAR for 15 minutes of continuous RF (a sequence or back to back series/scan without breaks) |
| MR Image Artifact | The presence of this implant may produce an image artifact of 10 mm. |

Note: MRI Conditional scan parameters are applicable to *rib plating only*. Sternum plating is MRI Unsafe due to temperature rise.

If more information about a specific parameter is not included, there are no conditions associated with that parameter.

Chest Wall Portfolio [continued]

**MR Unsafe**

A patient with the RibLoc U Plus wall plating system **implanted on the sternum** should not receive an MR exam. Failure to follow these conditions may result in injury to the patient.


**RibLoc U Plus Chest Wall Plating System**

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|------------------------------|
| N/A | RBL1221 | 6 mm x 2.7 mm Locking Screw |
| N/A | RBL1222 | 8 mm x 2.7 mm Locking Screw |
| N/A | RBL1223 | 10 mm x 2.7 mm Locking Screw |
| N/A | RBL1224 | 12 mm x 2.7 mm Locking Screw |
| N/A | RBL1225 | 14 mm x 2.7 mm Locking Screw |
| N/A | RBL1301 | 50 mm Rib Plate |
| N/A | RBL1302 | 75 mm Rib Plate |
| N/A | RBL1303 | 115 mm Rib Plate |
| N/A | RBL1304 | 155 mm Rib Plate |
| N/A | RBL1305 | 215 mm Rib Plate |
| N/A | RBL1401 | 126 mm Straight Plate |

Craniomaxillofacial (CMF) Portfolio

- ▶ 2.0 Orthognathic
- ▶ OSA™
- ▶ OsteoMed CFx™ Craniomaxillofacial Fixation System
- ▶ OsteoMed MFx™ Mandibular Fixation System
- ▶ OsteoMed ICON™ Facial Plating System
- ▶ 2.0 mm Auto Drive MMF Screws

MRI Restricted Zone Summary

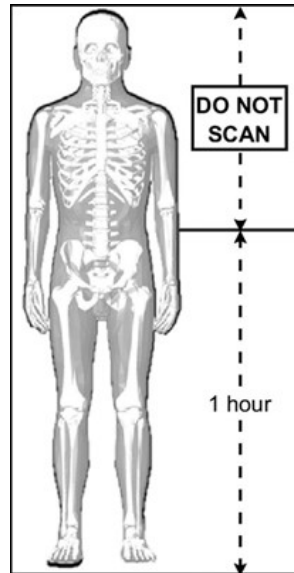

Conditions for Safe Scanning in the MRI



A patient with the screw and plate of 2.0 mm Auto Drive MMF Screws, 2.0 Orthognathic, OsteoMed CFx, OsteoMed ICON Facial Plating System, OsteoMed MFx, or OSA Systems may be safely scanned at 1.5 T or 3.0 T under the following conditions. Failure to follow these conditions may result in injury to the patient. See below for a complete list of part numbers included in these systems.

OsteoMed's Craniomaxillofacial Products

| Parameter | Condition |
|--|---|
| Device Name | OsteoMed's CMF Products |
| Static Magnetic Field Strength (B ₀) | 1.5T and 3T |
| MR Scanner Type | Cylindrical |
| B ₀ Field Orientation | Horizontal |
| Maximum Spatial Field Gradient | 23.5 T/m (2,350 G/cm) |
| Maximum Gradient Slew Rate | 200 T/m/s per axis |
| RF Excitation | Circularly Polarized (CP) |
| RF Transmit Coil Type | Integrated Whole Body Transmit Coil |
| Operating Mode | Normal Operating Mode |
| Scan Duration | Scan duration up to 1 hour without a cooling period when landmarked at or below the navel (as detailed in the MRI Restricted Zone Summary) |
| Scan Regions | Landmark locations at or below the navel are permissible (as detailed in the MRI Restricted Zone Summary) |
| Image Artifact | The presence of OsteoMed's CMF and Neuro Products may produce an image artifact of 3.0 cm. Some manipulation of scan parameters may be needed to compensate for the artifact. |
| Parameters may be needed to compensate for the artifact. | |

Craniomaxillofacial (CMF) Portfolio [continued]

2.0 Orthognathic

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|----------------------------|
| 202-2004-SP | 202-2004 | 2.0 mm x 4 mm Screw |
| 202-2005-SP | 202-2005 | 2.0 mm x 5 mm Screw |
| 202-2006-SP | 202-2006 | 2.0 mm x 6 mm Screw |
| 202-2008-SP | 202-2008 | 2.0 mm x 8 mm Screw |
| 202-2010-SP | 202-2010 | 2.0 mm x 10 mm Screw |
| 202-2012-SP | 202-2012 | 2.0 mm x 12 mm Screw |
| 202-2014-SP | 202-2014 | 2.0 mm x 14 mm Screw |
| 202-2016-SP | 202-2016 | 2.0 mm x 16 mm Screw |
| 202-2018-SP | 202-2018 | 2.0 mm x 18 mm Screw |
| 202-2304-SP | 202-2304 | 2.3 mm x 4 mm Safety Screw |
| N/A | 202-2305 | 2.3 mm x 5 mm Safety Screw |
| 202-2306-SP | 202-2306 | 2.3 mm x 6 mm Safety Screw |
| N/A | 202-2308 | 2.3 mm x 8 mm Safety Screw |
| 206-2404-SP | 206-2404 | 2.4 mm x 4 mm Screw |
| 206-2405-SP | 206-2405 | 2.4 mm x 5 mm Screw |
| 206-2406-SP | 206-2406 | 2.4 mm x 6 mm Screw |
| 206-2408-SP | 206-2408 | 2.4 mm x 8 mm Screw |
| 206-2410-SP | 206-2410 | 2.4 mm x 10 mm Screw |
| 206-2412-SP | 206-2412 | 2.4 mm x 12 mm Screw |
| 206-2414-SP | 206-2414 | 2.4 mm x 14 mm Screw |
| 206-2416-SP | 206-2416 | 2.4 mm x 16 mm Screw |
| 210-0040-SP | 210-0040 | 4 mm Chin Plate |
| 210-0041-SP | 210-0041 | 2 mm Chin Plate |
| N/A | 210-0042 | 6 mm Chin Plate |
| 210-0043-SP | 210-0043 | 3 mm Chin Plate |
| 210-0044-SP | 210-0044 | 8 mm Chin Plate |
| 210-0046-SP | 210-0046 | 10 mm Chin Plate |
| 210-0048-SP | 210-0048 | 12 mm Chin Plate |
| N/A | 210-0062 | 5 mm Chin Plate |

Craniomaxillofacial (CMF) Portfolio [continued]

2.0 Orthognathic

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|------------------------------------|
| N/A | 210-0063 | 7 mm Chin Plate |
| N/A | 210-0064 | 9 mm Chin Plate |
| N/A | 210-0065 | 11 mm Chin Plate |
| 210-0300-SP | 210-0300 | 5 mm Zygomatic Rigid Plate |
| 210-0301-SP | 210-0301 | 5 mm Pyriform Rigid Plate |
| 210-0302-SP | 210-0302 | 8 mm Zygomatic Rigid Plate |
| 210-0303-SP | 210-0303 | 8 mm Pyriform Rigid Plate |
| 210-0304-SP | 210-0304 | 12 mm Zygomatic Rigid Plate |
| 210-0305-SP | 210-0305 | 12 mm Pyriform Rigid Plate |
| 210-0306-SP | 210-0306 | 8 mm "Y" Rigid Plate |
| 210-0307-SP | 210-0307 | 13 mm "Y" Rigid Plate |
| 210-0308-SP | 210-0308 | 13 mm Segmental Rigid Plate |
| 210-0309-SP | 210-0309 | 16 mm, Zygomatic Rigid Plate |
| 210-0310-SP | 210-0310 | 16 mm, Pyriform Rigid Plate |
| 210-0320-SP | 210-0320 | 5 mm Zygomatic Plate |
| 210-0321-SP | 210-0321 | 5 mm Pyriform Plate |
| 211-2004-SP | 211-2004 | 2.0 mm x 4 mm Auto-Drive Screw |
| 211-2005-SP | 211-2005 | 2.0 mm x 5 mm Auto-Drive Screw |
| 211-2006-SP | 211-2006 | 2.0 mm x 6 mm Auto-Drive Screw |
| 211-2008-SP | 211-2008 | 2.0 mm x 8 mm Auto-Drive Screw |
| 214-0380-SP-01 | 214-0380 | BSSO Plate, 2.0 mm System, Medium |
| 214-0382-SP-01 | 214-0382 | BSSO Plate, 2.0 mm System, Long |
| 214-0386-SP-01 | 214-0386 | BSSO Plate, 2.0 mm System, X-Long |
| N/A | 214-0388 | BSSO Plate, 2.0 mm System, XX-Long |
| N/A | 202-2310 | 2.3 mm x 10 mm Safety Screw |

Craniomaxillofacial (CMF) Portfolio [continued]

OSA™

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|--|
| N/A | 202-2008 | 2.0 mm x 8 mm Screw |
| N/A | 202-2010 | 2.0 mm x 10 mm Screw |
| N/A | 202-2012 | 2.0 mm x 12 mm Screw |
| N/A | 202-2014 | 2.0 mm x 14 mm Screw |
| N/A | 202-2304 | 2.3 mm x 4 mm Safety Screw |
| N/A | 202-2305 | 2.3 mm x 5 mm Safety Screw |
| N/A | 202-2308 | 2.3 mm x 8 mm Safety Screw |
| N/A | 202-2310 | 2.3 mm x 10 mm Safety Screw |
| N/A | 208-1904 | 1.9 mm x 4 mm Safety Screw |
| N/A | 208-1905 | 1.9 mm x 5 mm Safety Screw |
| N/A | 208-1906 | 1.9 mm x 6 mm Safety Screw |
| N/A | 208-1908 | 1.9 mm x 8 mm Safety Screw |
| N/A | 210-0404 | ISO Chin Plate with Window, 4 mm Step |
| N/A | 210-0406 | ISO Chin Plate with Window, 6 mm Step |
| N/A | 210-0408 | ISO Chin Plate with Window, 8 mm Step |
| N/A | 210-0410 | ISO Chin Plate with Window, 10 mm Step |
| N/A | 210-0418 | Double Row BSSO Plate, 6 mm |
| N/A | 210-0419 | Double Row BSSO Plate, 8 mm |
| N/A | 210-0420 | Double Row BSSO Plate, 10 mm |
| N/A | 210-0421 | Double Row BSSO Plate, 12 mm |
| N/A | 210-0422 | Double Row BSSO Plate, 14 mm |
| N/A | 211-1604 | 1.6 mm x 4 mm Auto-Drive Screw |
| N/A | 211-1605 | 1.6 mm x 5 mm Auto-Drive Screw |
| N/A | 211-1606 | 1.6 mm x 6 mm Auto-Drive Screw |
| N/A | 211-1608 | 1.6 mm x 8 mm Auto-Drive Screw |
| N/A | 211-2004 | 2.0 mm x 4 mm Auto-Drive Screw |
| N/A | 211-2005 | 2.0 mm x 5 mm Auto-Drive Screw |
| N/A | 211-2006 | 2.0 mm x 6 mm Auto-Drive Screw |
| N/A | 211-2008 | 2.0 mm x 8 mm Auto-Drive Screw |
| N/A | 213-0048 | 1.6 mm Pre-bent Plt 2 mm Stp Lf |

Craniomaxillofacial (CMF) Portfolio [continued]

OSA™

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|---|
| N/A | 213-0049 | 1.6 mm Pre-bent Plt 2 mm Stp Rt |
| N/A | 213-0050 | 1.6mm Pre-bent Plt 4 mm Stp Lf |
| N/A | 213-0051 | 1.6 mm Pre-bent Plt 4 mm Stp Rt |
| N/A | 213-0052 | 1.6 mm Pre-bent Plt 6 mm Stp Lf |
| N/A | 213-0053 | 1.6 mm Pre-bent Plt 6 mm Stp Rt |
| N/A | 213-0054 | 1.6 mm Pre-bent Plt 8 mm Stp Lf |
| N/A | 213-0055 | 1.6 mm Pre-bent Plt 8 mm Stp Rt |
| N/A | 213-0056 | 1.6 mm Pre-bent Plt 10 mm Stp Lf |
| N/A | 213-0057 | 1.6 mm Pre-bent Plt 10 mm Stp Rt |
| N/A | 213-0058 | 1.6 mm Pre-bent Plt 12 mm Stp Lf |
| N/A | 213-0059 | 1.6 mm Pre-bent Plt 12 mm Stp Rt |
| N/A | 213-0062 | Buttress Plate Right |
| N/A | 213-0063 | Buttress Plate Left |
| N/A | 213-0064 | 1.6 mm Buttress Plate, Left, Modified |
| N/A | 213-0065 | 1.6 mm Buttress Plate, Right, Modified |
| N/A | 213-0080 | 1.6 mm Pre-bent Plt 3 mm Stp Lf |
| N/A | 213-0081 | 1.6 mm Pre-bent Plt 3 mm Stp Rt |
| N/A | 213-0082 | 1.6 mm Pre-bent Plt 5 mm Stp Lf |
| N/A | 213-0083 | 1.6 mm Pre-bent Plt 5 mm Stp Rt |
| N/A | 213-0084 | 1.6 mm Pre-bent Plt 7 mm Stp Lf |
| N/A | 213-0085 | 1.6 mm Pre-bent Plt 7 mm Stp Rt |
| N/A | 213-0200 | 1.6 mm Modified LeFort Plate, Left, 2 mm |
| N/A | 213-0201 | 1.6 mm Modified LeFort Plate, Right, 2 mm |
| N/A | 213-0202 | 1.6 mm Modified LeFort Plate, Left, 3 mm |
| N/A | 213-0203 | 1.6 mm Modified LeFort Plate, Right, 3 mm |
| N/A | 213-0204 | 1.6 mm Modified LeFort Plate, Left, 4 mm |
| N/A | 213-0205 | 1.6 mm Modified LeFort Plate, Right, 4 mm |
| N/A | 213-0206 | 1.6 mm Modified LeFort Plate, Left, 5 mm |
| N/A | 213-0207 | 1.6 mm Modified LeFort Plate, Right, 5 mm |
| N/A | 213-0208 | 1.6 mm Modified LeFort Plate, Left, 6 mm |

Craniomaxillofacial (CMF) Portfolio [continued]

OSA™

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|---|
| N/A | 213-0209 | 1.6 mm Modified LeFort Plate, Right, 6 mm |
| N/A | 213-0210 | 1.6 mm Modified LeFort Plate, Left, 7 mm |
| N/A | 213-0212 | 1.6 mm Modified LeFort Plate, Left, 8 mm |
| N/A | 213-0213 | 1.6 mm Modified LeFort Plate, Right, 8 mm |

OsteoMed CFx™ Craniomaxillofacial Fixation System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|---|
| N/A | 207-2412 | 2.4 mm x 12 mm MMF Screw |
| N/A | 207-2414 | 2.4 mm x 14 mm MMF Screw |
| N/A | 207-2416 | 2.4 mm x 16 mm MMF Screw |
| N/A | 207-2418 | 2.4 mm x 18 mm MMF Screw |
| N/A | 209-2008 | 2.0 mm x 8 mm MMF Auto-Drive Screw |
| N/A | 209-2011 | 2.0 mm x 11 mm MMF Auto-Drive Screw |
| N/A | 209-2014 | 2.0 mm x 14 mm MMF Auto-Drive Screw |
| N/A | 210-3018 | 4 Hole Straight Rigid Plate |
| N/A | 214-0111 | 2.0 mm 3 Hole Lozenge Straight Reduction Plate Short |
| N/A | 214-0112 | 2.0 mm 3 Hole Lozenge Straight Reduction Plate Medium |
| N/A | 214-0113 | 2.0 mm 3 Hole Lozenge Straight Reduction Plate Long |
| N/A | 214-0202 | Recon Plate, LT, 4 x 12 Hole |

Craniomaxillofacial (CMF) Portfolio [continued]

OsteoMed MFx™ Mandibular Fixation System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|------------------------------|
| N/A | 202-0206 | 2.0 mm x 6 mm Locking Screw |
| N/A | 202-0208 | 2.0 mm x 8 mm Locking Screw |
| N/A | 202-0210 | 2.0 mm x 10 mm Locking Screw |
| N/A | 202-0212 | 2.0 mm x 12 mm Locking Screw |
| N/A | 202-0214 | 2.0 mm x 14 mm Locking Screw |
| N/A | 202-0216 | 2.0 mm x 16 mm Locking Screw |
| N/A | 202-0218 | 2.0 mm x 18 mm Locking Screw |
| N/A | 202-2004 | 2.0 mm x 4 mm Screw |
| N/A | 202-2005 | 2.0 mm x 5 mm Screw |
| 202-2006-SP | 202-2006 | 2.0 mm x 6 mm Screw |
| 202-2008-SP | 202-2008 | 2.0 mm x 8 mm Screw |
| 202-2010-SP | 202-2010 | 2.0 mm x 10 mm Screw |
| 202-2012-SP | 202-2012 | 2.0 mm x 12 mm Screw |
| 202-2014-SP | 202-2014 | 2.0 mm x 14 mm Screw |
| 202-2016-SP | 202-2016 | 2.0 mm x 16 mm Screw |
| 202-2018-SP | 202-2018 | 2.0 mm x 18 mm Screw |
| 202-2020-SP | 202-2020 | 2.0 mm x 20 mm Screw |
| N/A | 202-2022 | 2.0 mm x 22 mm Screw |
| N/A | 202-2304 | 2.3 mm x 4 mm Safety Screw |
| N/A | 202-2305 | 2.3 mm x 5 mm Safety Screw |
| N/A | 202-2308 | 2.3 mm x 8 mm Safety Screw |
| N/A | 202-2310 | 2.3 mm x 10 mm Safety Screw |
| N/A | 204-1202 | 1.2 mm x 2 mm Screw |
| N/A | 204-1203 | 1.2 mm x 3 mm Screw |
| 204-1204-SP | 204-1204 | 1.2 mm x 4 mm Screw |
| N/A | 204-1206 | 1.2 mm x 6 mm Screw |
| 204-1208-SP | 204-1208 | 1.2 mm x 8 mm Screw |
| 204-1210-SP | 204-1210 | 1.2 mm x 10 mm Screw |
| 204-1212-SP | 204-1212 | 1.2 mm x 12 mm Screw |
| N/A | 204-1214 | 1.2 mm x 14 mm Screw |

Craniomaxillofacial (CMF) Portfolio [continued]

OsteoMed MFx™ Mandibular Fixation System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|-----------------------------------|
| N/A | 204-1503 | 1.5 mm x 3 mm Safety Screw |
| N/A | 204-1506 | 1.5 mm x 6 mm Safety Screw |
| N/A | 206-0206 | 2.4 mm x 6 mm Locking Screw |
| N/A | 206-0208 | 2.4 mm x 8 mm Locking Screw |
| N/A | 206-0210 | 2.4 mm x 10 mm Locking Screw |
| N/A | 206-0212 | 2.4 mm x 12 mm Locking Screw |
| N/A | 206-0214 | 2.4 mm x 14 mm Locking Screw |
| N/A | 206-0216 | 2.4 mm x 16 mm Locking Screw |
| N/A | 206-0218 | 2.4 mm x 18 mm Locking Screw |
| N/A | 206-2404 | 2.4 mm x 4 mm Screw |
| N/A | 206-2405 | 2.4 mm x 5 mm Screw |
| 206-2406-SP | 206-2406 | 2.4 mm x 6 mm Screw |
| 206-2408-SP | 206-2408 | 2.4 mm x 8 mm Screw |
| N/A | 206-2410 | 2.4 mm x 10 mm Screw |
| N/A | 206-2412 | 2.4 mm x 12 mm Screw |
| N/A | 206-2414 | 2.4 mm x 14 mm Screw |
| N/A | 206-2416 | 2.4 mm x 16 mm Screw |
| N/A | 206-2418 | 2.4 mm x 18 mm Screw |
| N/A | 206-2420 | 2.4 mm x 20 mm Screw |
| N/A | 206-2422 | 2.4 mm x 22 mm Screw |
| N/A | 206-2704 | Maxi Safety Screw, 2.7 mm x 4 mm |
| N/A | 206-2706 | Maxi Safety Screw, 2.7 mm x 6 mm |
| N/A | 206-2708 | 2.7 mm x 8 mm Safety Screw |
| N/A | 206-2710 | Maxi Safety Screw, 2.7 mm x 10 mm |
| N/A | 206-2712 | 2.7 mm x 12 mm Safety Screw |
| N/A | 206-2714 | Maxi Safety Screw, 2.7 mm x 14 mm |
| N/A | 206-2716 | 2.7 mm x 16 mm Safety Screw |
| N/A | 206-2718 | Maxi Safety Screw, 2.7 mm x 18 mm |
| N/A | 206-2720 | 2.7 mm x 20 mm Safety Screw |
| N/A | 208-1603 | 1.6 mm x 3 mm Screw |

Craniomaxillofacial (CMF) Portfolio [continued]

OsteoMed MFix™ Mandibular Fixation System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|--------------------------------------|
| 208-1604-SP | 208-1604 | 1.6 mm x 4 mm Screw |
| 208-1606-SP | 208-1606 | 1.6 mm x 6 mm Screw |
| 208-1608-SP | 208-1608 | 1.6 mm x 8 mm Screw |
| 208-1610-SP | 208-1610 | 1.6 mm x 10 mm Screw |
| 208-1612-SP | 208-1612 | 1.6 mm x 12 mm Screw |
| 208-1614-SP | 208-1614 | 1.6 mm x 14 mm Screw |
| N/A | 208-1904 | 1.9 mm x 4 mm Safety Screw |
| N/A | 208-1905 | 1.9 mm x 5 mm Safety Screw |
| N/A | 208-1906 | 1.9 mm x 6 mm Safety Screw |
| N/A | 208-1908 | 1.9 mm x 8 mm Safety Screw |
| N/A | 208-2004 | 2.0 mm x 4 mm Midface Screw |
| N/A | 208-2005 | 2.0 mm x 5 mm Midface Screw |
| N/A | 208-2006 | 2.0 mm x 6 mm Midface Screw |
| N/A | 208-2008 | 2.0 mm x 8 mm Midface Screw |
| N/A | 208-2010 | 2.0 mm x 10 mm Midface Screw |
| N/A | 208-2012 | 2.0 mm x 12 mm Midface Screw |
| N/A | 208-2014 | 2.0 mm x 14 mm Midface Screw |
| N/A | 208-2016 | 2.0 mm x 16 mm Midface Screw |
| N/A | 209-1904 | 1.9 mm x 4 mm Midface Safety Screw |
| N/A | 209-1905 | 1.9 mm x 5 mm Midface Safety Screw |
| N/A | 209-1906 | 1.9 mm x 6 mm Midface Safety Screw |
| N/A | 209-1908 | 1.9 mm x 8 mm Midface Safety Screw |
| N/A | 209-1910 | 1.9 mm x 10 mm Midface Safety Screw |
| N/A | 209-1912 | 1.9 mm x 12 mm Midface Safety Screw |
| N/A | 209-1914 | 1.9 mm x 14 mm Midface Safety Screw |
| N/A | 209-2304 | 2.3 mm x 4 mm Mandible Safety Screw |
| N/A | 209-2305 | 2.3 mm x 5 mm Mandible Safety Screw |
| N/A | 209-2306 | 2.3 mm x 6 mm Mandible Safety Screw |
| N/A | 209-2308 | 2.3 mm x 8 mm Mandible Safety Screw |
| N/A | 209-2310 | 2.3 mm x 10 mm Mandible Safety Screw |

Craniomaxillofacial (CMF) Portfolio [continued]

OsteoMed MFx™ Mandibular Fixation System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|--------------------------------------|
| N/A | 209-2312 | 2.3 mm x 12 mm Mandible Safety Screw |
| N/A | 209-2314 | 2.3 mm x 14 mm Mandible Safety Screw |
| N/A | 209-2316 | 2.3 mm x 16 mm Mandible Safety Screw |
| N/A | 209-2318 | 2.3 mm x 18 mm Mandible Safety Screw |
| N/A | 209-2320 | 2.3 mm x 20 mm Mandible Safety Screw |
| N/A | 209-2322 | 2.3 mm x 22 mm Mandible Safety Screw |
| N/A | 209-2324 | 2.3 mm x 24 mm Mandible Safety Screw |
| N/A | 210-0002 | 15 mm 90° Right "L" Plate |
| N/A | 210-0003 | 20 mm 90° Right "L" Plate |
| N/A | 210-0004 | 15 mm 90° Left "L" Plate |
| N/A | 210-0005 | 20 mm 90° Left "L" Plate |
| 210-0006-SP | 210-0006 | 15 mm Right "L" Plate |
| 210-0008-SP | 210-0008 | 15 mm Left "L" Plate |
| N/A | 210-0010 | 18 mm Right "L" Plate |
| N/A | 210-0012 | 18 mm Left "L" Plate |
| 210-0014-SP | 210-0014 | 23 mm Right "L" Plate |
| 210-0016-SP | 210-0016 | 23 mm Left "L" Plate |
| N/A | 210-0018 | 26 mm "Y" Plate |
| N/A | 210-0020 | 15 mm Right ZED Plate |
| N/A | 210-0022 | 15 mm Left ZED Plate |
| N/A | 210-0024 | 2.0 mm, ZED Plate, 18 mm, Right |
| N/A | 210-0026 | 2.0 mm, ZED Plate, 18 mm, Left |
| N/A | 210-0030 | 24 mm 4 Hole Straight Plate |
| N/A | 210-0032 | 10 Hole Straight Plate |
| N/A | 210-0034 | 16 Hole Straight Plate |
| N/A | 210-0035 | 16 mm "I" Plate |
| N/A | 210-0036 | 22 mm "T" Plate |
| N/A | 210-0038 | 18 Hole Straight Plate |
| N/A | 210-0039 | Large Burr Hole Plate |
| N/A | 210-0040 | 4 mm Chin Plate |

Craniomaxillofacial (CMF) Portfolio [continued]

OsteoMed MFx™ Mandibular Fixation System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|-------------------------------------|
| N/A | 210-0041 | 2 mm Chin Plate |
| N/A | 210-0042 | 6 mm Chin Plate |
| N/A | 210-0043 | 3 mm Chin Plate |
| N/A | 210-0044 | 8 mm Chin Plate |
| N/A | 210-0046 | 10 mm Chin Plate |
| N/A | 210-0048 | 12 mm Chin Plate |
| N/A | 210-0050 | 40 Hole Straight Plate |
| N/A | 210-0051 | 6 Hole Straight Plate |
| N/A | 210-0052 | 7 mm "X" Plate |
| N/A | 210-0054 | 10 mm Square Plate |
| N/A | 210-0062 | 5 mm Chin Plate |
| N/A | 210-0063 | 7 mm Chin Plate |
| N/A | 210-0064 | 9 mm Chin Plate |
| N/A | 210-0065 | 11 mm Chin Plate |
| N/A | 210-0076 | 127 mm x 127 mm Mesh Plate |
| N/A | 210-0078 | 90 mm x 127 mm Mesh Plate |
| N/A | 210-0080 | Extended "Y" Plate |
| N/A | 210-0081 | Extended Right "L" Plate |
| N/A | 210-0082 | Extended Left "L" Plate |
| N/A | 210-0090 | 17 mm Right "L" Plate, 1.0 mm Thick |
| N/A | 210-0091 | 17 mm Left "L" Plate, 1.0 mm Thick |
| N/A | 210-0092 | 19 mm Right "L" Plate, 1.0 mm Thick |
| N/A | 210-0093 | 19 mm Left "L" Plate, 1.0 mm Thick |
| N/A | 210-0102 | 4 Hole Compression Plate |
| N/A | 210-0104 | 5 Hole Compression Plate |
| N/A | 210-0105 | 4 Hole Curved Compression Plate |
| N/A | 210-0106 | 6 Hole Compression Plate |
| N/A | 210-0107 | 5 Hole Curved Compression Plate |
| N/A | 210-0108 | 5 Hole Curved Plate |
| N/A | 210-0109 | 7 Hole Curved Compression Plate |

Craniomaxillofacial (CMF) Portfolio [continued]

OsteoMed MFx™ Mandibular Fixation System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|---------------------------------------|
| N/A | 210-0110 | 7 Hole Curved Plate |
| N/A | 210-0111 | 2.0 mm, 8 Hole Curved Plate |
| N/A | 210-0112 | 9 Hole Curved Plate |
| N/A | 210-0115 | 150 mm x 150 mm Mesh Plate |
| N/A | 210-0116 | 150 mm x 300 mm Mesh Plate |
| N/A | 210-0120 | 2.0 mm, 2 Hole Straight Plate |
| N/A | 210-0206 | 15 mm Right "L" Rigid Plate |
| N/A | 210-0208 | 15 mm Left "L" Rigid Plate |
| N/A | 210-0210 | 18 mm Right "L" Rigid Plate |
| N/A | 210-0212 | 18 mm Left "L" Rigid Plate |
| N/A | 210-0214 | 23 mm Right "L" Rigid Plate |
| N/A | 210-0216 | 23 mm Left "L" Rigid Plate |
| N/A | 210-0218 | 26 mm "Y" Rigid Plate |
| N/A | 210-0220 | 15 mm Right ZED Rigid Plate |
| N/A | 210-0222 | 15 mm Left ZED Rigid Plate |
| N/A | 210-0224 | 2.0 mm, Rigid ZED Plate, 18 mm, Right |
| N/A | 210-0226 | 2.0 mm, Rigid ZED Plate, 18 mm, Left |
| N/A | 210-0230 | 24 mm 4 Hole Straight Rigid Plate |
| N/A | 210-0231 | 29 mm 4 Hole Straight Rigid Plate |
| N/A | 210-0232 | 10 Hole Straight Rigid Plate |
| N/A | 210-0235 | 16 mm "I" Rigid Plate |
| N/A | 210-0236 | "Y" Rigid Plate |
| N/A | 210-0237 | 30 mm 4 Hole Straight Rigid Plate |
| N/A | 210-0239 | Double "Y" Rigid Plate |
| N/A | 210-0250 | 2.0 mm Right L Plate, 15 mm, Rigid |
| N/A | 210-0251 | 2.0 mm Left L Plate, 15 mm, Rigid |
| N/A | 210-0252 | 2.0 mm Right L Plate, 18 mm, Rigid |
| N/A | 210-0253 | 2.0 mm Left L Plate, 18 mm, Rigid |
| N/A | 210-0254 | 2.0 mm Right L Plate, 23 mm, Rigid |
| N/A | 210-0255 | 2.0 mm Left L Plate, 23 mm, Rigid |

Craniomaxillofacial (CMF) Portfolio [continued]

OsteoMed MFx™ Mandibular Fixation System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|---|
| N/A | 210-0501 | Mini, Locking Strut Plate, Curved, 2 x 4 Hole |
| N/A | 210-0502 | Mini, Locking Strut Plate, Curved, 2 x 6 Hole |
| N/A | 210-0503 | Mini, Locking Strut Plate, Curved, 2 x 8 Hole |
| N/A | 210-1010 | 3 x 2 Hole Right "L" Plate |
| N/A | 210-1012 | 3 x 2 Hole Left "L" Plate |
| N/A | 210-1013 | 12 Hole Straight Plate, 1.0 mm Thick |
| N/A | 210-1014 | 48 Hole Double Row Straight Plate |
| N/A | 210-1020 | Mini, 34 Hole Strip, Rigid |
| N/A | 210-1021 | Mini 34 Hole Strip |
| N/A | 210-1022 | Mini 68 Hole Strip, Rigid |
| N/A | 210-1023 | Mini 68 Hole Strip |
| N/A | 210-1024 | Mini 124 mm x 124 mm Mesh, Rigid |
| N/A | 210-1025 | Mini 124 mm x 124 mm Mesh |
| N/A | 210-1026 | Mini, 68 Hole Strip, Low Profile |
| N/A | 210-1030 | Mini 4 Hole Str. Plate, 1.0 mm |
| N/A | 210-1031 | Mini 6 Hole Straight Plate, 1.0 mm |
| N/A | 210-1032 | Mini 8 Hole Straight Plate, 1.0 mm |
| 210-1036-SP | 210-1036 | Mini 16 Hole Straight Plate, 1.0 mm |
| N/A | 210-1037 | 2.0 mm, 4 Hole Straight Rigid Plate |
| N/A | 210-1038 | 2.0 mm, 12 Hole Double Row Straight Rigid Plate |
| N/A | 210-1039 | 2.0 mm, 10 Hole Double Row Straight Rigid Plate |
| N/A | 210-1210 | 8 Hole "U" Plate |
| N/A | 210-1213 | 7 Hole "U" Plate |
| N/A | 210-3018 | 4 Hole Straight Rigid Plate |
| N/A | 211-1203 | 1.2 mm x 3 mm Auto-Drive® Screw |
| N/A | 211-1204 | 1.2 mm x 4 mm Auto-Drive Screw |
| N/A | 211-1205 | 1.2 mm x 5 mm Auto-drive Screw |
| N/A | 211-1603 | 1.6 mm x 3.5 mm Auto-Drive Screw |
| N/A | 211-1604 | 1.6 mm x 4 mm Auto-Drive Screw |
| N/A | 211-1605 | 1.6 mm x 5 mm Auto-Drive Screw |

Craniomaxillofacial (CMF) Portfolio [continued]

OsteoMed MFx™ Mandibular Fixation System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|--------------------------------|
| N/A | 211-1606 | 1.6 mm x 6 mm Auto-Drive Screw |
| 211-1608-SP | 211-1608 | 1.6 mm x 8 mm Auto-Drive Screw |
| N/A | 211-2004 | 2.0 mm x 4 mm Auto-Drive Screw |
| 211-2005-SP | 211-2005 | 2.0 mm x 5 mm Auto-Drive Screw |
| N/A | 211-2006 | 2.0 mm x 6 mm Auto-Drive Screw |
| N/A | 211-2008 | 2.0 mm x 8 mm Auto-Drive Screw |
| N/A | 212-0030 | 7 Hole "Y" Plate |
| N/A | 212-0031 | 11 Hole "Y" Plate |
| N/A | 212-0032 | "H" Plate |
| N/A | 212-0034 | 18 Hole Straight Plate |
| N/A | 212-0036 | 12 Hole Curved Plate |
| N/A | 212-0038 | 5 Hole Straight Plate |
| N/A | 212-0039 | 10 Hole Straight Plate |
| N/A | 212-0040 | 2 x 4 Hole Quad Plate |
| N/A | 212-0041 | Short 2 x 4 Hole Quad Plate |
| N/A | 212-0042 | 37 x 50 mm Mesh Plate |
| N/A | 212-0043 | 74 x 100 mm Mesh Plate |
| N/A | 212-0044 | 13 mm Left "L" Plate |
| N/A | 212-0046 | 13 mm Right "L" Plate |
| N/A | 212-0048 | 16 mm Left "L" Plate |
| N/A | 212-0050 | 16 mm Right "L" Plate |
| 212-0052-SP | 212-0052 | 19 mm "T" Plate |
| N/A | 212-0062 | Large Burr Hole Plate |
| N/A | 212-0064 | Small Burr Hole Plate |
| N/A | 212-0068 | "X" Plate |
| 212-0069-SP | 212-0069 | Medium 2 x 4 Hole Quad Plate |
| N/A | 212-0080 | Orbital Plate |
| N/A | 212-0086 | Orbital Floor Plate, Large |
| N/A | 212-1030 | Micro 4 Hole Straight Plate |
| N/A | 212-1031 | Micro 6 Hole Straight Plate |

Craniomaxillofacial (CMF) Portfolio [continued]

OsteoMed MFx™ Mandibular Fixation System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|--------------------------------------|
| N/A | 212-1032 | Micro 8 Hole Straight Plate |
| N/A | 212-1040 | Micro 6 Hole Double "Y" Plate |
| N/A | 213-0000 | 13 mm Right "L" Plate |
| N/A | 213-0001 | 13 mm Left "L" Plate |
| N/A | 213-0002 | 15 mm Right "L" Plate |
| N/A | 213-0003 | 15 mm Left "L" Plate |
| N/A | 213-0004 | 22 mm Right "L" Plate |
| N/A | 213-0005 | 22 mm Left "L" Plate |
| N/A | 213-0006 | 9 mm Right ZED Plate |
| N/A | 213-0007 | 9 mm Left ZED Plate |
| N/A | 213-0008 | 11 mm Right ZED Plate |
| N/A | 213-0009 | 11 mm Left ZED Plate |
| N/A | 213-0010 | 15 mm 5 Hole "Y" Plate |
| 213-0011-SP | 213-0011 | 10 Hole Curved Plate |
| N/A | 213-0013 | 4 Hole "I" Plate |
| N/A | 213-0014 | 18 mm 4 Hole Straight Plate |
| N/A | 213-0016 | 12 Hole Straight Plate |
| N/A | 213-0017 | 11 Hole Extended "Y" Plate |
| 213-0018-SP | 213-0018 | 1.6 mm Rigid Curved Plate, 10 Hole |
| N/A | 213-0019 | 1.6 mm Rigid Straight Plate, 10 Hole |
| N/A | 213-0028 | Plate, Left "L", 8 Hole, 1.6 mm |
| N/A | 213-0029 | Plate, Right "L", 8 Hole, 1.6 mm |
| N/A | 213-0031 | 1.6 mm, Straight Plate, 18 Hole |
| N/A | 213-0040 | Right "L" Plate, 6 mm, Rigid |
| N/A | 213-0041 | Left "L" Plate, 6 mm, Rigid |
| N/A | 213-0042 | Right "L" Plate, 8 mm, Rigid |
| N/A | 213-0043 | Left "L" Plate, 8 mm, Rigid |
| N/A | 213-0044 | Right "L" Plate, 10 mm, Rigid |
| N/A | 213-0045 | Left "L" Plate, 10 mm, Rigid |
| N/A | 213-0046 | Right "L" Plate, 12 mm, Rigid |

Craniomaxillofacial (CMF) Portfolio [continued]

OsteoMed MFx™ Mandibular Fixation System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|---|
| N/A | 213-0047 | Left "L" Plate, 12 mm, Rigid |
| N/A | 213-0048 | 1.6 mm Pre-bent Plt 2 mm Stp Lf |
| N/A | 213-0049 | 1.6 mm Pre-bent Plt 2 mm Stp Rt |
| N/A | 213-0050 | 1.6 mm Pre-bent Plt 4 mm Stp Lf |
| N/A | 213-0051 | 1.6 mm Pre-bent Plt 4 mm Stp Rt |
| N/A | 213-0052 | 1.6 mm Pre-bent Plt 6 mm Stp Lf |
| N/A | 213-0053 | 1.6 mm Pre-bent Plt 6 mm Stp Rt |
| N/A | 213-0054 | 1.6 mm Pre-bent Plt 8 mm Stp Lf |
| N/A | 213-0055 | 1.6 mm Pre-bent Plt 8 mm Stp Rt |
| N/A | 213-0056 | 1.6 mm Pre-bent Plt 10 mm Stp Lf |
| N/A | 213-0057 | 1.6 mm Pre-bent Plt 10 mm Stp Rt |
| N/A | 213-0058 | 1.6 mm Pre-bent Plt 12 mm Stp Lf |
| N/A | 213-0059 | 1.6 mm Pre-bent Plt 12 mm Stp Rt |
| N/A | 213-0066 | 1.6 mm 6 Hole L Plate, Right |
| N/A | 213-0067 | 1.6 mm 6 Hole L Plate, Left |
| N/A | 213-0116 | 1.6 mm Inferior Orbital Rim Plate, Left |
| N/A | 213-0117 | 1.6 mm Inferior Orbital Rim Plate, Right |
| N/A | 213-0120 | 1.6 mm, Right "L" Plate, 8 mm, Rigid |
| N/A | 213-0121 | 1.6 mm, Left "L" Plate, 8 mm, Rigid |
| N/A | 213-0122 | 1.6 mm Right "L" Plate 10 mm, Rigid |
| N/A | 213-0123 | 1.6 mm Left "L" Plate 10 mm, Rigid |
| N/A | 213-0124 | 1.6 mm Right "L" Plate 12 mm, Rigid |
| N/A | 213-0125 | 1.6 mm Left "L" Plate 12 mm, Rigid |
| N/A | 213-0126 | 1.6 mm, Right "Zed" Plate, 11 mm |
| N/A | 213-0127 | 1.6 mm, Left "Zed" Plate, 11 mm |
| N/A | 213-0214 | 1.6 mm Modified LeFort Plate, Left, 0 mm |
| N/A | 213-0215 | 1.6 mm Modified LeFort Plate, Right, 0 mm |
| N/A | 213-1600 | Low Profile "Double Y" Plate |
| N/A | 213-1601 | Low Profile "T" Plate |
| N/A | 213-1602 | Low Profile "Y" Plate |

Craniomaxillofacial (CMF) Portfolio [continued]

OsteoMed MFX™ Mandibular Fixation System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|--|
| N/A | 213-1603 | Low Profile 4 Hole Straight Plate |
| N/A | 213-1604 | Low Profile "L" Plate Left |
| N/A | 213-1605 | Low Profile "L" Plate Right |
| N/A | 213-1606 | Low Profile 8 Hole Straight Plate |
| N/A | 213-1607 | Low Profile 2 x 2 Hole Large Plate |
| N/A | 213-1608 | Low Profile 5 x 2 Hole Large Plate |
| N/A | 213-1609 | Low Profile 4 Hole Arc Plate |
| N/A | 213-1610 | Low Profile 8 Hole Arc Plate |
| N/A | 213-1611 | Low Profile 24 Hole Straight Plate |
| N/A | 213-1612 | Low Profile 93 mm x 93 mm x 0.5 mm Mesh |
| N/A | 213-1613 | Low Profile 93 mm x 93 mm x 0.6 mm Mesh |
| N/A | 213-1614 | Low Profile 203 mm x 203 mm x 0.5 mm Mesh |
| N/A | 213-1615 | Low Profile 2 x 2 Hole Small Plate |
| N/A | 213-1616 | Low Profile 5 x 2 Hole Small Plate |
| N/A | 213-1620 | LPM Mesh, 93 mm x 93 mm x 0.6 mm, Rigid |
| N/A | 213-1621 | LPM Mesh, 203 mm x 203 mm x 0.6 mm, Rigid |
| N/A | 213-1622 | LPM Mesh, 187 mm x 161 mm x 0.6 mm, Rigid |
| N/A | 213-1623 | LPM Mesh, 187 mm x 161 mm x 0.6 mm |
| N/A | 213-1624 | LPM Mesh, 203 mm x 203 mm x 0.6 mm |
| N/A | 213-1634 | 2.0 mm Angled Locking L Plate, 6 x 3 Hole, Left |
| N/A | 213-1635 | 2.0 mm Angled Locking L Plate, 6 x 3 Hole, Right |
| N/A | 213-1640 | 2.0 mm Angled Locking Curved Plate, 8 Hole |
| N/A | 214-0050 | 33 mm 4 Hole Semi-Compression Plate |
| N/A | 214-0051 | 33 mm 4 Hole EDC Plate |
| N/A | 214-0052 | 48 mm 6 Hole EDC Plate |
| N/A | 214-0053 | 38 mm 4 Hole EDC Plate |
| N/A | 214-0055 | Short Angled EDC Plate |
| N/A | 214-0056 | 15 Hole Comminuted Fracture Plate |
| N/A | 214-0057 | Long Angled EDC Plate |
| N/A | 214-0058 | 28 mm 4 Hole Fracture Plate |

Craniomaxillofacial (CMF) Portfolio [continued]

OsteoMed MFx™ Mandibular Fixation System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|---|
| N/A | 214-0059 | 40 mm 6 Hole Fracture Plate |
| N/A | 214-0060 | Right Reconstruction Plate w/Angle |
| N/A | 214-0062 | Left Reconstruction Plate w/Angle |
| N/A | 214-0064 | Straight Reconstruction Plate |
| 214-0080-SP | 214-0080 | Medium 4 Hole BSSO Plate |
| 214-0082-SP | 214-0082 | Long 4 Hole BSSO Plate |
| N/A | 214-0083 | Long 6 Hole BSSO Plate |
| 214-0084-SP | 214-0084 | Short 4 Hole BSSO Plate |
| N/A | 214-0085 | TR BSSO Plate, 10 mm |
| N/A | 214-0086 | Short TR BSSO Plate |
| N/A | 214-0088 | Long TR BSSO Plate |
| 214-0092-SP | 214-0092 | Medium 8 Hole BSSO Plate |
| 214-0094-SP | 214-0094 | Long 8 Hole BSSO Plate |
| N/A | 214-0099 | Short 6 Hole BSSO Plate |
| N/A | 214-0100 | Medium 6 Hole BSSO Plate |
| N/A | 214-0101 | Short Right Reconstruction Plate w/Angle |
| N/A | 214-0102 | Short Left Reconstruction Plate w/Angle |
| N/A | 214-0103 | Short Straight Reconstruction Plate |
| N/A | 214-0104 | 6 Hole Angled EDC Plate |
| N/A | 214-0107 | 6 Hole, Small Curved BSSO Plate |
| N/A | 214-0108 | 6 Hole, Medium Curved BSSO Plate |
| N/A | 214-0109 | 6 Hole, Large Curved BSSO Plate |
| N/A | 214-0110 | 2.0 mm 5 Hole Lozenge Curved Reduction Plate |
| N/A | 214-0111 | 2.0 mm 3 Hole Lozenge Straight Reduction Plate Short |
| N/A | 214-0112 | 2.0 mm 3 Hole Lozenge Straight Reduction Plate Medium |
| N/A | 214-0113 | 2.0 mm 3 Hole Lozenge Straight Reduction Plate Long |
| 214-0114-SP | 214-0114 | 2.0 mm, 6 Hole Offset Curved Lozenge Reduction Plate |
| 214-0115-SP | 214-0115 | 2.0 mm, 7 Hole Offset Curved Lozenge Reduction Plate |
| N/A | 214-0116 | 2.0 mm 5 Hole Offset Curved Lozenge Reducion Plate |
| N/A | 214-0117 | 2.0 mm 3 Hole Lozenge STR PL, SHORT 0.8 mm |

Craniomaxillofacial (CMF) Portfolio [continued]

OsteoMed MFx™ Mandibular Fixation System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|--|
| N/A | 214-0118 | 2.0 mm 3 Hole Lozenge STR, MED, 0.8 mm |
| N/A | 214-0119 | 2.0 mm 3 Hole Lozenge STR, LONG, 0.8 mm |
| N/A | 214-0120 | 2.0 mm 5 Hole Lozenge CURVED PL, 0.8 mm |
| N/A | 214-0150 | 4 Hole Mandibular Fracture Compression Plate |
| N/A | 214-0154 | 6 Hole Curved Plate |
| N/A | 214-0201 | Recon Plate, RT, 4 x 12 Hole |
| N/A | 214-0202 | Recon Plate, LT, 4 x 12 Hole |
| N/A | 214-0203 | Recon Plate, RT, 5 x 16 Hole |
| N/A | 214-0204 | Recon Plate, LT 5 x 16 Hole |
| N/A | 214-0205 | Recon Plate, Full, 5 x 20 x 5 Hole |
| N/A | 214-0206 | Recon Plate, Full, 5 x 16 x 5 Hole |
| N/A | 214-0207 | Recon Plate, Strt, 14 Hole |
| N/A | 214-0208 | Recon Plate, Strt, 18 Hole |
| N/A | 214-0210 | Frac Plate, 4 Hole DCP, 32 mm |
| N/A | 214-0211 | Frac Plate, 4 Hole DCP 34 mm |
| N/A | 214-0212 | Frac Plate, 4 Hole DCP 37 mm |
| N/A | 214-0213 | Frac Plate, 6 Hole DCP 47 mm |
| N/A | 214-0214 | Frac Plate, 6 Hole DCP 49 mm |
| N/A | 214-0215 | Frac Plate, 6 Hole DCP 52 mm |
| N/A | 214-0216 | Frac Plate, 8 Hole DCP 59 mm |
| N/A | 214-0217 | Frac Plate, 4 Hole EDCP 32 mm |
| N/A | 214-0218 | Frac Plate, 4 Hole EDCP 34 mm |
| N/A | 214-0219 | 2 x 2 Hole Angled Fracture Plate |
| N/A | 214-0220 | 3 x 3 Hole Angled Fracture Plate |
| N/A | 214-0222 | Rigid Comminuted Fracture Plate, Locking |
| N/A | 214-0231 | MFx 2.0 mm Reduction Plate, 4 Hole |
| N/A | 214-0251 | MFx 2.4 mm Reduction Plate, 4 Hole |
| N/A | 214-0261 | MFx 125° Contoured Reconstruction Plate, 4 x 12, Right |
| N/A | 214-0262 | MFx 125° Contoured Reconstruction Plate, 4 x 12, Left |
| N/A | 214-0263 | MFx 125° Contoured Reconstruction Plate, 5 x 16, Right |

Craniomaxillofacial (CMF) Portfolio [continued]

OsteoMed MFX™ Mandibular Fixation System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|---|
| N/A | 214-0264 | MFX 125° Contoured Reconstruction Plate, 5 x 16, Left |
| N/A | 214-0266 | MFX 125° Contoured Reconstruction Plate, 5 x 16 x 5, Full |
| N/A | 214-0311 | 2.4 Locking Fracture Plate, 4 Hole, 34 mm |
| N/A | 214-0315 | 2.4 Locking Fracture Plate, 6 Hole, 52 mm |
| N/A | 214-0400 | 2.0 mm Locking Fracture Plate, 4 Hole, Short |
| N/A | 214-0401 | 2.0 mm Locking Fracture Plate, 4 Hole, Medium |
| N/A | 214-0402 | 2.0 mm Locking Fracture Plate, 6 Hole |
| N/A | 214-0403 | 2.0 mm Locking Plate, Comminution |
| N/A | 214-1080 | Ti Alloy Medium 4 Hole BSSO Plate |
| N/A | 214-1082 | Ti Alloy Long 4 Hole BSSO Plate |
| N/A | 214-1084 | Ti Alloy Short 4 Hole BSSO Plate |
| N/A | 214-2001 | MFX Angulated Locking BSSO Plate, 4 Hole, Medium |
| N/A | 214-2002 | MFX Angulated Locking BSSO Plate, 4 Hole, Long |
| N/A | 214-2003 | MFX Angulated Locking BSSO Plate, 4 Hole, Short |
| N/A | 214-2004 | MFX Angulated Locking Comminution Plate |
| N/A | 214-2005 | MFX Angulated Locking Straight Plate, 18 Hole |
| N/A | 214-2006 | MFX Angulated Locking Straight Plate, 4 Hole |
| N/A | 214-2007 | MFX Angulated Locking Curved Plate, 9 Hole |
| N/A | 214-2008 | MFX Angulated Locking Curved Plate, 7 Hole |
| N/A | 214-2009 | MFX Angulated Locking Straight Plate, 6 Hole |
| N/A | 214-2010 | MFX Angulated Locking Curved Plate, 5 Hole |
| N/A | 214-2011 | MFX Angulated Locking BSSO Plate, 6 Hole, Short |
| N/A | 214-2012 | MFX Angulated Locking BSSO Plate, 6 Hole, Medium |
| N/A | 214-2013 | MFX Angulated Locking Curved Strut Plate, 2 x 4 Hole |
| N/A | 214-2014 | MFX Angulated Locking Curved Strut Plate, 2 x 8 Hole |
| N/A | 214-2017 | MFX Angulated Locking Angle Plate, 2 x 2 Hole |
| N/A | 214-2018 | MFX Angulated Locking Angle Plate, 3 x 3 Hole |
| N/A | 214-2019 | 2.0 mm Champy Plate, Short, Right |
| N/A | 214-2020 | 2.0 mm Champy Plate, Medium, Right |
| N/A | 214-2021 | 2.0 mm Champy Plate, Long, Right |

Craniomaxillofacial (CMF) Portfolio [continued]

OsteoMed MFX™ Mandibular Fixation System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|--|
| N/A | 214-2022 | 2.0 mm Champy Plate, Short, Left |
| N/A | 214-2023 | 2.0 mm Champy Plate, Medium, Left |
| N/A | 214-2024 | 2.0 mm Champy Plate, Long, Left |
| N/A | 214-2025 | 2.0 mm, Locking Reconstruction Plate, 18 Hole, Straight |
| N/A | 214-2026 | 2.0 mm, Locking Reconstruction Plate, 5 x 16 Hole, Right |
| N/A | 214-2027 | 2.0 mm, Locking Reconstruction Plate, 5 x 16 Hole, Left |
| N/A | 214-2028 | 2.0 mm, Locking Reconstruction Plate, 5 x 20 x 5 Hole |
| N/A | 214-2030 | 2.0 mm Medium 4 Hole Locking Plate, 2.0 mm Thick |
| N/A | 214-2031 | 2.0 mm Long 4 Hole Locking Plate, 2.0 mm Thick |
| N/A | 214-2032 | 2.0 mm Short 4 Hole Locking Plate, 2.0 mm Thick |
| N/A | 214-2033 | 2.0 mm Comminution Locking Plate, 2.0 mm Thick |
| N/A | 214-2034 | 2.0 mm 18 Hole Straight Locking Plate, 2.0 mm Thick |
| N/A | 214-2035 | 2.0 mm 4 Hole Straight Locking Plate, 2.0 mm Thick |
| N/A | 214-2036 | 2.0 mm 9 Hole Curved Locking Plate, 2.0 mm Thick |
| N/A | 214-2037 | 2.0 mm 7 Hole Curved Locking Plate, 2.0 mm Thick |
| N/A | 214-2038 | 2.0 mm 6 Hole Straight Locking Plate, 2.0 mm Thick |
| N/A | 214-2039 | 2.0 mm 5 Hole Curved Locking Plate, 2.0 mm Thick |
| N/A | 214-2040 | 2.0 mm 6 Hole Short Locking Plate, 2.0 mm Thick |
| N/A | 214-2041 | 2.0 mm 6 Hole Medium Locking Plate, 2.0 mm Thick |
| N/A | 214-2042 | 2.0 mm 2 x 4 Strut, Locking Curved Plate, 2.0 mm Thick |
| N/A | 214-2043 | 2.0 mm 2 x 8 Strut, Locking Curved Plate, 2.0 mm Thick |
| N/A | 214-2044 | 2.0 mm Locking 2 x 2 Angle Plate, 2.0 mm Thick |
| N/A | 214-2045 | 2.0 mm Locking 3 x 3 Angle Plate, 2.0 mm Thick |
| N/A | 214-2046 | 2.0 mm Mandible Positional Plate |
| N/A | 214-2050 | MFx Angulated Locking 4 Hole Condyle Trauma Plate |
| N/A | 214-2051 | MFx Angulated Locking 5 Hole Condyle Trauma Plate |
| N/A | 214-2401 | MFx Angulated Locking Fracture Plate, 4 Hole, 32 mm |
| N/A | 214-2402 | MFx Angulated Locking Fracture Plate, 4 Hole, 34 mm |
| N/A | 214-2403 | MFx Angulated Locking Fracture Plate, 4 Hole, 37 mm |
| N/A | 214-2404 | MFx Angulated Locking Fracture Plate, 6 Hole, 47 mm |

Craniomaxillofacial (CMF) Portfolio [continued]

OsteoMed MFX™ Mandibular Fixation System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|---|
| N/A | 214-2405 | MFX Angulated Locking Fracture Plate, 6 Hole, 49 mm |
| N/A | 214-2406 | MFX Angulated Locking Fracture Plate, 6 Hole, 52 mm |
| N/A | 214-2407 | MFX Angulated Locking Fracture Plate, 8 Hole, 59 mm |
| N/A | 214-2408 | MFX Angulated Locking Fracture Plate, 2 x 2 Hole |
| N/A | 214-2409 | MFX Angulated Locking Fracture Plate, 3 x 3 Hole |
| N/A | 214-2410 | MFX Angulated Locking Comminuted Fracture Plate, Rigid |
| N/A | 214-2413 | MFX 125° Angulated Locking Contoured Reconstruction Plate, 4 x 12 Hole, Right |
| N/A | 214-2414 | MFX 125° Angulated Locking Contoured Reconstruction Plate, 4 x 12 Hole, Left |
| N/A | 214-2415 | MFX 125° Angulated Locking Contoured Reconstruction Plate, Left 5 x 20 x 5 Hole |
| N/A | 214-2416 | MFX 125° Angulated Locking Contoured Reconstruction Plate, Left 5 x 16 x 5 Hole |
| N/A | 214-2417 | MFX Angulated Locking Contoured Reconstruction Plate, 14 Hole |
| N/A | 214-2418 | MFX Angulated Locking Contoured Reconstruction Plate, 18 Hole |
| N/A | 214-2419 | MFX 125° Angulated Locking Contoured Reconstruction Plate, 5 x 16 Hole, Right |
| N/A | 214-2420 | MFX 125° Angulated Locking Contoured Reconstruction Plate, 5 x 16 Hole, Left |
| N/A | 216-1604 | 1.6 mm X 4 mm Screw |
| N/A | 216-1605 | 1.6 mm X 5 mm Screw |
| N/A | 216-1606 | 1.6 mm X 6 mm Screw |
| N/A | 216-1608 | 1.6 mm X 8 mm Screw |
| N/A | 216-1610 | 1.6 mm X 10 mm Screw |
| N/A | 216-1612 | 1.6 mm X 12 mm Screw |
| N/A | 216-1614 | 1.6 mm X 14 mm Screw |
| N/A | 218-0028 | Double "Y" Plate |
| N/A | 218-0038 | 10 Hole Straight Plate |
| N/A | 218-1603 | 1.6 mm x 3.5 mm Profile Zero Auto-Drive Screw |
| N/A | 218-1903 | 1.9 mm x 3 mm Profile Zero Safety Screw |
| N/A | 223-2004 | MFX 2.0 mm x 4 mm Angulated Locking Auto-Drive® Screw |
| N/A | 223-2005 | MFX 2.0 mm x 5 mm Angulated Locking Auto-Drive® Screw |

Craniomaxillofacial (CMF) Portfolio [continued]

OsteoMed MFx™ Mandibular Fixation System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|---|
| N/A | 223-2006 | MFx 2.0 mm x 6 mm Angulated Locking Auto-Drive® Screw |
| N/A | 223-2008 | MFx 2.0 mm x 8 mm Angulated Locking Auto-Drive® Screw |
| N/A | 223-2406 | MFx 2.4 mm x 6 mm Angulated Locking Auto-Drive® Screw |
| N/A | 223-2408 | MFx 2.4 mm x 8 mm Angulated Locking Auto-Drive® Screw |
| N/A | 225-2004 | MFx 2.0 mm x 4 mm Angulated Locking Standard Screw |
| N/A | 225-2005 | MFx 2.0 mm x 5 mm Angulated Locking Standard Screw |
| N/A | 225-2006 | MFx 2.0 mm x 6 mm Angulated Locking Standard Screw |
| N/A | 225-2008 | MFx 2.0 mm x 8 mm Angulated Locking Standard Screw |
| N/A | 225-2010 | MFx 2.0 mm x 10 mm Angulated Locking Standard Screw |
| N/A | 225-2012 | MFx 2.0 mm x 12 mm Angulated Locking Standard Screw |
| N/A | 225-2014 | MFx 2.0 mm x 14 mm Angulated Locking Standard Screw |
| N/A | 225-2016 | MFx 2.0 mm x 16 mm Angulated Locking Standard Screw |
| N/A | 225-2018 | MFx 2.0 mm x 18 mm Angulated Locking Standard Screw |
| N/A | 225-2406 | MFx 2.4 mm x 6 mm Angulated Locking Standard Screw |
| N/A | 225-2408 | MFx 2.4 mm x 8 mm Angulated Locking Standard Screw |
| N/A | 225-2410 | MFx 2.4 mm x 10 mm Angulated Locking Standard Screw |
| N/A | 225-2412 | MFx 2.4 mm x 12 mm Angulated Locking Standard Screw |
| N/A | 225-2414 | MFx 2.4 mm x 14 mm Angulated Locking Standard Screw |
| N/A | 225-2416 | MFx 2.4 mm x 16 mm Angulated Locking Standard Screw |
| N/A | 225-2418 | MFx 2.4 mm x 18 mm Angulated Locking Standard Screw |
| N/A | 225-2420 | MFx 2.4 mm x 20 mm Angulated Locking Standard Screw |
| N/A | 225-2422 | MFx 2.4 mm x 22 mm Angulated Locking Standard Screw |
| N/A | 225-2706 | MFx 2.7 mm x 6 mm Angulated Locking Safety Screw |
| N/A | 225-2708 | MFx 2.7 mm x 8 mm Angulated Locking Safety Screw |
| N/A | 225-2710 | MFx 2.7 mm x 10 mm Angulated Locking Safety Screw |
| N/A | 225-2712 | MFx 2.7 mm x 12 mm Angulated Locking Safety Screw |
| N/A | 225-2714 | MFx 2.7 mm x 14 mm Angulated Locking Safety Screw |
| N/A | 225-2716 | MFx 2.7 mm x 16 mm Angulated Locking Safety Screw |
| N/A | 225-2718 | MFx 2.7 mm x 18 mm Angulated Locking Safety Screw |

Craniomaxillofacial (CMF) Portfolio [continued]

OsteoMed ICON™ Facial Plating System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|---|
| N/A | 212-1661 | 1.2 mm Orbital Floor Plate, Fan, .2 mm Thick |
| N/A | 212-1662 | 1.2 mm Orbital Floor Plate, Fan, .3 mm Thick |
| N/A | 212-1663 | 1.2 mm Orbital Floor Plate, Fan, .4 mm Thick |
| N/A | 212-1664 | 1.2 mm Orbital Floor Plate, Fan, .5 mm Thick |
| N/A | 212-1671 | 1.2 mm Orbital Floor Plate, Pick, .2 mm Thick |
| N/A | 212-1672 | 1.2 mm Orbital Floor Plate, Pick, .3 mm Thick |
| N/A | 212-1673 | 1.2 mm Orbital Floor Plate, Pick, .4 mm Thick |
| N/A | 212-1674 | 1.2 mm Orbital Floor Plate, Pick, .5 mm Thick |
| N/A | 212-1681 | 1.2 mm Orbital Floor Plate, Preformed, Left, .2 mm Thick |
| N/A | 212-1682 | 1.2 mm Orbital Floor Plate, Preformed, Left, .3 mm Thick |
| N/A | 212-1683 | 1.2 mm Orbital Floor Plate, Preformed, Left, .4 mm Thick |
| N/A | 212-1684 | 1.2 mm Orbital Floor Plate, Preformed, Left, .5 mm Thick |
| N/A | 212-1691 | 1.2 mm Orbital Floor Plate, Preformed, Right, .2 mm Thick |
| N/A | 212-1692 | 1.2 mm Orbital Floor Plate, Preformed, Right, .3 mm Thick |
| N/A | 212-1693 | 1.2 mm Orbital Floor Plate, Preformed, Right, .4 mm Thick |
| N/A | 212-1694 | 1.2 mm Orbital Floor Plate, Preformed, Right, .5 mm Thick |
| N/A | 213-1661 | 1.6 mm Orbital Floor Plate, Fan, .2 mm Thick |
| N/A | 213-1663 | 1.6 mm Orbital Floor Plate, Fan, .3 mm Thick |
| N/A | 213-1665 | 1.6 mm Orbital Floor Plate, Fan, .4 mm Thick |
| N/A | 213-1667 | 1.6 mm Orbital Floor Plate, Fan, .5 mm Thick |
| N/A | 213-1671 | 1.6 mm Orbital Floor Plate, Pick, .2 mm Thick |
| N/A | 213-1673 | 1.6 mm Orbital Floor Plate, Pick, .3 mm Thick |
| N/A | 213-1675 | 1.6 mm Orbital Floor Plate, Pick, .4 mm Thick |
| N/A | 213-1677 | 1.6 mm Orbital Floor Plate, Pick, .5 mm Thick |
| N/A | 213-1681 | 1.6 mm Orbital Floor Plate, Preformed, Left, .2 mm Thick |
| N/A | 213-1683 | 1.6 mm Orbital Floor Plate, Preformed, Left, .3 mm Thick |
| N/A | 213-1685 | 1.6 mm Orbital Floor Plate, Preformed, Left, .4 mm Thick |
| N/A | 213-1687 | 1.6 mm Orbital Floor Plate, Preformed, Left, .5 mm Thick |
| N/A | 213-1691 | 1.6 mm Orbital Floor Plate, Preformed, Right, .2 mm Thick |
| N/A | 213-1693 | 1.6 mm Orbital Floor Plate, Preformed, Right, .3 mm Thick |

Craniomaxillofacial (CMF) Portfolio [continued]

OsteoMed ICON™ Facial Plating System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|---|
| N/A | 213-1695 | 1.6 mm Orbital Floor Plate, Preformed, Right, .4 mm Thick |
| N/A | 213-1697 | 1.6 mm Orbital Floor Plate, Preformed, Right, .5 mm Thick |


2.0 mm Auto Drive MMF Screws

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|-------------------------------------|
| N/A | 209-2008 | 2.0 mm x 8 mm MMF Auto-Drive Screw |
| N/A | 209-2011 | 2.0 mm x 11 mm MMF Auto-Drive Screw |
| N/A | 209-2014 | 2.0 mm x 14 mm MMF Auto-Drive Screw |

Neurosurgical Portfolio

- ▶ OsteoMed Fast Flap™ Neurosurgical Fixation System
- ▶ OsteoMed Profile Zero™ Low Profile Neuro Plating System

MRI Restricted Zone Summary

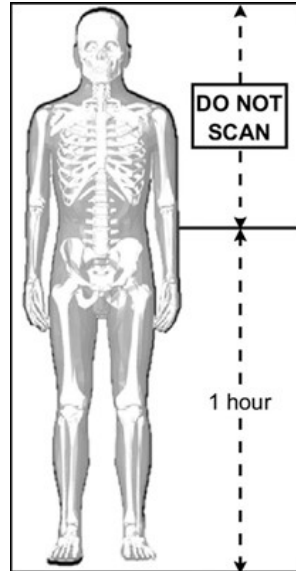

Conditions for Safe Scanning in the MRI

A patient with the screw and plate of OsteoMed Fast Flap Neurosurgical Fixation System, OsteoMed Profile Zero Low Profile Neuro Plating System or Acumed Cranial Mesh may be safely scanned at 1.5T or 3T under the following conditions. Failure to follow these conditions may result in injury to the patient. See below for a complete list of part numbers included in these systems.

OsteoMed's Neurosurgical Products

| Parameter | Condition |
|--|---|
| Device Name | OsteoMed's Neuro Products |
| Static Magnetic Field Strength (B0) | 1.5T and 3T |
| MR Scanner Type | Cylindrical |
| B0 Field Orientation | Horizontal |
| Maximum Spatial Field Gradient | 32 T/m (3,200 G/cm) |
| Maximum Gradient Slew Rate | 200 T/m/s per axis |
| RF Excitation | Circularly Polarized (CP) |
| RF Transmit Coil Type | Integrated Whole Body Transmit Coil |
| Operating Mode | Normal Operating Mode |
| Scan Duration | Scan duration up to 1 hour without a cooling period when landmarked at or below the navel (as detailed in the MRI Restricted Zone Summary) |
| Scan Regions | Landmark locations at or below the navel are permissible (as detailed in the MRI Restricted Zone Summary) |
| Image Artifact | The presence of OsteoMed's Neuro Products may produce an image artifact of 3.0 cm. Some manipulation of scan parameters may be needed to compensate for the artifact. |
| Parameters may be needed to compensate for the artifact. | |

Neurosurgical Portfolio [continued]

OsteoMed Fast Flap™ Neurosurgical Fixation System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|-----------------------------------|
| N/A | 204-1503 | 1.5 mm x 3 mm Safety Screw |
| N/A | 208-1904 | 1.9 mm x 4 mm Safety Screw |
| N/A | 211-1203 | 1.2 mm x 3 mm Auto-drive Screw |
| N/A | 211-1204 | 1.2 mm x 4 mm Auto-drive Screw |
| N/A | 211-1205 | 1.2 mm x 5 mm Auto-drive Screw |
| N/A | 211-1603 | 1.6 mm x 3.5 mm Auto-Drive Screw |
| N/A | 211-1604 | 1.6 mm x 4 mm Auto-Drive Screw |
| N/A | 211-1605 | 1.6 mm x 5 mm Auto-Drive Screw |
| N/A | 211-1606 | 1.6 mm x 6 mm Auto-Drive Screw |
| N/A | 211-2006 | 2.0 mm x 6 mm Auto-Drive Screw |
| N/A | 212-0034 | 18 Hole Straight Plate |
| N/A | 212-0036 | 12 Hole Curved Plate |
| N/A | 212-0038 | 5 Hole Straight Plate |
| N/A | 212-0040 | 2 x 4 Hole Quad Plate |
| N/A | 212-0130 | Medium 2 Hole Straight Plate |
| N/A | 212-0132 | Long 2 Hole Straight Plate |
| N/A | 212-0134 | Square Plate |
| N/A | 212-0136 | Extended "X" Plate |
| N/A | 212-0138 | Extended Double "Y" Plate |
| N/A | 212-0140 | Small Burr Hole Plate |
| N/A | 212-0142 | Large Burr Hole Plate |
| N/A | 212-0144 | X-Large Burr Hole Plate |
| N/A | 212-0146 | 40 x 40 mm Osteoform Mesh, 1.2 mm |
| N/A | 212-0148 | 90 x 90 mm Osteoform Mesh, 1.2 mm |
| N/A | 212-0150 | Micro Mesh, 74 mm x 50 mm |
| 218-0020-SP | 218-0020 | 9 mm Square Plate |
| N/A | 218-0022 | Short 2 Hole Straight Plate |
| 218-0024-SP | 218-0024 | Medium 2 Hole Straight Plate |
| 218-0026-SP | 218-0026 | Long 2 Hole Straight Plate |
| N/A | 218-0028 | Double "Y" Plate |

Neurosurgical Portfolio [continued]

OsteoMed Fast Flap™ Neurosurgical Fixation System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|---|
| N/A | 218-0030 | 5 Hole Straight Plate |
| 218-0032-SP | 218-0032 | Small Burr Hole Plate |
| N/A | 218-0033 | Small Burr Hole Plate, Long Arm |
| N/A | 218-0034 | Large Burr Hole Plate |
| N/A | 218-0036 | 36 x 50 mm Mesh Plate |
| N/A | 218-0038 | 10 Hole Straight Plate |
| N/A | 218-0039 | 74 x 50 mm Mesh Plate |
| N/A | 218-0040 | 111 x 191 mm Mesh Plate |
| N/A | 218-0041 | 150 x 300 mm Mesh Plate |
| N/A | 218-0042 | 25 mm Gap Plate |
| N/A | 218-0044 | 35 mm Gap Plate |
| N/A | 218-0046 | 45 mm Gap Plate |
| N/A | 218-0048 | 21 x 33 mm Sub-Temporal Plate |
| N/A | 218-0050 | 30 x 48 mm Sub-Temporal Plate |
| N/A | 218-0052 | 74 x 50 mm Mesh Plate, Rigid |
| N/A | 218-0053 | 125 mm x 180 mm Mesh Plate, Rigid |
| 218-0059-SP | 218-0059 | 2 x 3 Hole Quad Plate |
| 218-0060-SP | 218-0060 | 2 x 4 Hole Quad Plate |
| N/A | 218-0061 | Bur Hole Shunt Plate |
| N/A | 218-0062 | X-Large Burr Hole Plate |
| N/A | 218-0064 | Burr Hole Shunt Plate, Extended Hole |
| N/A | 218-0065 | X Plate, X-Small |
| N/A | 218-0066 | Burr Hole Plate, Medium |
| N/A | 218-0067 | Burr Hole Plate, Small |
| N/A | 218-0069 | Burr Hole Groove Shunt Plate |
| N/A | 218-0080 | 9 mm Square Plate, Rigid |
| N/A | 218-0081 | Medium 2 Hole Straight Plate, Rigid |
| N/A | 218-0082 | Long 2 Hole Straight Plate, Rigid |
| N/A | 218-0090 | 1.6 mm, Low Profile 9 mm Square Plate |
| N/A | 218-0091 | 1.6 mm, Low Profile Short 2 Hole Straight Plate |

Neurosurgical Portfolio [continued]

OsteoMed Fast Flap™ Neurosurgical Fixation System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|---|
| N/A | 218-0092 | 1.6 mm, Low Profile Medium 2 Hole Straight Plate |
| N/A | 218-0093 | 1.6 mm, Low Profile Long 2 Hole Straight Plate |
| N/A | 218-0094 | 1.6 mm, Low Profile Small Bur Hole Plate |
| N/A | 218-0095 | 1.6 mm, Low Profile Large Bur Hole Plate |
| N/A | 218-0096 | 1.6 mm, Low Profile X-Large Bur Hole Plate |
| N/A | 218-0097 | Low Profile Cover Shunt Plate with Groove |
| N/A | 218-1603 | 1.6 mm x 3.5 mm Profile Zero Auto-Drive Screw |
| N/A | 218-1903 | 1.9 mm x 3 mm Profile Zero Safety Screw |
| 218-2001-SP | 218-2001 | 40 mm x 40 mm Malleable Osteoform Mesh Plate |
| 218-2002-SP | 218-2002 | 40 mm x 40 mm Rigid Osteoform Mesh Plate |
| 218-2003-SP | 218-2003 | 90 mm x 90 mm Malleable Osteoform Mesh Plate |
| 218-2004-SP | 218-2004 | 90 mm x 90 mm Rigid Osteoform Mesh Plate |
| 218-2005-SP | 218-2005 | 165 mm x 95 mm Pre-Formed Osteoform Mesh Plate |
| 218-2006-SP | 218-2006 | 90 mm x 90 mm Pre-Formed Osteoform Mesh Plate |
| N/A | 218-2007 | 165 mm x 95 mm Rigid Osteoform Mesh Plate |
| N/A | 218-2008 | 165 mm x 95 mm Pre-Formed Osteoform Mesh, Extra-Rigid |
| N/A | 218-2014 | 165 x 95 mm OsteoForm Mesh, Extra-Rigid |
| N/A | 218-4500 | 135 mm X 135 mm Rigid Osteoform Mesh |
| N/A | 218-4501 | 157 mm X 157 mm Rigid Osteoform Mesh |
| N/A | 218-4502 | 203 mm X 203 mm Rigid Osteoform Mesh |

Neurosurgical Portfolio [continued]

OsteoMed Profile Zero™ Low Profile Neuro Plating System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|--|
| N/A | 218-0201 | Profile Zero 2 Hole Plate |
| N/A | 218-0202 | Profile Zero 4 Hole Straight Plate, Long |
| N/A | 218-0203 | Profile Zero 4 Hole Straight Plate |
| N/A | 218-0204 | Profile Zero 16 Hole Straight Plate |
| N/A | 218-0205 | Profile Zero Y Plate |
| N/A | 218-0206 | Profile Zero Double Y Plate |
| N/A | 218-0207 | Profile Zero 9 mm Square Plate |
| N/A | 218-0208 | Profile Zero 2 X 3 Hole Quad Plate |
| N/A | 218-0209 | Profile Zero 2 X 4 Hole Quad Plate |
| N/A | 218-0210 | Profile Zero Burr Hole Plate Small |
| N/A | 218-0211 | Profile Zero Large Burr Hole Plate |
| N/A | 218-0212 | Profile Zero X-Large Burr Hole Plate |
| N/A | 218-0213 | Profile Zero Shunt Plate, Groove |
| N/A | 218-0214 | Profile Zero 25 mm Gap Plate |
| N/A | 218-0215 | Profile Zero 35 mm Gap Plate |
| N/A | 218-0216 | Profile Zero 45 mm Gap Plate |
| N/A | 218-0217 | Profile Zero 2 Hole Plate, Long |
| N/A | 218-0218 | Profile Zero Medium Burr Hole Plate |
| N/A | 218-0219 | Profile Zero, 25 mm Gap Plate, Preformed |
| N/A | 218-0220 | Profile Zero, 35 mm Gap Plate, Preformed |
| N/A | 218-0221 | Profile Zero, 45 mm Gap Plate, Preformed |
| N/A | 218-0223 | Profile Zero Temporal Plate |
| N/A | 218-0224 | Profile Zero MVD Plate |
| N/A | 218-0225 | Profile Zero Translabyrinthine Plate |
| N/A | 218-0226 | Profile Zero Retrosigmoid Plate |
| N/A | 218-0227 | Profile Zero Suboccipital Plate |
| N/A | 218-0228 | Profile Zero Shunt Plate, Wide Groove |
| N/A | 218-0229 | Neuro Temporal Plate, Small |
| N/A | 218-0230 | Neuro Arc Mesh Plate, Small |
| N/A | 218-0231 | Neuro Mastoid Plate, Small |


Neurosurgical Portfolio [continued]

OsteoMed Profile Zero™ Low Profile Neuro Plating System

| Sterile Part Number | Non-Sterile Part Number | Description |
|---------------------|-------------------------|---|
| N/A | 218-0232 | Neuro Mastoid Plate, Medium |
| N/A | 218-0233 | Neuro Mastoid Plate, Large |
| 218-0234-SP | 218-0234 | Translabyrinthine Plate, S |
| 218-0235-SP | 218-0235 | Translabyrinthine Plate, L |
| N/A | 218-0236 | Profile Zero , XXL Burr Hole Plate |
| N/A | 218-0285 | Profile Zero 80 mm x 50 mm x .36 mm Mesh |
| N/A | 218-0290 | Profile Zero 65 mm x 90 mm x .36 mm Mesh |
| N/A | 218-1603 | 1.6 mm x 3.5 mm Profile Zero Auto-Drive Screw |
| 218-1604-SP | 218-1604 | 1.6 mm x 4 mm Profile Zero Auto-Drive Screw |
| 218-1605-SP | 218-1605 | 1.6 mm X 5 mm Profile Zero Auto-Drive Screw |
| N/A | 218-1606 | 1.6 mm x 6 mm Profile Zero Auto-Drive Screw |
| 218-1904-SP | 218-1904 | 1.9 mm x 4 mm Low Profile Emergency Screw |
| N/A | 218-1905 | 1.9 mm x 5 mm Low Profile Emergency Screw |
| N/A | 218-2015 | Pterional Mesh Plate, Right |
| N/A | 218-2016 | Pterional Mesh Plate, Left |

OsteoMed® Upper Extremities Portfolio

- ▶ OsteoMed Hand Plating System (hps™)
- ▶ OsteoMed hps™ Hand Fusion System
- ▶ OsteoMed M3-X™ Extremity Fixation System
- ▶ OsteoMed InstaFix™ Shape Memory Fixation System
- ▶ OsteoMed ExtremiLock™ Wrist Plating System

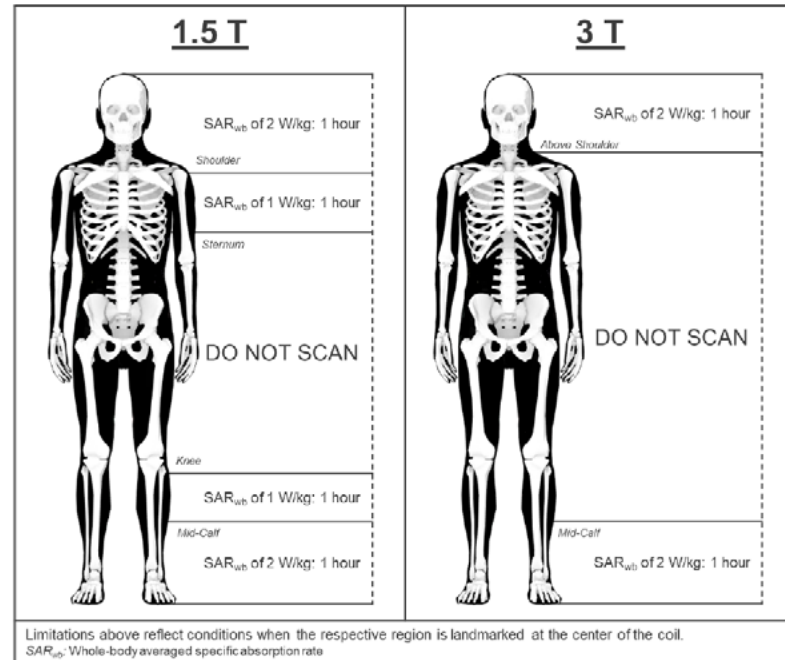

Conditions for Safe Scanning in the MRI

A patient with OsteoMed's Upper Extremity Products, **excluding any stand-alone screws or K-wires**, may be safely scanned at 1.5 T or 3.0 T under the following conditions. Failure to follow these conditions may result in injury.

OsteoMed's Upper Extremity Portfolio

| Parameter | Condition |
|---|---|
| Device Name | OsteoMed's Upper Extremity Products, excluding stand-alone screws and K-wires |
| Nominal value(s) of Static Magnetic Field | 1.5 T or 3 T |
| MR Scanner Type | Cylindrical |
| BO Field Orientation | Horizontal |
| Maximum Spatial Field Gradient | 40 T/m (4,000 G/cm) |
| RF Excitation | Circularly Polarized (CP) |
| RF Transmit Coil Type | Integrated Whole Body Transmit Coil |

OsteoMed® Upper Extremities Portfolio [continued]

OsteoMed's Upper Extremity Portfolio

| Parameter | Condition |
|----------------|---|
| Operating Mode | <p><u>1.5T:</u> Normal Operating Mode from the top of head to the shoulder and from mid-calf to bottom of the foot</p> <p>SARwb≤1W/kg from the shoulder to the sternum and from the knee to mid-calf</p> <p>Do Not Scan between the sternum and the knee</p> |
| | <p><u>3 T:</u> Normal Operating Mode from the top of head to directly above the shoulder and from mid-calf to the bottom of the foot</p> <p>Do Not Scan from the top of shoulder to the mid-calf</p> |
| | <p>Scan duration up to 1 hour without a cooling period with landmark restrictions according to Operating Mode section above and as detailed in the Upper Extremity MRI Restricted Zone Summary</p> |
| Scan Duration | <p>Do not scan between sternum and knee for 1.5 T and from the top of shoulder to the mid-calf in 3 T. All other landmark locations are permissible with restrictions according to Operating Mode section above and as detailed in the Upper Extremity MRI Restricted Zone Summary</p> |
| Scan Regions | <p>The presence of OsteoMed's Upper Extremity products may produce an image artifact of 2.3 cm. Some manipulation of scan parameters may be needed to compensate for the artifact.</p> |
| Image Artifact | |

OsteoMed® Upper Extremities Portfolio [continued]

OsteoMed® Hand Plating System (hps™)

| Non-Sterile Part Number | Descriptor |
|-------------------------|-------------------------------|
| 302-2006 | 2.0 mm X 6 mm Lag Screw |
| 302-2007 | 2.0 mm X 7 mm Lag Screw |
| 302-2009 | 2.0 mm X 9 mm Lag Screw |
| 302-2011 | 2.0 mm X 11 mm Lag Screw |
| 302-2013 | 2.0 mm X 13 mm Lag Screw |
| 302-2024 | 2.0 mm X 24 mm Mini Lag Screw |
| 302-2026 | 2.0 mm X 26 mm Mini Lag Screw |
| 302-2028 | 2.0 mm X 28 mm Mini Lag Screw |
| 302-2030 | 2.0 mm X 30 mm Mini Lag Screw |
| 302-2032 | 2.0 mm X 32 mm Mini Lag Screw |
| 302-2034 | 2.0 mm X 34 mm Mini Lag Screw |
| 302-2036 | 2.0 mm X 36 mm Lag Screw |
| 306-2406 | 2.4 mm X 6 mm Lag Screw |
| 306-2407 | 2.4 mm X 7 mm Lag Screw |
| 306-2409 | 2.4 mm X 9 mm Lag Screw |
| 306-2411 | 2.4 mm X 11 mm Lag Screw |
| 306-2413 | 2.4 mm X 13 mm Lag Screw |
| 306-2434 | 2.4 mm X 34 mm Lag Screw |
| 306-2436 | 2.4 mm X 36 mm Lag Screw |
| 308-1606 | 1.6 mm X 6 mm Lag Screw |
| 308-1607 | 1.6 mm X 7 mm Lag Screw |
| 308-1608 | 1.6 mm X 8 mm Lag Screw |
| 308-1609 | 1.6 mm X 9 mm Lag Screw |
| 308-1610 | 1.6 mm X 10 mm Lag Screw |
| 308-1611 | 1.6 mm X 11 mm Lag Screw |
| 308-1612 | 1.6 mm X 12 mm Lag Screw |
| 308-1613 | 1.6 mm X 13 mm Lag Screw |
| 308-1614 | 1.6 mm X 14 mm Lag Screw |
| 308-1616 | 1.6 mm X 16 mm Lag Screw |
| 308-1618 | 1.6 mm X 18 mm Lag Screw |

OsteoMed® Upper Extremities Portfolio [continued]

OsteoMed® Hand Plating System (hps™)

| Non-Sterile Part Number | Descriptor |
|-------------------------|--|
| 308-1620 | 1.6 mm X 20 mm Lag Screw |
| 308-1622 | 1.6 mm X 22 mm Lag Screw |
| 308-1624 | 1.6 mm X 24 mm Lag Screw |
| 317-2011 | 2.0 mm X 11 mm Headless Screw |
| 317-2013 | 2.0 mm X 13 mm Headless Screw |
| 317-2411 | 2.4 mm X 11 mm Headless Screw |
| 317-2413 | 2.4 mm X 13 mm Headless Screw |
| 317-3011 | 3.0 mm X 11 mm Headless Screw |
| 317-3013 | 3.0 mm X 13 mm Headless Screw |
| 319-2007 | 2.0 mm X 7 mm Cannulated Lag Screw |
| 319-2009 | 2.0 mm X 9 mm Cannulated Lag Screw |
| 319-2011 | 2.0 mm X 11 mm Cannulated Lag Screw |
| 319-2013 | 2.0 mm X 13 mm Cannulated Lag Screw |
| 319-2407 | 2.4 mm X 7 mm Cannulated Lag Screw |
| 319-2409 | 2.4 mm X 9 mm Cannulated Lag Screw |
| 319-2411 | 2.4 mm X 11 mm Cannulated Lag Screw |
| 319-2413 | 2.4 mm X 13 mm Cannulated Lag Screw |
| 330-1606 | 1.6 mm X 6 mm Fully Threaded Screw, Locking |
| 330-1607 | 1.6 mm X 7 mm Fully Threaded Screw, Locking |
| 330-1608 | 1.6 mm X 8 mm Fully Threaded Screw, Locking |
| 330-1609 | 1.6 mm X 9 mm Fully Threaded Screw, Locking |
| 330-1610 | 1.6 mm X 10 mm Fully Threaded Screw, Locking |
| 330-1611 | 1.6 mm X 11 mm Fully Threaded Screw, Locking |
| 330-1612 | 1.6 mm X 12 mm Fully Threaded Screw, Locking |
| 330-1613 | 1.6 mm X 13 mm Fully Threaded Screw, Locking |
| 330-1614 | 1.6 mm X 14 mm Fully Threaded Screw, Locking |
| 330-1616 | 1.6 mm X 16 mm Fully Threaded Screw, Locking |
| 330-1618 | 1.6 mm X 18 mm Fully Threaded Screw, Locking |
| 330-1620 | 1.6 mm X 20 mm Fully Threaded Screw, Locking |
| 330-1622 | 1.6 mm X 22 mm Fully Threaded Screw, Locking |
| 330-1624 | 1.6 mm X 24 mm Fully Threaded Screw, Locking |
| 330-2006 | 2.0 mm X 6 mm Fully Threaded Screw, Locking |
| 330-2007 | 2.0 mm X 7 mm Fully Threaded Screw, Locking |

OsteoMed® Upper Extremities Portfolio [continued]

OsteoMed® Hand Plating System (hps™)

| Non-Sterile Part Number | Descriptor |
|-------------------------|--|
| 330-2008 | 2.0 mm X 8 mm Fully Threaded Screw, Locking |
| 330-2009 | 2.0 mm X 9 mm Fully Threaded Screw, Locking |
| 330-2010 | 2.0 mm X 10 mm Fully Threaded Screw, Locking |
| 330-2011 | 2.0 mm X 11 mm Fully Threaded Screw, Locking |
| 330-2012 | 2.0 mm X 12 mm Fully Threaded Screw, Locking |
| 330-2013 | 2.0 mm X 13 mm Fully Threaded Screw, Locking |
| 330-2014 | 2.0 mm X 14 mm Fully Threaded Screw, Locking |
| 330-2016 | 2.0 mm X 16 mm Fully Threaded Screw, Locking |
| 330-2018 | 2.0 mm X 18 mm Fully Threaded Screw, Locking |
| 330-2020 | 2.0 mm X 20 mm Fully Threaded Screw, Locking |
| 330-2022 | 2.0 mm X 22 mm Fully Threaded Screw, Locking |
| 330-2024 | 2.0 mm X 24 mm Fully Threaded Screw, Locking |
| 330-2026 | 2.0 mm X 26 mm Fully Threaded Screw, Locking |
| 330-2028 | 2.0 mm X 28 mm Fully Threaded Screw, Locking |
| 330-2030 | 2.0 mm X 30 mm Fully Threaded Screw, Locking |
| 330-2032 | 2.0 mm X 32 mm Fully Threaded Screw, Locking |
| 330-2034 | 2.0 mm X 34 mm Fully Threaded Screw, Locking |
| 330-2036 | 2.0 mm X 36 mm Fully Threaded Screw, Locking |
| 330-2406 | 2.4 mm X 6 mm Fully Threaded Screw, Locking |
| 330-2407 | 2.4 mm X 7 mm Fully Threaded Screw, Locking |
| 330-2408 | 2.4 mm X 8 mm Fully Threaded Screw, Locking |
| 330-2409 | 2.4 mm X 9 mm Fully Threaded Screw, Locking |
| 330-2410 | 2.4 mm X 10 mm Fully Threaded Screw, Locking |
| 330-2411 | 2.4 mm X 11 mm Fully Threaded Screw, Locking |
| 330-2412 | 2.4 mm X 12 mm Fully Threaded Screw, Locking |
| 330-2413 | 2.4 mm X 13 mm Fully Threaded Screw, Locking |
| 330-2414 | 2.4 mm X 14 mm Fully Threaded Screw, Locking |
| 330-2416 | 2.4 mm X 16 mm Fully Threaded Screw, Locking |
| 330-2418 | 2.4 mm X 18 mm Fully Threaded Screw, Locking |
| 330-2420 | 2.4 mm X 20 mm Fully Threaded Screw, Locking |
| 330-2422 | 2.4 mm X 22 mm Fully Threaded Screw, Locking |

OsteoMed® Upper Extremities Portfolio [continued]

OsteoMed® Hand Plating System (hps™)

| Non-Sterile Part Number | Descriptor |
|-------------------------|--|
| 330-2424 | 2.4 mm X 24 mm Fully Threaded Screw, Locking |
| 330-2426 | 2.4 mm X 26 mm Fully Threaded Screw, Locking |
| 330-2428 | 2.4 mm X 28 mm Fully Threaded Screw, Locking |
| 330-2430 | 2.4 mm X 30 mm Fully Threaded Screw, Locking |
| 330-2432 | 2.4 mm X 32 mm Fully Threaded Screw, Locking |
| 330-2434 | 2.4 mm X 34 mm Fully Threaded Screw, Locking |
| 330-2436 | 2.4 mm X 36 mm Fully Threaded Screw, Locking |
| 331-1204 | 1.2 mm X 4 mm Fully Threaded Screw |
| 331-1205 | 1.2 mm X 5 mm Fully Threaded Screw |
| 331-1206 | 1.2 mm X 6 mm Fully Threaded Screw |
| 331-1207 | 1.2 mm X 7 mm Fully Threaded Screw |
| 331-1208 | 1.2 mm X 8 mm Fully Threaded Screw |
| 331-1209 | 1.2 mm X 9 mm Fully Threaded Screw |
| 331-1210 | 1.2 mm X 10 mm Fully Threaded Screw |
| 331-1211 | 1.2 mm X 11 mm Fully Threaded Screw |
| 331-1212 | 1.2 mm X 12 mm Fully Threaded Screw |
| 331-1213 | 1.2 mm X 13 mm Fully Threaded Screw |
| 331-1214 | 1.2 mm X 14 mm Fully Threaded Screw |
| 331-1216 | 1.2 mm X 16 mm Fully Threaded Screw |
| 331-1218 | 1.2 mm X 18 mm Fully Threaded Screw |
| 331-1604 | 1.6 mm X 4 mm Fully Threaded Screw |
| 331-1605 | 1.6 mm X 5 mm Fully Threaded Screw |
| 331-1606 | 1.6 mm X 6 mm Fully Threaded Screw |
| 331-1607 | 1.6 mm X 7 mm Fully Threaded Screw |
| 331-1608 | 1.6 mm X 8 mm Fully Threaded Screw |
| 331-1609 | 1.6 mm X 9 mm Fully Threaded Screw |
| 331-1610 | 1.6 mm X 10 mm Fully Threaded Screw |
| 331-1611 | 1.6 mm X 11 mm Fully Threaded Screw |
| 331-1612 | 1.6 mm X 12 mm Fully Threaded Screw |
| 331-1613 | 1.6 mm X 13 mm Fully Threaded Screw |

OsteoMed® Upper Extremities Portfolio [continued]

OsteoMed® Hand Plating System (hps™)

| Non-Sterile Part Number | Descriptor |
|-------------------------|-------------------------------------|
| 331-1614 | 1.6 mm X 14 mm Fully Threaded Screw |
| 331-1616 | 1.6 mm X 16 mm Fully Threaded Screw |
| 331-1618 | 1.6 mm X 18 mm Fully Threaded Screw |
| 331-1620 | 1.6 mm X 20 mm Fully Threaded Screw |
| 331-1622 | 1.6 mm X 22 mm Fully Threaded Screw |
| 331-1624 | 1.6 mm X 24 mm Fully Threaded Screw |
| 331-2006 | 2.0 mm X 6 mm Fully Threaded Screw |
| 331-2007 | 2.0 mm X 7 mm Fully Threaded Screw |
| 331-2008 | 2.0 mm X 8 mm Fully Threaded Screw |
| 331-2009 | 2.0 mm X 9 mm Fully Threaded Screw |
| 331-2010 | 2.0 mm X 10 mm Fully Threaded Screw |
| 331-2011 | 2.0 mm X 11 mm Fully Threaded Screw |
| 331-2012 | 2.0 mm X 12 mm Fully Threaded Screw |
| 331-2013 | 2.0 mm X 13 mm Fully Threaded Screw |
| 331-2014 | 2.0 mm X 14 mm Fully Threaded Screw |
| 331-2016 | 2.0 mm X 16 mm Fully Threaded Screw |
| 331-2018 | 2.0 mm X 18 mm Fully Threaded Screw |
| 331-2020 | 2.0 mm X 20 mm Fully Threaded Screw |
| 331-2022 | 2.0 mm X 22 mm Fully Threaded Screw |
| 331-2024 | 2.0 mm X 24 mm Fully Threaded Screw |
| 331-2026 | 2.0 mm X 26 mm Fully Threaded Screw |
| 331-2028 | 2.0 mm X 28 mm Fully Threaded Screw |
| 331-2030 | 2.0 mm X 30 mm Fully Threaded Screw |
| 331-2032 | 2.0 mm X 32 mm Fully Threaded Screw |
| 331-2034 | 2.0 mm X 34 mm Fully Threaded Screw |
| 331-2036 | 2.0 mm X 36 mm Fully Threaded Screw |
| 331-2406 | 2.4 mm X 6 mm Fully Threaded Screw |
| 331-2407 | 2.4 mm X 7 mm Fully Threaded Screw |
| 331-2408 | 2.4 mm X 8 mm Fully Threaded Screw |
| 331-2409 | 2.4 mm X 9 mm Fully Threaded Screw |

OsteoMed® Upper Extremities Portfolio [continued]

OsteoMed® Hand Plating System (hps™)

| Non-Sterile Part Number | Descriptor |
|-------------------------|-------------------------------------|
| 331-2410 | 2.4 mm X 10 mm Fully Threaded Screw |
| 331-2411 | 2.4 mm X 11 mm Fully Threaded Screw |
| 331-2412 | 2.4 mm X 12 mm Fully Threaded Screw |
| 331-2413 | 2.4 mm X 13 mm Fully Threaded Screw |
| 331-2414 | 2.4 mm X 14 mm Fully Threaded Screw |
| 331-2416 | 2.4 mm X 16 mm Fully Threaded Screw |
| 331-2418 | 2.4 mm X 18 mm Fully Threaded Screw |
| 331-2420 | 2.4 mm X 20 mm Fully Threaded Screw |
| 331-2422 | 2.4 mm X 22 mm Fully Threaded Screw |
| 331-2424 | 2.4 mm X 24 mm Fully Threaded Screw |
| 331-2426 | 2.4 mm X 26 mm Fully Threaded Screw |
| 331-2428 | 2.4 mm X 28 mm Fully Threaded Screw |
| 331-2430 | 2.4 mm X 30 mm Fully Threaded Screw |
| 331-2432 | 2.4 mm X 32 mm Fully Threaded Screw |
| 331-2434 | 2.4 mm X 34 mm Fully Threaded Screw |
| 331-2436 | 2.4 mm X 36 mm Fully Threaded Screw |
| 332-1206 | 1.2 mm X 6 mm Lag Screw |
| 332-1207 | 1.2 mm X 7 mm Lag Screw |
| 332-1208 | 1.2 mm X 8 mm Lag Screw |
| 332-1209 | 1.2 mm X 9 mm Lag Screw |
| 332-1210 | 1.2 mm X 10 mm Lag Screw |
| 332-1211 | 1.2 mm X 11 mm Lag Screw |
| 332-1212 | 1.2 mm X 12 mm Lag Screw |
| 332-1213 | 1.2 mm X 13 mm Lag Screw |
| 332-1214 | 1.2 mm X 14 mm Lag Screw |
| 332-1216 | 1.2 mm X 16 mm Lag Screw |
| 332-1218 | 1.2 mm X 18 mm Lag Screw |
| 333-1201 | 1.2 mm 6 Hole Straight Plate |
| 333-1202 | 1.2 mm 12 Hole Straight Plate |
| 333-1203 | 1.2 mm Y Plate |

OsteoMed® Upper Extremities Portfolio [continued]

OsteoMed® Hand Plating System (hps™)

| Non-Sterile Part Number | Descriptor |
|-------------------------|--|
| 333-1204 | 1.2 mm L Plate, Left |
| 333-1205 | 1.2 mm L Plate, Right |
| 333-1206 | 1.2 mm 3 X 8 Hole T Plate |
| 333-1207 | 1.2 mm 4 X 8 Hole T Plate |
| 333-1208 | 1.2 mm Grid Plate, Offset Left |
| 333-1209 | 1.2 mm Grid Plate, Offset Right |
| 333-1601 | 1.6 mm 6 Hole Straight Plate, Locking |
| 333-1602 | 1.6 mm 12 Hole Straight Plate, Locking |
| 333-1603 | 1.6 mm Y Plate, Locking |
| 333-1604 | 1.6 mm L Plate, Left, Locking |
| 333-1605 | 1.6 mm L Plate, Right, Locking |
| 333-1606 | 1.6 mm 3 X 8 Hole T Plate, Locking |
| 333-1607 | 1.6 mm 4 X 8 Hole T Plate, Locking |
| 333-1608 | 1.6 mm Grid Plate, Offset Left, Locking |
| 333-1609 | 1.6 mm Grid Plate, Offset Right, Locking |
| 333-1611 | 1.6 mm Condylar Plate, Left, Locking |
| 333-1612 | 1.6 mm Condylar Plate, Right, Locking |
| 333-1620 | 1.6 mm Subcondylar Plate, Locking |
| 333-1651 | 1.6 mm 6 Hole Straight Plate, Tia |
| 333-1652 | 1.6 mm Subcondylar Plate, Tia |
| 333-1653 | 1.6 mm 8 Hole Straight Plate, Tia |
| 333-2001 | 2.0 mm 6 Hole Straight Plate, Locking |
| 333-2002 | 2.0 mm 12 Hole Straight Plate, Locking |
| 333-2003 | 2.0 mm Y Plate, Locking |
| 333-2004 | 2.0 mm L Plate, Left, Locking |
| 333-2005 | 2.0 mm L Plate, Right, Locking |
| 333-2006 | 2.0 mm 2 X 8 Hole T Plate, Locking |
| 333-2007 | 2.0 mm 3 X 8 Hole T Plate, Locking |
| 333-2008 | 2.0 mm Grid Plate, Offset Left, Locking |
| 333-2009 | 2.0 mm Grid Plate, Offset Right, Locking |

OsteoMed® Upper Extremities Portfolio [continued]

OsteoMed® Hand Plating System (hps™)

| Non-Sterile Part Number | Descriptor |
|-------------------------|--|
| 333-2010 | 2.0 mm Subcondylar Plate, Locking |
| 333-2011 | 2.0 mm Condylar Plate, Left, Locking |
| 333-2012 | 2.0 mm Condylar Plate, Right, Locking |
| 333-2013 | 2.0 mm Z Plate, Locking |
| 333-2020 | 2.0 mm 4 Hole Straight Lcdcp |
| 333-2021 | 2.0 mm 6 Hole Straight Lcdcp |
| 333-2022 | 2.0 mm 8 Hole Straight Lcdcp |
| 333-2051 | 2.0 mm 6 Hole Straight Plate, Tia |
| 333-2052 | 2.0 mm Subcondylar Plate, Tia |
| 333-2053 | 2.0 mm 8 Hole Straight Plate, Tia |
| 333-2054 | 2.0 mm Cmc Plate, Left |
| 333-2055 | 2.0 mm Cmc Plate, Right |
| 333-2401 | 2.4 mm 6 Hole Straight Plate, Locking |
| 333-2402 | 2.4 mm 12 Hole Straight Plate, Locking |
| 333-2403 | 2.4 mm Y Plate, Locking |
| 333-2404 | 2.4 mm L Plate, Locking, Left |
| 333-2405 | 2.4 mm L Plate, Locking, Right |
| 333-2406 | 2.4 mm 2 X 8 T Plate, Locking |
| 333-2407 | 2.4 mm 3 X 8 T Plate, Locking |
| 333-2408 | 2.4 mm Offset Grid Plate, Left, Locking |
| 333-2409 | 2.4 mm Offset Grid Plate, Right, Locking |
| 333-2410 | 2.4 mm Sybcondylar Plate, Locking |
| 333-2411 | 2.4 mm Condylar Plate, Left, Locking |
| 333-2412 | 2.4 mm Condylar Plate, Right, Locking |
| 333-2413 | 2.4 mm Z Plate, Locking |
| 333-2420 | 2.4 mm 4 Hole Straight Lcdcp |
| 333-2421 | 2.4 mm 6 Hole Straight Lcdcp |
| 333-2422 | 2.4 mm 8 Hole Straight Lcdcp |
| 333-2451 | 2.4 mm 6 Hole Straight Plate, Tia |
| 333-2452 | 2.4 mm Subcondylar Plate, Tia |

OsteoMed® Upper Extremities Portfolio [continued]

OsteoMed® Hand Plating System (hps™)

| Non-Sterile Part Number | Descriptor |
|-------------------------|---|
| 333-2453 | 2.4 mm 8 Hole Straight Plate, Tia |
| 333-2454 | 2.4 mm Cmc Plate, Left |
| 333-2455 | 2.4 mm Cmc Plate, Right |
| 333-2456 | 2.4 mm 6 Hole Straight Plate, Tia, Gold |
| 333-2457 | 2.4 mm Subcondylar Plate, Tia, Gold |
| 333-2458 | 2.4 mm 8 Hole Straight Plate, Tia, Gold |
| 316-0125 | 2.0 mm Washer |
| 316-0126 | 2.4 mm Washer |
| 316-0225 | 3.0 mm Washer |
| 316-0226 | 4.0 mm Washer |
| 319-2006 | 2.0 mm X 6 mm Cannulated Lag Screw |
| 319-2008 | 2.0 mm X 8 mm Cannulated Lag Screw |
| 319-2010 | 2.0 mm X 10 mm Cannulated Lag Screw |
| 319-2012 | 2.0 mm X 12 mm Cannulated Lag Screw |
| 319-2014 | 2.0 mm X 14 mm Cannulated Lag Screw |
| 319-2016 | 2.0 mm X 16 mm Cannulated Lag Screw |
| 319-2018 | 2.0 mm X 18 mm Cannulated Lag Screw |
| 319-2020 | 2.0 mm X 20 mm Cannulated Lag Screw |
| 319-2022 | 2.0 mm X 22 mm Cannulated Lag Screw |
| 319-2024 | 2.0 mm X 24 mm Cannulated Lag Screw |
| 319-2026 | 2.0 mm X 26 mm Cannulated Lag Screw |
| 319-2028 | 2.0 mm X 28 mm Cannulated Lag Screw |
| 319-2030 | 2.0 mm X 30 mm Cannulated Lag Screw |
| 319-2032 | 2.0 mm X 32 mm Cannulated Lag Screw |
| 319-2034 | 2.0 mm X 34 mm Cannulated Lag Screw |
| 319-2036 | 2.0 mm X 36 mm Cannulated Lag Screw |
| 319-2038 | 2.0 mm X 38 mm Cannulated Lag Screw |
| 319-2040 | 2.0 mm X 40 mm Cannulated Lag Screw |
| 319-2042 | 2.0 mm X 42 mm Cannulated Lag Screw |
| 319-2406 | 2.4 mm X 6 mm Cannulated Lag Screw |

OsteoMed® Upper Extremities Portfolio [continued]

OsteoMed® Hand Plating System (hps™)

| Non-Sterile Part Number | Descriptor |
|-------------------------|-------------------------------------|
| 319-2408 | 2.4 mm X 8 mm Cannulated Lag Screw |
| 319-2410 | 2.4 mm X 10 mm Cannulated Lag Screw |
| 319-2412 | 2.4 mm X 12 mm Cannulated Lag Screw |
| 319-2414 | 2.4 mm X 14 mm Cannulated Lag Screw |
| 319-2416 | 2.4 mm X 16 mm Cannulated Lag Screw |
| 319-2418 | 2.4 mm X 18 mm Cannulated Lag Screw |
| 319-2420 | 2.4 mm X 20 mm Cannulated Lag Screw |
| 319-2422 | 2.4 mm X 22 mm Cannulated Lag Screw |
| 319-2424 | 2.4 mm X 24 mm Cannulated Lag Screw |
| 319-2426 | 2.4 mm X 26 mm Cannulated Lag Screw |
| 319-2428 | 2.4 mm X 28 mm Cannulated Lag Screw |
| 319-2430 | 2.4 mm X 30 mm Cannulated Lag Screw |
| 319-2432 | 2.4 mm X 32 mm Cannulated Lag Screw |
| 319-2434 | 2.4 mm X 34 mm Cannulated Lag Screw |
| 319-2436 | 2.4 mm X 36 mm Cannulated Lag Screw |
| 319-2438 | 2.4 mm X 38 mm Cannulated Lag Screw |
| 319-2440 | 2.4 mm X 40 mm Cannulated Lag Screw |
| 319-2442 | 2.4 mm X 42 mm Cannulated Lag Screw |
| 319-2444 | 2.4 mm X 44 mm Cannulated Lag Screw |
| 319-2446 | 2.4 mm X 46 mm Cannulated Lag Screw |
| 319-2448 | 2.4 mm X 48 mm Cannulated Lag Screw |
| 319-2450 | 2.4 mm X 50 mm Cannulated Lag Screw |
| 319-3010 | 3.0 mm X 10 mm Cannulated Lag Screw |
| 319-3012 | 3.0 mm X 12 mm Cannulated Lag Screw |
| 319-3014 | 3.0 mm X 14 mm Cannulated Lag Screw |
| 319-3016 | 3.0 mm X 16 mm Cannulated Lag Screw |
| 319-3018 | 3.0 mm X 18 mm Cannulated Lag Screw |
| 319-3020 | 3.0 mm X 20 mm Cannulated Lag Screw |
| 319-3022 | 3.0 mm X 22 mm Cannulated Lag Screw |
| 319-3024 | 3.0 mm X 24 mm Cannulated Lag Screw |

OsteoMed® Upper Extremities Portfolio [continued]

OsteoMed® Hand Plating System (hps™)

| Non-Sterile Part Number | Descriptor |
|-------------------------|-------------------------------------|
| 319-3026 | 3.0 mm X 26 mm Cannulated Lag Screw |
| 319-3028 | 3.0 mm X 28 mm Cannulated Lag Screw |
| 319-3030 | 3.0 mm X 30 mm Cannulated Lag Screw |
| 319-3032 | 3.0 mm X 32 mm Cannulated Lag Screw |
| 319-3034 | 3.0 mm X 34 mm Cannulated Lag Screw |
| 319-3036 | 3.0 mm X 36 mm Cannulated Lag Screw |
| 319-3038 | 3.0 mm X 38 mm Cannulated Lag Screw |
| 319-3040 | 3.0 mm X 40 mm Cannulated Lag Screw |
| 319-4012 | 4.0 mm X 12 mm Cannulated Lag Screw |
| 319-4014 | 4.0 mm X 14 mm Cannulated Lag Screw |
| 319-4016 | 4.0 mm X 16 mm Cannulated Lag Screw |
| 319-4018 | 4.0 mm X 18 mm Cannulated Lag Screw |
| 319-4020 | 4.0 mm X 20 mm Cannulated Lag Screw |
| 319-4022 | 4.0 mm X 22 mm Cannulated Lag Screw |
| 319-4024 | 4.0 mm X 24 mm Cannulated Lag Screw |
| 319-4026 | 4.0 mm X 26 mm Cannulated Lag Screw |
| 319-4028 | 4.0 mm X 28 mm Cannulated Lag Screw |
| 319-4030 | 4.0 mm X 30 mm Cannulated Lag Screw |
| 319-4032 | 4.0 mm X 32 mm Cannulated Lag Screw |
| 319-4034 | 4.0 mm X 34 mm Cannulated Lag Screw |
| 319-4036 | 4.0 mm X 36 mm Cannulated Lag Screw |
| 319-4038 | 4.0 mm X 38 mm Cannulated Lag Screw |
| 319-4040 | 4.0 mm X 40 mm Cannulated Lag Screw |
| 319-4044 | 4.0 mm X 44 mm Cannulated Lag Screw |
| 319-4048 | 4.0 mm X 48 mm Cannulated Lag Screw |
| 319-4052 | 4.0 mm X 52 mm Cannulated Lag Screw |

OsteoMed® Upper Extremities Portfolio [continued]

OsteoMed® hps™ Hand Fusion System

| Non-Sterile Part Number | Descriptor |
|-------------------------|-----------------------------|
| 333-1614 | 1.6 mm Hand Fusion Plate |
| 333-2014 | 2.0 mm Hand Fusion Plate |
| 334-2016 | 2.0 mm X 16 mm Fusion Screw |
| 334-2018 | 2.0 mm X 18 mm Fusion Screw |
| 334-2020 | 2.0 mm X 20 mm Fusion Screw |
| 334-2022 | 2.0 mm X 22 mm Fusion Screw |
| 334-2024 | 2.0 mm X 24 mm Fusion Screw |
| 334-2026 | 2.0 mm X 26 mm Fusion Screw |
| 334-2028 | 2.0 mm X 28 mm Fusion Screw |
| 334-2030 | 2.0 mm X 30 mm Fusion Screw |
| 334-2032 | 2.0 mm X 32 mm Fusion Screw |
| 334-2420 | 2.4 mm X 20 mm Fusion Screw |
| 334-2422 | 2.4 mm X 22 mm Fusion Screw |
| 334-2424 | 2.4 mm X 24 mm Fusion Screw |
| 334-2426 | 2.4 mm X 26 mm Fusion Screw |
| 334-2428 | 2.4 mm X 28 mm Fusion Screw |
| 334-2430 | 2.4 mm X 30 mm Fusion Screw |
| 334-2432 | 2.4 mm X 32 mm Fusion Screw |
| 334-2434 | 2.4 mm X 34 mm Fusion Screw |
| 334-2436 | 2.4 mm X 36 mm Fusion Screw |

OsteoMed® Upper Extremities Portfolio [continued]

OsteoMed® M3-X™ Extremity Fixation System

| Non-Sterile Part Number | Descriptor |
|-------------------------|-------------------------------------|
| 301-2008 | 2.0 mm X 8 mm Fully Threaded Screw |
| 301-2010 | 2.0 mm X 10 mm Fully Threaded Screw |
| 301-2012 | 2.0 mm X 12 mm Fully Threaded Screw |
| 301-2014 | 2.0 mm X 14 mm Fully Threaded Screw |
| 301-2016 | 2.0 mm X 16 mm Fully Threaded Screw |
| 301-2018 | 2.0 mm X 18 mm Fully Threaded Screw |
| 301-2020 | 2.0 mm X 20 mm Fully Threaded Screw |
| 301-2022 | 2.0 mm X 22 mm Fully Threaded Screw |
| 302-2008 | 2.0 mm X 8 mm Lag Screw |
| 302-2010 | 2.0 mm X 10 mm Lag Screw |
| 302-2012 | 2.0 mm X 12 mm Lag Screw |
| 302-2014 | 2.0 mm X 14 mm Lag Screw |
| 302-2016 | 2.0 mm X 16 mm Lag Screw |
| 302-2018 | 2.0 mm X 18 mm Lag Screw |
| 302-2020 | 2.0 mm X 20 mm Lag Screw |
| 302-2022 | 2.0 mm X 22 mm Lag Screw |
| 303-1204 | 1.2 mm X 4 mm Fully Threaded Screw |
| 303-1206 | 1.2 mm X 6 mm Fully Threaded Screw |
| 303-1208 | 1.2 mm X 8 mm Fully Threaded Screw |
| 303-1210 | 1.2 mm X 10 mm Fully Threaded Screw |
| 303-1212 | 1.2 mm X 12 mm Fully Threaded Screw |
| 304-1206 | 1.2 mm X 6 mm Lag Screw |
| 304-1208 | 1.2 mm X 8 mm Lag Screw |
| 304-1210 | 1.2 mm X 10 mm Lag Screw |
| 304-1212 | 1.2 mm X 12 mm Lag Screw |
| 304-1214 | 1.2 mm X 14 mm Lag Screw |
| 304-1216 | 1.2 mm X 16 mm Lag Screw |
| 305-2408 | 2.4 mm X 8 mm Fully Threaded Screw |
| 305-2410 | 2.4 mm X 10 mm Fully Threaded Screw |
| 305-2412 | 2.4 mm X 12 mm Fully Threaded Screw |

OsteoMed® Upper Extremities Portfolio [continued]

OsteoMed® M3-X™ Extremity Fixation System

| Non-Sterile Part Number | Descriptor |
|-------------------------|-------------------------------------|
| 305-2414 | 2.4 mm X 14 mm Fully Threaded Screw |
| 305-2416 | 2.4 mm X 16 mm Fully Threaded Screw |
| 305-2418 | 2.4 mm X 18 mm Fully Threaded Screw |
| 305-2420 | 2.4 mm X 20 mm Fully Threaded Screw |
| 305-2422 | 2.4 mm X 22 mm Fully Threaded Screw |
| 306-2408 | 2.4 mm X 8 mm Lag Screw |
| 306-2410 | 2.4 mm X 10 mm Lag Screw |
| 306-2412 | 2.4 mm X 12 mm Lag Screw |
| 306-2414 | 2.4 mm X 14 mm Lag Screw |
| 306-2416 | 2.4 mm X 16 mm Lag Screw |
| 306-2418 | 2.4 mm X 18 mm Lag Screw |
| 306-2420 | 2.4 mm X 20 mm Lag Screw |
| 306-2422 | 2.4 mm X 22 mm Lag Screw |
| 306-2424 | 2.4 mm X 24 mm Lag Screw |
| 306-2426 | 2.4 mm X 26 mm Lag Screw |
| 306-2428 | 2.4 mm X 28 mm Lag Screw |
| 306-2430 | 2.4 mm X 30 mm Lag Screw |
| 306-2432 | 2.4 mm X 32 mm Lag Screw |
| 308-3012 | 3.0 mm X 12 mm Lag Screw |
| 308-3014 | 3.0 mm X 14 mm Lag Screw |
| 308-3016 | 3.0 mm X 16 mm Lag Screw |
| 308-3018 | 3.0 mm X 18 mm Lag Screw |
| 308-3020 | 3.0 mm X 20 mm Lag Screw |
| 308-3022 | 3.0 mm X 22 mm Lag Screw |
| 308-3024 | 3.0 mm X 24 mm Lag Screw |
| 308-3026 | 3.0 mm X 26 mm Lag Screw |
| 308-3028 | 3.0 mm X 28 mm Lag Screw |
| 308-3030 | 3.0 mm X 30 mm Lag Screw |
| 308-3032 | 3.0 mm X 32 mm Lag Screw |
| 308-3034 | 3.0 mm X 34 mm Lag Screw |

OsteoMed® Upper Extremities Portfolio [continued]

OsteoMed® M3-X™ Extremity Fixation System

| Non-Sterile Part Number | Descriptor |
|-------------------------|--|
| 308-3036 | 3.0 mm X 36 mm Lag Screw |
| 310-2114 | 2 X 4 Hole Compression Left "L" Plate |
| 310-2124 | 2 X 4 Hole Compression Right "L" Plate |
| 310-2134 | 2 X 4 Hole Compression Left Oblique "L" Plate |
| 310-2144 | 2 X 4 Hole Compression Right Oblique "L" Plate |
| 310-2204 | 2 X 4 Hole Compression "T" Plate |
| 310-2234 | 2 X 4 Hole Compression Left Ob |
| 310-2244 | 2 X 4 Hole Compression Right Oblique "T" Plate |
| 310-2304 | 3 X 4 Hole Compression "T" Plate |
| 310-2404 | 3 X 4 Hole Compression "Y" Plate |
| 310-3004 | 4 Hole Straight Compression Plate |
| 310-3005 | 5 Hole Straight Compression Plate |
| 310-3006 | Mini 6 Hole Compression Plate |
| 310-3304 | 2 X 1 X 4 Compression "T" Plate |
| 310-4005 | 5 Hole, Straight Compression Plate Rigid |
| 310-4006 | 6 Hole Straight Compression Plate Rigid |
| 310-4007 | 7 Hole Straight Compression Plate Rigid |
| 312-1003 | 2 X 3 Hole Quad Plate |
| 312-2203 | 2 X 3 Hole Compression "T" Plate |
| 312-2233 | 2 X 3 Hole Oblique Left Compression "T" Plate |
| 312-2243 | 2 X 3 Hole Oblique Right Compression "T" Plate |
| 312-3004 | 4 Hole Straight Compression Plate |
| 314-0040 | D-Maxi 4 Hole Straight Plate |
| 314-2114 | 2 X 4 Hole Compression Left "L" Plate |
| 314-2124 | 2 X 4 Hole Compression Right "L" Plate |
| 314-2204 | 2 X 4 Hole Compression "T" Plate |
| 314-2234 | 2 X 4 Hole Oblique Left Compression "T" Plate |
| 314-2244 | 2 X 4 Hole Oblique Right Compression "T" Plate |
| 314-3004 | 4 Hole Straight Compression Plate |
| 314-3005 | 5 Hole Straight Compression Plate |

OsteoMed® Upper Extremities Portfolio [continued]

OsteoMed® M3-X™ Extremity Fixation System

| Non-Sterile Part Number | Descriptor |
|----------------------------|-----------------------------------|
| 314-3006 | 6 Hole Straight Compression Plate |
| 306-2438 | 2.4 mm X 38 mm Lag Screw |
| 306-2440 | 2.4 mm X 40 mm Lag Screw |
| 306-2442 | 2.4 mm X 42 mm Lag Screw |

OsteoMed® Upper Extremities Portfolio [continued]

OsteoMed® InstaFix™ Shape Memory Fixation System

| Non-Sterile Part Number | Descriptor |
|-------------------------|-------------------------------|
| 500-0020-00 | Implant, 2 Arm, 25c, 8x8 |
| 500-0020-01 | Implant, 2 Arm, 25c, 10x10 |
| 500-0020-02 | Implant, 2 Arm, 25c, 15x10 |
| 500-0020-03 | Implant, 2 Arm, 25c, 10x15/17 |
| 500-0020-04 | Implant, 2 Arm, 25c, 15x15/17 |
| 500-0020-05 | Implant, 2 Arm, 25c, 20x15 |
| 500-0020-06 | Implant, 2 Arm, 25c, 20x20 |
| 500-0020-07 | Implant, 2 Arm, 25c, 25x20 |
| 500-0020-08 | Implant, 2 Arm, 25c, 30x20 |

OsteoMed® Upper Extremities Portfolio [continued]

OsteoMed® ExtremiLock™ Wrist Plating System

| Non-Sterile Part Number | Descriptor |
|-------------------------|--|
| 349-2412 | 2.4 mm X 12 mm Cann Lag Va Locking Screw |
| 349-2413 | 2.4 mm X 13 mm Cann Lag Va Locking Screw |
| 349-2414 | 2.4 mm X 14 mm Cann Lag Va Locking Screw |
| 349-2416 | 2.4 mm X 16 mm Cann Lag Va Locking Screw |
| 349-2418 | 2.4 mm X 18 mm Cann Lag Va Locking Screw |
| 349-2419 | 2.4 mm X 19 mm Cann Lag Va Locking Screw |
| 349-2420 | 2.4 mm X 20 mm Cann Lag Va Locking Screw |
| 349-2421 | 2.4 mm X 21 mm Cann Lag Va Locking Screw |
| 349-2422 | 2.4 mm X 22 mm Cann Lag Va Locking Screw |
| 349-2423 | 2.4 mm X 23 mm Cann Lag Va Locking Screw |
| 349-2424 | 2.4 mm X 24 mm Cann Lag Va Locking Screw |
| 349-2426 | 2.4 mm X 26 mm Cann Lag Va Locking Screw |
| 349-2428 | 2.4 mm X 28 mm Cann Lag Va Locking Screw |
| 349-2430 | 2.4 mm X 30 mm Cann Lag Va Locking Screw |
| 349-2432 | 2.4 mm X 32 mm Cann Lag Va Locking Screw |
| 349-2434 | 2.4 mm X 34 mm Cann Lag Va Locking Screw |
| 349-2436 | 2.4 mm X 36 mm Cann Lag Va Locking Screw |
| 349-2438 | 2.4 mm X 38 mm Cann Lag Va Locking Screw |
| 349-2440 | 2.4 mm X 40 mm Cann Lag Va Locking Screw |
| 354-1808 | 2.4 mm X 8 mm Partial Thread Locking Screw |
| 354-1809 | 2.4 mm X 9 mm Partial Thread Locking Screw |
| 354-1810 | 2.4 mm X 10 mm Partial Thread Lockingscrew |
| 354-1811 | 2.4 mm X 11 mm Partial Thread Lockingscrew |
| 354-1812 | 2.4 mm X 12 mm Partial Thread Lockingscrew |
| 354-1813 | 2.4 mm X 13 mm Partial Thread Lockingscrew |
| 354-1814 | 2.4 mm X 14 mm Partial Thread Lockingscrew |
| 354-1815 | 2.4 mm X 15 mm Partial Thread Lockingscrew |
| 354-1816 | 2.4 mm X 16 mm Partial Thread Lockingscrew |
| 354-1818 | 2.4 mm X 18 mm Partial Thread Lockingscrew |
| 354-1819 | 2.4 mm X 19 mm Partial Thread Lockingscrew |

OsteoMed® Upper Extremities Portfolio [continued]

OsteoMed® ExtremiLock™ Wrist Plating System

| Non-Sterile Part Number | Descriptor |
|-------------------------|--|
| 354-1820 | 2.4 mm X 20 mm Partial Thread Lockingscrew |
| 354-1821 | 2.4 mm X 21 mm Partial Thread Lockingscrew |
| 354-1822 | 2.4 mm X 22 mm Partial Thread Lockingscrew |
| 354-1823 | 2.4 mm X 23 mm Partial Thread Lockingscrew |
| 354-1824 | 2.4 mm X 24 mm Partial Thread Lockingscrew |
| 354-1826 | 2.4 mm X 26 mm Partial Thread Lockingscrew |
| 354-1828 | 2.4 mm X 28 mm Partial Thread Lockingscrew |
| 354-1830 | 2.4 mm X 30 mm Partial Thread Lockingscrew |
| 354-1832 | 2.4 mm X 32 mm Partial Thread Lockingscrew |
| 354-1834 | 2.4 mm X 34 mm Partial Thread Lockingscrew |
| 354-1836 | 2.4 mm X 36 mm Partial Thread Lockingscrew |
| 354-1838 | 2.4 mm X 38 mm Partial Thread Lockingscrew |
| 354-1840 | 2.4 mm X 40 mm Partial Thread Lockingscrew |
| 355-0101 | 4 Hole Vdr Narrow Plate, Left |
| 355-0102 | Vdr Narrow/Left, Long |
| 355-0106 | 4 Hole Vdr Narrow Plate, Right |
| 355-0107 | Vdr Narrow/Right, Long |
| 355-0111 | 4 Hole Vdr Standard Plate, Left |
| 355-0112 | 7 Hole Vdr Standard Plate, Left |
| 355-0116 | 4 Hole Vdr Standard Plate, Right |
| 355-0117 | 7 Hole Vdr Standard Plate, Right |
| 355-0200 | 1 Hole Vdr Narrow Plate, Left |
| 355-0201 | 1 Hole Vdr Standard Plate, Left |
| 355-0205 | 1 Hole Vdr Narrow Plate, Right |
| 355-0206 | 1 Hole Vdr Standard Plate, Right |
| 355-0300 | 4 Hole Vdr Narrow Pq Sparing Left |
| 355-0301 | 4 Hole Vdr Standard Pq Sparing, Left |
| 355-0305 | 4 Hole Vdr Narrow Pq Sparing, Right |
| 355-0306 | 4 Hole Vdr Standard Pq Sparing, Right |
| 355-0400 | Vdr Radial Column Standard, Left |

OsteoMed® Upper Extremities Portfolio [continued]

OsteoMed® ExtremiLock™ Wrist Plating System

| Non-Sterile Part Number | Descriptor |
|-------------------------|---|
| 355-0401 | Vdr Radial Column Standard, Right |
| 355-0402 | Vdr Intermediate Column Standard, Left |
| 355-0403 | Vdr Intermediate Column Standard, Right |
| 355-0601 | 6 Hole Radial Styloid Plate |
| 358-0099 | 2.4/2.7 Washer |
| 358-2408 | 2.4 mm X 8 mm Cortex Screw |
| 358-2409 | 2.4 mm X 9 mm Cortex Screw |
| 358-2410 | 2.4 mm X 10 mm Cortex Screw |
| 358-2411 | 2.4 mm X 11 mm Cortex Screw |
| 358-2412 | 2.4 mm X 12 mm Cortex Screw |
| 358-2413 | 2.4 mm X 13 mm Cortex Screw |
| 358-2414 | 2.4 mm X 14 mm Cortex Screw |
| 358-2415 | 2.4 mm X 15 mm Cortex Screw |
| 358-2416 | 2.4 mm X 16 mm Cortex Screw |
| 358-2418 | 2.4 mm X 18 mm Cortex Screw |
| 358-2419 | 2.4 mm X 19 mm Cortex Screw |
| 358-2420 | 2.4 mm X 20 mm Cortex Screw |
| 358-2421 | 2.4 mm X 21 mm Cortex Screw |
| 358-2422 | 2.4 mm X 22 mm Cortex Screw |
| 358-2423 | 2.4 mm X 23 mm Cortex Screw |
| 358-2424 | 2.4 mm X 24 mm Cortex Screw |
| 358-2426 | 2.4 mm X 26 mm Cortex Screw |
| 358-2428 | 2.4 mm X 28 mm Cortex Screw |
| 358-2430 | 2.4 mm X 30 mm Cortex Screw |
| 358-2432 | 2.4 mm X 32 mm Cortex Screw |
| 358-2434 | 2.4 mm X 34 mm Cortex Screw |
| 358-2436 | 2.4 mm X 36 mm Cortex Screw |
| 358-2438 | 2.4 mm X 38 mm Cortex Screw |
| 358-2440 | 2.4 mm X 40 mm Cortex Screw |
| 358-2708 | 2.7 mm X 8 mm Cortex Screw |

OsteoMed® Upper Extremities Portfolio [continued]

OsteoMed® ExtremiLock™ Wrist Plating System

| Non-Sterile Part Number | Descriptor |
|-------------------------|--|
| 358-2709 | 2.7 mm X 9 mm Cortex Screw |
| 358-2710 | 2.7 mm X 10 mm Cortex Screw |
| 358-2711 | 2.7 mm X 11 mm Cortex Screw |
| 358-2712 | 2.7 mm X 12 mm Cortex Screw |
| 358-2713 | 2.7 mm X 13 mm Cortex Screw |
| 358-2714 | 2.7 mm X 14 mm Cortex Screw |
| 358-2715 | 2.7 mm X 15 mm Cortex Screw |
| 358-2716 | 2.7 mm X 16 mm Cortex Screw |
| 358-2718 | 2.7 mm X 18 mm Cortex Screw |
| 358-2720 | 2.7 mm X 20 mm Cortex Screw |
| 358-2722 | 2.7 mm X 22 mm Cortex Screw |
| 358-2724 | 2.7 mm X 24 mm Cortex Screw |
| 358-2726 | 2.7 mm X 26 mm Cortex Screw |
| 358-2728 | 2.7 mm X 28 mm Cortex Screw |
| 358-2730 | 2.7 mm X 30 mm Cortex Screw |
| 358-2732 | 2.7 mm X 32 mm Cortex Screw |
| 358-2734 | 2.7 mm X 34 mm Cortex Screw |
| 358-2736 | 2.7 mm X 36 mm Cortex Screw |
| 358-2738 | 2.7 mm X 38 mm Cortex Screw |
| 358-2740 | 2.7 mm X 40 mm Cortex Screw |
| 359-2408 | 2.4 mm X 8 mm Variable Angle Locking Screw |
| 359-2409 | 2.4 mm X 9 mm Variable Angle Locking Screw |
| 359-2410 | 2.4 mm X 10 mm Variable Angle Lockingscrew |
| 359-2411 | 2.4 mm X 11 mm Variable Angle Lockingscrew |
| 359-2412 | 2.4 mm X 12 mm Variable Angle Lockingscrew |
| 359-2413 | 2.4 mm X 13 mm Variable Angle Lockingscrew |
| 359-2414 | 2.4 mm X 14 mm Variable Angle Lockingscrew |
| 359-2415 | 2.4 mm X 15 mm Variable Angle Lockingscrew |
| 359-2416 | 2.4 mm X 16 mm Variable Angle Lockingscrew |
| 359-2418 | 2.4 mm X 18 mm Variable Angle Lockingscrew |

OsteoMed® Upper Extremities Portfolio [continued]

OsteoMed® ExtremiLock™ Wrist Plating System

| Non-Sterile Part Number | Descriptor |
|-------------------------|--|
| 359-2419 | 2.4 mm X 19 mm Variable Angle Lockingscrew |
| 359-2420 | 2.4 mm X 20 mm Variable Angle Lockingscrew |
| 359-2421 | 2.4 mm X 21 mm Variable Angle Lockingscrew |
| 359-2422 | 2.4 mm X 22 mm Variable Angle Lockingscrew |
| 359-2423 | 2.4 mm X 23 mm Variable Angle Lockingscrew |
| 359-2424 | 2.4 mm X 24 mm Variable Angle Lockingscrew |
| 359-2426 | 2.4 mm X 26 mm Variable Angle Lockingscrew |
| 359-2428 | 2.4 mm X 28 mm Variable Angle Lockingscrew |
| 359-2430 | 2.4 mm X 30 mm Variable Angle Lockingscrew |
| 359-2432 | 2.4 mm X 32 mm Variable Angle Lockingscrew |
| 359-2434 | 2.4 mm X 34 mm Variable Angle Lockingscrew |
| 359-2436 | 2.4 mm X 36 mm Variable Angle Lockingscrew |
| 359-2438 | 2.4 mm X 38 mm Variable Angle Lockingscrew |
| 359-2440 | 2.4 mm X 40 mm Variable Angle Lockingscrew |
| 359-2708 | 2.7 mm X 8 mm Variable Angle Locking Screw |
| 359-2709 | 2.7 mm X 9 mm Variable Angle Locking Screw |
| 359-2710 | 2.7 mm X 10 mm Variable Angle Lockingscrew |
| 359-2711 | 2.7 mm X 11 mm Variable Angle Lockingscrew |
| 359-2712 | 2.7 mm X 12 mm Variable Angle Lockingscrew |
| 359-2713 | 2.7 mm X 13 mm Variable Angle Lockingscrew |
| 359-2714 | 2.7 mm X 14 mm Variable Angle Lockingscrew |
| 359-2715 | 2.7 mm X 15 mm Variable Angle Lockingscrew |
| 359-2716 | 2.7 mm X 16 mm Variable Angle Lockingscrew |
| 359-2718 | 2.7 mm X 18 mm Variable Angle Lockingscrew |
| 359-2720 | 2.7 mm X 20 mm Variable Angle Lockingscrew |
| 359-2722 | 2.7 mm X 22 mm Variable Angle Lockingscrew |
| 359-2724 | 2.7 mm X 24 mm Variable Angle Lockingscrew |
| 359-2726 | 2.7 mm X 26 mm Variable Angle Lockingscrew |
| 359-2728 | 2.7 mm X 28 mm Variable Angle Lockingscrew |
| 359-2730 | 2.7 mm X 30 mm Variable Angle Lockingscrew |


OsteoMed® Upper Extremities Portfolio [continued]

OsteoMed® ExtremiLock™ Wrist Plating System

| Non-Sterile Part Number | Descriptor |
|-------------------------|--|
| 359-2732 | 2.7 mm X 32 mm Variable Angle Lockingscrew |
| 359-2734 | 2.7 mm X 34 mm Variable Angle Lockingscrew |
| 359-2736 | 2.7 mm X 36 mm Variable Angle Lockingscrew |
| 359-2738 | 2.7 mm X 38 mm Variable Angle Lockingscrew |
| 359-2740 | 2.7 mm X 40 mm Variable Angle Lockingscrew |

OsteoMed® Lower Extremities Portfolio

- ▶ OsteoMed Talar-Fit™ Subtalar Arthroereisis Implant
- ▶ OsteoMed Hemi™ Great Toe Implant
- ▶ OsteoMed FPS™ Foot Plating System
- ▶ OsteoMed ExtremiFuse™ Hammertoe Fixation System
- ▶ OsteoMed ExtremiLock™ Ankle Plating System
- ▶ OsteoMed EnCompass™ Metatarsal Resurfacing Implant
- ▶ OsteoMed ReFlexion™ 1st MTP Implant System

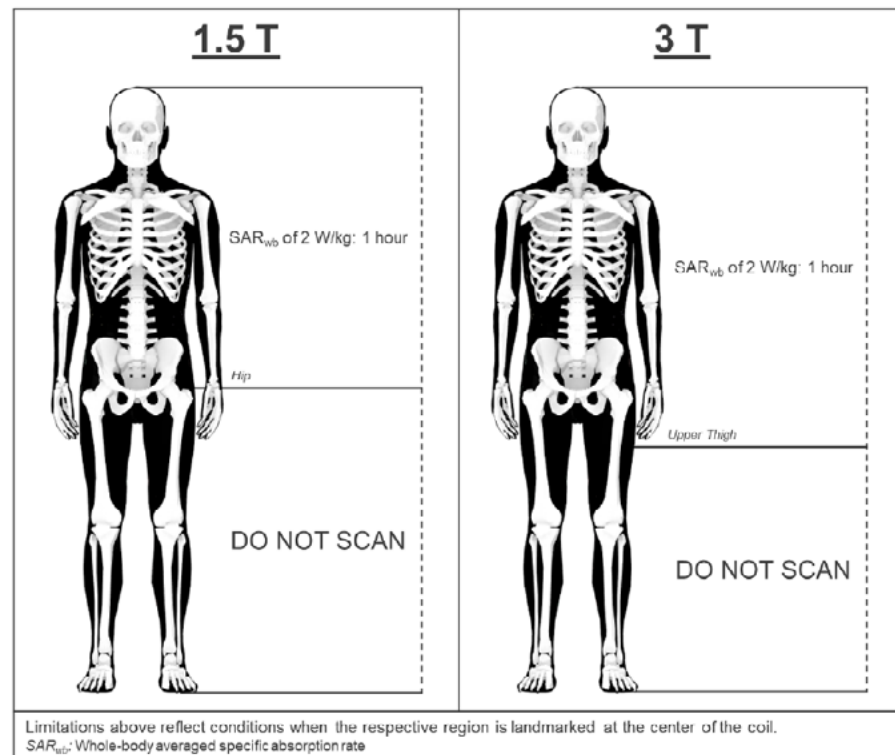

Conditions for Safe Scanning in the MRI

A patient with OsteoMed's Lower Extremity Products, **excluding any stand-alone screws or K-wires**, may be safely scanned at 1.5 T or 3.0 T under the following conditions. Failure to follow these conditions may result in injury.

OsteoMed® Lower Extremities Portfolio [continued]

OsteoMed's Lower Extremity Portfolio

| Parameter | Condition |
|---|---|
| Device Name | OsteoMed's Lower Extremity Products, excluding stand-alone screws and K-wires |
| Nominal value(s) of Static Magnetic Field | 1.5 T or 3 T |
| MR Scanner Type | Cylindrical |
| BO Field Orientation | Horizontal |
| Maximum Spacial Field Gradient | 40 T/m (4,000 G/cm) |
| RF Excitation | Circularly Polarized (CP) |
| RF Transmit Coil Type | Integrated Whole Body Transmit Coil |
| Operating Mode | <u>1.5 T:</u> Normal Operating Mode from the top of head to the hip |
| | Do Not Scan below the hip |
| | <u>3 T:</u> Normal Operating Mode from the top of head to the upper thigh |
| Scan Duration | Do Not Scan below the upper thigh |
| | Scan duration up to 1 hour without a cooling period with landmark restrictions according to Operating Mode section above and as detailed in the Lower Extremity MRI Restricted Zone Summary |
| Scan Regions | Do not scan below the hip for 1.5 T. All other landmark locations are permissible with restrictions according to Operating Mode section above and as detailed in the Lower Extremity MRI Restricted Zone Summary |
| Image Artifact | The presence of OsteoMed's Lower Extremity products may produce an image artifact of 3.0 cm. Some manipulation of scan parameters may be needed to compensate for the artifact. |

OsteoMed® Lower Extremities Portfolio [continued]

OsteoMed® Talar-Fit™ Subtalar Arthroereisis Implant

| Non-Sterile Part Number | Descriptor |
|----------------------------|--------------------------|
| 340-0001 | Subtalar Implant, Ø8 mm |
| 340-0002 | Subtalar Implant, Ø9 mm |
| 340-0003 | Subtalar Implant, Ø10 mm |
| 340-0004 | Subtalar Implant, Ø11 mm |
| 340-0005 | Subtalar Implant, Ø12 mm |

OsteoMed® Lower Extremities Portfolio [continued]

OsteoMed® Hemi™ Great Toe Implant

| Non-Sterile Part Number | Descriptor |
|----------------------------|-------------------------------------|
| 375-0001 | Hemi Great Toe Implant, Extra Small |
| 375-0002 | Hemi Great Toe Implant, Small |
| 375-0003 | Hemi Great Toe Implant, Medium |
| 375-0004 | Hemi Great Toe Implant, Large |

OsteoMed® Lower Extremities Portfolio [continued]

OsteoMed® FPS™ Foot Plating System

| Non-Sterile Part Number | Descriptor |
|-------------------------|------------------------------------|
| 321-2004 | Fps™ 2.0 mm X 4 mm Standard Screw |
| 321-2005 | Fps™ 2.0 mm X 5 mm Standard Screw |
| 321-2006 | Fps™ 2.0 mm X 6 mm Standard Screw |
| 321-2008 | Fps™ 2.0 mm X 8 mm Standard Screw |
| 321-2010 | Fps™ 2.0 mm X 10 mm Standard Screw |
| 321-2012 | Fps™ 2.0 mm X 12 mm Standard Screw |
| 321-2014 | Fps™ 2.0 mm X 14 mm Standard Screw |
| 321-2016 | Fps™ 2.0 mm X 16 mm Standard Screw |
| 321-2018 | Fps™ 2.0 mm X 18 mm Standard Screw |
| 321-2020 | Fps™ 2.0 mm X 20 mm Standard Screw |
| 321-2022 | Fps™ 2.0 mm X 22 mm Standard Screw |
| 321-2024 | Fps™ 2.0 mm X 24 mm Standard Screw |
| 321-2026 | Fps™ 2.0 mm X 26 mm Standard Screw |
| 321-2028 | Fps™ 2.0 mm X 28 mm Standard Screw |
| 321-2030 | Fps™ 2.0 mm X 30 mm Standard Screw |
| 321-2032 | Fps™ 2.0 mm X 32 mm Standard Screw |
| 321-2099 | Fps™ 2.0 mm Washer |
| 321-2499 | Fps™ 2.7 mm Washer |
| 321-2710 | Fps™ 2.7 mm X 10 mm Standard Screw |
| 321-2712 | Fps™ 2.7 mm X 12 mm Standard Screw |
| 321-2714 | Fps™ 2.7 mm X 14 mm Standard Screw |
| 321-2716 | Fps™ 2.7 mm X 16 mm Standard Screw |
| 321-2718 | Fps™ 2.7 mm X 18 mm Standard Screw |
| 321-2720 | Fps™ 2.7 mm X 20 mm Standard Screw |
| 321-2722 | Fps™ 2.7 mm X 22 mm Standard Screw |
| 321-2724 | Fps™ 2.7 mm X 24 mm Standard Screw |
| 321-2726 | Fps™ 2.7 mm X 26 mm Standard Screw |
| 321-2728 | Fps™ 2.7 mm X 28 mm Standard Screw |
| 321-2730 | Fps™ 2.7 mm X 30 mm Standard Screw |
| 321-2732 | Fps™ 2.7 mm X 32 mm Standard Screw |

OsteoMed® Lower Extremities Portfolio [continued]

OsteoMed® FPS™ Foot Plating System

| Non-Sterile Part Number | Descriptor |
|-------------------------|------------------------------------|
| 321-2734 | Fps™ 2.7 mm X 34 mm Standard Screw |
| 321-2736 | Fps™ 2.7 mm X 36 mm Standard Screw |
| 321-2738 | Fps™ 2.7 mm X 38 mm Standard Screw |
| 321-2740 | Fps™ 2.7 mm X 40 mm Standard Screw |
| 321-2745 | Fps™ 2.7 mm X 45 mm Standard Screw |
| 321-2750 | Fps™ 2.7 mm X 50 mm Standard Screw |
| 321-2755 | Fps™ 2.7 mm X 50 mm Standard Screw |
| 321-2760 | Fps™ 2.7 mm X 60 mm Standard Screw |
| 321-2765 | Fps™ 2.7 mm X 65 mm Standard Screw |
| 321-2770 | Fps™ 2.7 mm X 70 mm Standard Screw |
| 321-2775 | Fps™ 2.7 mm X 75 mm Standard Screw |
| 321-3510 | Fps™ 3.5 mm X 10 mm Standard Screw |
| 321-3512 | Fps™ 3.5 mm X 12 mm Standard Screw |
| 321-3514 | Fps™ 3.5 mm X 14 mm Standard Screw |
| 321-3516 | Fps™ 3.5 mm X 16 mm Standard Screw |
| 321-3518 | Fps™ 3.5 mm X 18 mm Standard Screw |
| 321-3520 | Fps™ 3.5 mm X 20 mm Standard Screw |
| 321-3522 | Fps™ 3.5 mm X 22 mm Standard Screw |
| 321-3524 | Fps™ 3.5 mm X 24 mm Standard Screw |
| 321-3526 | Fps™ 3.5 mm X 26 mm Standard Screw |
| 321-3528 | Fps™ 3.5 mm X 28 mm Standard Screw |
| 321-3530 | Fps™ 3.5 mm X 30 mm Standard Screw |
| 321-3532 | Fps™ 3.5 mm X 32 mm Standard Screw |
| 321-3534 | Fps™ 3.5 mm X 34 mm Standard Screw |
| 321-3536 | Fps™ 3.5 mm X 36 mm Standard Screw |
| 321-3538 | Fps™ 3.5 mm X 38 mm Standard Screw |
| 321-3540 | Fps™ 3.5 mm X 40 mm Standard Screw |
| 321-3545 | Fps™ 3.5 mm X 45 mm Standard Screw |
| 321-3550 | Fps™ 3.5 mm X 50 mm Standard Screw |
| 321-3555 | Fps™ 3.5 mm X 55 mm Standard Screw |

OsteoMed® Lower Extremities Portfolio [continued]

OsteoMed® FPS™ Foot Plating System

| Non-Sterile Part Number | Descriptor |
|-------------------------|------------------------------------|
| 321-3560 | Fps™ 3.5 mm X 60 mm Standard Screw |
| 321-3565 | Fps™ 3.5 mm X 65 mm Standard Screw |
| 321-3570 | Fps™ 3.5 mm X 70 mm Standard Screw |
| 321-3575 | Fps™ 3.5 mm X 75 mm Standard Screw |
| 321-3599 | Fps™ 3.5/4.0 mm Washer |
| 321-4010 | Fps™ 4.0 mm X 10 mm Standard Screw |
| 321-4012 | Fps™ 4.0 mm X 12 mm Standard Screw |
| 321-4014 | Fps™ 4.0 mm X 14 mm Standard Screw |
| 321-4016 | Fps™ 4.0 mm X 16 mm Standard Screw |
| 321-4018 | Fps™ 4.0 mm X 18 mm Standard Screw |
| 321-4020 | Fps™ 4.0 mm X 20 mm Standard Screw |
| 321-4022 | Fps™ 4.0 mm X 22 mm Standard Screw |
| 321-4024 | Fps™ 4.0 mm X 24 mm Standard Screw |
| 321-4026 | Fps™ 4.0 mm X 26 mm Standard Screw |
| 321-4028 | Fps™ 4.0 mm X 28 mm Standard Screw |
| 321-4030 | Fps™ 4.0 mm X 30 mm Standard Screw |
| 321-4032 | Fps™ 4.0 mm X 32 mm Standard Screw |
| 321-4034 | Fps™ 4.0 mm X 34 mm Standard Screw |
| 321-4036 | Fps™ 4.0 mm X 36 mm Standard Screw |
| 321-4038 | Fps™ 4.0 mm X 38 mm Standard Screw |
| 321-4040 | Fps™ 4.0 mm X 40 mm Standard Screw |
| 321-4045 | Fps™ 4.0 mm X 45 mm Standard Screw |
| 321-4050 | Fps™ 4.0 mm X 50 mm Standard Screw |
| 321-4055 | Fps™ 4.0 mm X 55 mm Standard Screw |
| 321-4060 | Fps™ 4.0 mm X 60 mm Standard Screw |
| 321-4065 | Fps™ 4.0 mm X 65 mm Standard Screw |
| 321-4070 | Fps™ 4.0 mm X 70 mm Standard Screw |
| 321-4075 | Fps™ 4.0 mm X 75 mm Standard Screw |
| 322-2006 | Fps™ 2.0 mm X 6 mm Locking Screw |
| 322-2007 | Fps™ 2.0 mm X 7 mm Locking Screw |

OsteoMed® Lower Extremities Portfolio [continued]

OsteoMed® FPS™ Foot Plating System

| Non-Sterile Part Number | Descriptor |
|-------------------------|-----------------------------------|
| 322-2008 | Fps™ 2.0 mm X 8 mm Locking Screw |
| 322-2009 | Fps™ 2.0 mm X 9 mm Locking Screw |
| 322-2010 | Fps™ 2.0 mm X 10 mm Locking Screw |
| 322-2011 | Fps™ 2.0 mm X 11 mm Locking Screw |
| 322-2012 | Fps™ 2.0 mm X 12 mm Locking Screw |
| 322-2013 | Fps™ 2.0 mm X 13 mm Locking Screw |
| 322-2014 | Fps™ 2.0 mm X 14 mm Locking Screw |
| 322-2016 | Fps™ 2.0 mm X 16 mm Locking Screw |
| 322-2018 | Fps™ 2.0 mm X 18 mm Locking Screw |
| 322-2020 | Fps™ 2.0 mm X 20 mm Locking Screw |
| 322-2022 | Fps™ 2.0 mm X 22 mm Locking Screw |
| 322-2024 | Fps™ 2.0 mm X 24 mm Locking Screw |
| 322-2026 | Fps™ 2.0 mm X 26 mm Locking Screw |
| 322-2028 | Fps™ 2.0 mm X 28 mm Locking Screw |
| 322-2030 | Fps™ 2.0 mm X 30 mm Locking Screw |
| 322-2032 | Fps™ 2.0 mm X 32 mm Locking Screw |
| 322-2710 | Fps™ 2.7 mm X 10 mm Locking Screw |
| 322-2712 | Fps™ 2.7 mm X 12 mm Locking Screw |
| 322-2714 | Fps™ 2.7 mm X 14 mm Locking Screw |
| 322-2716 | Fps™ 2.7 mm X 16 mm Locking Screw |
| 322-2718 | Fps™ 2.7 mm X 18 mm Locking Screw |
| 322-2720 | Fps™ 2.7 mm X 20 mm Locking Screw |
| 322-2722 | Fps™ 2.7 mm X 22 mm Locking Screw |
| 322-2724 | Fps™ 2.7 mm X 24 mm Locking Screw |
| 322-2726 | Fps™ 2.7 mm X 26 mm Locking Screw |
| 322-2728 | Fps™ 2.7 mm X 28 mm Locking Screw |
| 322-2730 | Fps™ 2.7 mm X 30 mm Locking Screw |
| 322-2732 | Fps™ 2.7 mm X 32 mm Locking Screw |
| 322-2734 | Fps™ 2.7 mm X 34 mm Locking Screw |
| 322-2736 | Fps™ 2.7 mm X 36 mm Locking Screw |

OsteoMed® Lower Extremities Portfolio [continued]

OsteoMed® FPS™ Foot Plating System

| Non-Sterile Part Number | Descriptor |
|-------------------------|--|
| 322-2738 | Fps™ 2.7 mm X 38 mm Locking Screw |
| 322-2740 | Fps™ 2.7 mm X 40 mm Locking Screw |
| 322-2745 | Fps™ 2.7 mm X 45 mm Locking Screw |
| 322-2750 | Fps™ 2.7 mm X 50 mm Locking Screw |
| 322-2755 | Screw, Cortical, Locking, 2.7 Diameter, 55 L |
| 322-2760 | Screw, Cortical, Locking, 2.7 Diameter, 60 L |
| 322-2765 | Screw, Cortical, Locking, 2.7 Diameter, 65 L |
| 322-2770 | Screw, Cortical, Locking, 2.7 Diameter, 70 L |
| 322-2775 | Fps™ 2.7 mm X 75 mm Locking Screw |
| 322-3510 | Fps™ 3.5 mm X 10 mm Locking Screw |
| 322-3512 | Fps™ 3.5 mm X 12 mm Locking Screw |
| 322-3514 | Fps™ 3.5 mm X 14 mm Locking Screw |
| 322-3516 | Fps™ 3.5 mm X 16 mm Locking Screw |
| 322-3518 | Fps™ 3.5 mm X 18 mm Locking Screw |
| 322-3520 | Fps™ 3.5 mm X 20 mm Locking Screw |
| 322-3522 | Fps™ 3.5 mm X 22 mm Locking Screw |
| 322-3524 | Fps™ 3.5 mm X 24 mm Locking Screw |
| 322-3526 | Fps™ 3.5 mm X 26 mm Locking Screw |
| 322-3528 | Fps™ 3.5 mm X 28 mm Locking Screw |
| 322-3530 | Fps™ 3.5 mm X 30 mm Locking Screw |
| 322-3532 | Fps™ 3.5 mm X 32 mm Locking Screw |
| 322-3534 | Fps™ 3.5 mm X 34 mm Locking Screw |
| 322-3536 | Fps™ 3.5 mm X 36 mm Locking Screw |
| 322-3538 | Fps™ 3.5 mm X 38 mm Locking Screw |
| 322-3540 | Fps™ 3.5 mm X 40 mm Locking Screw |
| 322-3545 | Fps™ 3.5 mm X 45 mm Locking Screw |
| 322-3550 | Fps™ 3.5 mm X 50 mm Locking Screw |
| 322-3555 | Fps™ 3.5 mm X 55 mm Locking Screw |
| 322-3560 | Fps™ 3.5 mm X 60 mm Locking Screw |
| 322-3565 | Fps™ 3.5 mm X 65 mm Locking Screw |

OsteoMed® Lower Extremities Portfolio [continued]

OsteoMed® FPS™ Foot Plating System

| Non-Sterile Part Number | Descriptor |
|-------------------------|------------------------------------|
| 322-3570 | Fps™ 3.5 mm X 70 mm Locking Screw |
| 322-3575 | Fps™ 3.5 mm X 75 mm Locking Screw |
| 322-4010 | Fps™ 4.0 mm X 10 mm Locking Screw |
| 322-4012 | Fps™ 4.0 mm X 12 mm Locking Screw |
| 322-4014 | Fps™ 4.0 mm X 14 mm Locking Screw |
| 322-4016 | Fps™ 4.0 mm X 16 mm Locking Screw |
| 322-4018 | Fps™ 4.0 mm X 18 mm Locking Screw |
| 322-4020 | Fps™ 4.0 mm X 20 mm Locking Screw |
| 322-4022 | Fps™ 4.0 mm X 22 mm Locking Screw |
| 322-4024 | Fps™ 4.0 mm X 24 mm Locking Screw |
| 322-4026 | Fps™ 4.0 mm X 26 mm Locking Screw |
| 322-4028 | Fps™ 4.0 mm X 28 mm Locking Screw |
| 322-4030 | Fps™ 4.0 mm X 30 mm Locking Screw |
| 322-4032 | Fps™ 4.0 mm X 32 mm Locking Screw |
| 322-4034 | Fps™ 4.0 mm X 34 mm Locking Screw |
| 322-4036 | Fps™ 4.0 mm X 36 mm Locking Screw |
| 322-4038 | Fps™ 4.0 mm X 38 mm Locking Screw |
| 322-4040 | Fps™ 4.0 mm X 40 mm Locking Screw |
| 322-4045 | Fps™ 4.0 mm X 45 mm Locking Screw |
| 322-4050 | Fps™ 4.0 mm X 50 mm Locking Screw |
| 322-4055 | Fps™ 4.0 mm X 55 mm Locking Screw |
| 322-4060 | Fps™ 4.0 mm X 60 mm Locking Screw |
| 322-4065 | Fps™ 4.0 mm X 65 mm Locking Screw |
| 322-4070 | Fps™ 4.0 mm X 70 mm Locking Screw |
| 322-4075 | Fps™ 4.0 mm X 75 mm Locking Screw |
| 324-1000 | Fps 2.0 mm X 4 Hole Straight Plate |
| 324-1001 | Fps 2.0 mm X 5 Hole Straight Plate |
| 324-1002 | Fps 2.0 mm X 6 Hole Straight Plate |
| 324-1010 | Fps 2.0 mm X 6 Hole Right L Plate |
| 324-1011 | Fps 2.0 mm X 6 Hole Left L Plate |

OsteoMed® Lower Extremities Portfolio [continued]

OsteoMed® FPS™ Foot Plating System

| Non-Sterile Part Number | Descriptor |
|-------------------------|--|
| 324-1012 | Fps 2.0 mm X 6 Hole Oblique Right L Plate |
| 324-1013 | Fps 2.0 mm X 6 Hole Oblique Left L Plate |
| 324-1020 | Fps 2.0 mm X 7 Hole T Plate, Locking |
| 324-1021 | Fps 2.0 mm X 7 Hole T Plate, Compression |
| 324-1022 | Fps 2.0 mm X 7 Hole Oblique Right T Plate |
| 324-1023 | Fps 2.0 mm X 7 Hole Oblique Left T Plate |
| 324-1030 | Fps 2.0 mm X 6 Hole Y Plate |
| 324-1031 | Fps 2.0 mm X 7 Hole Y Plate |
| 324-1090 | Fps 2.0 mm Subcondylar Plate |
| 324-1100 | Fps 2.7 mm 4 Hole Straight Compression Pl |
| 324-1101 | Fps 2.7 mm 5 Hole Straight Compression Pl |
| 324-1102 | Fps 2.7 mm 6 Hole Straight Compression Pl |
| 324-1103 | Fps 2.7 mm 7 Hole Straight Compression Pl |
| 324-1110 | Fps 2.7 mm 6 Hole L Plate Right |
| 324-1111 | Fps 2.7 mm 6 Hole L Plate Left |
| 324-1112 | Fps 2.7 mm 6 Hole Oblique L Compression P, R |
| 324-1113 | Fps 2.7 mm 6 Hole Oblique L Compression P, L |
| 324-1190 | Fps 2.7 mm Subcondylar Plate |
| 324-1220 | Fps 2.7 mm 6 Hole T Plate |
| 324-1221 | Fps 2.7 mm 6 Hole Oblique T Plate Right |
| 324-1222 | Fps 2.7 mm 6 Hole Oblique T Plate Left |
| 324-1250 | Fps 2.7 0 mm Step Plate |
| 324-1251 | Fps 2.7 1 mm Step Plate |
| 324-1252 | Fps 2.7 2 mm Step Plate |
| 324-1253 | Fps 2.7 3 mm Step Plate |
| 324-1254 | Fps 2.7 4 mm Step Plate |
| 324-1255 | Fps 2.7 5 mm Step Plate |
| 324-1260 | Fps 2.7 mm 0 mm Wedge Plate |
| 324-1261 | Fps 2.7 mm 2.0 mm Wedge Plate |
| 324-1262 | Fps 2.7 mm 2.5 mm Wedge Plate |

OsteoMed® Lower Extremities Portfolio [continued]

OsteoMed® FPS™ Foot Plating System

| Non-Sterile Part Number | Descriptor |
|-------------------------|--|
| 324-1263 | Fps 2.7 mm 3.0 mm Wedge Plate |
| 324-1264 | Fps 2.7 mm 3.5 mm Wedge Plate |
| 324-1265 | Fps 2.7 mm 4.0 mm Wedge Plate |
| 324-1266 | Fps 2.7 mm 4.5 mm Wedge Plate |
| 324-1267 | Fps 2.7 mm 5.0 mm Wedge Plate |
| 324-1270 | Fps 2.7 mm Hook Plate Right |
| 324-1271 | Fps 2.7 mm Hook Plate Left |
| 324-1272 | Fps 2.7 mm Hook Plate, Short |
| 324-1280 | Fps 2.7 mm 1st Mpt/Mpj Plate Small Right |
| 324-1281 | Fps 2.7 mm 1st Mpt/Mpj Plate Small Left |
| 324-1284 | Fps 2.7 mm 1st Mpt/Mpj Plate Large Right |
| 324-1285 | Fps 2.7 mm 1st Mpt/Mpj Plate Large Left |
| 324-1286 | Fps 2.7 mm 1st Mpt/Mpj Primary Plate Right |
| 324-1287 | Fps 2.7 mm 1st Mpt/Mpj Primary Plate Left |
| 324-1290 | Fps 2.7 mm Lapidus Plate Right |
| 324-1291 | Fps 2.7 mm Lapidus Plate Left |
| 324-1292 | Fps 2.7 Nc Fusion Plate Large |
| 324-1293 | Fps 2.7 Nc Fusion Plate Small |
| 324-1296 | Fps 3.5/4.0 mm Medial Column Plate Small |
| 324-1297 | Fps 3.5/4.0 mm Medial Column Plate Large |
| 324-1386 | 1st Mpj Plate, Primary, Right, Modified |
| 324-1387 | 1st Mpj Plate, Primary, Left, Modified |
| 324-2204 | Fps 3.5/4.0 mm 4 Hole Straight Plate |
| 324-2206 | Fps 3.5/4.0 mm 6 Hole Straight Plate |
| 324-2208 | Fps 3.5/4.0 mm 8 Hole Straight Plate |
| 324-2210 | Fps 3.5/4.0 mm 10 Hole Straight Plate |
| 324-2212 | Fps 3.5/4.0 mm 12 Hole Straight Plate |
| 324-2220 | Fps 3.5/4.0 mm 6 Hole T Plate |
| 324-2221 | Fps 3.5/4.0 mm 8 Hole T Plate |
| 324-2222 | Fps 3.5/4.0 mm 10 Hole T Plate |

OsteoMed® Lower Extremities Portfolio [continued]

OsteoMed® FPS™ Foot Plating System

| Non-Sterile Part Number | Descriptor |
|-------------------------|---|
| 324-2223 | Fps 3.5/4.0 mm 6 Hole Oblique T Plate Right |
| 324-2224 | Fps 3.5/4.0 mm 6 Hole Oblique T Plate Left |
| 324-2225 | Plate, Oblique T, Right,3.5/4.0 mm,8 Holes |
| 324-2226 | Plate, Oblique T, Left,3.5/4.0 mm,8 Holes |
| 324-2227 | Fps 3.5/4.0 mm 10 Hole Oblique T Right |
| 324-2228 | Fps 3.5/4.0 mm 10 Hole Oblique T Left |
| 324-2240 | Fps 3.5/4.0 mm X 15 H Plate |
| 324-2241 | Fps 3.5/4.0 mm X 20 H Plate |
| 324-2242 | Fps 3.5/4.0 mm X 30 H Plate |
| 324-2250 | Fps 3.5/4.0 mm 6 mm Step Plate |
| 324-2251 | Fps 3.5/4.0 mm 8 mm Step Plate |
| 324-2252 | Fps 3.5/4.0 mm 10 mm Step Plate |
| 324-2260 | Fps 3.5 mm 2 mm Wedge Plate |
| 324-2261 | Fps 3.5 mm 4 mm Wedge Plate |
| 324-2262 | Fps 3.5 mm 6 mm Wedge Plate |
| 324-2263 | Fps 3.5 mm 8 mm Wedge Plate |
| 324-2264 | Fps 3.5 mm 10 mm Wedge Plate |
| 324-3004 | 3.5/4.0 mm 4 Hole Lcdc Fibula Plate |
| 324-3006 | 3.5/4.0 mm 6 Hole Lcdc Fibula Plate |
| 324-3008 | 3.5/4.0 mm 8 Hole Lcdc Fibula Plate |
| 324-3010 | 3.5/4.0 mm 10 Hole Lcdc Fibula Plate |
| 324-3012 | 3.5/4.0 mm 12 Hole Lcdc Fibula Plate |
| 324-3456 | Calfix Calcaneal Plate, 56 mm |
| 324-3465 | Calfix Calcaneal Plate, 65 mm |
| 324-3474 | Calfix Calcaneal Plate, 74 mm |
| 305-3518 | 3.5 mm X 18 mm Screw |
| 305-3520 | 3.5 mm X 20 mm Screw |
| 305-3522 | 3.5 mm X 22 mm Screw |
| 305-3524 | 3.5 mm X 24 mm Screw |
| 305-3526 | 3.5 mm X 26 mm Screw |

OsteoMed® Lower Extremities Portfolio [continued]

OsteoMed® FPS™ Foot Plating System

| Non-Sterile Part Number | Descriptor |
|-------------------------|---|
| 305-3528 | 3.5 mm X 28 mm Screw |
| 305-3530 | 3.5 mm X 30 mm Screw |
| 305-3532 | 3.5 mm X 32 mm Screw |
| 305-3534 | 3.5 mm X 34 mm Screw |
| 305-3536 | 3.5 mm X 36 mm Screw |
| 305-3538 | 3.5 mm X 38 mm Screw |
| 305-3540 | 3.5 mm X 40 mm Screw |
| 305-3542 | 3.5 mm X 42 mm Screw |
| 305-3544 | 3.5 mm X 44 mm Screw |
| 305-3546 | 3.5 mm X 46 mm Screw |
| 305-3548 | 3.5 mm X 48 mm Screw |
| 305-3550 | 3.5 mm X 50 mm Screw |
| 305-3555 | 3.5 mm X 55 mm Screw |
| 313-3518 | 3.5 mm X 18 mm Locking Screw, Calcaneal Plt |
| 313-3520 | 3.5 mm X 20 mm Locking Screw, Calcaneal Plt |
| 313-3522 | 3.5 mm X 22 mm Locking Screw, Calcaneal Plt |
| 313-3524 | 3.5 mm X 24 mm Locking Screw, Calcaneal Plt |
| 313-3526 | 3.5 mm X 26 mm Locking Screw, Calcaneal Plt |
| 313-3528 | 3.5 mm X 28 mm Locking Screw, Calcaneal Plt |
| 313-3530 | 3.5 mm X 30 mm Locking Screw, Calcaneal Plt |
| 313-3532 | 3.5 mm X 32 mm Locking Screw, Calcaneal Plt |
| 313-3534 | 3.5 mm X 34 mm Locking Screw, Calcaneal Plt |
| 313-3536 | 3.5 mm X 36 mm Locking Screw, Calcaneal Plt |
| 313-3538 | 3.5 mm X 38 mm Locking Screw, Calcaneal Plt |
| 313-3540 | 3.5 mm X 40 mm Locking Screw, Calcaneal Plt |
| 313-3542 | 3.5 mm X 42 mm Locking Screw, Calcaneal Plt |
| 313-3544 | 3.5 mm X 44 mm Locking Screw, Calcaneal Plt |
| 313-3546 | 3.5 mm X 46 mm Locking Screw, Calcaneal Plt |
| 313-3548 | 3.5 mm X 48 mm Locking Screw, Calcaneal Plt |
| 313-3550 | 3.5 mm X 50 mm Locking Screw, Calcaneal Plt |

OsteoMed® Lower Extremities Portfolio [continued]

OsteoMed® FPS™ Foot Plating System

| Non-Sterile Part Number | Descriptor |
|-------------------------|--|
| 313-3555 | 3.5 mm X 55 mm Locking Screw, Calcaneal Plt |
| 314-1456 | Calfix Plate, 56 mm |
| 314-1465 | Calfix Plate, 65 mm |
| 314-1474 | Calfix Plate, 74 mm |
| 336-2700 | 2.7 mm - 5 Hole Straight Plate |
| 336-2701 | 2.7 mm - 10 Hole Y Compression Plate |
| 336-2702 | 2.7 mm - 7 Hole T Compression Plate |
| 336-2703 | 2.7 mm - 4 Hole Straight Compression Plate |
| 336-2704 | 2.7 mm - 2 Hole Straight Compression Plate |
| 336-2710 | 2.7 mm - 10 Hole Subcondylar Plate |
| 336-2720 | 2.7 mm - 0 mm Step Plate (H-Plate) |
| 336-2721 | 2.7 mm - 1 mm Step Plate |
| 336-2722 | 2.7 mm - 2 mm Step Plate |
| 336-2723 | 2.7 mm - 3 mm Step Plate |
| 336-2724 | 2.7 mm - 4 mm Step Plate |
| 336-2725 | 2.7 mm - 5 mm Step Plate |
| 336-2730 | 2.7 mm - 0 mm Wedge Plate (Closing Wedge) |
| 336-2733 | 2.7 mm - 3 mm Wedge Plate |
| 336-2734 | 2.7 mm - 4 mm Wedge Plate |
| 336-2735 | 2.7 mm - 5 mm Wedge Plate |
| 336-2736 | 2.7 mm - 6 mm Wedge Plate |
| 336-2740 | 2.7 mm - 6 Hole Hook Plate, Left |
| 336-2744 | 2.7 mm - 4 Hole Hook Plate |
| 336-2745 | 2.7 mm - 6 Hole Hook Plate, Right |
| 336-2750 | 2.7 mm - 1st Mtp Primary Plate, Left 0° |
| 336-2751 | 2.7 mm - 1st Mtp Primary Plate, Left 5° |
| 336-2752 | 2.7 mm - 1st Mtp Primary Plate Left - 10° |
| 336-2753 | 2.7 mm - 1st Mtp Primary Plate, Left 15° |
| 336-2754 | 2.7 mm - 1st Mtp Primary Plate, 6h, 0°/5°, L |
| 336-2755 | 2.7 mm - 1st Mtp Primary Plate, Right 0° |

OsteoMed® Lower Extremities Portfolio [continued]

OsteoMed® FPS™ Foot Plating System

| Non-Sterile Part Number | Descriptor |
|-------------------------|--|
| 336-2756 | 2.7 mm - 1st Mtp Primary Plate, Right 5° |
| 336-2757 | 2.7 mm - 1st Mtp Primary Plate Right, 10° |
| 336-2758 | 2.7 mm -1st Mtp Primary Plate, Right 15° |
| 336-2759 | 2.7 mm 1st Mtp Primary Plate, 6h, 0°/5°, R |
| 336-2760 | 2.7 mm - Lapidus Plate |
| 336-2770 | 2.7 mm 1st Transfixation Mtp Plate, L,0°/10° |
| 336-2771 | 2.7 mm 1st Transfixation Mtp Plate, L,5°/10° |
| 336-2772 | 2.7 mm 1st Transfixation Mtp Plate, L10°/10° |
| 336-2773 | 2.7 mm 1st Transfixation Mtp Plate, L15°/10° |
| 336-2775 | 2.7 mm 1st Transfixation Mtp Plate, R,0°/10° |
| 336-2776 | 2.7 mm 1st Transfixation Mtp Plate, R,5°/10° |
| 336-2777 | 2.7 mm 1st Transfixation Mtp Plate, R10°/10° |
| 336-2778 | 2.7 mm 1st Transfixation Mtp Plate, R15°/10° |
| 336-2779 | 2.7 mm 1st Mtp Fusion Plate, 5 Hole Narrow |
| 336-2780 | 2.7 mm 1st Mtp Fusion Plate, 5 Hole Narrow |
| 336-3500 | 3.5 mm - 15 mm H Plate |
| 336-3501 | 3.5 mm - 20 mm H Plate |
| 336-3502 | 3.5 mm - 30 mm H Plate |
| 336-3506 | 3.5 mm - 6 mm Step Plate |
| 336-3508 | 3.5 mm - 8 mm Step Plate |
| 336-3510 | 3.5 mm - 10 mm Step Plate |
| 336-3524 | 3.5 mm - 4 mm Wedge Plate |
| 336-3526 | 3.5 mm - 6 mm Wedge Plate |
| 336-3528 | 3.5 mm - 8 mm Wedge Plate |
| 336-3530 | 3.5 mm - 56 mm Calcaneal Mesh Plate |
| 336-3531 | 3.5 mm - 65 mm Calcaneal Mesh Plate |
| 336-3540 | 3.5/4.0 mm - 58 mm Mi Calcaneal Plate Left |
| 336-3541 | 3.5 mm - 65 mm Mi Calcaneal Plate Left |
| 336-3545 | 3.5/4.0 mm - 58 mm Mi Calcaneal Plate Right |
| 336-3546 | 3.5 mm - 65 mm Mi Calcaneal Plate Right |

OsteoMed® Lower Extremities Portfolio [continued]

OsteoMed® FPS™ Foot Plating System

| Non-Sterile Part Number | Descriptor |
|-------------------------|---|
| 336-3554 | 3.5 mm 1st Mtp Primary Plate, 6 Hole, 0°/5° |
| 336-3559 | 3.5 mm 1st Mtp Primary Plate, 6 Hole, 0°/5° |
| 336-3560 | 3.5 mm - Lapidus Plate |
| 336-3561 | 3.5/4.0 mm - Lapidus Plate, Dorsalmedial |
| 336-3562 | 3.5/4.0 Medial Column Fusion Plate Short |
| 336-3563 | 3.5/4.0 Medial Column Fusion Plate Mediu |
| 336-3564 | 3.5/4.0 Medial Column Fusion Plate Long |
| 336-3570 | 3.5/4.0 mm - 3 Hole Straight Compression PI |
| 336-3571 | 3.5/4.0 mm - 2 Hole Straight Compression PI |
| 336-3572 | 3.5/4.0 mm - 2 Hole Straight Compression PI |
| 337-2710 | 2.7 X 10 mm Fully Threaded Screw, T8 |
| 337-2712 | 2.7 X 12 mm Fully Threaded Screw, T8 |
| 337-2714 | 2.7 X 14 mm Fully Threaded Screw, T8 |
| 337-2716 | 2.7 X 16 mm Fully Threaded Screw, T8 |
| 337-2718 | 2.7 X 18 mm Fully Threaded Screw, T8 |
| 337-2720 | 2.7 X 20 mm Fully Threaded Screw, T8 |
| 337-2722 | 2.7 X 22 mm Fully Threaded Screw, T8 |
| 337-2724 | 2.7 X 24 mm Fully Threaded Screw, T8 |
| 337-2726 | 2.7 X 26 mm Fully Threaded Screw, T8 |
| 337-2728 | 2.7 X 28 mm Fully Threaded Screw, T8 |
| 337-2730 | 2.7 X 30 mm Fully Threaded Screw, T8 |
| 337-2799 | 2.7 mm Washer |
| 337-3010 | 3.0 X 10 mm Double-Lead Non-Locking Screw |
| 337-3012 | 3.0 X 12 mm Double-Lead NI Screw |
| 337-3014 | 3.0 X 14 mm Double-Lead Non-Locking Screw |
| 337-3016 | 3.0 X 16 mm Double-Lead NI Screw |
| 337-3018 | 3.0 X 18 mm Double-Lead Non-Locking Screw |
| 337-3020 | 3.0 X 20 mm Double-Lead NI Screw |
| 337-3022 | 3.0 X 22 mm Double-Lead Non-Locking Screw |
| 337-3026 | 3.0 X 26 mm Double-Lead Non-Locking Screw |

OsteoMed® Lower Extremities Portfolio [continued]

OsteoMed® FPS™ Foot Plating System

| Non-Sterile Part Number | Descriptor |
|-------------------------|---|
| 337-3030 | 3.0 X 30 mm Double-Lead Non-Locking Screw |
| 337-3510 | 3.5 X 10 mm Double-Lead Non-Locking Screw |
| 337-3512 | 3.5 X 12 mm Double-Lead Non-Locking Screw |
| 337-3514 | 3.5 X 14 mm Double-Lead Non-Locking Screw |
| 337-3516 | 3.5 X 16 mm Double-Lead Non-Locking Screw |
| 337-3518 | 3.5 X 18 mm Double-Lead Non-Locking Screw |
| 337-3520 | 3.5 X 20 mm Double Lead Non-Locking Screw |
| 337-3522 | 3.5 X 22 mm Double-Lead Non-Locking Screw |
| 337-3524 | 3.5 X 24 mm Double-Lead Non-Locking Screw |
| 337-3526 | 3.5 X 26 mm Double-Lead Non-Locking Screw |
| 337-3528 | 3.5 X 28 mm Double-Lead Non-Locking Screw |
| 337-3530 | 3.5 X 30 mm Double-Lead Non-Locking Screw |
| 337-3599 | 3.5/4.0 mm Washer |
| 337-4020 | 4.0 X 20 mm Double-Lead Non-Locking Screw |
| 337-4022 | 4.0 X 22 mm Double-Lead Non-Locking Screw |
| 337-4024 | 4.0 X 24 mm Double-Lead Non-Locking Screw |
| 337-4026 | 4.0 X 26 mm Double-Lead Non-Locking Screw |
| 337-4028 | 4.0 X 28 mm Double-Lead Non-Locking Screw |
| 337-4030 | 4.0 X 30 mm Double-Lead Non-Locking Screw |
| 337-4032 | 4.0 X 32 mm Double-Lead Non-Locking Screw |
| 337-4034 | 4.0 X 34 mm Double-Lead Non-Locking Screw |
| 337-4036 | 4.0 X 36 mm Double-Lead Non-Locking Screw |
| 337-4038 | 4.0 X 38 mm Double-Lead Non-Locking Screw |
| 337-4040 | 4.0 X 40 mm Double-Lead, Non-Locking Scr |
| 337-4042 | 4.0 X 42 mm Double-Lead, Non-Locking Scr |
| 337-4044 | 4.0 X 44 mm Double-Lead, Non-Locking Scr |
| 337-4045 | 4.0 X 45 mm Double-Lead, Non-Locking Scr |
| 337-4046 | 4.0 X 46 mm Double-Lead, Non-Locking Scr |
| 337-4048 | 4.0 X 48 mm Double-Lead, Non-Locking Scr |
| 337-4050 | 4.0 X 50 mm Double-Lead Non-Locking Screw |

OsteoMed® Lower Extremities Portfolio [continued]

OsteoMed® FPS™ Foot Plating System

| Non-Sterile Part Number | Descriptor |
|-------------------------|---|
| 337-4055 | 4.0 X 55 mm Double-Lead Non-Locking Screw |
| 337-4060 | 4.0 X 60 mm Double-Lead Non-Locking Screw |
| 338-2710 | 2.7 X 10 mm Fully Threaded Angled Locking |
| 338-2712 | 2.7 X 12 mm Fully Threaded Angled Locking |
| 338-2714 | 2.7 X 14 mm Fully Threaded Angled Locking |
| 338-2716 | 2.7 X 16 mm Fully Threaded Angled Locking |
| 338-2718 | 2.7 X 18 mm Fully Threaded Angled Locking |
| 338-2720 | 2.7 X 20 mm Fully Threaded Angled Locking |
| 338-2722 | 2.7 X 22 mm Fully Threaded Angled Locking |
| 338-2724 | 2.7 X 24 mm Fully Threaded Angled Locking |
| 338-2726 | 2.7 X 26 mm Fully Threaded Angled Locking |
| 338-2728 | 2.7 X 28 mm Fully Threaded Angled Locking |
| 338-2730 | 2.7 X 30 mm Fully Threaded Angled Locking |
| 338-3510 | 3.5 X 10 mm Double-Lead Locking T15 |
| 338-3512 | 3.5 X 12 mm Double-Lead Locking T15 |
| 338-3514 | 3.5 X 14 mm Double-Lead Locking T15 |
| 338-3516 | 3.5 X 16 mm Double-Lead Locking T15 |
| 338-3518 | 3.5 X 18 mm Double-Lead Locking T15 |
| 338-3520 | 3.5 X 20 mm Double-Lead Locking T15 |
| 338-3522 | 3.5 X 22 mm Double-Lead Locking T15 |
| 338-3524 | 3.5 X 24 mm Double-Lead Locking T15 |
| 338-3526 | 3.5 X 26 mm Double-Lead Locking T15 |
| 338-3528 | 3.5 X 28 mm Double-Lead Locking T15 |
| 338-3530 | 3.5 X 30 mm Double-Lead Locking T15 |
| 338-4020 | 4.0 X 20 mm Double-Lead Locking T15 |
| 338-4022 | 4.0 X 22 mm Double-Lead Locking T15 |
| 338-4024 | 4.0 X 24 mm Double-Lead Locking T15 |
| 338-4026 | 4.0 X 26 mm Double-Lead Locking T15 |
| 338-4028 | 4.0 X 28 mm Double-Lead Locking T15 |
| 338-4030 | 4.0 X 30 mm Double-Lead Locking T15 |

OsteoMed® Lower Extremities Portfolio [continued]

OsteoMed® FPS™ Foot Plating System

| Non-Sterile Part Number | Descriptor |
|-------------------------|--|
| 338-4032 | 4.0 X 32 mm Double-Lead Locking T15 |
| 338-4034 | 4.0 X 34 mm Double-Lead Locking T15 |
| 338-4036 | 4.0 X 36 mm Double-Lead Locking T15 |
| 338-4038 | 4.0 X 38 mm Double-Lead Locking T15 |
| 338-4040 | 4.0 X 40 mm Double-Lead,Locking Screw, T |
| 338-4042 | 4.0 X 42 mm Double-Lead,Locking Screw, T |
| 338-4044 | 4.0 X 44 mm Double-Lead,Locking Screw, T |
| 338-4045 | 4.0 X 45 mm Double-Lead,Locking Screw, T |

OsteoMed® Lower Extremities Portfolio [continued]

OsteoMed® ExtremiFuse™ Hammertoe Fixation System

| Non-Sterile Part Number | Descriptor |
|----------------------------|----------------------------------|
| 319-0020 | 2.0 mm X 0° ExtremiFuse Implant |
| 319-0024 | 2.4 mm X 0° ExtremiFuse Implant |
| 319-0030 | 3.0 mm X 0° ExtremiFuse Implant |
| 319-1020 | 2.0 mm X 10° ExtremiFuse Implant |
| 319-1024 | 2.4 mm X 10° ExtremiFuse Implant |
| 319-1030 | 3.0 mm X 10° ExtremiFuse Implant |

OsteoMed® Lower Extremities Portfolio [continued]

OsteoMed® ExtremiLock™ Ankle Plating System

| Non-Sterile Part Number | Descriptor |
|-------------------------|--------------------------------------|
| 336-0002 | 2h 1/3 Tubular Plate |
| 336-0003 | 3h 1/3 Tubular Plate |
| 336-0004 | 4h 1/3 Tubular Plate |
| 336-0005 | 5h 1/3 Tubular Plate |
| 336-0006 | 6h 1/3 Tubular Plate |
| 336-0007 | 7h 1/3 Tubular Plate |
| 336-0008 | 8h 1/3 Tubular Plate |
| 336-0009 | 9h 1/3 Tubular Plate |
| 336-0010 | 10h 1/3 Tubular Plate |
| 336-0011 | 11h 1/3 Tubular Plate |
| 336-0012 | 12h 1/3 Tubular Plate |
| 336-0013 | 13h 1/3 Tubular Plate |
| 336-0014 | 14h 1/3 Tubular Plate |
| 336-1003 | 3h Anterior Tibia Plate |
| 336-1006 | 6h Anterior Tibia Plate |
| 336-2107 | 7h Anterolateral Tibia Plate, Left |
| 336-2110 | 10h Anterolateral Tibia Plate, Left |
| 336-2113 | 13h Anterolateral Tibia Plate, Left |
| 336-2115 | 15h Anterolateral Tibia Plate, Left |
| 336-2117 | 17h Anterolateral Tibia Plate, Left |
| 336-2119 | 19h Anterolateral Tibia Plate, Left |
| 336-2121 | 21h Anterolateral Tibia Plate, Left |
| 336-2123 | 23h Anterolateral Tibia Plate, Left |
| 336-2125 | 25h Anterolateral Tibia Plate, Left |
| 336-2127 | 27h Anterolateral Tibia Plate, Left |
| 336-2207 | 7h Anterolateral Tibia Plate, Right |
| 336-2210 | 10h Anterolateral Tibia Plate, Right |
| 336-2213 | 13h Anterolateral Tibia Plate, Right |
| 336-2215 | 15h Anterolateral Tibia Plate, Right |
| 336-2217 | 17h Anterolateral Tibia Plate, Right |

OsteoMed® Lower Extremities Portfolio [continued]

OsteoMed® ExtremiLock™ Ankle Plating System

| Non-Sterile Part Number | Descriptor |
|-------------------------|--------------------------------------|
| 336-2219 | 19h Anterolateral Tibia Plate, Right |
| 336-2221 | 21h Anterolateral Tibia Plate, Right |
| 336-2223 | 23h Anterolateral Tibia Plate, Right |
| 336-2225 | 25h Anterolateral Tibia Plate, Right |
| 336-2227 | 27h Anterolateral Tibia Plate, Right |
| 336-4002 | 2h Medial Tibia Plate |
| 336-4004 | 4h Medial Tibia Plate |
| 336-4006 | 6h Medial Tibia Plate |
| 336-4008 | 8h Medial Tibia Plate |
| 336-4010 | 10h Medial Tibia Plate |
| 336-4012 | 12h Medial Tibia Plate |
| 336-5102 | 2h Fibula Plate, Left |
| 336-5103 | 3h Fibula Plate, Left |
| 336-5104 | 4h Fibula Plate, Left |
| 336-5105 | 5h Fibula Plate, Left |
| 336-5106 | 6h Fibula Plate, Left |
| 336-5107 | 7h Fibula Plate, Left |
| 336-5108 | 8h Fibula Plate, Left |
| 336-5109 | 9h Fibula Plate, Left |
| 336-5110 | 10h Fibula Plate, Left |
| 336-5111 | 11h Fibula Plate, Left |
| 336-5112 | 12h Fibula Plate, Left |
| 336-5113 | 13h Fibula Plate, Left |
| 336-5114 | 14h Fibula Plate, Left |
| 336-5202 | 2h Fibula Plate, Right |
| 336-5203 | 3h Fibula Plate, Right |
| 336-5204 | 4h Fibula Plate, Right |
| 336-5205 | 5h Fibula Plate, Right |
| 336-5206 | 6h Fibula Plate, Right |
| 336-5207 | 7h Fibula Plate, Right |

OsteoMed® Lower Extremities Portfolio [continued]

OsteoMed® ExtremiLock™ Ankle Plating System

| Non-Sterile Part Number | Descriptor |
|-------------------------|-------------------------------|
| 336-5208 | 8h Fibula Plate, Right |
| 336-5209 | 9h Fibula Plate, Right |
| 336-5210 | 10h Fibula Plate, Right |
| 336-5211 | 11h Fibula Plate, Right |
| 336-5212 | 12h Fibula Plate, Right |
| 336-5213 | 13h Fibula Plate, Right |
| 336-5214 | 14h Fibula Plate, Right |
| 336-6005 | 5h Hook Plate, Universal |
| 336-6007 | 7h Hook Plate, Universal |
| 336-6102 | 2h Distal Fibula Plate, Left |
| 336-6103 | 3h Distal Fibula Plate, Left |
| 336-6104 | 4h Distal Fibula Plate, Left |
| 336-6105 | 5h Distal Fibula Plate, Left |
| 336-6106 | 6h Distal Fibula Plate, Left |
| 336-6107 | 7h Distal Fibula Plate, Left |
| 336-6108 | 8h Distal Fibula Plate, Left |
| 336-6109 | 9h Distal Fibula Plate, Left |
| 336-6110 | 10h Distal Fibula Plate, Left |
| 336-6111 | 11h Distal Fibula Plate, Left |
| 336-6112 | 12h Distal Fibula Plate, Left |
| 336-6113 | 13h Distal Fibula Plate, Left |
| 336-6114 | 14h Distal Fibula Plate, Left |
| 336-6202 | 2h Distal Fibula Plate, Right |
| 336-6203 | 3h Distal Fibula Plate, Right |
| 336-6204 | 4h Distal Fibula Plate, Right |
| 336-6205 | 5h Distal Fibula Plate, Right |
| 336-6206 | 6h Distal Fibula Plate, Right |
| 336-6207 | 7h Distal Fibula Plate, Right |
| 336-6208 | 8h Distal Fibula Plate, Right |
| 336-6209 | 9h Distal Fibula Plate, Right |

OsteoMed® Lower Extremities Portfolio [continued]

OsteoMed® ExtremiLock™ Ankle Plating System

| Non-Sterile Part Number | Descriptor |
|-------------------------|--|
| 336-6210 | 10h Distal Fibula Plate, Right |
| 336-6211 | 11h Distal Fibula Plate, Right |
| 336-6212 | 12h Distal Fibula Plate, Right |
| 336-6213 | 13h Distal Fibula Plate, Right |
| 336-6214 | 14h Distal Fibula Plate, Right |
| 336-7103 | 3h Anterior Fusion Plate, Left |
| 336-7105 | 5h Anterior Fusion Plate, Left |
| 336-7107 | 7h Anterior Fusion Plate, Left |
| 336-7203 | 3h Anterior Fusion Plate, Right |
| 336-7205 | 5h Anterior Fusion Plate, Right |
| 336-7207 | 7h Ankle Fusion Plate, Right |
| 336-8003 | 3h Posterior Fusion Plate, Universal |
| 336-8006 | 6h Posterior Fusion Plate, Universal |
| 336-9001 | Lateral Tt Fusion Plate, Ankle |
| 336-9002 | Lateral Ttc Fusion Plate, Ankle |
| 337-3532 | 3.5 mm X 32 mm, Double Lead, Non-Locking Scr |
| 337-3534 | 3.5 mm X 34 mm, Double Lead, Non-Locking Scr |
| 337-3536 | 3.5 mm X 36 mm, Double Lead, Non-Locking Scr |
| 337-3538 | 3.5 mm X 38 mm, Double Lead, Non-Locking Scr |
| 337-3540 | 3.5 mm X 40 mm, Double Lead, Non-Locking Scr |
| 337-3542 | 3.5 mm X 42 mm, Double Lead, Non-Locking Scr |
| 337-3544 | 3.5 mm X 44 mm, Double Lead, Non-Locking Scr |
| 337-3546 | 3.5 mm X 46 mm, Double Lead, Non-Locking Scr |
| 337-3548 | 3.5 mm X 48 mm, Double Lead, Non-Locking Scr |
| 337-3550 | 3.5 mm X 50 mm, Double Lead, Non-Locking Scr |
| 337-3552 | 3.5 mm X 52 mm, Double Lead, Non-Locking Scr |
| 337-3554 | 3.5 mm X 54 mm, Double Lead, Non-Locking Scr |
| 337-3555 | 3.5 mm X 55 mm, Double Lead, Non-Locking Scr |
| 337-3556 | 3.5 mm X 56 mm, Double Lead, Non-Locking Scr |
| 337-3558 | 3.5 mm X 58 mm, Double Lead, Non-Locking Scr |

OsteoMed® Lower Extremities Portfolio [continued]

OsteoMed® ExtremiLock™ Ankle Plating System

| Non-Sterile Part Number | Descriptor |
|-------------------------|--|
| 337-3560 | 3.5 mm X 60 mm, Double Lead, Non-Locking Scr |
| 337-3565 | 3.5 mm X 65 mm, Double Lead, Non-Locking Scr |
| 337-3570 | 3.5 mm X 70 mm, Double Lead, Non-Locking Scr |
| 337-3575 | 3.5 mm X 75 mm, Double Lead, Non-Locking Scr |
| 337-3580 | 3.5 mm X 80 mm, Double Lead, Non-Locking Scr |
| 337-3585 | 3.5 mm X 85 mm, Double Lead, Non-Locking Scr |
| 337-3590 | 3.5 mm X 90 mm, Double Lead, Non-Locking Scr |
| 337-3595 | 3.5 mm X 95 mm, Double Lead, Non-Locking Scr |
| 337-3600 | 3.5 mm X 100 mm, Double Lead, Non-Locking Sc |
| 337-3605 | 3.5 mm X 105 mm, Double Lead, Non-Locking Sc |
| 337-3610 | 3.5 mm X 110 mm, Double Lead, Non-Locking Sc |
| 337-4010 | 4.0 mm X 10 mm, Double Lead, Non-Locking Scr |
| 337-4012 | 4.0 mm X 12 mm, Double Lead, Non-Locking Scr |
| 337-4014 | 4.0 mm X 14 mm, Double Lead, Non-Locking Scr |
| 337-4016 | 4.0 mm X 16 mm, Double Lead, Non-Locking Scr |
| 337-4018 | 4.0 mm X 18 mm, Double Lead, Non-Locking Scr |
| 337-4052 | 4.0 mm X 52 mm, Double-Lead Non-Locking Scr |
| 337-4054 | 4.0 mm X 54 mm, Double-Lead Non-Locking Scr |
| 337-4056 | 4.0 mm X 56 mm, Double-Lead Non-Locking Scr |
| 337-4058 | 4.0 mm X 58 mm, Double-Lead Non-Locking Scr |
| 337-4065 | 4.0 mm X 65 mm, Double Lead, Non-Locking Sc |
| 337-4070 | 4.0 mm X 70 mm, Double Lead, Non-Locking Scr |
| 337-4075 | 4.0 mm X 75 mm, Double Lead, Non-Locking Scr |
| 337-4080 | 4.0 mm X 80 mm, Double Lead, Non-Locking Scr |
| 337-4085 | 4.0 mm X 85 mm, Double Lead, Non-Locking Scr |
| 337-4090 | 4.0 mm X 90 mm, Double Lead, Non-Locking Scr |
| 337-4095 | 4.0 mm X 95 mm, Double Lead, Non-Locking Scr |
| 337-4100 | 4.0 mm X 100 mm, Double Lead, Non-Locking Sc |
| 337-4105 | 4.0 mm X 105 mm, Double Lead, Non-Locking Sc |
| 337-4110 | 4.0 mm X 110 mm, Double Lead, Non-Locking Sc |

OsteoMed® Lower Extremities Portfolio [continued]

OsteoMed® ExtremiLock™ Ankle Plating System

| Non-Sterile Part Number | Descriptor |
|-------------------------|---|
| 338-3532 | 3.5 mm X 32 mm, Double Lead, Locking Screw, T |
| 338-3534 | 3.5 mm X 34 mm, Double Lead, Locking Screw, T |
| 338-3536 | 3.5 mm X 36 mm, Double Lead, Locking Screw, T |
| 338-3538 | 3.5 mm X 38 mm, Double Lead, Locking Screw, T |
| 338-3540 | 3.5 mm X 40 mm, Double Lead, Locking Screw, T |
| 338-3542 | 3.5 mm X 42 mm, Double Lead, Locking Screw, T |
| 338-3544 | 3.5 mm X 44 mm, Double Lead, Locking Screw, T |
| 338-3546 | 3.5 mm X 46 mm, Double Lead, Locking Screw, T |
| 338-3548 | 3.5 mm X 48 mm, Double Lead, Locking Screw, T |
| 338-3550 | 3.5 mm X 50 mm, Double Lead, Locking Screw, T |
| 338-3552 | 3.5 mm X 52 mm, Double Lead, Locking Screw, T |
| 338-3555 | 3.5 mm X 55 mm, Double Lead, Locking Screw, T |
| 338-3556 | 3.5 mm X 56 mm, Double Lead, Locking Screw, T |
| 338-3560 | 3.5 mm X 60 mm, Double Lead, Locking Screw, T |
| 338-3565 | 3.5 mm X 65 mm, Double Lead, Locking Screw, T |
| 338-3570 | 3.5 mm X 70 mm, Double Lead, Locking Screw, T |
| 338-3575 | 3.5 mm X 75 mm, Double Lead, Locking Screw, T |
| 338-3580 | 3.5 mm X 80 mm, Double Lead, Locking Screw, T |
| 338-3585 | 3.5 mm X 85 mm, Double Lead, Locking Screw, T |
| 338-3590 | 3.5 mm X 90 mm, Double Lead, Locking Screw, T |
| 338-3595 | 3.5 mm X 95 mm, Double Lead, Locking Screw, T |
| 338-3600 | 3.5 mm X 100 mm, Double Lead, Locking Screw, |
| 338-3605 | 3.5 mm X 105 mm, Double Lead, Locking Screw, |
| 338-3610 | 3.5 mm X 110 mm, Double Lead, Locking Screw, |
| 338-4010 | 4.0 mm X 10 mm, Double Lead, Locking Screw, T |
| 338-4012 | 4.0 mm X 12 mm, Double Lead, Locking Screw, T |
| 338-4014 | 4.0 mm X 14 mm, Double Lead, Locking Screw, T |
| 338-4016 | 4.0 mm X 16 mm, Double Lead, Locking Screw, T |
| 338-4018 | 4.0 mm X 18 mm, Double Lead, Locking Screw, T |
| 338-4046 | 4.0 X 46 mm, Double Lead, Locking Screw, T |

OsteoMed® Lower Extremities Portfolio [continued]

OsteoMed® ExtremiLock™ Ankle Plating System

| Non-Sterile Part Number | Descriptor |
|-------------------------|---|
| 338-4048 | 4.0 X 48 mm Double Lead, Locking Screw, T |
| 338-4050 | 4.0 mm X 50 mm, Double Lead, Locking Screw, T |
| 338-4055 | 4.0 mm X 55 mm, Double Lead, Locking Screw, T |
| 338-4060 | 4.0 mm X 60 mm, Double Lead, Locking Screw, T |
| 338-4065 | 4.0 mm X 65 mm, Double Lead, Locking Screw, T |
| 338-4070 | 4.0 mm X 70 mm, Double Lead, Locking Screw, T |
| 338-4075 | 4.0 mm X 75 mm, Double Lead, Locking Screw, T |
| 338-4080 | 4.0 mm X 80 mm, Double Lead, Locking Screw, T |
| 338-4085 | 4.0 mm X 85 mm, Double Lead, Locking Screw, T |
| 338-4090 | 4.0 mm X 90 mm, Double Lead, Locking Screw, T |
| 338-4095 | 4.0 mm X 95 mm, Double Lead, Locking Screw, T |
| 338-4100 | 4.0 mm X 100 mm, Double Lead, Locking Screw, |
| 338-4105 | 4.0 mm X 105 mm, Double Lead, Locking Screw, |
| 338-4110 | 4.0 mm X 110 mm, Double Lead, Locking Screw, |
| 339-4010 | 4.0 mm X 10 mm, Double Lead, Non-Locking Scr |
| 339-4012 | 4.0 mm X 12 mm, Double Lead, Non-Locking Scr |
| 339-4014 | 4.0 mm X 14 mm, Double Lead, Non-Locking Scr |
| 339-4016 | 4.0 mm X 16 mm, Double Lead, Non-Locking Scr |
| 339-4018 | 4.0 mm X 18 mm, Double Lead, Non-Locking Scr |
| 339-4020 | 4.0 mm X 20 mm, Double Lead, Non-Locking Scr |
| 339-4022 | 4.0 mm X 22 mm, Double Lead, Non-Locking Scr |
| 339-4024 | 4.0 mm X 24 mm, Double Lead, Non-Locking Scr |
| 339-4026 | 4.0 mm X 26 mm, Double Lead, Non-Locking Scr |
| 339-4028 | 4.0 mm X 28 mm, Double Lead, Non-Locking Scr |
| 339-4030 | 4.0 mm X 30 mm, Double Lead, Non-Locking Scr |
| 339-4032 | 4.0 mm X 32 mm, Double Lead, Non-Locking Scr |
| 339-4034 | 4.0 mm X 34 mm, Double Lead, Non-Locking Scr |
| 339-4036 | 4.0 mm X 36 mm, Double Lead, Non-Locking Scr |
| 339-4038 | 4.0 mm X 38 mm, Double Lead, Non-Locking Scr |
| 339-4040 | 4.0 mm X 40 mm, Double Lead, Non-Locking Scr |

OsteoMed® Lower Extremities Portfolio [continued]

OsteoMed® ExtremiLock™ Ankle Plating System

| Non-Sterile Part Number | Descriptor |
|-------------------------|---|
| 339-4042 | 4.0 mm X 42 mm, Double Lead, Non-Locking Scr |
| 339-4044 | 4.0 mm X 44 mm, Double Lead, Non-Locking Scr |
| 339-4046 | 4.0 mm X 46 mm, Double Lead, Non-Locking Scr |
| 339-4048 | 4.0 mm X 48 mm, Double Lead, Non-Locking Scr |
| 339-4052 | 4.0 mm X 52 mm, Double Lead, Non-Locking Scr |
| 339-4056 | 4.0 mm X 56 mm, Double Lead, Non-Locking Scr |
| 339-4060 | 4.0 mm X 60 mm, Double Lead, Non-Locking Scr |
| 339-4064 | 4.0 mm X 64 mm Non-Locking Ft Cann T15 |
| 339-4068 | 4.0 mm X 68 mm Non-Locking Ft Cann T15 |
| 339-4072 | 4.0 mm X 72 mm Non-Locking Ft Cann T15 |
| 339-4076 | 4.0 mm X 76 mm Non-Locking Ft Cann T15 |
| 339-4080 | 4.0 mm X 80 mm Non-Locking Ft Cann T15 |
| 342-2710 | 2.7 mm X 10 mm, Double Lead, Non-Locking Scr |
| 342-2712 | 2.7 mm X 12 mm, Double Lead, Non-Locking Scr |
| 342-2714 | 2.7 mm X 14 mm, Double Lead, Non-Locking Scr |
| 342-2716 | 2.7 mm X 16 mm, Double Lead, Non-Locking Scr |
| 342-2718 | 2.7 mm X 18 mm, Double Lead, Non-Locking Scr |
| 342-2720 | 2.7 mm X 20 mm, Double Lead, Non-Locking Scr |
| 342-2722 | 2.7 mm X 22 mm, Double Lead, Non-Locking Scr |
| 342-2724 | 2.7 mm X 24 mm, Double Lead, Non-Locking Scr |
| 342-2726 | 2.7 mm X 26 mm, Double Lead, Non-Locking Scr |
| 342-2728 | 2.7 mm X 28 mm, Double Lead, Non-Locking Scr |
| 342-2730 | 2.7 mm X 30 mm, Double Lead, Non-Locking Scr |
| 342-2732 | 2.7 mm X 32 mm, Double Lead, Non-Locking Scr |
| 342-2734 | 2.7 mm X 34 mm, Double Lead, Non-Locking Scr |
| 342-2736 | 2.7 mm X 36 mm, Double Lead, Non-Locking Scr |
| 342-2738 | 2.7 mm X 38 mm, Double Lead, Non-Locking Scr |
| 342-2740 | 2.7 mm X 40 mm, Double Lead, Non-Locking Scr |
| 343-2710 | 2.7 mm X 10 mm, Double Lead, Locking Screw, T |
| 343-2712 | 2.7 mm X 12 mm, Double Lead, Locking Screw, T |

OsteoMed® Lower Extremities Portfolio [continued]

OsteoMed® ExtremiLock™ Ankle Plating System

| Non-Sterile Part Number | Descriptor |
|-------------------------|---|
| 343-2714 | 2.7 mm X 14 mm, Double Lead, Locking Screw, T |
| 343-2716 | 2.7 mm X 16 mm, Double Lead, Locking Screw, T |
| 343-2718 | 2.7 mm X 18 mm, Double Lead, Locking Screw, T |
| 343-2720 | 2.7 mm X 20 mm, Double Lead, Locking Screw, T |
| 343-2722 | 2.7 mm X 22 mm, Double Lead, Locking Screw, T |
| 343-2724 | 2.7 mm X 24 mm, Double Lead, Locking Screw, T |
| 343-2726 | 2.7 mm X 26 mm, Double Lead, Locking Screw, T |
| 343-2728 | 2.7 mm X 28 mm, Double Lead, Locking Screw, T |
| 343-2730 | 2.7 mm X 30 mm, Double Lead, Locking Screw, T |
| 343-2732 | 2.7 mm X 32 mm, Double Lead, Locking Screw, T |
| 343-2734 | 2.7 mm X 34 mm, Double Lead, Locking Screw, T |
| 343-2736 | 2.7 mm X 36 mm, Double Lead, Locking Screw, T |
| 343-2738 | 2.7 mm X 38 mm, Double Lead, Locking Screw, T |
| 343-2740 | 2.7 mm X 40 mm, Double Lead, Locking Screw, T |

OsteoMed® Lower Extremities Portfolio [continued]

OsteoMed® EnCompass™ Metatarsal Resurfacing Implant

| Non-Sterile Part Number | Descriptor |
|-------------------------|---|
| 385-0008 | Encompass 8 mm Ha Coated Implant |
| 385-0010 | Encompass 10 mm Ha Coated Implant |
| 385-0012 | Encompass 12 mm Ha Coated Implant |
| 385-0014 | Encompass 14 mm Ha Coated Implant |
| 385-0015 | Encompass 15 mm Met Resurfacing Implant |
| 385-0016 | Encompass 16 mm Met Resurfacing Implant |
| 385-0017 | Encompass 17 mm Met Resurfacing Implant |
| 385-0018 | Encompass 18 mm Met Resurfacing Implant |
| 385-0019 | Encompass 19 mm Met Resurfacing Implant |
| 385-0021 | Encompass 21 mm Ha Coated Implant |

OsteoMed® Lower Extremities Portfolio [continued]

OsteoMed® ReFlexion™ 1st MTP Implant System

| Non-Sterile Part Number | Descriptor |
|----------------------------|-----------------------------|
| 370-0203-01 | Metatarsal Head, Standard |
| 370-0205-01 | Metatarsal Head, Long |
| 370-0207-01 | Metatarsal Head, Extra-Long |
| 370-0223 | Metatarsal Stem, Small |
| 370-0225 | Metatarsal Stem, Medium |
| 370-0227 | Metatarsal Stem, Large |
| 370-0423 | Phalangeal Base, Small |
| 370-0425 | Phalangeal Base, Medium |
| 370-0427 | Phalangeal Base, Large |



www.acumed.net

Acumed Oregon Campus
5885 NE Cornelius Pass Road
Hillsboro, OR 97124
+1.888.627.9957

Acumed Texas Campus
3885 Arapaho Road
Addison, TX 75001
+1.800.456.7779

Acumed Iberica Campus
C. Proción, 1
Edificio Oficor
28023 Madrid, Spain
+34.913.51.63.57

GEN70-56-W

| Effective: 2025/04 | © 2025 Acumed® LLC

Reference

1. Woods T, Delfino J, Rajan S. Assessment of magnetically induced displacement force and torque on metal alloys used in medical devices. *J Test Eval* 49. Published online July 3, 2019. <https://doi.org/10.1520/JTE20190096>.

These materials contain information about products that may or may not be available in any particular country or may be available under different trademarks in different countries. The products may be approved or cleared by governmental regulatory organizations for sale or use with different indications or restrictions in different countries. Products may not be approved for use in all countries. Nothing contained in these materials should be construed as a promotion or solicitation for any product or for the use of any product in a particular way that is not authorized under the laws and regulations of the country where the reader is located. Nothing in these materials should be construed as a representation or warranty as to the efficacy or quality of any product, nor the appropriateness of any product to treat any specific condition. Physicians may direct questions about the availability and use of the products described in these materials to their authorized Acumed distributor. Specific questions patients may have about the use of the products described in these materials or the appropriateness for their own conditions should be directed to their own physician.

Acumed®, Polarus®, Acu-Loc®, Acutrak®, Acutrak 2®, AcuTwist®, Biotrak®, Acu-Sinch®, FlexBraid®, Hub Cap®, Tap-Loc®, Frag-Loc®, and RibLoc® are registered trademarks of Acumed LLC

ExsoMed™ Corporation is a wholly owned subsidiary of Acumed LLC. ArcPhix™, InFrame™, INnate™, NanoPhix™ are registered trademarks of ExsoMed Corporation

OsteoMed® LLC is a wholly owned subsidiary of Acumed LLC. OsteoMed®, OSA™, MFx™, CFx™, ICON™, Profile Zero™, and Fast Flap™ are registered trademarks of OsteoMed® LLC